

This PDF is for our eIFU site prototype only.
Do not share or release this file!

SELECT LANGUAGE:

ENGLISH

FRENCH • FRANÇAIS

SPANISH • ESPAÑOL

PORTUGUESE • PORTUGUÊS

ITALIAN • ITALIANO

GERMAN • DEUTSCH

FINNISH • SUOMI

SWEDISH • SVENSKA

GREEK • Νέα Ελληνικά

KOREAN • 한국인

ARABIC • اَلْعَرَبِيَّةُ

This PDF is for our eIFU site prototype only.
Do not share or release this file!

This PDF is for our eIFU site prototype only.
Do not share or release this file!

BONEFOAM™

BLOCK & WEDGE

UPPER EXTREMITY POSITIONERS



Upper Extremity Blocks & Wedges are reusable freestanding positioners designed to position the patient's upper extremities for upper extremity procedures in the lateral decubitus position.




BoneFoam, Inc.
20175 County Road 50
Corcoran, MN 55340
USA

Email: customerservice@bonefoam.com
Phone: (763) 559-1830
Toll-free: 1 (877) 861-2663
Fax: (763) 559-1822


bonefoam.com

UEBW-eIFU Rev.0 2025-07-01


INSTRUCTIONS FOR USE

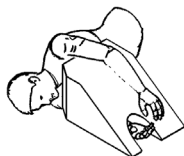


925.LUEWL

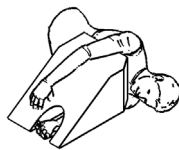

925.RUEWL

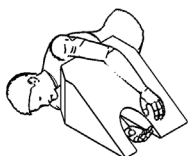

925.LUEWS

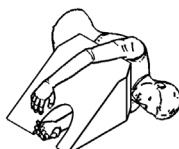


925.RUEWS



925.UUELW



925.UEBL

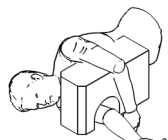

925.UEBS

CLEANING AND DISINFECTING INSTRUCTIONS

- Do not sterilize the product. This product is not intended to be sterilized.
- For coated products, wipe the product with a clean wipe and disinfectants after every use. The following disinfecting materials are known to be compatible with coated BoneFoam products:
 - Rubbing alcohol (70% isopropyl alcohol)
 - General purpose/all purpose cleaner
 - Read and follow the cleaning product's instructions.
- Make sure that the coated product is dry before you store it or use it again.

WARNINGS/PRECAUTIONS/CONTRAINDICATIONS

- Intended to contact intact skin only.
- Incorrect positioning of the patient with this product can increase the risk of:
 - Pressure injury resulting in markings or ulcers
 - Peripheral nerve injury, which can lead to temporary or permanent paralysis
 - Venous pressure, which can lead to ischemia
 - Instability, which can lead to patient falls
- Product is to be used by licensed medical professionals only.
- Not intended for patient transport.
- This product has a vinyl coating. While this coating is intended to enhance product durability and usability and is safe for skin contact, it may cause allergic reactions in individuals sensitive to vinyl or its derivatives.

- Possibility of artifacts in imaging applications depending on setup and imaging settings.
- Risk of fire. Below are flammability ratings of the materials:
 - Coating: UL 94 HBF
 - Foam: Rated either UL 94 HBF or CAL 117

CARE, HANDLING AND USEFUL LIFE

- Avoid using the product if package is damaged during transportation.
- Allow coated products to reach room temperature after storage or transportation in cold temperatures.
- Transport and store in dry conditions and out of direct light where possible.
- To prevent patient and/or user injury and/or equipment damage, examine the product for potential damage or wear prior to use. Do not use the product if damage is visible, if parts are missing, or if it does not function as expected.
- Avoid storing where sharp objects could pierce the coating of the product.
- Avoid stacking items on top of product when storing.
- Do not submerge the product in water. Equipment damage can occur.
- Prior to using this product, read the Instructions for Use.
- Discontinue use if foam has degraded or lost resilience, or if coating is torn, cracked, or loose and wrinkles, as this may lead to acute pressure points.

SAFE DISPOSAL INSTRUCTIONS

Customers should adhere to all federal, state, regional, and/or local laws and regulations as it pertains to the safe disposal of medical devices and accessories.










PACKAGING









The product is packaged in a clear polythene bag inside a corrugated cardboard box.

SAFETY INFORMATION

The techniques detailed in this manual are only manufacturer's suggestions. The final responsibility for patient care with respect to this product remains with the attending physician. Any serious incident that has occurred in relation to this product should be reported to the manufacturer and regulatory authority, where required. Keep this manual available for future reference.

SKU	DESCRIPTION	INTENDED USE	INDICATIONS
925.LUEWL	Wedge Left Upper Extremity Large	Nonsterile, reusable upper extremity positioner intended to support, stabilize, and offset the arm in a bent position in the lateral decubitus position.	Indicated for use by healthcare professionals to position the patient’s left arm for left arm procedures.
925.LUEWS	Wedge Left Upper Extremity Small		
925.RUEWL	Wedge Right Upper Extremity Large		Indicated for use by healthcare professionals to position the patient’s right arm for right arm procedures.
925.RUEWS	Wedge Right Upper Extremity Small		
925.UUELW	Wedge Universal Upper Extremity	Nonsterile, reusable upper extremity positioner intended to support and stabilize the arm in a bent position in the lateral decubitus position.	Indicated for use by healthcare professionals to position the patient’s arm for arm procedures.
925.UEBL	Block Upper Extremity Large	Nonsterile, reusable upper extremity positioner intended to support, stabilize, and offset the arm in a bent position in the lateral decubitus position.	
925.UEBS	Block Upper Extremity Small		

SYMBOL	DESCRIPTION
	Indicates the device is a medical device
	Indicates the medical device manufacturer
	Indicates the country where the medical device was manufactured and date of manufacture
	Indicates the manufacturer's lot number
	Indicates the medical device Global Trade Item Number
	Indicates the manufacturer's catalogue number
	Indicates the unique product identifier information for the device
	Indicates the need for the user to consult the instruction for use
	Indicates the device does not contain natural rubber or dry natural rubber latex

SYMBOL	DESCRIPTION
	Indicates the temperature limits to which the medical device can be stored, transported or used
	Indicates a medical device should not be used if the package has been damaged or opened
	Indicates quantity in package
	Indicates the European Authorized Representative name and address in the European Community
	Signifies European technical conformity
	Signifies Great Britain technical conformity
	Indicates a UK Responsible Person on medical device information supplied by the manufacturer
	Indicates an Authorized Representative in Switzerland on medical information supplied by the manufacturer