## SELECT LANGUAGE:

ENGLISH

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SPANISH • ESPAÑOL

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GERMAN • DEUTSCH

FINNISH • SUOMI

SWEDISH • SVENSKA

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# HEAD CRADLE

## SUPINE HEAD PROTECTION & SUPPORT

Head Cradle is a freestanding patient positioner designed to position the patient's head and neck anatomy for medical procedures in the supine position, and comes in both single-use and reusable variations.

To reorder the following Head Cradle disposable products, please contact [sales@bonefoam.com](mailto:sales@bonefoam.com)

- 10-910-02 - Head Cradle - Disposable - Box of 4

BoneFoam, Inc.  
20175 County Road 50  
Corcoran, MN 55340  
USA

Email: [customerservice@bonefoam.com](mailto:customerservice@bonefoam.com)  
Phone: (763) 559-1830  
Toll-free: 1 (877) 861-2663  
Fax: (763) 559-1822


BoneFoam B.V.  
P.O. Box 10110, NL-5060 GA Oisterwijk  
Kantstraat 19, 5076 NP Haaren  
[EAR@qnetbv.eu](mailto:EAR@qnetbv.eu)

[bonefoam.com](https://bonefoam.com)

HC-eIFU Rev.0 2025-10-01

## INSTRUCTIONS FOR USE
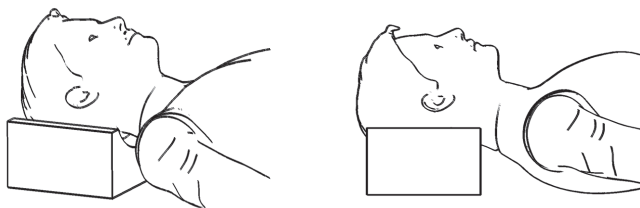
## CLEANING AND DISINFECTING INSTRUCTIONS

- Do not sterilize the product. This product is not intended to be sterilized.
- For coated products, wipe the product with a clean wipe and disinfectants after every use. The following disinfecting materials are known to be compatible with coated BoneFoam products:
  - Rubbing alcohol (70% isopropyl alcohol)
  - General purpose/all purpose cleaner
    - Read and follow the cleaning product's instructions.
- Make sure that the coated product is dry before you store it or use it again.

## WARNINGS/PRECAUTIONS/CONTRAINDICATIONS

- Intended to contact intact skin only.
- For single-use product, risk of cross-contamination if not disposed of immediately after use.
- Incorrect positioning of the patient with this product can increase the risk of:
  - Pressure injury resulting in markings or ulcers
  - Peripheral nerve injury, which can lead to temporary or permanent paralysis
  - Venous pressure, which can lead to ischemia
  - Instability, which can lead to patient falls
- Product is to be used by licensed medical professionals only.
- Not intended for patient transport.
- Ensure device is compatible with work surface, procedural instruments, and equipment.
- Possibility of artifacts in imaging applications depending on setup and imaging settings.
- Risk of fire. Below are flammability ratings of the materials:
  - Foam: Rated either UL 94 HBF or CAL 117

## CARE, HANDLING AND USEFUL LIFE

- Avoid using the product if package is damaged during transportation.
- Allow coated products to reach room temperature after storage or transportation in cold temperatures.
- Transport and store in dry conditions and out of direct light where possible.
- To prevent patient and/or user injury and/or equipment damage, examine the product for potential damage or wear prior to use. Do not use the product if damage is visible, if parts are missing, or if it does not function as expected.
- Avoid storing where sharp objects could pierce the coating of the product.
- Avoid stacking items on top of product when storing.
- Do not submerge the product in water. Equipment damage can occur.
- Prior to using this product, read the Instructions for Use.
- Discontinue use if foam has degraded or lost resilience, or if coating is torn, cracked, or loose and wrinkles, as this may lead to acute pressure points.

## SAFE DISPOSAL INSTRUCTIONS

Customers should adhere to all federal, state, regional, and/or local laws and regulations as it pertains to the safe disposal of medical devices and accessories.

## PACKAGING

The product is packaged in a clear polythene bag inside a corrugated cardboard box.

## SAFETY INFORMATION

The techniques detailed in this manual are only manufacturer's suggestions. The final responsibility for patient care with respect to this product remains with the attending physician. Any serious incident that has occurred in relation to this product should be reported to the manufacturer and regulatory authority, where required. Keep this manual available for future reference.

| SKU                           | DESCRIPTION                          | INTENDED USE                                                                                         | INDICATIONS                                                                                                                                                                                                                                     |
|-------------------------------|--------------------------------------|------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <b>REUSABLE POSITIONERS</b>   |                                      |                                                                                                      |                                                                                                                                                                                                                                                 |
| 10-910-02-C                   | Head Cradle - Coated<br>Box of 4     | Nonsterile, reusable positioner intended to support and stabilize the head in the supine position.   | Indicated for use by healthcare professionals to position the patient's head during medical procedures. Can be used as a standalone positioner or with the Troop Airway Management Elevation Complete System, Standard System, 3-Set, or 2-Set. |
| <b>SINGLE-USE POSITIONERS</b> |                                      |                                                                                                      |                                                                                                                                                                                                                                                 |
| 10-910-02*                    | Head Cradle - Disposable<br>Box of 4 | Nonsterile, single-use positioner intended to support and stabilize the head in the supine position. | Indicated for use by healthcare professionals to position the patient's head during medical procedures.                                                                                                                                         |

\*  applies only to these products.  does not apply to these products.

| SYMBOL                                                                                                   | DESCRIPTION                                                                             |
|----------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------|
|                         | Indicates the device is a medical device                                                |
|                        | Indicates the medical device manufacturer                                               |
|                       | Indicates the country where the medical device was manufactured and date of manufacture |
|                       | Indicates the manufacturer's lot number                                                 |
|                       | Indicates the medical device Global Trade Item Number                                   |
|                       | Indicates the manufacturer's catalogue number                                           |
|                       | Indicates the unique product identifier information for the device                      |
| <br>bonefoam.com/eIFU | Indicates the need for the user to consult the instruction for use                      |
|                       | Indicates the device does not contain natural rubber or dry natural rubber latex        |

| SYMBOL                                                                              | DESCRIPTION                                                                                               |
|-------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------|
|    | Indicates the device is intended for one single-use only                                                  |
|   | Indicates the temperature limits to which the medical device can be stored, transported or used           |
|  | Indicates a medical device should not be used if the package has been damaged or opened                   |
|  | Indicates quantity in package                                                                             |
|  | Indicates the European Authorized Representative name and address in the European Community               |
|  | Signifies European technical conformity                                                                   |
|  | Signifies Great Britain technical conformity                                                              |
|  | Indicates a UK Responsible Person on medical device information supplied by the manufacturer              |
|  | Indicates an Authorized Representative in Switzerland on medical information supplied by the manufacturer |

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