


Please take time to read the instructions before use.

NOTE: Gloves should be worn when using this device. For blood collection using a syringe, or for short-term IV administration, refer to manufacturer's instruction included in the box.

STEP 1

Open package and ensure the multi luer adapter is securely attached.

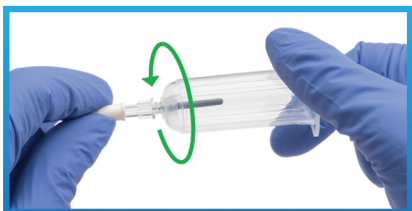

STEP 2

Remove needle sheath and attach the tubeholder




STEP 3

Prepare venipuncture site per your facility's protocol. Remove the protective needle sheath.



STEP 4

Perform venipuncture by holding the wings and per your facility's protocol. (Insert needle in the center of the tube when collecting the sample.)

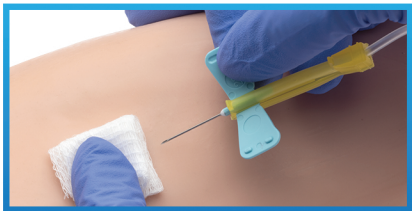

STEP 5

Withdraw blood collection set from site by grasping the translucent yellow safety shield grip area with thumb and index finger while grasping the tubing. Apply pressure to site.

STEP 6

A. One-hand Safety Activation - Holding the tubing in hand and advance the translucent yellow safety shield with thumb and index finger until an audible click is heard indicating the needle is completely retracted and locked in place.

B1. Two-hand Safety Activation - Using dominant hand, place the thumb behind the tab below the wings with index finger underneath push forward.

B2. Using the other hand, grasp the tubing and pull back until an audible click is heard indicating the needle is completely retracted and locked in place.

STEP 7

Dispose of device in the nearest Sharps container.

Blood Collection Using a Syringe

STEP 1 Open package and remove the multi luer adapter and attach the syringe. Proceed with steps 3-5. (Upon completion, use a blood transfer device to transfer blood from the syringe to the tube.)

Short Term Infusion Set (up to 2 hours)

STEP 1 Remove multi luer adapter.

STEP 2 Attach IV line or syringe to the Female Luer Adapter. Remove needle shield.

STEP 3 Prime set in accordance with facility procedure.

STEP 4 Perform venipuncture by holding wings.

Administer medication

It is recommended that this device be changed per the institution's procedures, but no longer than 2 hours.

STEP 5 See Steps 5-7 under Blood Collection



MYCO
MEDICAL
TheSharpChoice

NMSDC Certified Diversity Supplier

919 460 2535

2015 Production Drive
Apex, NC 27539 USA

www.mycomedical.com

info@mycomedical.com

Some medical devices are not licensed for sale in Canada.

© 2025 MYCO Medical, Inc. All rights reserved.

MYCO® is a registered trademark of MYCO Medical, Inc.

LIT-01037 Rev 0 Eff 12.11.2025