


PUSH BUTTON Safety Blood Collection Set

The **RELI[®] Push Button Safety Blood Collection Set** is a sterile, single-use, winged blood collection device designed for safety and efficiency. Depressing the button on the top of the device activates in-vein needle retraction for simple and effective prevention of needle-stick injury.

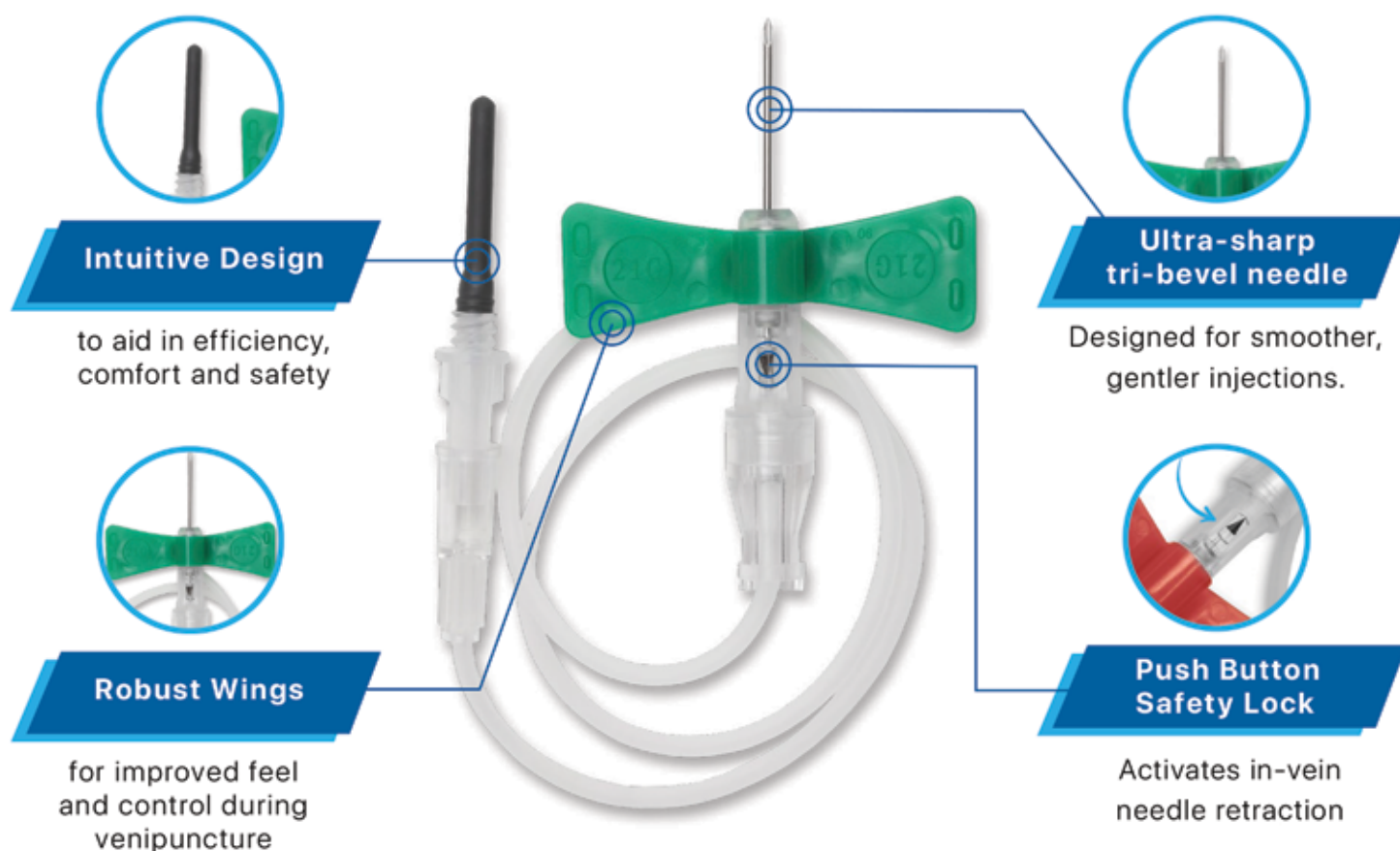

Pictured above: RELI[®] Push Button Safety Blood Collection Set 21G x 3/4" (PBSBCS21G-12)

Specifications

- Approved for short-term infusion (<2 hours)
- High flow rate design
- One-handed safety activation with audible "click"
- Not made with natural rubber latex
- User-friendly packaging designed to reduce coil memory in tubing
- Compatible with most brands of tube holders
- Visible indicator arrow for activation

Instructions for Use

NOTE: Wear gloves when using. For syringe blood collection or short-term IV use, see instructions in box.

STEP 1



Twist sheath to remove from Luer adapter and attach tube holder, transfer device, or syringe.

STEP 2



Twist and pull off needle cover and perform phlebotomy per your facility's procedure.

STEP 3



Hold the wings and insert needle into the vein. Observe for the presence of blood in the chamber.

STEP 4

Collect blood specimen(s) per your facility's procedure.

STEP 5



Cover the puncture site with a sterile gauze.

STEP 6



Press button to retract needle and check visually. Avoid blocking retraction. Apply pressure or bandage.

After Use: Dispose of all used materials in the nearest sharps container per your facility's procedure.



**Scan to watch
Instructions for Use**

Order Information

MYCO CODE	NEEDLE GAUGE	TUBE LENGTH	COLOR	UNIT OF MEASURE
PBSBCS21G-12	21G x ¾"	30cm (12")	Green	50 ea/bx; 4 bx/cs
PBSBCS23G-12	23G x ¾"	30cm (12")	Blue	50 ea/bx; 4 bx/cs
PBSBCS25G-12	25G x ¾"	30cm (12")	Orange	50 ea/bx; 4 bx/cs



MYCO
MEDICAL

TheSharpChoice

NMSDC Certified Diversity Supplier



1.919.460.2535



2015 Production Drive
Apex, NC 27539 USA



www.mycomedical.com



info@mycomedical.com



Scan to learn more

or receive complimentary
evaluation samples



<https://mycomed.co/reli-push-button>

Some medical devices are not licensed for sale in Canada.

© 2026 MYCO Medical, Inc. All rights reserved.

MYCO® is a registered trademark of MYCO Medical, Inc.

LIT-01035 Rev 3 Eff 08.19.2026