

# Sterile Glove

## Instructions for Use | Donning

### PRODUCT DESCRIPTION & INTENDED USE

Inspect the outer package for damage before use. Do not use if the package is opened or damaged in a way that could compromise the sterile integrity of the product.




#### Perform Hand Hygiene

1. Perform hand hygiene and allow hands to dry completely before handling sterile gloves.



#### Open Sterile Glove Packaging

2. Open the glove package carefully, ensuring the sterile contents remain uncontaminated.

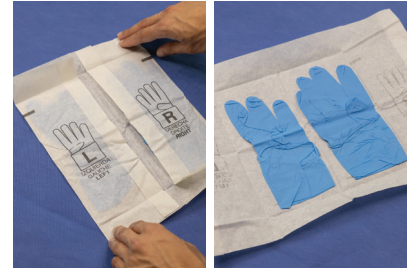

#### Expose Sterile Gloves

3. Open the inner sterile wrapper by touching only the outer edges.




#### Grasp First Glove at Folded Cuff

4. Pick up the first glove by its folded cuff, touching only the inner surface.




#### Don First Sterile Glove

5. Insert the dominant hand into the glove without contacting the outer sterile surface.

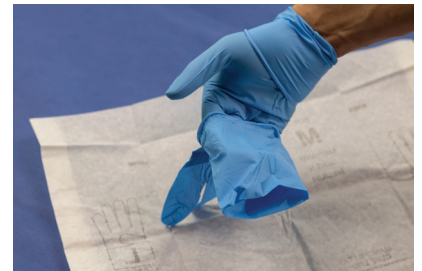

#### Insert Gloved Fingers Below Next Cuff

6. Slide the gloved fingers under the folded cuff of the second glove.

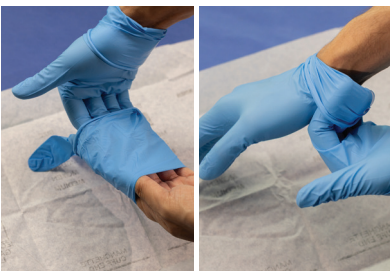

#### Don Second Sterile Glove

7. Insert the ungloved hand into the second glove while maintaining sterility.



#### Adjust Gloves Using Sterile-to-Sterile Contact

8. Adjust the gloves by touching only sterile glove surfaces together.



**PRIMED Medical Products Inc.**  
200, 2003-91 Street SW, Edmonton, AB, Canada, T6X 0W8  
TF +1.877.877.4633 | F +1.780.497.7670 | W PRIMED.com | E info@PRIMED.com  
\*Registered Trademark or Trademark of PRIMED Medical Products Inc. or its affiliates.

# Sterile Glove

## Instructions for Use | Doffing

---



### Grasp Exterior of First Glove

1. Pinch the outside of one glove near the wrist without touching exposed skin.



### Remove First Glove

2. Peel the glove away from the hand, turning it inside out as it is removed.



### Retain Removed Glove

3. Hold the removed glove in the palm of the remaining gloved hand.



### Slide Fingers Under Second Cuff

4. Insert unglved fingers beneath the cuff of the remaining glove, contacting only the clean inner surface.



### Remove Second Glove

5. Peel the second glove off from the hand, turning it inside out over the first glove.



### Encapsulate First Glove

6. Ensure the first glove is contained inside the second glove during removal.



### Discard Gloves

7. Dispose of the removed gloves immediately in the appropriate waste receptacle according to organizational policy.



### Perform Hand Hygiene

8. Perform hand hygiene immediately after glove removal using soap and water or an alcohol-based hand rub.



PRIMED Medical Products Inc.

200, 2003-91 Street SW, Edmonton, AB, Canada, T6X 0W8

TF +1.877.877.4633 | F +1.780.497.7670 | W PRIMED.com | E info@PRIMED.com

\*Registered Trademark or Trademark of PRIMED Medical Products Inc. or its affiliates.