



BEFORE & AFTER

PICOHI™ 300

Pigmentation Treatment



PICOHI™ 300

PICOHI™ 300
High Power High Effect High Convenience



Lentigine

1 time treatment (1 month later)



Before

After

532nm Zoom : 3mm · 0.17J/cm² · 1pass

PICOHI™ 300
High Power High Effect High Convenience



Lentigine

1 time treatment



532nm Zoom : 3mm · 0.3J/cm² · 1pass

PICOHI™ 300
High Power High Effect High Convenience



Solar Lentigo

1 time treatment (2months later)



Before



After

532nm Zoom : 3mm · 0.18J/cm² · 1pass

PICOHI™ 300
High Power High Effect High Convenience



ABNOM

1 time treatment




1064nm Zoom : 5mm · 1.60 ~ 1.90J/cm² · 1pass

PICOHI™ 300
High Power High Effect High Convenience




ABNOM

1 time treatment

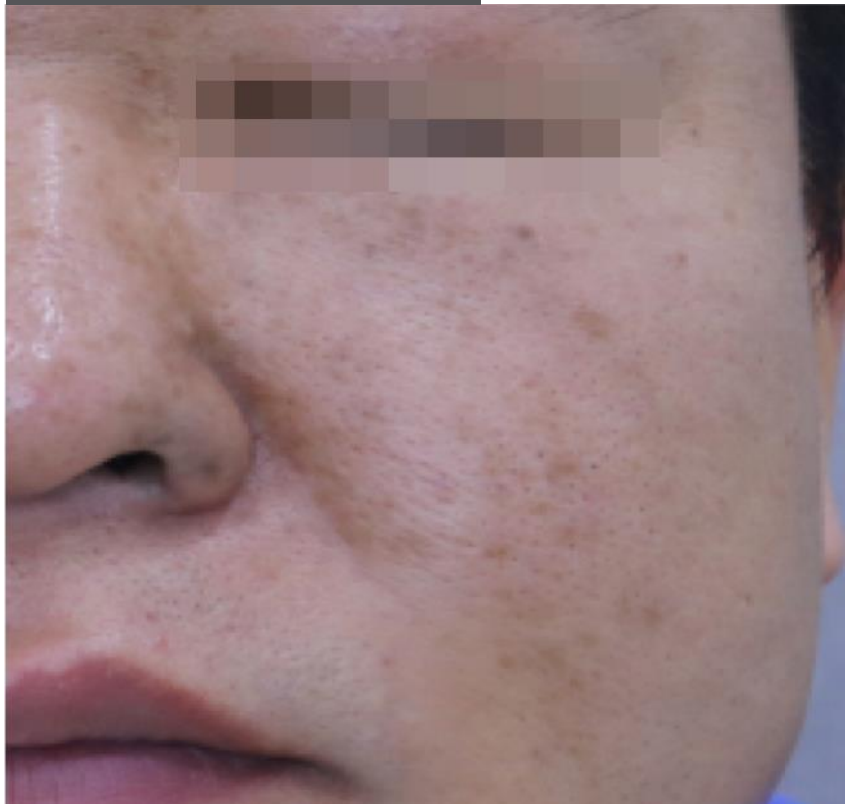

Before

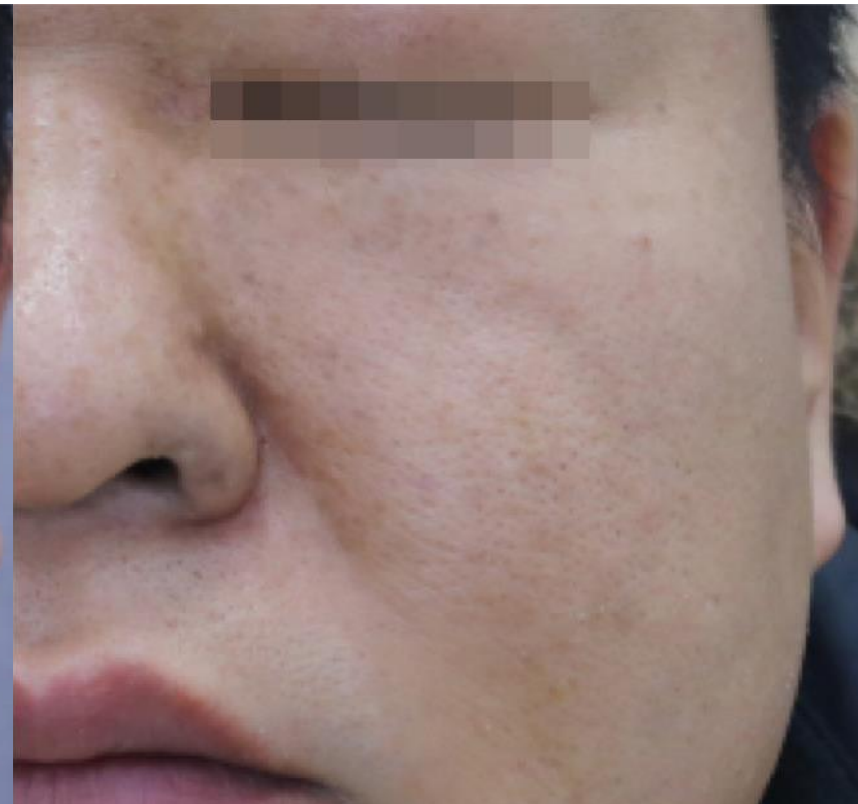

After

1064nm Zoom : 5mm · 1.60 ~ 1.90J/cm² · 1pass

PICOHI™ 300
High Power High Effect High Convenience



ABNOM + Lentigo

12 time treatment



Before



After

[ABNOM]
1064nm Collimated
0.6J/cm² · 10Hz · Stacking

[Lentigo]
1064nm Collimated
0.3J/cm² · 10Hz · 2~3passes

PICOHI™ 300
High Power High Effect High Convenience



PIH

8 time treatment



Before

After

1064nm Collimated : $0.3\text{J}/\text{cm}^2 \cdot 10\text{Hz} \cdot 1\sim 2\text{passes}$

PICOHI™ 300
High Power High Effect High Convenience



Melasma

10 time treatment



Before



After

1064nm Collimated : $0.3\sim0.4\text{J}/\text{cm}^2$ · 10Hz · 1~2passes

PICOHI™ 300
High Power High Effect High Convenience



Melasma

After 10 days



① Zoom HP

1064nm 10mm 0.3J 10Hz 1883 shots (Full Face) / 1064nm 7mm 0.7J 10Hz 275 shots (pigment lesion)
532nm **2-3mm** 0.18-0.25J 5Hz 1542 shots

② Collimation HP

1064nm 10mm 0.3J 10Hz 1377 shots (Full Face)

③ DOE 532

10x10mm 0.03J 10Hz 267 shots (pigment lesion)

Toning Treatment


PICOHI™ 300

PICOHI™ 300
High Power High Effect High Convenience




Freckles + Toning

10 times treatment

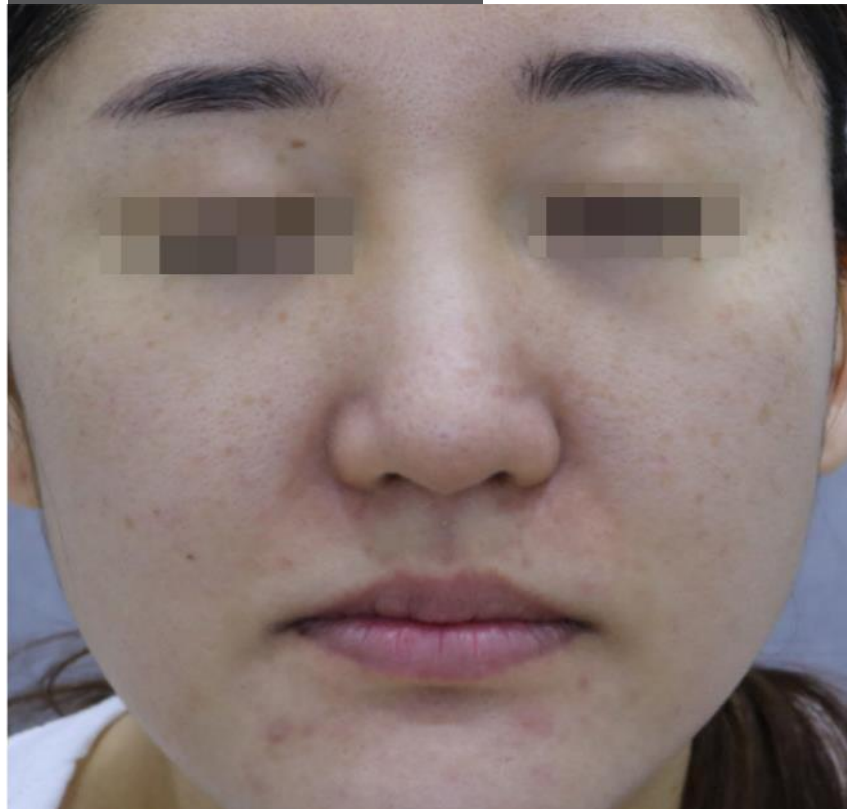

Before

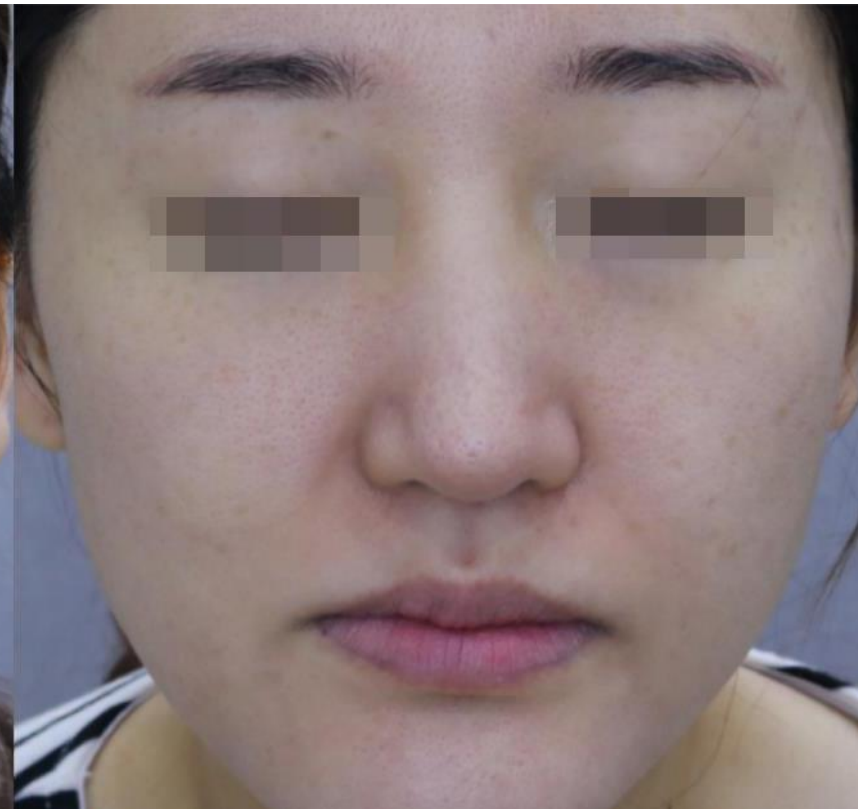

After

532nm Zoom : 3mm · 0.18J/cm² · 1pass · 1session +
1064nm Zoom : 10mm · 0.30J/cm² · 1500shots

PICOHI™ 300
High Power High Effect High Convenience



Lentigines + Toning

10 times treatment



Before



After

532nm Zoom : 3mm · 0.17J/cm² · 1pass · 1session +
1064nm Zoom : 10mm · 0.30J/cm² · 1500shots

PICOHI™ 300
High Power High Effect High Convenience



Freckles + Toning

10 times treatment



Before



After


532nm Zoom : 4mm · 0.17J/cm² · 1pass · 1session +
1064nm Zoom : 10mm · 0.30J/cm² · 1500shots

PICOHI™ 300
High Power High Effect High Convenience




Toning

10 times treatment

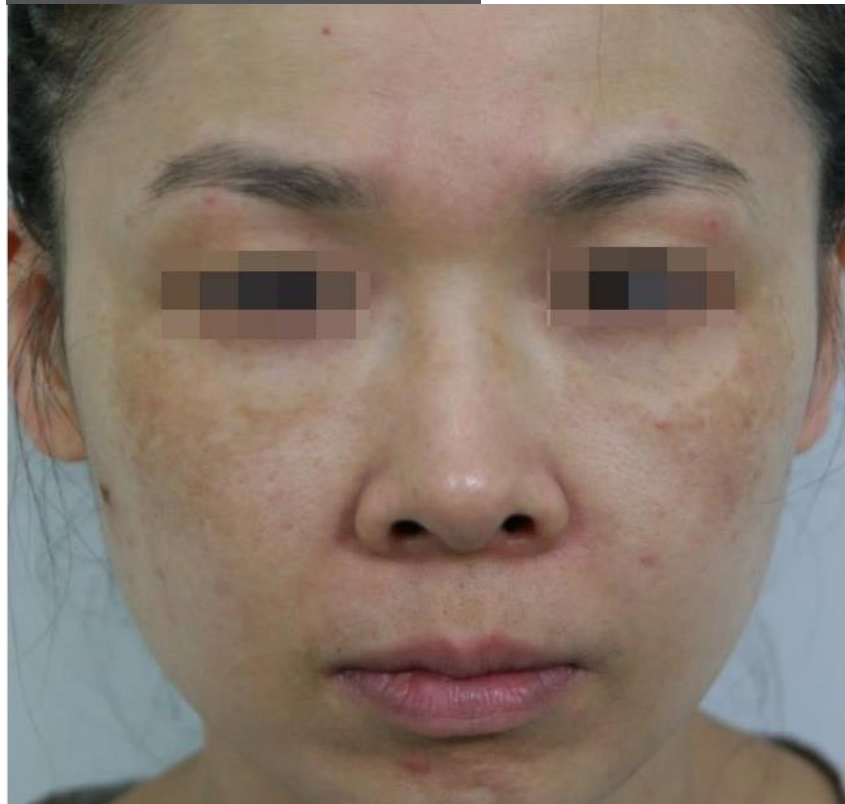

Before

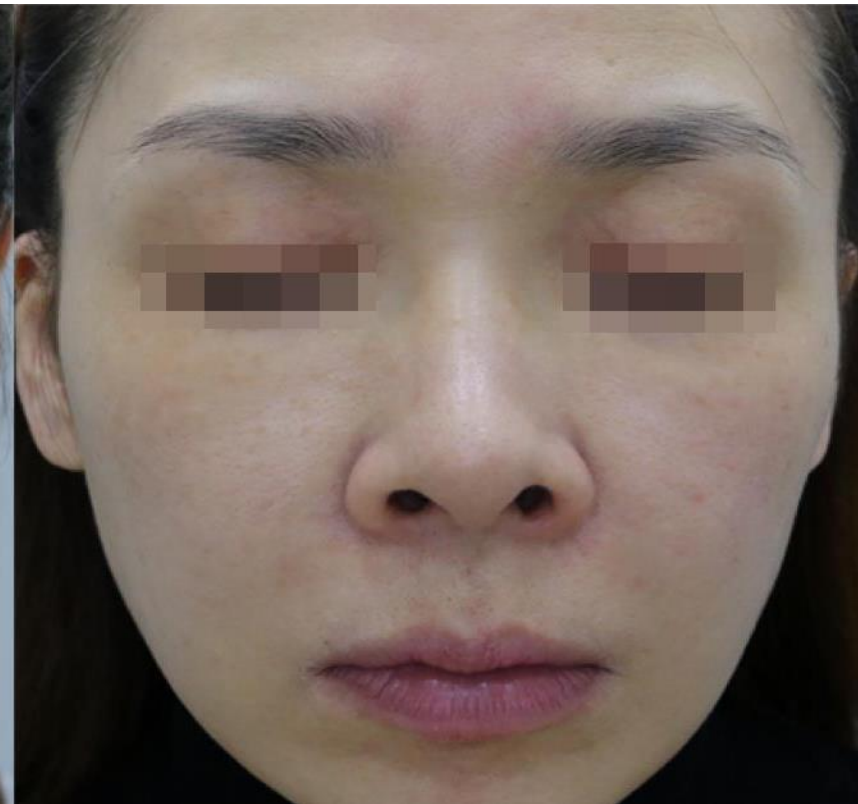

After


1064nm Zoom : 10mm · 0.20J/cm² · 1500~2000shots

PICOHI™ 300
High Power High Effect High Convenience




Toning

10 times treatment

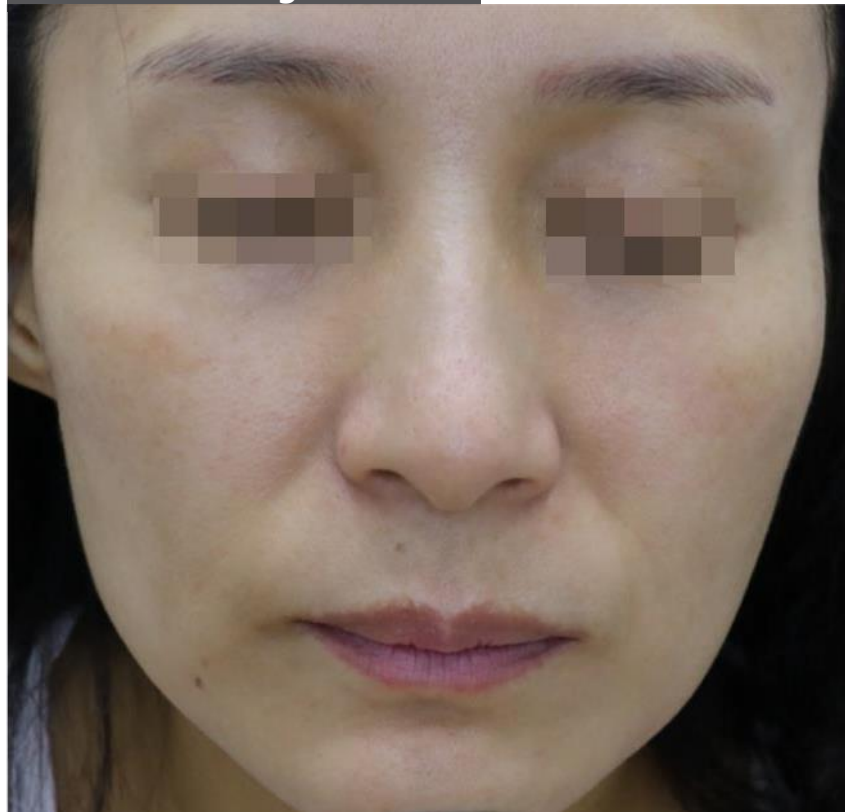

Before

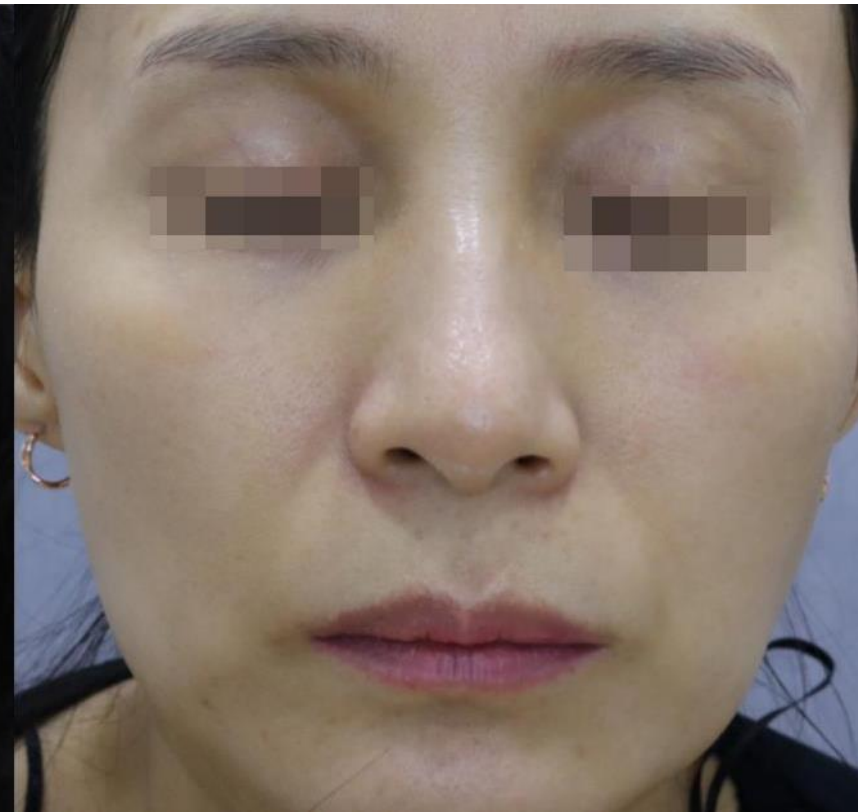

After

1064nm Zoom : 10mm · 0.20J/cm² · 1500~2000shots

Tattoo Removal




PICOHI™ 300

PICOHI™ 300
High Power High Effect High Convenience




Tattoo removal

1 time treatment

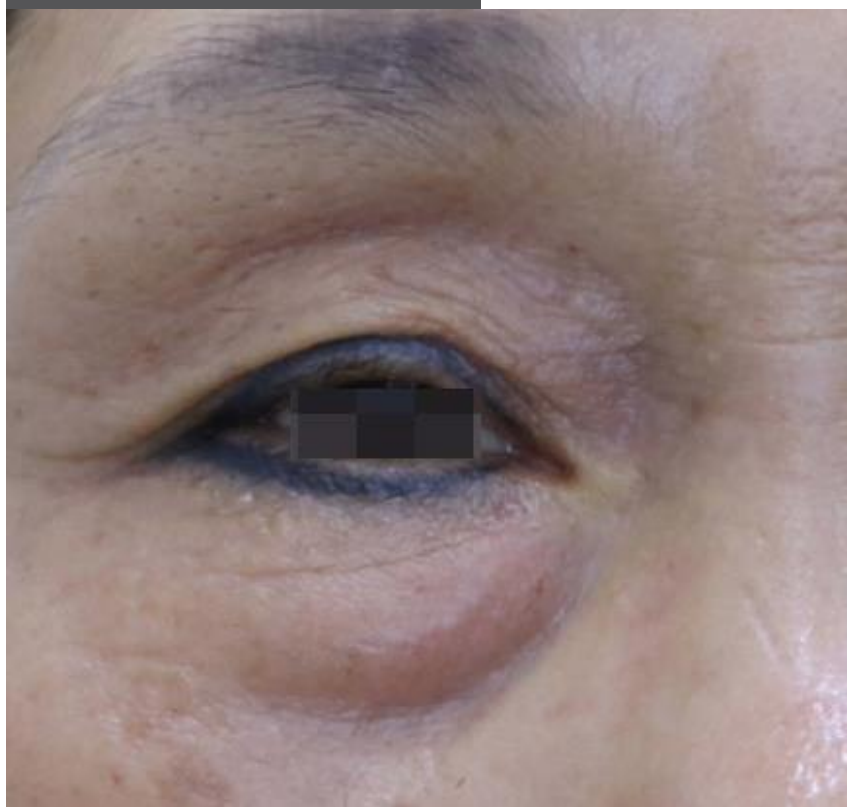

Before

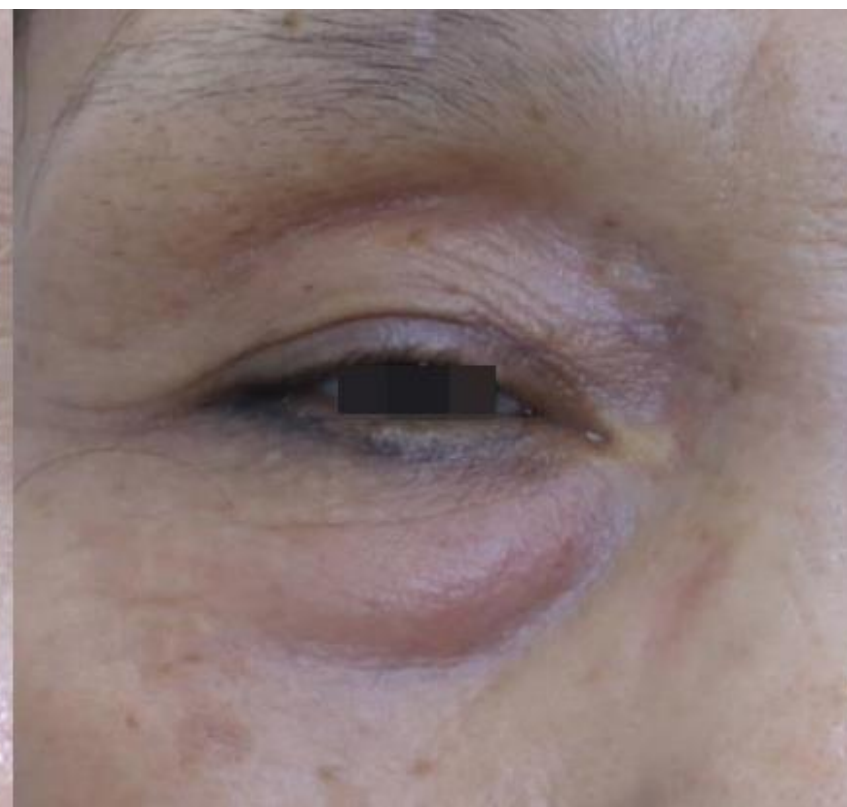

After

1064nm Zoom : 3mm · 1.90J/cm² · 1pass

PICOHI™ 300
High Power High Effect High Convenience



Tattoo removal

1 time treatment



Before



After

1064nm Zoom : 4mm · 1.20J/cm² · 1pass

PICOHI™ 300
High Power High Effect High Convenience



Tattoo removal

3 times treatment



Before



After

1064nm Zoom : 4mm · 1.20J/cm² · 1pass

PICOHI™ 300
High Power High Effect High Convenience



Tattoo removal

3 times treatment



Before



After

1064nm Zoom: 4mm · 1.20J/cm² · 1pass
532nm Zoom: 4mm · 0.12~0.20J/cm² · 1pass

Rejuvenation Treatment

PICOHI™ 300

PICOHI™ 300
High Power High Effect High Convenience



Scratched Scar

1 time treatment



ZMLA 1064nm : 5mm · 1.30J/cm² · Stacking

PICOHI™ 300
High Power High Effect High Convenience




Traumatic tattoo

3 times treatment



Before

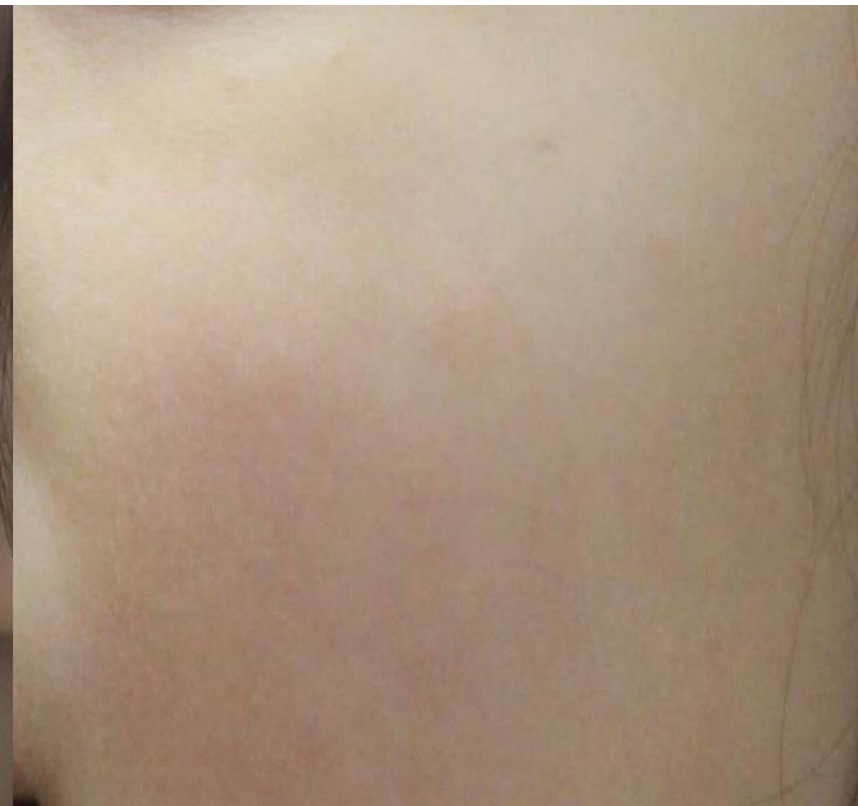

After

1064nm Zoom : 4mm · 1.10J/cm² · 1pass

PICOHI™ 300
High Power High Effect High Convenience



Toning + Rejuvenation

10 times treatment



Before

After

1064nm VMLA : level 3 · 0.05~0.10J/cm² · 3pass

PICOHI™ 300
High Power High Effect High Convenience



Rejuvenation

5 times treatment (1 week interval)



Before



After

1064nm DOE : 0.06J/cm² · 1500shots

PICOHI™ 300
High Power High Effect High Convenience




Rejuvenation

3 times treatment (1 week interval)



Before

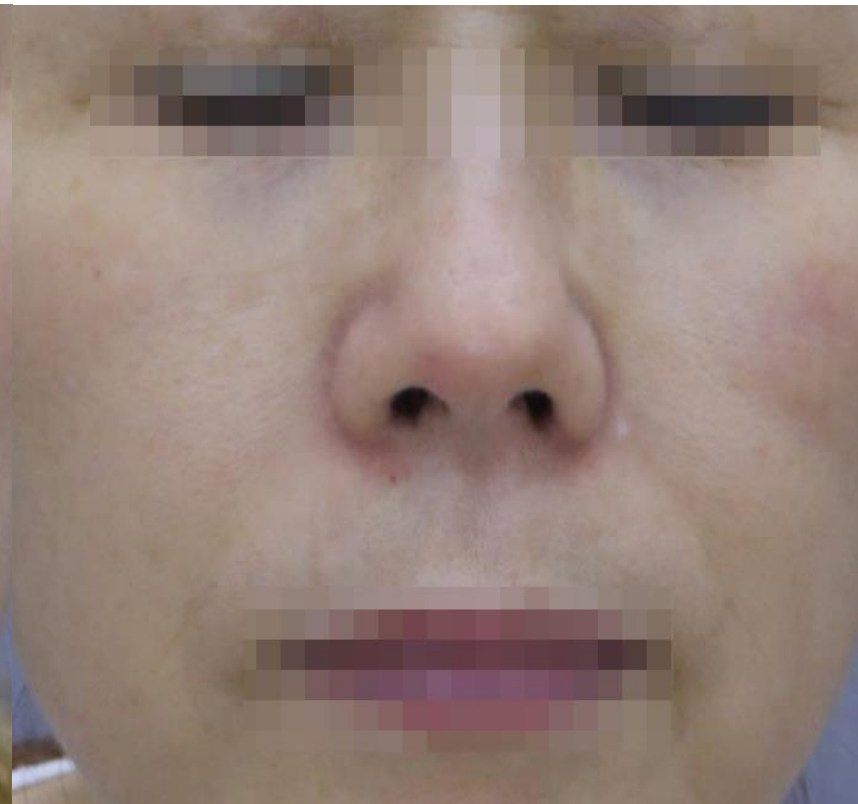

After

1064nm DOE : 0.06J/cm² · 3~4pass

PICOHI™ 300
High Power High Effect High Convenience



Large pores

5 times treatment



Before

After

1064nm VMLA : level 3 · 0.05J/cm² · 3pass

PICOHI™ 300
High Power High Effect High Convenience



Large pores & Acne scar

5 times treatment



Before



After

1064nm ZMLA : 4mm · 1.20~1.30J/cm² · Stacking

PICOHI™ 300
High Power High Effect High Convenience



Scratched Scar

4 times treatment



Before

After

1064nm ZMLA : 4mm · 1.00~1.10J/cm² · Stacking

PICOHI™ 300
High Power High Effect High Convenience



Acne Scar

15 times treatment



Before



After

1064nm VMLA : level 3 · 0.15J/cm² · 3pass
(+ 1540nm FRX + NRF)

PICOHI™ 300
High Power High Effect High Convenience



Acne Scar

15 times treatment



Before



After

1064nm VMLA : level 3 · 0.15J/cm² · 3pass
(+ 1540nm FRX + NRF)

PICOHI™ 300
High Power High Effect High Convenience



Acne Scar

5 times treatment



Before



After

1064nm VMLA : level 3 · 0.15J/cm² · 3pass
(+ 1540nm FRX)

PICOHI™ 300
High Power High Effect High Convenience



Acne Scar

14 times treatment



Before

After

1064nm VMLA : level 3 · 0.20~0.25J/cm² · 10Hz · 3~4passes