

DetechBeads Influenza B RT-LAMP LyoMix



Rapid Detection of Influenza B

A major challenge of diagnostics is time. PCR being the gold standard, often utilizes expensive equipment, labor and reagents, and reactions taking more than an hour. Detection of Influenza B virus is heavily reliant on PCR for diagnosis.

Here we present Detechgene's alternative to PCR and rapid antigen tests, offering levels of specificity and sensitivity rivalling that of PCR in the form of lyophilized LAMP beads. Loop mediated Isothermal Amplification (LAMP) is an isothermal nucleic acid amplification technique which offers specificity and sensitivity matching PCR. The results are obtained under 30 minutes and as the lyophilized bead could be stored in room temperature, it is a perfect point-of-care (POC) option for quick detection.

Specificity and Sensitivity


The DetechBeads Influenza B RT-LAMP LyoMix contains a special LAMP primer-set designed for the detection of Influenza B viruses (Victoria and Yamagata lineages) and can detect as low as 100 Influenza B copies of RNA per reaction at 65°C within 30 minutes.

Moreover, the primer-set is compatible with the RNA of the Influenza B vaccine strains recommended by the World Health Organization (WHO) and the Paul Ehrlich Institute (PEI) among other relevant strains*.

**tested in-silico*



Product Features

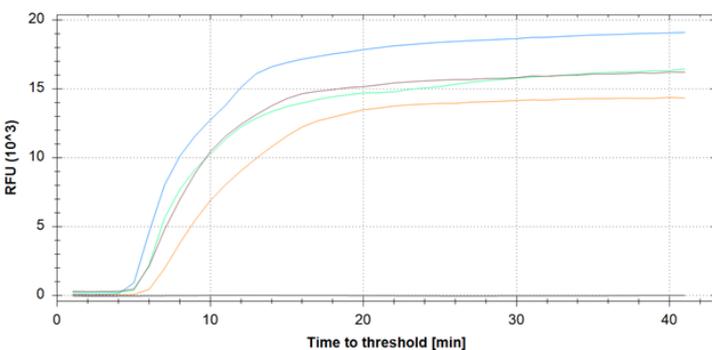

- 2500 copies/μl (Blue line)
- 125 copies/μl (Green line)
- 500 copies/μl (Brown line)
- 50 copies/μl (Orange line)

Negative control:

- Nuclease-free water: 23μl;
- Template RNA: 2μl;
- Final volume of LAMP reaction: 25μl




Sensitivity tests performed using Influenza B RT-LAMP Beads with fluorescent dye and synthetic Influenza B RNA Control.

Colorimetric change is noticed when reaction turns from pink (no-reaction/negative) to yellow (positive).

- Utilizes robust Bst DNA polymerase for LAMP
- Thermostable reverse transcriptase for RNA samples
- Sensitive and specific to detect various strains of Influenza B (Victoria and Yamagata lineages)
- Ready-to-use format with accurate results in 30 minutes
- Compatible with fluorescent dye for real-time readouts
- Suitable for use in common PCR and RT-PCR thermocyclers

Applications Overview

The use of the lyophilized LAMP bead is as simple as 1-2-3. First, 23µl of water is resuspended to the PCR tube containing the bead. Second, 2µl of RNA extracted from saliva or nasal mucus is added. Finally, the tube is incubated at 65°C for 30 minutes. The change in color, from pink to yellow, indicates the presence of Influenza B in the sample. The colors are very bright and saturated, and moreover, the beads' low foaming property makes it easy to read the results.

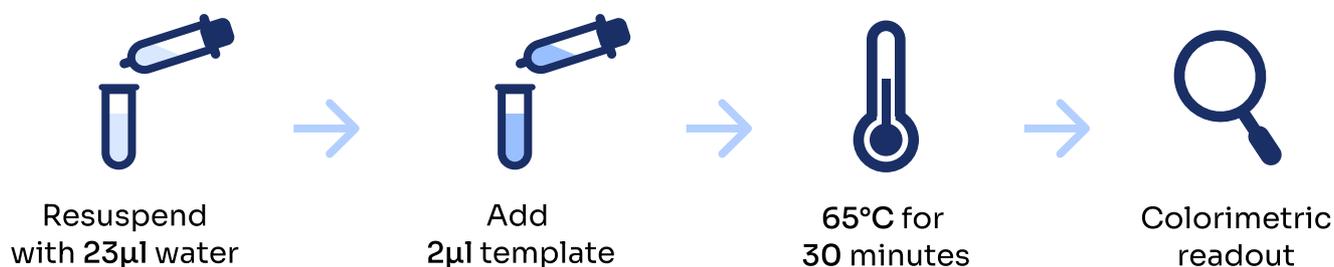

Product Specifications

The DetechBeads Influenza B RT-LAMP LyoMix are available in various formats to suit various needs. A singular Influenza B RT-LAMP bead consists of a ready-made pre-mix consisting of *Bst* polymerase, reverse transcriptase, Influenza B specific LAMP primers, dNTPs, dUTPs and the necessary salts required to carry out the reaction.

Product	Quantity	Product-ID	Price
DetechBeads Influenza B RT-LAMP LyoMix, 1 Pack	8 PCR-Strips with 8 beads each	DG-DB-FLUB-PK01	249€
DetechBeads Influenza B RT-LAMP LyoMix, 5 Pack	40 PCR-Strips with 8 beads each	DG-DB-FLUB-PK05	1229€
DetechBeads Influenza B RT-LAMP LyoMix, 10 Pack	80 PCR-Strips with 8 beads each	DG-DB-FLUB-PK010	2419€

All RT-LAMP beads are intended for Research Use Only (RUO) and not for use in diagnostic procedures

