

SINUS LATERAL APPROACH KIT



1/ MINIMALLY INVASIVE

- Minimal flap size
- Smaller window size than conventional techniques
- Instruments provide method of opening the lateral wall near the crest

2/ SPEED

- Reamers designed to be used with conventional surgical motor handpiece at 2,000rpm
- Tapered trunk of LS reamer designed to control drilling depth up to 3.5mm without use of drill stops

3/ COMPLETE


- Kit provides a complete set of tools for performing a lateral window sinus lift

Sinus Lateral Approach Kit Technique Guide




1 Flap

- Open the soft tissue.

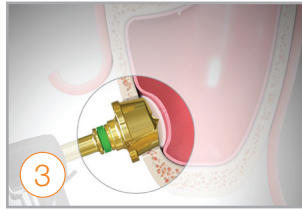


3 Bone Disk

- The LS reamer will leave a thin bone disk.

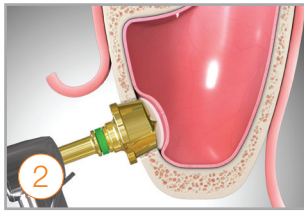


2 LS-Reamer

- Make a lateral hole with the LS-reamer.
- Keep the drill upright and maintain a steady drilling motion; Both sides of the LS-reamer blade should be in contact with the bone while drilling.
- Continue drilling.
- Speed: 2,000rpm
- NOTE: C-Reamer can be used instead of the LS-Reamer if you wish to save a circular bone core disk from the lateral wall.

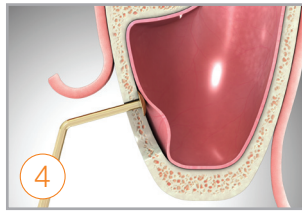


4 Detach sinus membrane

- #1 elevator is used to initiate detachment of the membrane in the mesial-distal area.

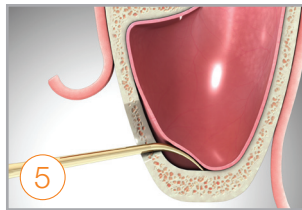


5 Elevate sinus membrane

- #2 elevator is used for detaching the membrane in the inferior area.
- #3 elevator is used for the anterior and posterior areas.

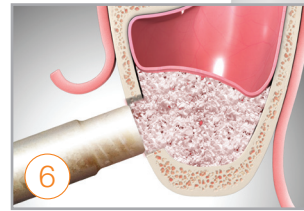

6 Bone graft

- Place bone graft material.
- Cover the lateral hole with a collagen membrane.
- Suture the soft tissue.

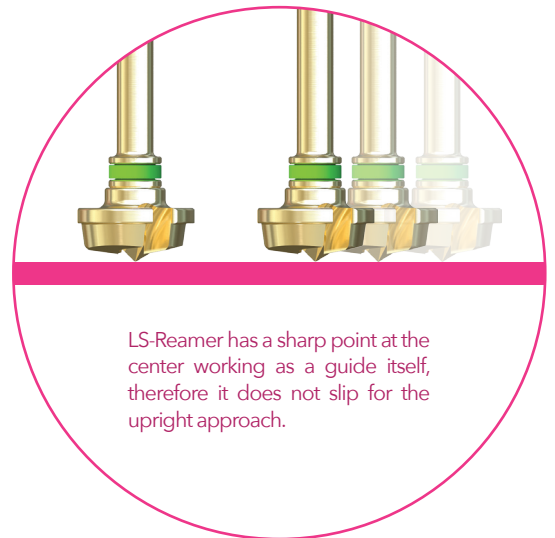

Ordering information

Catalog Number	Description
SLAKIT	Sinus Lateral Approach Instrument Kit

Note: Use copious irrigation during drilling!

Note: This guide is for educational use only.

See instructions for use for additional information on this product.

Neobiotech Co.,Ltd

70 Sinpyeongro, Jijeongmyeon, Wonjusi, Gangwondo, 26351, Republic of Korea

Tel. 82-2-582-2885, Fax. 82-2-582-2883 www.neobiotech.com / Made in Korea

**Neo
Biotech**