

New York **KERATOR**


Overdenture Attachment Systems




Your Best Choice in Overdenture



ADVANTAGE

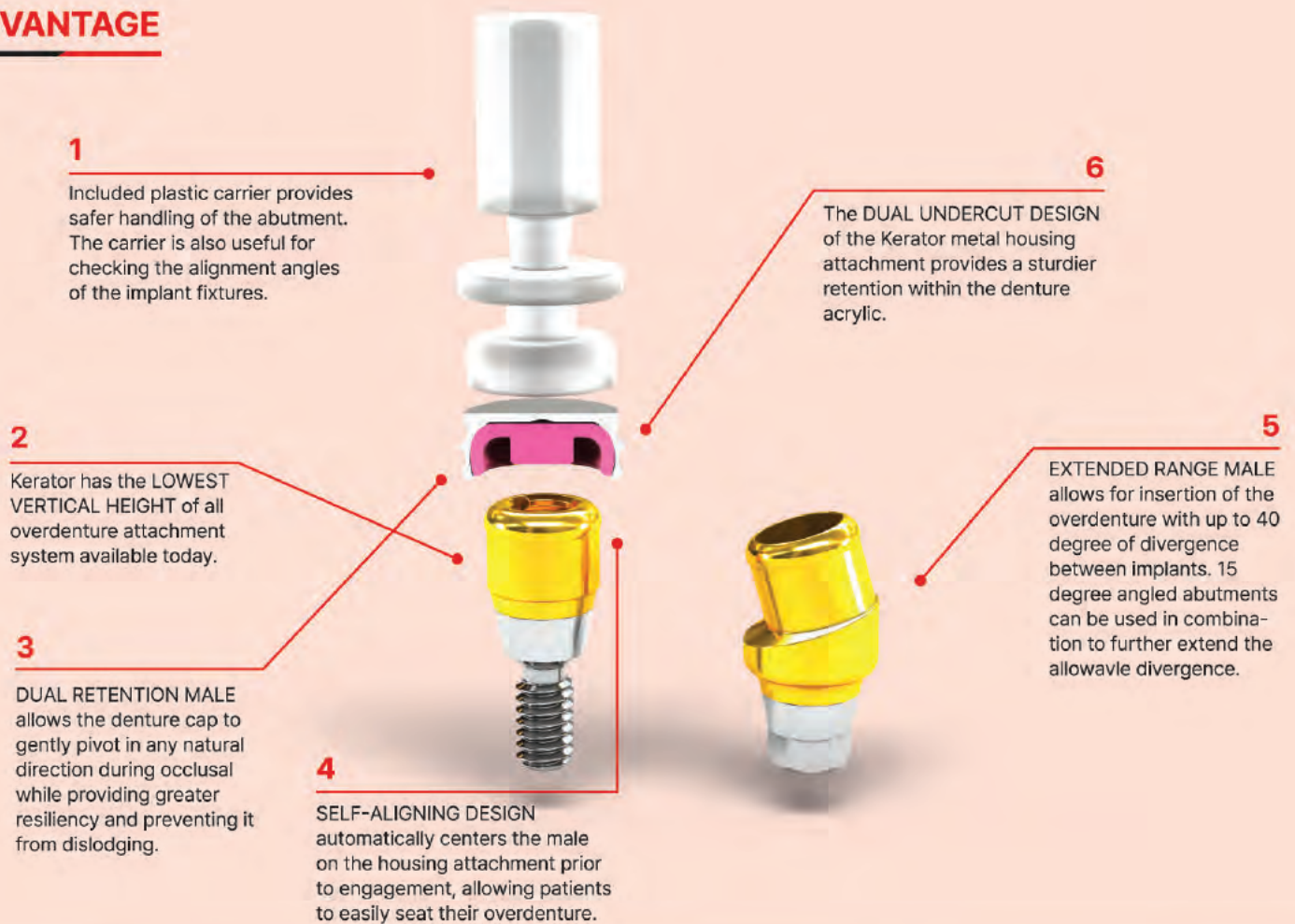

The titanium self-aligning housing that corrects divergency up to 40°

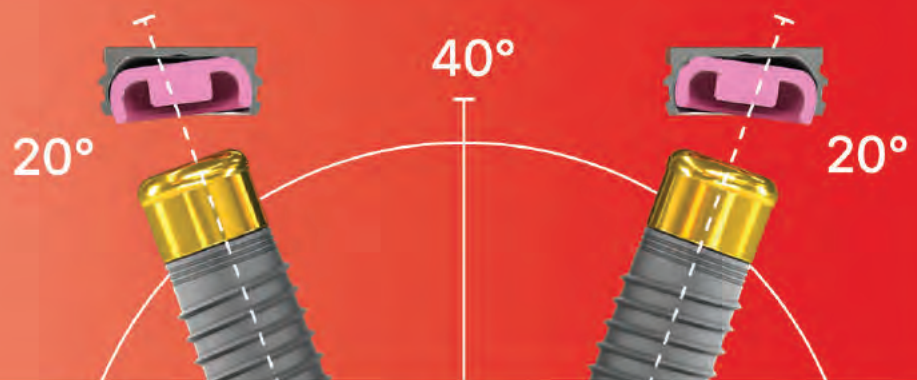

KERATOR ATTACHMENTS

For Implant



- | | | | | | | |
|-----------------|--------------|------------------|---------------|-----------|---------------|----------------|
| • Nobel Biocare | • ITI | • Bicon | • Zimmer | • Neoss | • Neo Biotech | • Megagen |
| • 3i | • Ankylos | • Camlog | • Uniti | • Thommen | • Dentium | • Cowell Medi |
| • Astra | • Biohorizon | • Frialit-2&Xive | • Euroteknika | • ICX | • Dio | • Adin, AB |
| • Alpha-Bio | • Bego | • Lifecore | • Tekka | • Osstem | • Dentis | • All Implants |

The others can be made within 10 days!

THE KERATOR COMPONENT LIST

KERATOR ABUTMENTS

IMPLANT ABUTMENTS



Kerator Implant
Abutment

BAR ABUTMENTS



Kerator bar
Thread



Kerator bar
Casting

ROOT ABUTMENTS



Root Post
Abutment
LRA1



10° Root Post
Abutment
LRA3



20° Root Post
Abutment
LRA5



Root Casting
Type
OA002

LAB PROCESSING



Impression
Coping
DK14845



Lab
Analog
DKA3854



White
Spacer Ring
WS



Denture
Spacer
KDSP



Special
Coping Carrier
SCC



Metal
Housing
(Octa)
CDPH



Metal
Housing
(Round)
CDPH(R)



Male Processing
Package
(Standard Range)
CDPH SET



Male Processing
Package
(Angled Range)
CDPH SET(A)



Black
Processing
Male
CP2

KERATOR MALE CAPS

STANDARD RANGE



0.8 lbs
CPR2



1.2 lbs
CPB2



2.4 lbs
CPP2

EXTENDED RANGE



0 lbs
CPY2



1.6 lbs
CPO2



3.2 lbs
CPG2



4 lbs
CPW2

TORQUE DRIVERS & TOOLS



Magic Tool
KMT002



Torque Wrench Tip
(Square)
KMD619



Torque Wrench Tip
(Round)
KMD719



Torque Wrench Tip
(Hexagonal)
KMD819




Housing
KMH001



Hex Driver 0.05
KMHD01


There are two types for root retained overdenture. The Post type is used in direct technique and there are straight, 10 degree, 20 degree angled females. The female part is made of Titanium alloy (Ti6AL4V) and is coated with a Titanium Nitride coating. And the Cast type has a pattern that is incorporated into wax-pattern for post coping. The Root Kit includes 4 attachments (Male&Female), 3 kinds of burs (initial drill, main pilot drill, countersink diamond bur).

The Diamond bur is used for the shoulder on the top of the root accept the KERATOR Female.

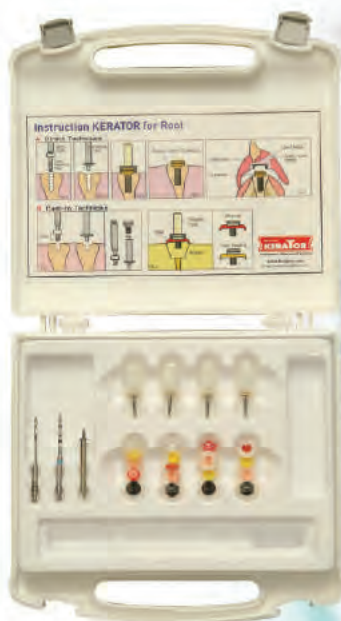


A. Direct Technique

■ NON-DIVERGENT ROOT

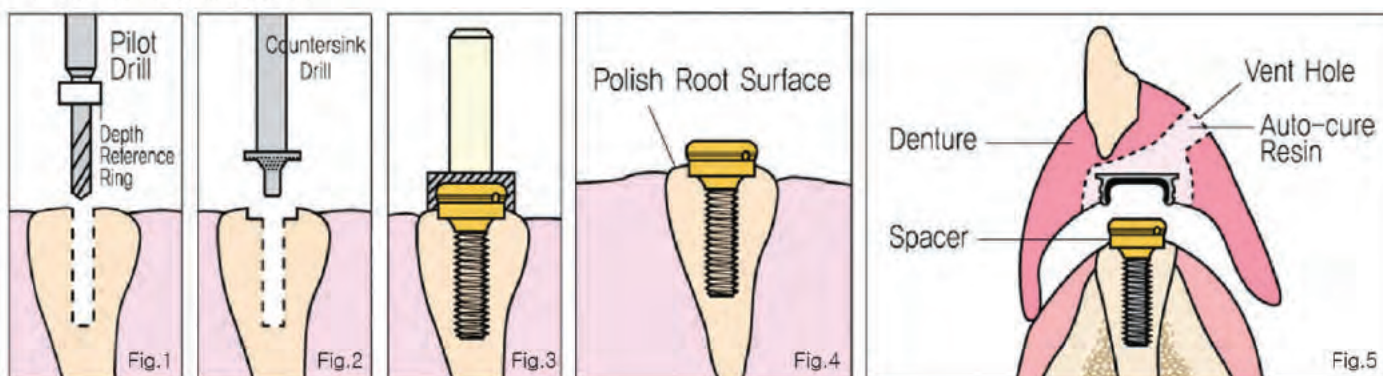


■ DIVERGENT ROOT

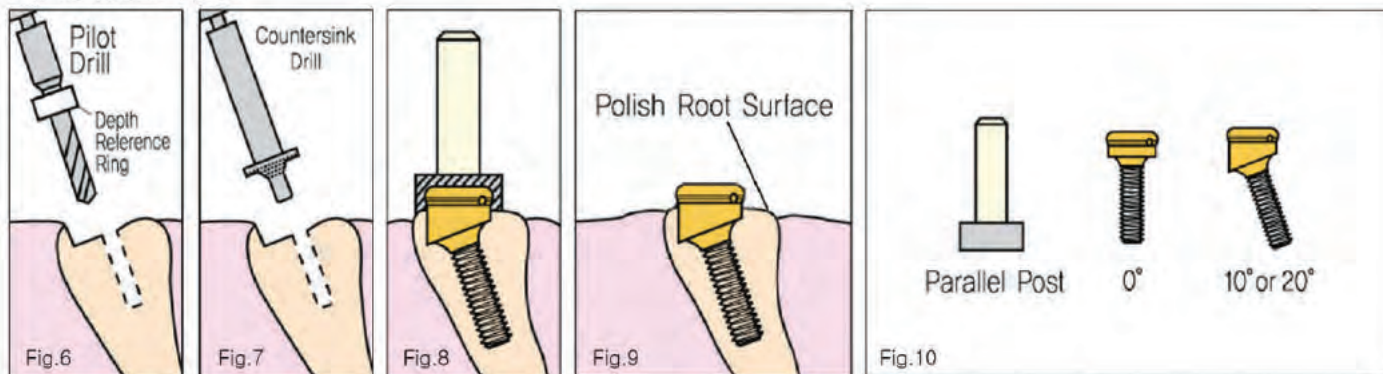

B. Cast-to Technique

■ Transfer Impression

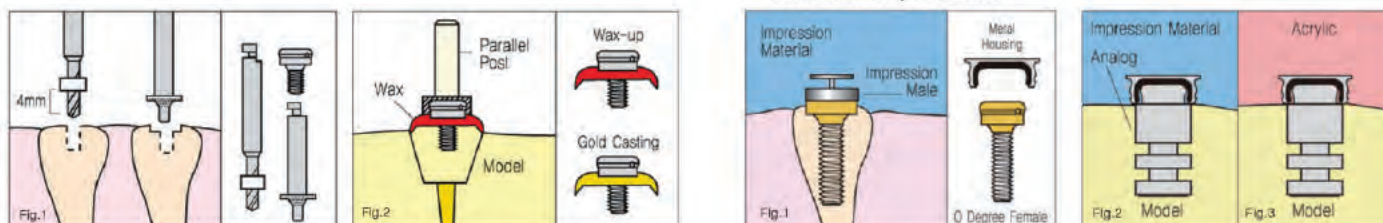

FABRICATION INSTRUCTIONS



1 Remove the healing abutment using hex-driver



2 Place the kerator abutment on to the implant body using the carrier



3 Tighten each Kerator abutment to 30Ncm using Kerator driver tip



4 Place a white block-out spacer around each abutment Snap a metal housing with a black cap onto each abutment



5 Place a white block-out spacer around each abutment Snap a metal housing with a black cap onto each abutment




6 Prepare a recess in the denture Make a vent hole on the lingual aspect to excess the resin



7 Place a small amount of resin into the recess in the denture




8 Insert the denture in position

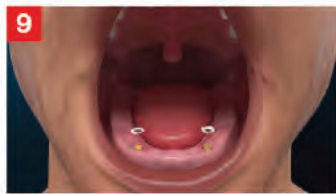

9 Remove the white spacer ring



10 Remove the black cap from the metal housing using Kerator Magic Tool

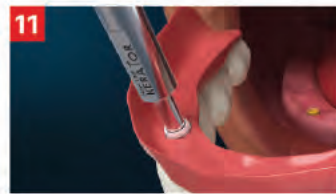

11 Remove the black cap from the metal housing in the denture and insert a final cap using Kerator Magic Tool



12 Seat the overdenture on the Kerator abutment

GLOBAL EXPORTING AREA



Major market

Europe

Germany, England, Portugal, Netherlands, Italy, Israel, Greece, Czech Republic and most of Europe.

North America

The United States, Canada and other North American regions.

Asia

China, Japan, Cambodia, Thailand, Singapore, and other Asian regions.

Middle East

Iran, Saudi Arabia, Lebanon, and other Middle East regions.

ETC

Egypt, South Africa, etc.

**HOW SOMETHING
SMALL, CAN
MAKE SUCH A
BIG DIFFERENCE!**



**Kerator has the
world's most
powerful com-
petitiveness
in overdenture
attachment
system**

Why?

KERATOR was
chosen over the
other overdenture
attachment, for the
comfortable retention
power of caps and
excellent mobility in
the denture metal
housing, in the world
market.



DAEKWANG I.D.M

Address : #802 Woori
Building, 42, Chilpae-ro,
Jung-gu, Seoul, Republic
of Korea

Tel : 82-2-3783-7041~3
Fax : 82-2-774-7350
overdenture@kerator.com
www.kerator.com

**NewYork
KERATOR**
Overdenture Attachment Systems