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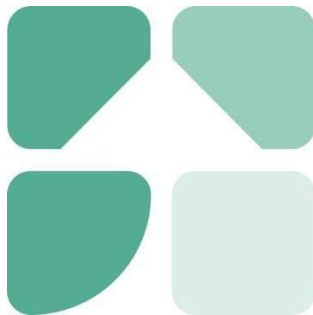
The screenshot displays the 'Exercises' section of the Healthcare Platform. It includes a 'Create exercise' button, a list of exercises (Glutes contraction, Icepacks Hip, Side-heel taps, Quadriceps contraction and extension of the knee), and a 'Daily message' section. Below this is a 'Manage patient scheme' table with columns for 'Active exercises', 'Day', and various exercise types (Icepacks Hip, Circulation, Heel raises, Glutes contraction). The table shows a grid of exercises with blue checkmarks indicating active status. A 'Plan exercise' dialog box is open, showing settings for 'Start day / End day' (5 / 5), 'Repeat' (Daily), 'Intensity start' (10), 'Intensity target' (10), 'Frequency start' (3), and 'Frequency target' (3). The dialog also includes 'Cancel', 'Apply', and 'Apply & save' buttons.

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Exercises Step goal **Recurring questionnaires** Medical questionnaires Video analysis Show in Patient language: No Advanced + Stop active follow-up Complete

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Name	Start day	End day	Frequency (Week)	Active from	Reminders	Order	Preview	Remove
Discharge	9	10	1	10:00	Reminders not set		Preview	Remove
Postoperative recovery	0	3	1	10:00	Reminders not set		Preview	Remove
Activities	15	100	1	10:00	Reminders not set		Preview	Remove
Exercises	5	100	1	10:00	Reminders not set		Preview	Remove

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Videoanalyse:

Exercises Step goal **Recurring questionnaires** **Medical questionnaires** **Video analysis** Show in Patient language: No Advanced + Stop active follow-up Complete

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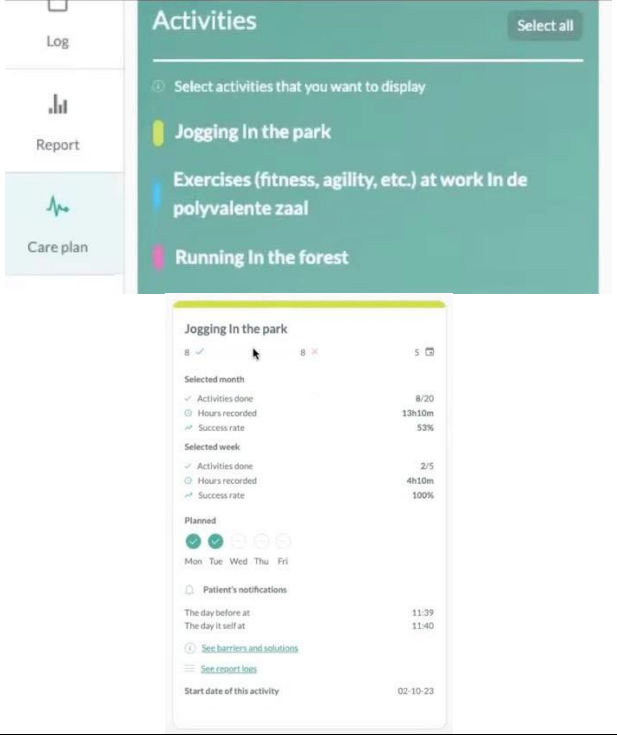
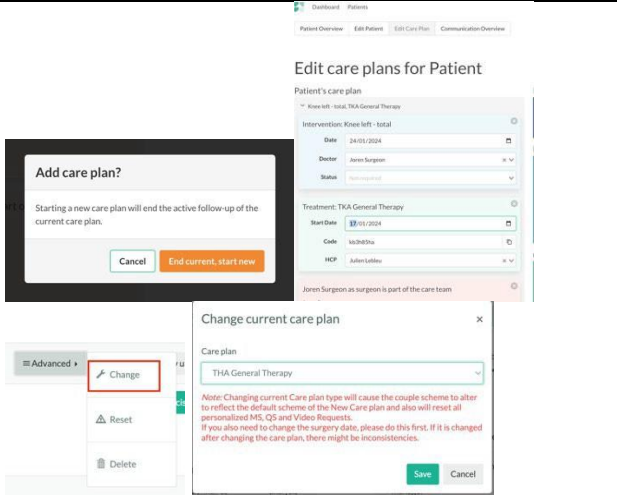
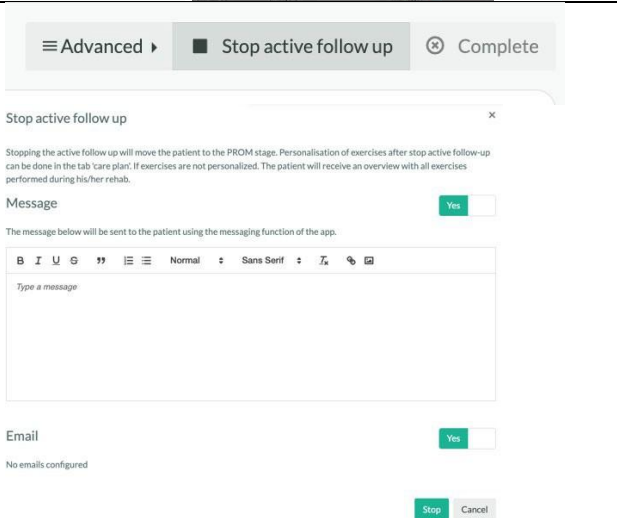
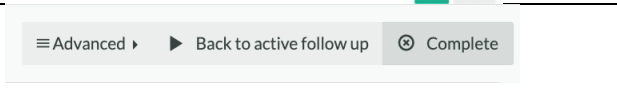
The screenshot shows the 'Edit Laurens Baes' page in the moveUP system. The page has a top navigation bar with 'moveUP', 'Patients', and 'Onboarding'. On the right, there are links for 'Config', 'Admin', and a user profile for 'Laurens'. The main content area is titled 'Edit Laurens Baes' and contains a form with the following fields:

- First name: Laurens
- Surname: Baes
- Username: laurens@moveup.care
- Email: laurens@moveup.care
- Telephone number: +32496123456
- Social security number: Social security number
- Address: Oudeplein 1
- City: Gent
- Postal code: 9000
- Country: Belgium
- Title: Title
- Image: Image
- Language: English

At the bottom of the form, there are two toggle switches:

- Authentication:** Two factor authentication (turned on)
- Wellness module:** Enable wellness module (turned on, highlighted with a green circle)

Below the toggles, there are three buttons: 'Save', 'Set new password', and 'Manage consent/contracts'.

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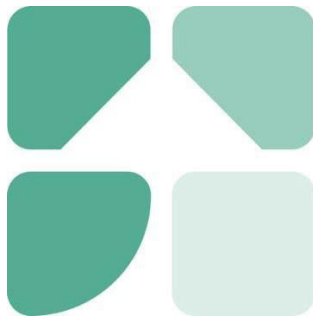
User manual for the Health Care Provider (HCP)

A personalized treatment, with more insight and involvement in the care pathway.

App version : 2.0.0 and later

Manual version : 1.5

Data of release manual : 13/11/2025

Link to e-learning: <https://moveup.talentlms.com/index>

Table of Contents

Table of Contents.....	2
PRODUCT DESCRIPTION.....	3
<i>Intended use.....</i>	<i>3</i>
<i>Intended users.....</i>	<i>3</i>
<i>Target population / indication for use</i>	<i>4</i>
<i>Claims</i>	<i>4</i>
<i>moveUP app & devices</i>	<i>4</i>
<i>Warnings and precautions</i>	<i>5</i>
<i>Contra-indications and potential side effects</i>	<i>5</i>
<i>Notification of serious incidents</i>	<i>5</i>
<i>Human body contact</i>	<i>6</i>
<i>Accessories/products used in combination</i>	<i>6</i>
<i>Device lifetime</i>	<i>6</i>
<i>Company details</i>	<i>7</i>
PATIENT APP.....	8
MEDICAL MODULE	8
HEALTHCARE PLATFORM ACCOUNT CREATION.....	8
ONBOARDING OF PATIENTS	8
HEALTHCARE PLATFORM GENERAL FUNCTIONS	11
HEALTHCARE PLATFORM CARE PLAN FUNCTIONS	18
WELLNESS MODULE.....	21
WELLNESS MODULE ACTIVATION	21
WELLNESS MODULE OVERVIEW	22
Wellness module overview	23
Wellness report	23

PRODUCT DESCRIPTION

Intended use

For healthcare providers, moveUP provides efficient clinical management of pathologies and treatments. With the use of moveUP, valuable insights are given to the HCP/care team about the status of the patient condition and its progression. The HCP/care team can provide targeted therapy advice & adapt treatment (exercise) schemes based on the needs of the patient through the use of the moveUP app.

For patients, moveUP gives insights about their progress, provides instructions related to their disease or treatment and helps them manage their symptoms. Based on the patients' needs and timepoint in the patient pathway, patients can be followed-up intensively and receive active therapy advice through the moveUP app.

moveUP is intended for the clinical follow-up on adult patients who underwent (postoperative) or will undergo (preoperative) a Total Knee Arthroplasty (TKA), Total Hip Arthroplasty (THA), Chronic obstructive pulmonary disease (COPD) or Bariatric surgery.

Intended users

moveUP is intended to be used by patients and healthcare providers.

Patient

The main user of the mobile app and patient website is the patient. Patient users should be alert, mentally competent and capable of operating a tablet or smartphone.

Inclusion:

- Age:
 - minimum: 18 years
 - maximum: no limit
- Weight: NA
- Health & condition: capable of performing basic activities of daily living
- Language: understanding one of the available languages of the app (Dutch, French, German, English)

Exclusion:

- Patients who are mentally incompetent or having trouble expressing what they are feeling (for instance, mentally diseased people, people staying in elderly care centres.) are excluded.
- Patients who are not capable of operating a tablet or smartphone.
- Patients who can't understand one of the available languages of the app (Dutch, French, German, English)

Healthcare provider

The main user of the medical web interface is a healthcare practitioner (group) or clinical researcher (group), named the care team. The care team is able to operate a web interface via web browser on PC/tablet/smartphone. The healthcare practitioner needs to understand one of the available languages of the web interface (today only available in English).

Only qualified HCP with a diploma or certificate in the healthcare domain or qualified clinical researchers can access the platform and complete the moveUP e-learning trainings.

Target population / indication for use

moveUP is used by patients:

- Who underwent or will undergo a Knee (Resurfacing) arthroplasty
- Who underwent or will undergo a Hip (Resurfacing) arthroplasty
- Who underwent or will undergo a gastric sleeve or gastric bypass surgery
- Patients with COPD.

Claims

- Patients using moveUP are highly satisfied.
- Patients using moveUP show a high adherence/compliance level to their therapy.
- The use of moveUP is more cost-effective than standard of care for hip & knee arthroplasty
- The length of stay for Knee & hip arthroplasty can be reduced with the use of moveUP
- Patients can rehabilitate without the need to leave the home environment & without receiving F2F physical therapy session for TKA and THA
- Using moveUP is as effective as standard of care
- Using moveUP is as safe as standard of care

moveUP app & devices

The moveUP app runs on:

- iPhones/iPads (iOS 11.0 or newer)
- Android systems (Android 7.0.0 and newer).

If the patient has no compatible tablet/smartphone, moveUP can lease these devices to the patient.

- Lenovo Tab 3 - Android 6.0.1 or Higher
- Huawei MediaPad T3 - Android 7.0.0 or higher

The moveUP application needs an internet connection when in use.

The moveUP application is compatible with the following browsers:

- Chrome
- Safari
- Edge

- Firefox

We strongly recommend to use the latest version of the browser.

For your data safety please enable the 2-factor authentication(2FA) and use a safe network

Warnings and precautions

	• moveUP is not an emergency tool. In case of emergency please contact your doctor or call 112.
	• moveUP is not intended to monitor/treat vital parameters of critical diseases.
	• Patients need to be aware if questionnaires are not filled out on frequent basis or not filled out in a trustworthy, the care team has more difficulty to know their health situation & evolution.
	• moveUP is not a replacement of any treatment you need to follow, but is used as an addition.
	• The moveUP App offers personalized healthcare based on the patient's profile and therefore can not be shared, borrowed or traded between users. The information and advice of the moveUP application do not apply to a patient other than the one whose profile is stored in the application. The patient may not grant access to the moveUP App to third parties. Any other person can, even accidentally, send wrong information or change the profile of the patient, with the result that the moveUP-guided advice no longer corresponds to the needs of the patient.
	• The platform includes modules displaying wellness data. These data are not part of the medical device and are not intended to be used for clinical decision making

* This data originates from a non-medical device product & may not be validated as such and the data shown in the wellness module cannot be used for medical purposes. The wellness module information is not part of the CE certificate of moveUP.

Contra-indications and potential side effects

There are no contra-indications or known side effects

Notification of serious incidents

A serious incident means any incident that directly or indirectly led, might have led or might lead to any of the following:

- the death of a patient, user or another person,
- the temporary or permanent serious deterioration of a patient's, user's or other person's state of health that resulted in any of the following:
 - life-threatening illness or injury,
 - permanent impairment of a body structure or a body function, hospitalisation
 - prolongation of patient hospitalisation, medical or surgical intervention to prevent life-threatening illness or injury
 - permanent impairment to a body structure or a body function,
 - chronic disease,
- a serious public health threat;

If you detect any of these cases please send an email to info@moveUP.care and your competent authority:

- Netherlands: meldpunt@igj.nl
- Belgium: meddev@fagg.be
- France: dedim.ugsv@ansm.sante.fr
- Germany: zlg@zlg.nrw.de

Human body contact

No human body contact with the patient or user, due to the nature of the product (software).

Accessories/products used in combination

There are no accessories.

Device lifetime

2-years, depending on the willingness of the user to update the app. We sent out a notification with the recommendation to update the app, to make sure the app will keep the performance & recently included features.

Company details

 moveUP NV/SA	Tel: +32 800 88 008	UDI: ++G166MOVEUPV200AU  (01)G166MOVEUPV200AU (10)+\$+200MDR (11) /16D20251031 (21) HI
Head office Brussel: Cantersteen 47, B-1000 Brussels, Belgium Secondary office: Kliniekstraat 27A 9050; Ghent, Belgium	Email: info@moveUP.care www.moveup.care	

If you have any questions about the App or the tablet, you can contact moveUP via the following channels - in order of priority: Via the message system of the App. Via email to support@moveup.care - Via phone during office hours: +32 800 88 008

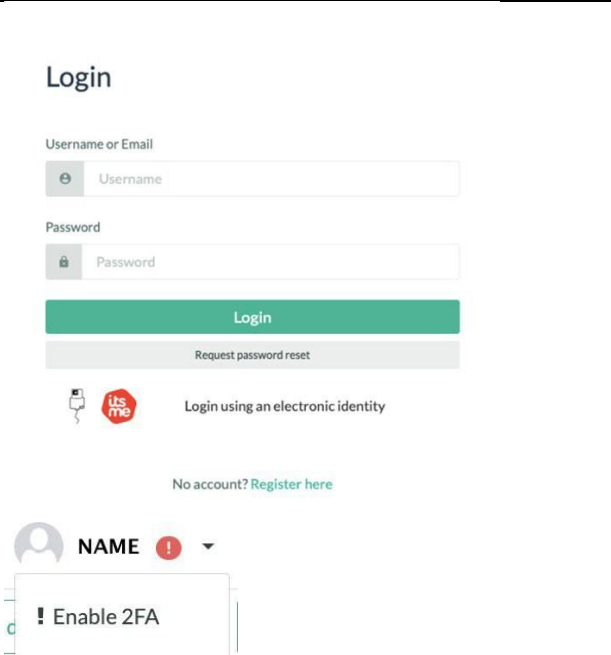
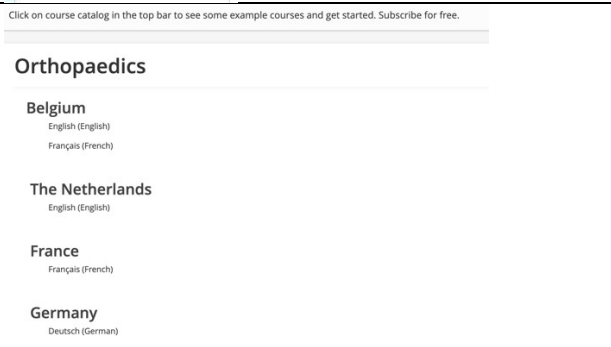
MoveUP guarantees delivery of the eIFU within seven (7) days upon request. Please contact info@moveUP.care if you want to obtain a free copy.

PATIENT APP

Please read the user manual for patients on www.moveup.care/appinfo to view the functionalities of the patient app.

MEDICAL MODULE

HEALTHCARE PLATFORM ACCOUNT CREATION

Step 1: The HCP registers on the healthcare platform on https://v1.medical.moveup.care/	
Step 2: moveUP will contact the HCP to add his/her profile to the correct organization.	
Step 3: Security measures: It is important to enable 2FA as you have access to sensitive health data.	
E-learning: For HCPs who will provide care through the use of the healthcare platform it is recommended to follow the e-learning courses for the indications listed on https://moveup.talentlms.com . We ask the HCPs to register for the e-learning course in their applicable language.	

ONBOARDING OF PATIENTS

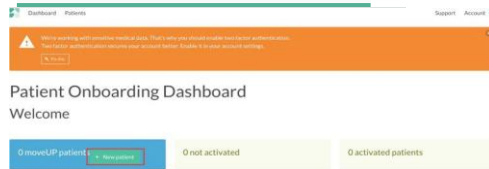
Onboarding is the process of registering a new patient in the healthcare platform.

A crucial part is adding a care plan to the patient's profile. A care plan determines what will be visible for the patient in the application, e.g.: type of exercises, information modules, message functionality, etc. A care plan usually consists out of an intervention, a treatment and (a) care

team(s).

Onboarding patients

Patients Onboarding

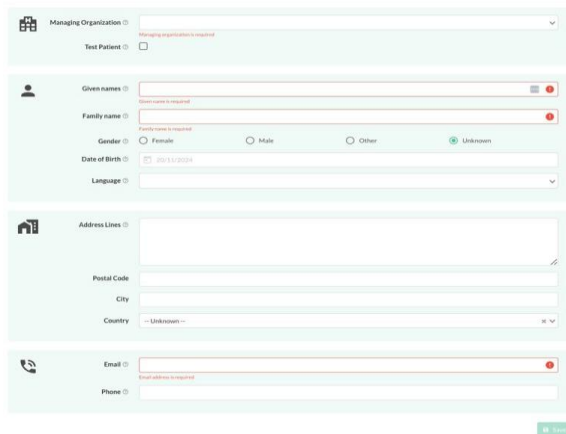
The onboarding dashboard is easy retrievable in the health care platform in the upper left corner when clicking 'Onboarding'. It is also accessible through <https://patient-onboarding-admin-moveup.care> using the HCP credentials of the health care platform.

Click on 'new patient' to onboard a new patient and fill out the patient information.

After the patient information is filled out, press 'Save'.

New patient

Important note: Only one account can be coupled to an e-mail address. It is not possible to reuse an e-mail address that is already known within our database, e.g.: reusing the same e-mail address for the account of a partner or parent.

Patient's care plans

Next up, the care plan needs to be selected: press 'Start new care plan'. Drag and drop an intervention and treatment to the care plan and fill out the dates, name of the treating doctor and care team.

Drag and drop a study if applicable. Click Save.

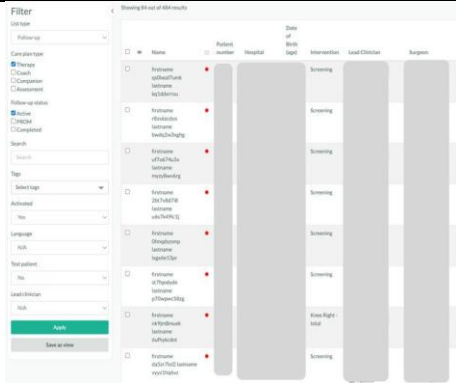
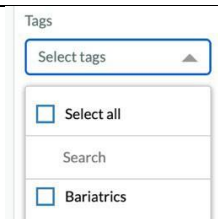
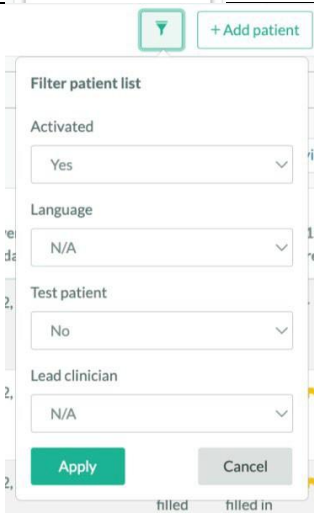
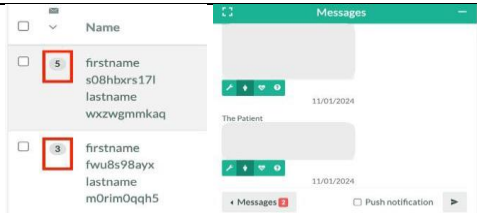
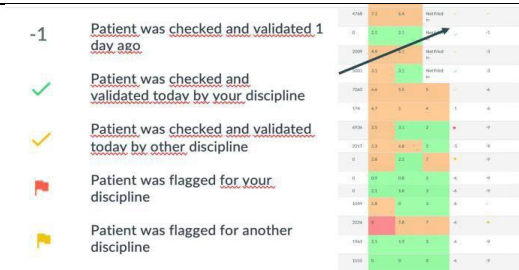
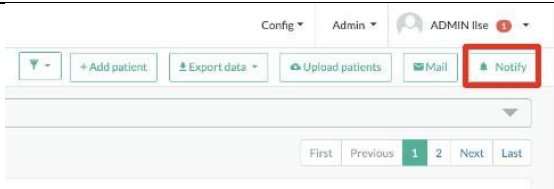
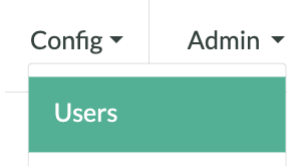
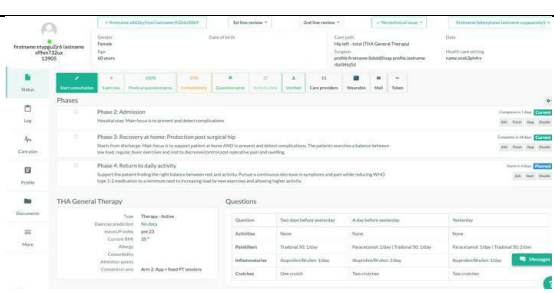
A mail template is shown. The HCP reviews the e-mail variables and clicks on 'send selected mails'.

An invitation mail will be sent to the patient with all the needed information to install and activate the application.

HEALTHCARE PLATFORM GENERAL FUNCTIONS

Below the HCP can retrieve all functions of the healthcare platform. The applicable functions will be shown based on the HCPs subscription. The extent of the available patient information is based on the patients subscription & his/her compliancy rate.

Patient overview		The health care platform gives a clear overview of the HCPs patients. In the top bar the relevant categories are shown that identify the patient and show relevant health parameters, unread messages, the lead clinician and patient compliance. By default, when opening the healthcare platform, the patient overview is sorted based on intervention date.
Tags		Select patients related to a specific pathology by using tags.
Filters		Select patients based on language, test patient, lead clinician or being activated or not. Activated can mean 2 things: either the patient registered him/herself with a code through the app or he has been marked as immediately activated in onboarding. Non activated: patients who haven't registered themselves through the app yet and haven't been marked as immediately activated.
Color codes		Color codes indicate which patients need closer attention.

Messaging		The number of unread messages are displayed in the overview.
Review / Flags		With the use of review flags it is known when the patient was checked or who needs to check the patient. See more info in Patient details – Review / Flags.
Notifications		Notifications can be sent directly to the patient's device home screen.
User management		Only limited HCPs of an organization have the permission to add users. In case a HCP has this permission, they can add users by clicking the 'New' button, fill out the credentials, select
		the applicable role and press 'Register'. A mail will be sent to the newly created user to follow the account creation steps, see above .
Patient details – Overview		In the patient details overview all relevant (health) information of the individual patient can be retrieved. With the use of the previous and next patient buttons the HCP can easily switch to another patient. In the left side bar, the HCP can navigate between status, log, care plan, profile, documents and more.

Patient details – Status

Task 1

← firstname y@bank 1ref lastname bankSource → → Flagged (no test) ←

	Gender	Female	Date of birth	
	Age	76 years	Language	Dutch
			Firstname	Shakeel
			Lastname	26719

Status	Start consultation	Exercises	75% Medical questionnaires	54.55% Consultations	Questionnaires	Reported events	Not verified
--------	--------------------	-----------	-------------------------------	-------------------------	----------------	-----------------	--------------

TKA General Therapy

- Tyot
- Exercise prediction
- moveUP index
- Current BMI
- Convention arm
- Patient variable
- Patient reported allergies
- Patient reported intolerance
- Patient reported GI issues
- Patient reported stationary bike
- Patient reported history of thrombosis
- Patient reported attention point
- HCP reported attention point
- HCP reported rehabilitation material

Questions

- Question
- Activities
- Walking aids
- Other medication
- Nonsmoking
- Optical

● Compliance OK for today or yesterday
● Non compliant > 2 days
● Non compliant > 7 days

The reported healthcare data are retrievable in the 'Status' section. From there it is possible to navigate to the recurring and medical questionnaire results, exercises and many more.

A **GREEN**, **YELLOW** or **RED** status dot shows how compliant the patient is to the care plan, for example filling out all the suggested recurring and medical questionnaires and performing the suggested exercises will result in a **GREEN** status dot.

A **GREEN**, **YELLOW** or **RED** status dot shows how compliant the patient is to the care plan, for example filling out all the suggested recurring and medical questionnaires and performing the suggested exercises will result in a **GREEN** status dot.

Patient details – Medical questionnaires

Status

Start consultation

7

30%

Exercises
Medical questionnaires

THA General Therapy

- ✓ Which medical questionnaires are active for a patient.
- ✓ How many of these are completed or not.
- ✓ What the reported results are on the individual questions within the medical questionnaire.

What the sub scores and total scores are to be able to see a quick indication on how the patient is progressing.

- ✓ Which medical questionnaires are active for a patient.
 - ✓ How many of these are completed or not.
 - ✓ What the reported results are on the individual questions within the medical questionnaire.
- What the sub scores and total scores are to be able to see a quick indication on how the patient is progressing.

What the sub scores and total scores are to be able to see a quick indication on how the patient is progressing.

Patient details – Recurring questionnaires
Patient details - Reported events

[illegible]

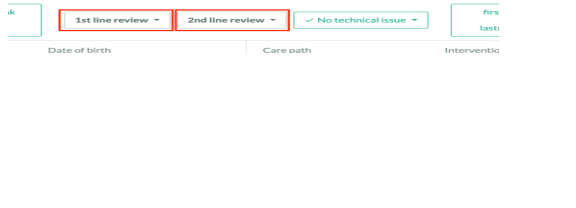
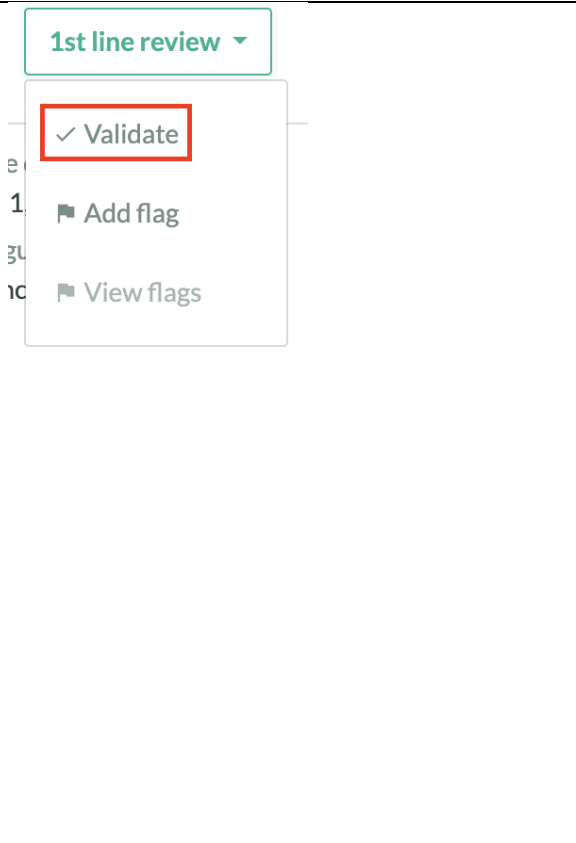
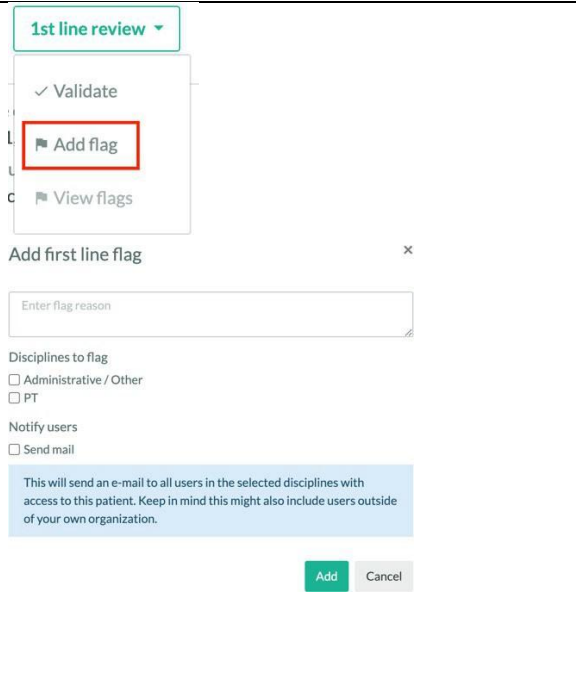
With these patient reported answers it is possible to follow the progress of important health parameters frequently.
The reported events are categorized per event. For each reported event you can view: (1) The date when the event was reported (+day in care plan) (2) The answers to all questions asked to the patient.

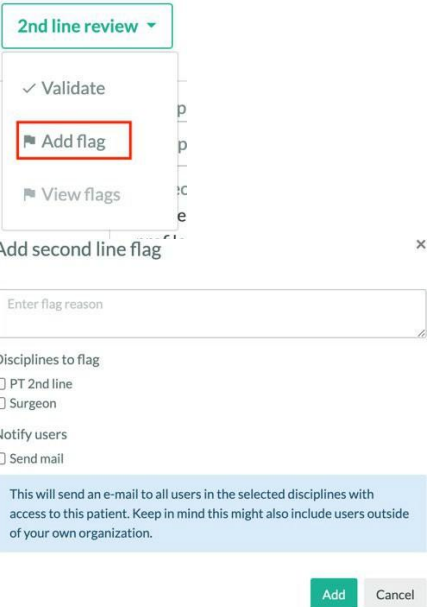
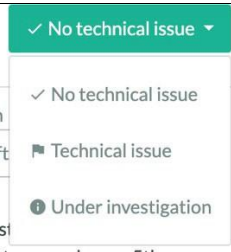
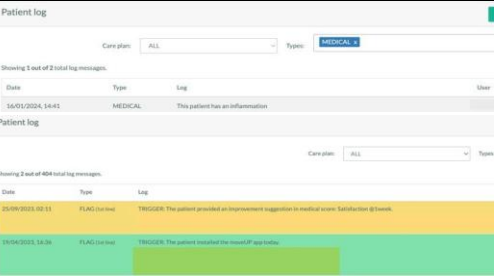
- Patient details
- Reported events

Reported events

Day	Reported at	Answers
242	Jan 23, 2015, 9:24:45 AM	0.52 0.81
		Yes Night : Back Exercises : Call No pain : Good

The reported events are categorized per event. For each reported event you can view: (1) The date when the event was reported (+day in care plan) (2) The answers to all questions asked to the patient.

Patient details - Review / Flags		1st line review is for the HCP that is not the lead clinician (e.g. physical therapist), 2nd line review is for the lead clinician (e.g. surgeon)
		Click validate when a patient has been monitored. Monitoring can mean: ✓ Questions answered - Advice given ✓ Adjusted or checked exercises ✓ Adjusted or checked step goal ✓ Checked daily parameters: pain, swelling, activity,... ✓ Checked postoperative protocol: use of crutches/walker, medication,... Same principle for 2nd line review.
		Flag someone from first line ✓ Patient needs to be checked by someone from first line (e.g. a physical therapist) ✓ Give a reason for flag ✓ Select the discipline to flag ✓ Optional: add an e-mail notification to flag. Flagged discipline will receive an email about the flag. A flag will add a new line in the patient log with the reason for flagging(see further).

		Flag someone from second line ✓ Patient needs to be checked by someone from second line (e.g. a surgeon) ✓ Give a reason to flag. ✓ Select the discipline to flag. Optional: add an e-mail notification to flag. The flagged discipline will receive an email about that flag.
		Report a technical issue.
Patient details - Log		Flags and technical updates automatically appear in the log of the patient. In case of flagging: ✓ GREEN : solved flag ✓ RED : Unsolved flag for discipline of logged in user (e.g. PT) ✓ YELLOW : Unsolved flag for other discipline than logged in user (e.g. Surgeon) Other important information can be added to the log manually. The applicable tags are added to the log in view of easy log searches

Patient details – Report

This report is used to visualize the main milestones and medical questionnaires of the patient, as well as the progression of the patient.

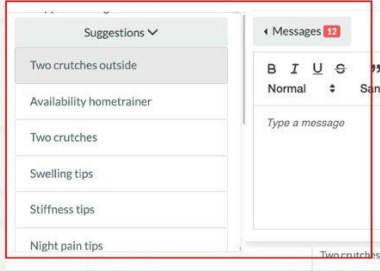
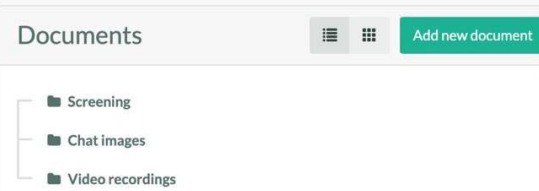
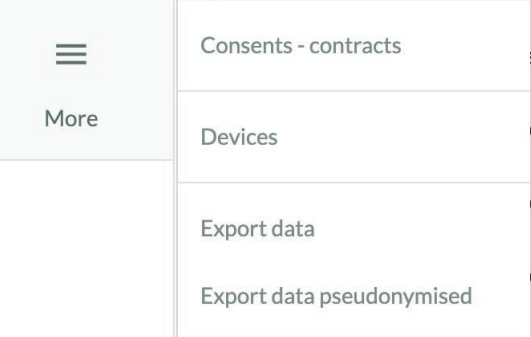
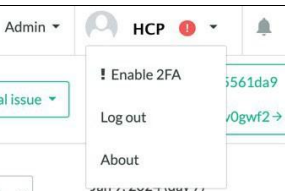
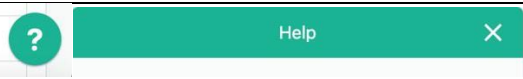
The report can be downloaded in pdf form by clicking the 'Download report' button. This way the HCP can send the report to the patient. It's also possible to share the link directly to the report, using the button 'Copy report link'.

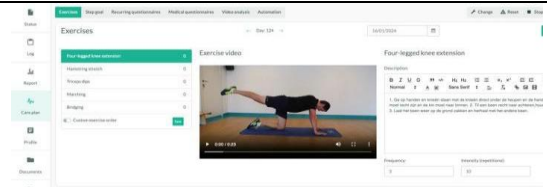
Patient details – Profile

The HCP can view, edit or add all relevant basic, clinical, administrative or care plan specific data in the patient profile

Patient details – Messaging

The message popup (which can be enlarged) is closed by default and the number of new messages is visible in the floating button. This functionality is used to give advice to patients and answer their questions. It can also be used to ask for pictures, for example to follow-up on wound healing, swelling, redness etc...

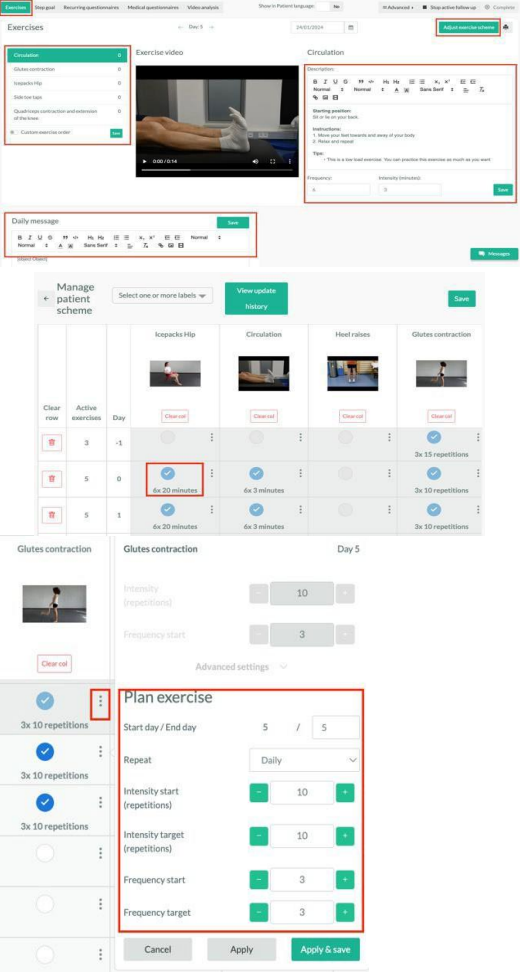
		When clicking 'messages', a set of prescribed personal messages is shown wherefrom the HCP can select an applicable message for that particular patient.
		Click on the "1st line symbol" to mark the question as read.
		Messages that are submitted by a HCP can be edited or removed. When hovering over a message, a button to edit or delete the message will become visible.
Patient details -- Documents		Documents can be uploaded in the 'Documents' section. Chat images and video recordings are automatically collected in this section in view of easy access to these materials.
Patient details - More		See which documents / contracts the patient has signed. See on which device the patient is using the application and when the application was first and last seen. Export data is only possible for limited HCPs due to security reasons and to protect data integrity.
Log out		In the upper right corner the HCP can log out. Know that the automatic log out time is 14 days, therefore it is recommended to log out more frequently to avoid unauthorized access.
Support		Reach out via support@moveup.care for questions or use the '?' button

		on the platform.
Care plan management		In care plan management it is possible to navigate between multiple functions that are applicable to the indication of the patient, such as exercises, step goal, recurring questionnaires, video analysis, automation, etc... For more information see further.

HEALTHCARE PLATFORM CARE PLAN FUNCTIONS

It is possible to provide care through the use of the healthcare platform. With the functions below the HCP can adapt the default treatment plan of their individual patient.

Only the applicable functions will be displayed conform the care plan.

Exercise management		The exercise management module includes the following possibilities of adaptation: ✓ Adapt description, frequency, and intensity of the exercises of the day. ✓ Adapt the daily message. ✓ Adapt the exercise scheme: ○ Blue check marks indicate active exercises, accompanied with frequency and intensity of these active exercises. ○ Schedule exercises over a longer period of time. Extend the standard exercise scheme with a set of exercises that are not part of the exercises defined in the care plan
---------------------	--	--

Recurring questionnaire

Exercises Step goal **Recurring questionnaires** Medical questionnaires Video analysis

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Using treatment scheme **Personalize**

Exercises Step goal **Recurring questionnaires** Medical questionnaires Video analysis Show in Patient language No Advanced Stop active follow-up Complete

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Name	Start day	End day	Frequency (times)	Active from	Reminders	Order	Preview	Remove
Discharge	9	10	1	10:00	No reminder set last		Preview	Remove
Postoperative rescue	0	3	3	10:00	No reminder set last		Preview	Remove
Activities	15	100	1	10:00	No reminder set last		Preview	Remove
Exercises	1	100	1	10:00	No reminder set last		Preview	Remove

Add questionnaire to personalized care plan

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Start date

End date

Frequency

Daily

Activation hour

15

Activation minutes

0

Save Cancel

Recurring questionnaires can be removed, reordered or a new recurring questionnaire can be added. The default list of recurring questionnaires, including their start day, end day, frequency and reminders is predetermined in the care plan.

Medical questionnaire

Exercises Step goal **Recurring questionnaires** **Medical questionnaires** Video analysis Change Reset Stop active follow-up Complete

Medical scores scheme

Name	Start day	End day	Preview	Remove
HIP PROM collection month 6	170	210	Preview	Remove
HIP PROM collection Year 1	340	500	Preview	Remove
Discharge Questionnaire HIP	0	15	Preview	Remove

Add medical questionnaire to care plan

Medical score

Select a medical score

Activation type

Day in treatment

Start Date

0

End Date

0

Save Cancel

Medical questionnaires can also be removed or added to the default scheme. The default list of medical questionnaires, including their start and end day is predetermined in the care plan.

Video analysis

Exercises Step goal **Recurring questionnaires** **Medical questionnaires** **Video analysis** Show in Patient language No Advanced Stop active follow-up Complete

Video analysis

Analysis enabled for this patient

Results NOT shown to this patient

Add request

Jan 12, 2024

Gait analysis

not recorded

Video requests

Name	Start day	End day	Preview	Remove
Gait analysis	-7	-1	Preview	Remove
Gait analysis	21	30	Preview	Remove
Gait analysis	40	45	Preview	Remove

Add video request to personalized care plan

Video request

Select a video type

Start date

0

End date

0

Save Cancel

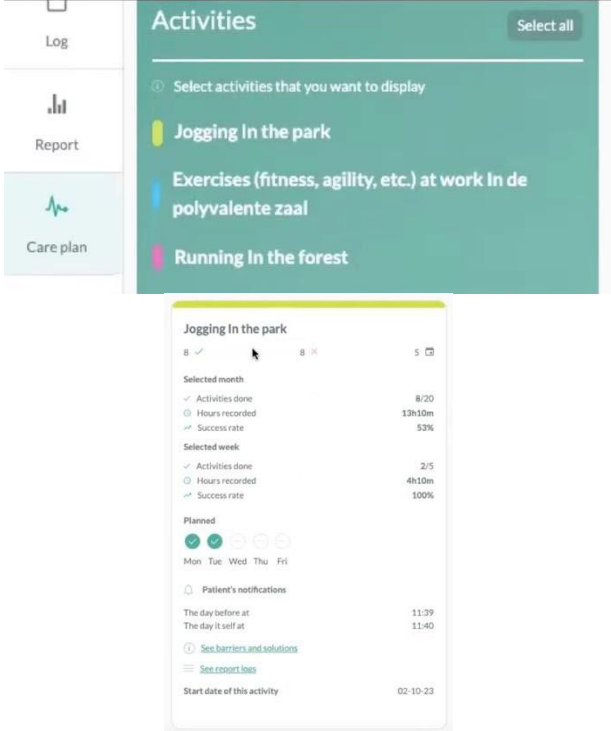
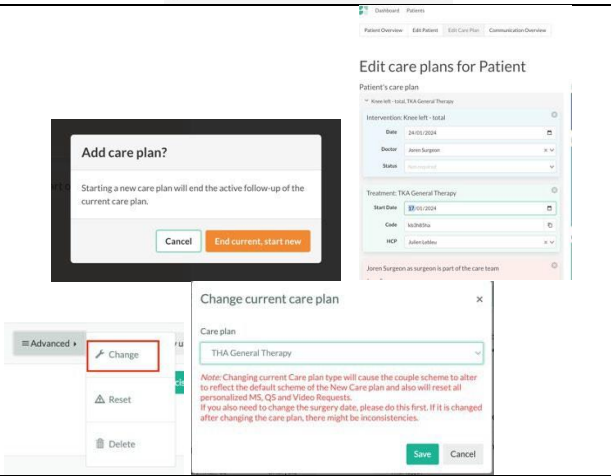
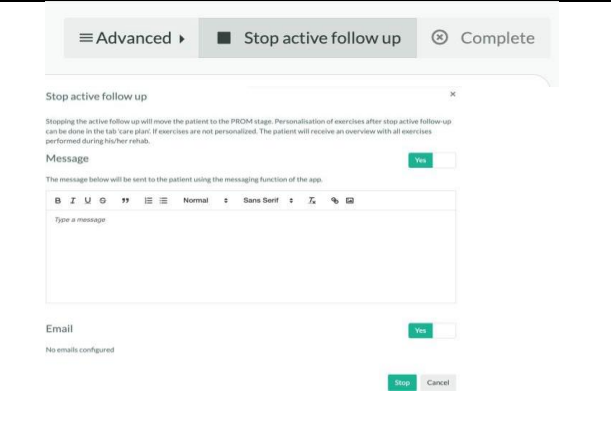
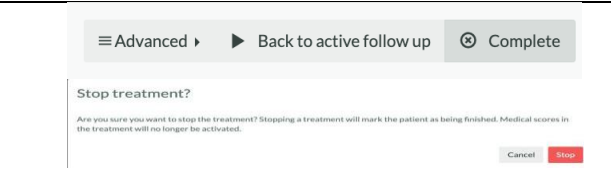
Request the patient to record a video for analysis

Video type

Select a Video type

Save Cancel

With the use of video requests, it is possible to gain more objective information, for example about the gait analysis or range of motion of an operated joint. The default list of video requests is predetermined in the care plan. A video request can be added directly by using the 'Add request' button or a video request can be scheduled for a longer period of time by using the 'Add' button and by adding a start and end day to put the video request active.

Goal manager		The HCP can enable and disable a behavioral change tool which enables patients to set their personal goal, define barriers and solutions in order to reach a healthy lifestyle. The HCP can see the progress on the patient's goals that are defined by the patient.
Change of care plan		The care plan can be adapted based on the patient's interventions, needs and wishes. This can be done via the onboarding tool or via the healthcare platform via the button 'Advanced'. Instructions how to change a care plan are clearly indicated within the tools.
Stop active follow-up		When the active follow-up is finished (what falls under active follow-up is predetermined in the care plan), press 'Stop active follow-up'. This will redirect the patient to a passive long term follow-up if applicable (also predetermined in the care plan if long term follow-up is needed).
		When no actions are required anymore from the HCP or the patient, the patient can be redirected to the 'Completed' status.

WELLNESS MODULE

The wellness module is an optional non-medical device module and therefore the data shown is not intended to be used for clinical decision making.

This data gives an overview of the general wellbeing of the patient.

The wellness module by default is blurred and needs to be activated to check the patient's wellness data.

WELLNESS MODULE ACTIVATION

To activate the wellness module you need to enter your profile settings and activate the option:

- Enable wellness data

moveUP Patients Onboarding Config Admin Laurens EN

Edit Laurens Baes

First name Laurens	Address Oktrooi plein 1
Surname Baes	City Gent
Username laurens@moveup.care	Postal code 9000
Email laurens@moveup.care	Country Belgium
Telephone number +32496123456	Title Title
Social security number Rizin nr	Image Image
Rizin nr Rizin nr	Language English
Authentication ☒ Two factor authentication	Wellness module ☒ Enable wellness module

Save Set new password Manage consent/contracts

Upon activation, a mandatory acknowledgment pop-up is displayed. It informs you that the data from the Wellness Module must not be used to support clinical decision-making.

Platform notice

Wellness module terms of use

Dear User,
You are about to activate the Wellness Module.

What is the Wellness Module?

The wellness module provides you with additional data about the general wellbeing and activity levels of your patients. The data in the Wellness Module can come from:

- imported data from general-purpose devices such as e.g. activity trackers
- manual entry by patients

This data does not come from medical devices. This means that neither the devices from which the data is imported, nor the manual entries by the patients have been verified or validated by moveUP, and therefore, it should not be used for clinical decision-making.

What is the purpose of this data?

We provide the option to see wellness data to give you a broader idea about the general wellbeing of the patients and to give you the opportunity to discuss and give further guidance regarding wellbeing.

This module is an optional module and is not linked to the medical data presented in the "Medical data" tab of the dashboard.

The data in the wellness module is not intended to be used for clinical decisions.

How to recognize the wellness data?

Wellness data is displayed in a separate tab in the dashboard or in a separate report.

Whenever wellness data is displayed in a column, together with other data it has a grey background to visually indicate that this does not concern medical data. Furthermore, it has a warning sign next to it.

Enabling/Disabling the Wellness Module

At any time, a user can enable/disable the Wellness Module in the user settings.

Acknowledgement

By clicking on the button below, you acknowledge that you have read and understood that the Wellness Module is an optional module that provides you with additional data regarding the general wellbeing of the patient. This data is not validated by moveUP and is not intended for clinical decision making.

☐ I agree

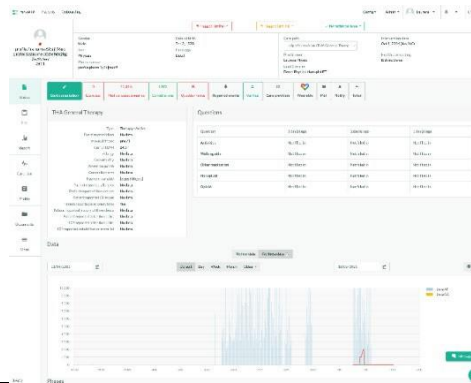
Sign

Cancel

Once acknowledged, you gain access to the wellness data. This activation allows you to view the wellness data and export the wellness reports.

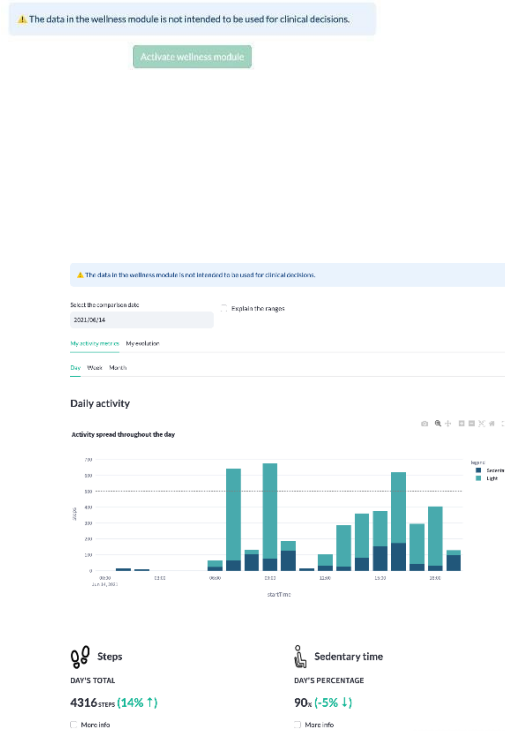
WELLNESS MODULE OVERVIEW

Wellness module overview

In this tab you'll be able to see the step data from the patient and see their evolution.

Wellness report

You can see the wellness reports when needed.

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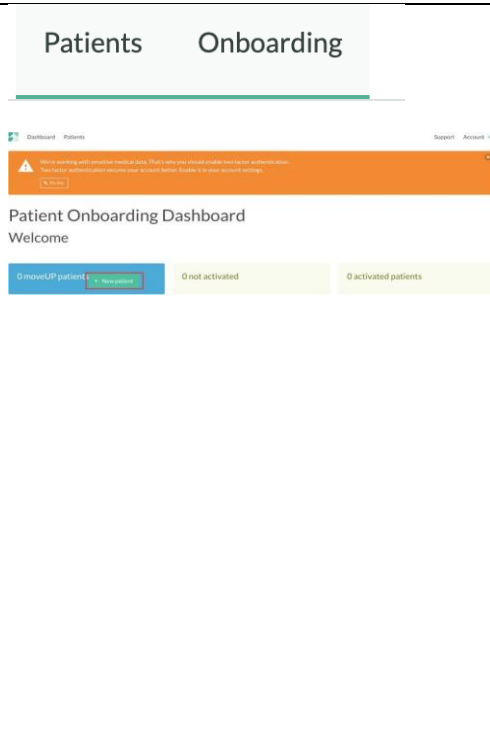
Accueil des patients

Patients Onboarding

Dashboard Patients Support Account

Warning: you are not logged in as a user. Please log in to your account to access the dashboard.

Patient Onboarding Dashboard

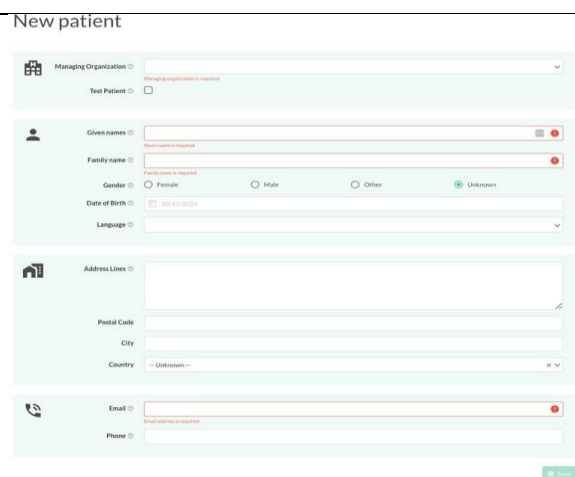
Welcome

0 moveup patient **New patient** 0 not activated 0 activated patients

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New patient

Managing Organization Managing organization is required

Test Patient ☐

Given names Given names is required

Family name Family name is required

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Patient's care plans

Start a care plan to add intervention and treatment for this patient

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The screenshot displays the 'Exercices' (Exercises) module interface. It includes a sidebar with a list of exercises, a main area with an exercise video and a 'Circulation' section, and a 'Manage patient scheme' section. A 'Plan exercise' dialog box is open, showing settings for 'Start day / End day', 'Repeat', 'Intensity start (repetitions)', 'Intensity target (repetitions)', 'Frequency start', and 'Frequency target'. The 'Plan exercise' dialog box is highlighted with a red border.

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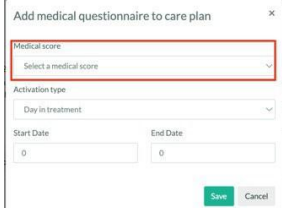
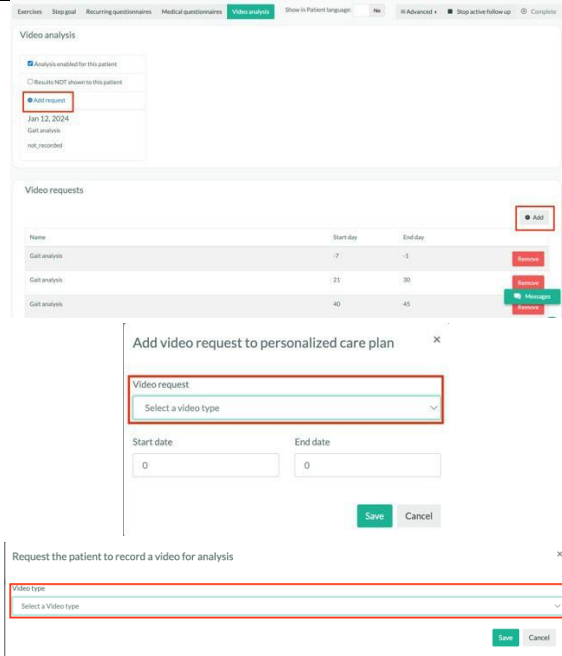
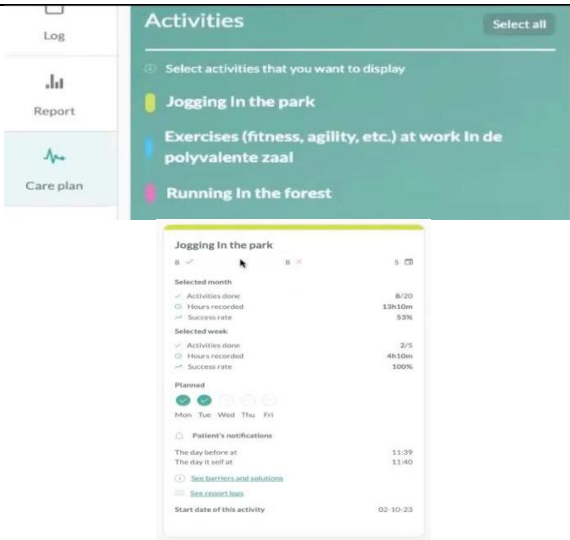
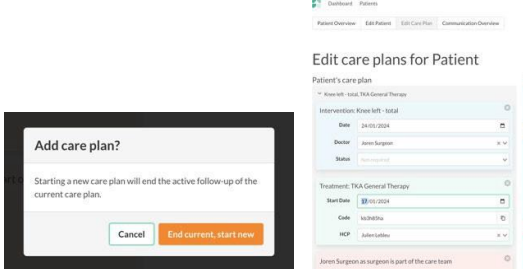
The screenshot displays the 'Questionnaires' module interface. It includes a sidebar with a list of questionnaires, a main area with a 'Questionnaires' section, and a 'Personalize' button. A 'Add questionnaire to personalized care plan' dialog box is open, showing fields for 'Questionnaire', 'Start date', 'End date', 'Frequency', 'Activation hour', and 'Activation minutes'. The 'Add questionnaire to personalized care plan' dialog box is highlighted with a red border.

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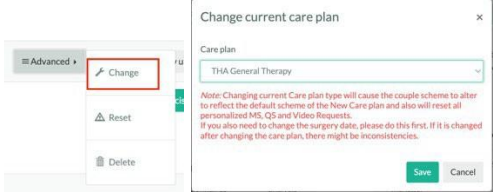
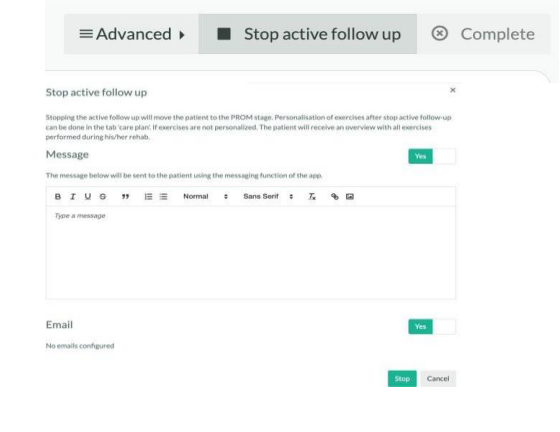
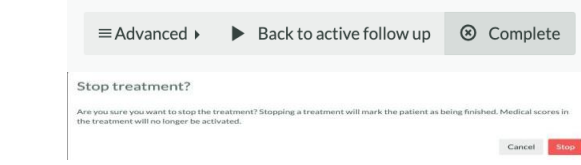
The screenshot displays the 'Medical scores scheme' module interface. It includes a sidebar with a list of medical scores, a main area with a 'Medical scores scheme' section, and an 'Add' button. A table lists medical scores with columns for 'Name', 'Start day', 'End day', 'Provider', and 'Remove'.

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The screenshot shows the 'Edit Laurens Baes' profile page. The page has a header with 'moveUP Patients Onboarding' and a user menu for 'Laurens'. The main content area is titled 'Edit Laurens Baes' and contains a form with the following fields:

- First name: Laurens
- Surname: Baes
- Username: laurens@moveup.care
- Email: laurens@moveup.care
- Telephone number: +32496123456
- Social security number: Social security number
- RIZIV nr: RIZIV nr
- Address: Oltroolplein 1
- City: Gent
- Postal code: 9000
- Country: Belgium
- Title: Title
- Image: Image
- Language: English
- Authentication: Two factor authentication (checked)
- Wellness module: Enable wellness module (unchecked)

At the bottom of the form, there are three buttons: 'Save', 'Set new password', and 'Manage consent/contracts'.

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Oefenbeheer

The screenshot displays the 'Oefenbeheer' (Exercise Management) module. It includes a 'Manage patient scheme' section with a grid of exercises. A 'Plan exercise' modal is open, showing configuration options for start/end day, repeat frequency, intensity, and frequency. Red boxes highlight the 'Add exercise' button, the 'Plan exercise' modal, and the 'Personalize' button in the 'Recurring questionnaires' section.

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The screenshot displays the 'Terugkerende vragenlijst' (Recurring Questionnaire) module. It includes a 'Questionnaires' section with a 'Personalize' button. Below, a table lists questionnaires with columns for Name, Start date, End date, Frequency, and Activation. A modal 'Add questionnaire to personalized care plan' is open, showing fields for Name, Start date, End date, Frequency, and Activation minutes. Red boxes highlight the 'Personalize' button and the 'Add' button.

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Medical score: Select a medical score

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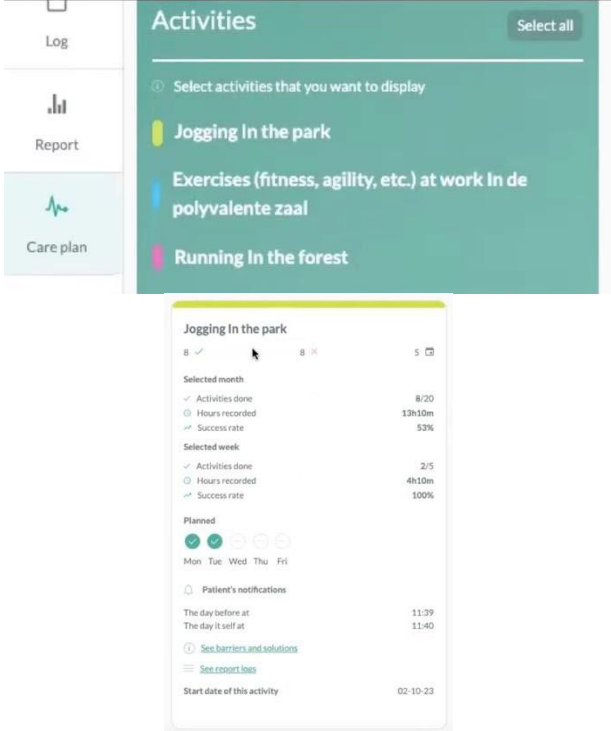
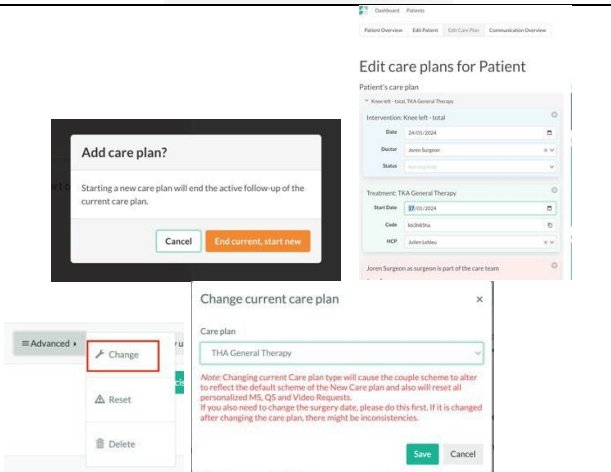
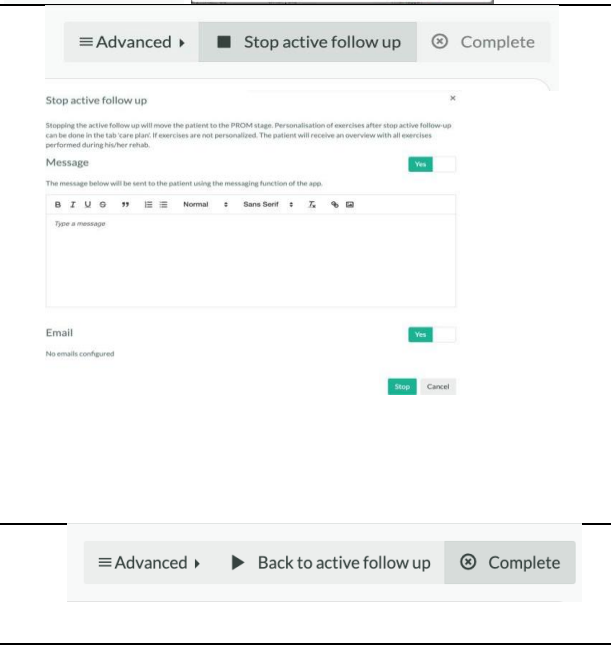
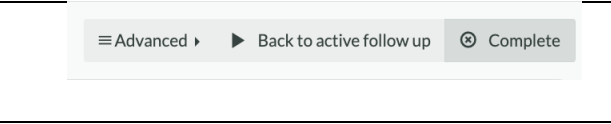
	Stop treatment? Are you sure you want to stop the treatment? Stopping a treatment will mark the patient as being finished. Medical scores in the treatment will no longer be activated. Cancel Stop	patiënt worden doorgestuurd naar de status 'Voltooid'.
--	---	---

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- Wellnessgegevens inschakelen

The screenshot shows the 'Edit Laurens Baes' form in the moveUP system. The form is divided into two columns. The left column contains fields for First name (Laurens), Surname (Baes), Username (laurens@moveup.care), Email (laurens@moveup.care), Telephone number (+32496123456), Social security number, and Rizin nr. The right column contains fields for Address (Oktrooi plein 1), City (Gent), Postal code (9000), Country (Belgium), Title, Image, and Language (English). At the bottom, there are two sections: 'Authentication' with a radio button for 'Two factor authentication' (selected) and 'Wellness module' with a radio button for 'Enable wellness module' (selected). The form has a 'Save' button and links for 'Set new password' and 'Manage consent/contracts'.

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