



Avive Solutions, Inc.
185 Valley Drive Brisbane, CA, 94005 USA
www.avive.life | 1 (800) 489-4428
Patents: www.avive.life/patents

© 2024 Avive Solutions, Inc. All Rights Reserved.




Avive® Power Case



MS-03104 Rev E

140x110mm


Get to know the Avive Power Case

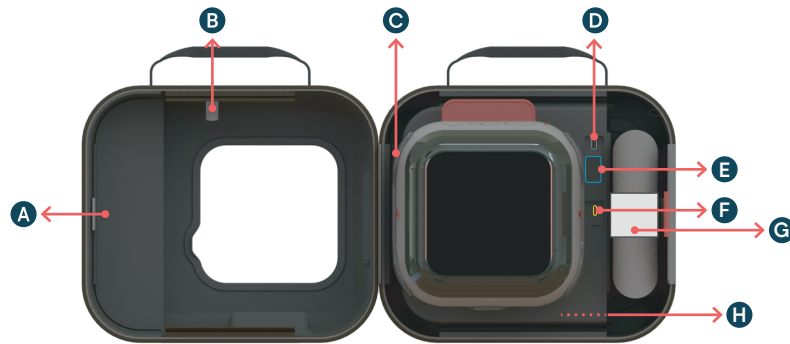

- A Storage Pocket** for the Avive USB Charging Accessories.
- B Status Lightpipe** shows the current state of the AED Status Light.
- C Holding Clamp** secures the Avive Connect AED® inside the Avive Power Case.
- D Power Bank Button** turns on the Display Screen on the Avive Power Bank.
- E Power Bank Display Screen** shows the charge level for the Power Bank. Will never show the battery percentage of the AED.
- F USB-C port** used to charge the AED and the Power Bank.
- G Elastic Strap** for the Avive Rescue Kit.
- H Avive Power Bank** for extending the standby time for the Avive Connect AED®.

Set up the Avive Power Case

Let's get started!

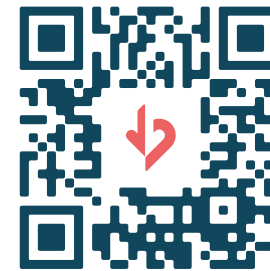

Scan for setup video!

Step 1: Set up the Avive Power Case


A. Verify the Power Bank's battery level by pressing the Power Bank Button (see page 2 for location). If the battery is not fully charged, follow the instructions on page 5 to charge the Power Bank.

B. Fully Charge the Avive Connect AED® before placing it in the Power Case. Use an Avive USB Charging Cable, an Avive USB Power Adapter, and a compatible wall outlet to charge.

C. Secure the Avive Connect AED® in the Power Case's Holding Clamp (see page 2 for location).

D. Locate the device charging cable inside of the Power Case and plug it into the Avive Connect AED®.

Important: The Power Case's device charging cable **MUST** be plugged into the Avive Connect AED® for the Power Case to function as intended.

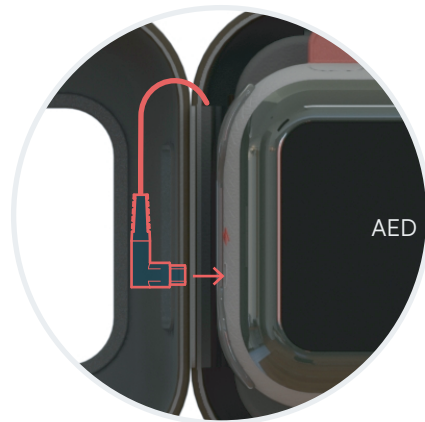


Step 2: Turn on Power Bank Mode

A. If using the Avive Connect AED®, navigate to Settings > Power on the device screen to view device power options.



B. Toggle "Power Bank Mode" **ON**.

NOTE: This will automatically configure your Avive Connect AED® to maximize Power Case standby time.


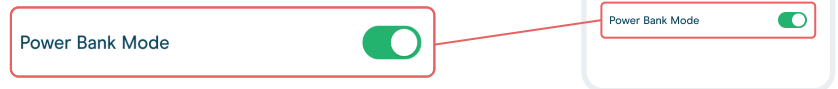


Done! The Avive Power Case will automatically turn on to extend the standby time of your Avive Connect AED®.

Recharging the Avive Power Case

A. When it's time to recharge the Power Bank in your Avive Power Case, first unplug the Power Bank Cable from the USB-C port on the Power Bank. →

B. Connect the USB-C end of a compatible charging cable to the USB-C port on the Power Bank. Connect the other end of the cable to a compatible adapter and wall outlet and let charge for 24 hrs or until the battery level reaches 100%.

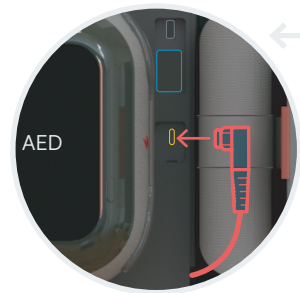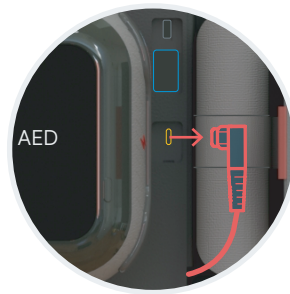

← **C.** Once the Power Bank is fully charged, unplug the cable and reconnect the Power Bank Cable into the USB-C port on the Power Bank.

Important: The Power Bank Cable **MUST** be plugged into the Power Bank so that the Avive Power Case can charge your device.

Operating & storage conditions

Temperature	32°F – 95°F (0°C - 35°C)
Humidity	5 - 95% non-condensing
Altitude	< 9,842 ft. (3,000 m)
Water & Ingress	Do not expose the Avive Power Case to any liquids.



Keep away from liquids



Do not drop



Do not expose to extreme temperatures



Keep away from direct sunlight



Do not disassemble



Keep away from fire

Avive AED and Power Case Standby Time: Up to 5 years **USB-C output:** 5V/3A
Battery Rated Capacity: 42000mAh (5V/3A) **USB-C input:** 5V/2A