

Assessment of Structural and Copy Number Variants in Constitutional Samples using Electronic Genome Mapping

Adam C. Smith^{1,2}, Michael Howard¹, Amanda Sussman¹, Karen K. Phillips¹, Patrick A. Lennon¹, Venkateswara R. Potluri¹, Robert G. Best¹, Lily Nasanovsky³, Anna Lozar³, Yang Zhang³, Hanae Sugiura⁴, Shuk Shukor⁴, and Eiichi Araki⁴

Introduction

The detection and interpretation of structural variants (SVs), including copy number variants (CNVs), is the key function of modern cytogenetic laboratories and is important for diagnosis, clinical management, and genetic counseling. Current cytogenetic laboratory methods such as chromosomal microarray analysis (CMA), fluorescence *in situ* hybridization (FISH), and chromosome banding analysis (CBA) can each provide structural genomic information, but each technique has specific limitations in their ability to resolve certain types or sizes of SVs. Electronic genome mapping (EGM), a solid-state nanochannel-based platform, is a high-resolution technology capable of detecting a broad range of structural genomic variations. We present the concordance of EGM for the detection of SVs in a cohort of constitutional and prenatal samples with findings previously characterized by conventional cytogenetic testing and Optical Genome Mapping (OGM).

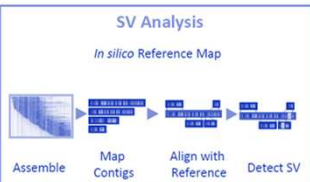
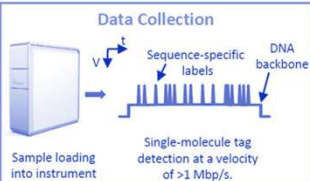
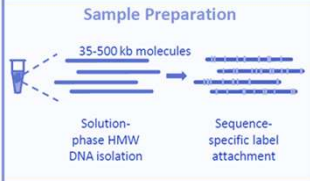


Figure 1. EGM sample processing workflow. Extraction of HMW DNA is labeled at specific sequences (top). Sample is loaded on detector and data collected in OhmX instrument (middle). Molecule data is assembled and SV and CNV events detected (bottom).

Reference: Kosyakova et al. Heteromorphic variants of chromosome 9. *Mol Cytogenet.* 2013;6(1):14.

Affiliations:
¹ Labcorp, Center for Molecular Biology and Pathology (Cytogenetics), Research Triangle Park, NC.
& ² Dept of Laboratory Medicine and Pathobiology, University of Toronto, Toronto, Canada,
³ Nabsys 2.0, LLC, Providence, RI, ⁴ Hitachi High-Tech America, San Diego, CA

American Society of Human Genetics Annual Meeting 2026; Oct 20-24; Montreal, Canada
© 2026 Labcorp. All rights reserved.

Methods

Cohort
Blood specimens referred for follow up of known chromosomal abnormality were analyzed in this study and included samples with genomic findings previously characterized using orthogonal methods such as CMA, quantitative PCR, CBA and OGM. Institutional IRB approval under # SCMM-RND-403.

EGM
DNA was extracted using NEB Monarch[®] high molecular weight DNA Extraction Kit and processed using a dual-enzyme workflow and the OhmX[™] Platform. SVs were identified using Human Chromosome Explorer[®] (HCE) *de novo* assembly, while copy number alterations were assessed using a CNV analysis pipeline.

OGM
Ultra high molecular weight DNA was extracted and labelled using the Bionano Blood Extraction kit. Labelled DNA was loaded on a nanochannel chip and imaged using a Stratys[™] instrument. Genome assembly was performed using the Guided Assembly pipeline and visualized in Access and/or VIA[™].

Array
DNA was run on Thermo-Fisher Cytoscan[™] HD Accel assay, processed using Genechip System 3000 and analyzed using Chromosome Analysis Suite software.

Table 1. Concordance of Results Between Array, EGM and OGM

Sample	Event	Region	EGM [GRCh38]	Size in kb	OGM [GRCh38]	Size in kb	Array [GRCh38]	Size in kb	Comment
EGM361	Inversion	6q14.1	6-8(q14.1;q14.3)(77,948,869-86,856,919)	8,908	6-8(q14.1;q14.3)(77,622,617-77,624,243-87,207,318-87,209,933)	9,583	Not Detected	-	Array cannot detect balanced rearrangements (Fig 2C)
	Duplication	6q14.1	6q14.1(77,641,282-77,881,326)x3	240	6q14.1(77,642,688-77,881,674)x3	239	6q14.1(78337721-78655840)x3	318	Deletions and duplications flanking inversion is a common mechanism
	Duplication	6q14.3	6q14.3(86,825,273-87,230,273)x3	405	6q14.3(86,829,785-87,229,762)x3	400	6q14.3(87565439-87920197)x3	355	Deletions and duplications flanking inversion is a common mechanism
	Loss	15q11.2	15q11.2(22,630,042-23,489,042)x1	850	15q11.2(22,522,095-23,392,094)x1	870	15q11.2(22595659-23102645)x1	507	Size discrepancy due to probe placement in repetitive region
EGM366	Duplication	22q11.1	22q11.1(16,844,635-17,244,635)x3	400	dup(22)(q11.1)(16,948,763-17,204,010)	255	22q11.1(16947158-17235989)x3	289	22q11.1-q11.2 region challenging due to presence of multiple low copy repeats
EGM378	Deletion	13q31.2q31.3	13q31.1q31.3(82,870,000-90,720,000)x1	7,850	del(13)(q31.1q31.3)(82,950,924-82,954,146-90,499,601-90,502,823)	7,545	13q31.1q31.3(82,947,255-90,503,450)x1	7,556	All techniques highly concordant
EGM349	Gain	Xp22.33	9p12p11.1(38,810,000-43,213,392)x4	4,400	dup(X)(p22.33)(1,653,916-2,464,217)	810	Xp22.33(1,630,017-2,468,877)x3	839	Chromosome 9 heteromorphism (Fig 3D).
	Gain	9q13q21.11	9q13q21.11(61,686,937-66,338,937)x4	6,600	9p13.1p11.2(36,780,001-41,246,000)x4	2,466	ND - last probe at ~38.9Mb	4,626	EGM/OGM has better coverage in 9p centromere than array (Fig 3A,B,C)
EGM371	Deletion	10q21.1	del(10)(q21.1)(54,850,341-54,939,277)	89	del(10)(q21.1)(54,846,835-54,943,344)	96	10q21.1(54,849,462-54,936,124)x1	87	All techniques highly concordant
EGM381	Deletion	9p21.3	del 9p21.3(22,224,822-22,749,272)	515	Sample unavailable	-	9p21.3(22,219,776-22,739,865)x1	520	EGM and array highly concordant


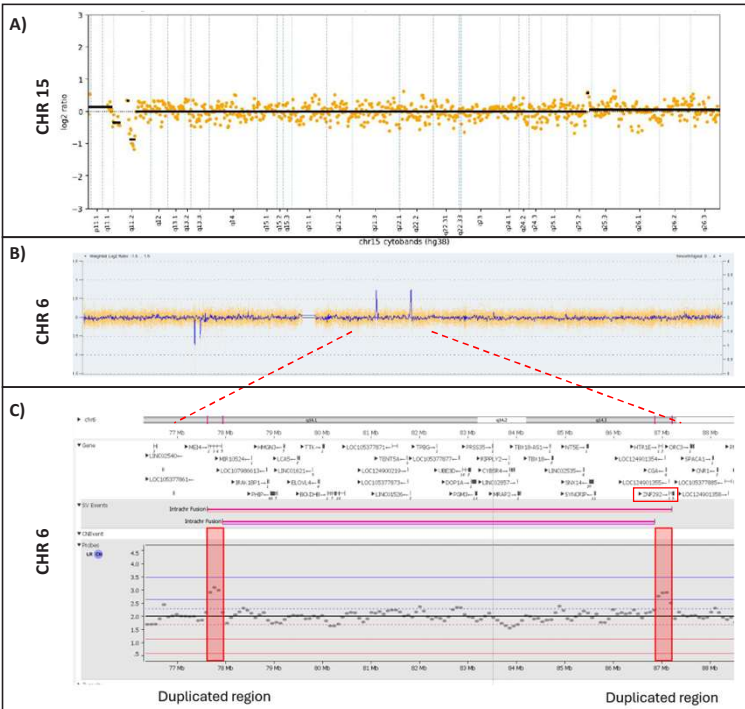


Figure 2. Data for EGM361. A) EGM copy number analysis confirms a deletion on proximal chromosome 15 in the BP1-BP2 region. B) An observation on the array showed two small duplications on chromosome 6. These regions did not meet positive reporting criteria for array. The mechanism (gain) did not match any genes in the region susceptible to triplosensitivity. C) Genome mapping shows that the duplications are actually flanking an intrachromosomal inversion. The distal duplication on chromosome 6 disrupts the coding region of the gene ZNF292 (red box), which has a Pathogenic/Likely Pathogenic ClinGen Score for haploinsufficiency.

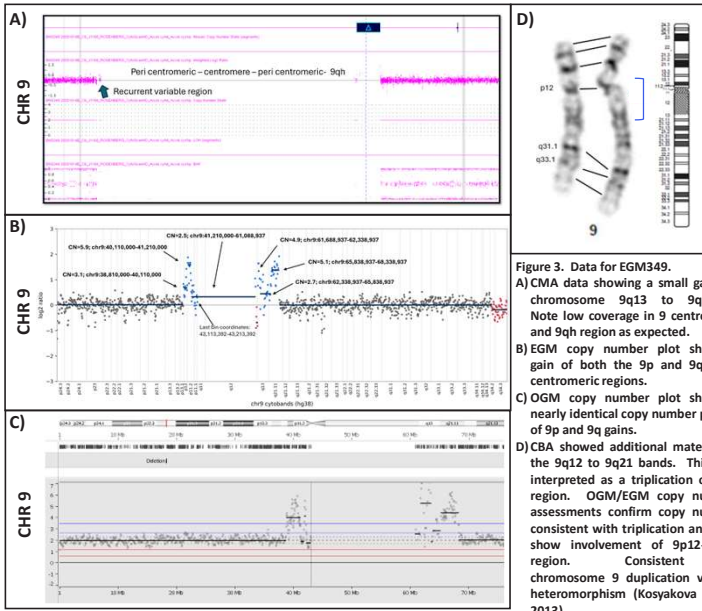

Results & Conclusions

EGM361 shows a 15q deletion concordantly detected by all platforms. An inversion on 6q was detected only by EGM and OGM flanked by two duplications (confirmed by array). The inversion disrupts a gene, ZNF292, that is relevant to the reason for referral showing benefit of comprehensive high resolution structural analysis. EGM349 shows a 9q variant with gain of the p- and q-arm. Array shows the triplication but has sparser probe coverage (especially in the p-arm) due to the highly repetitive nature of the sequence, but with more coverage in this region for EGM and OGM that provides additional clarity to the karyotype results. EGM, OGM and array results were highly concordant and show EGM to be a promising, emerging technology for integrated genome-wide assessment of SVs in constitutional and prenatal samples.

Poster #: 2025F