



B/S[®]
SPANGE
Classic

Instruction

A step-by-step guide

Instruction - B/S Spange Classic

A step-by-step guide

Points to note before treatment

The B/S brace must not be applied to diseased nails, e.g., fungal disease, psoriasis, onycholysis, etc.

- No footbaths
- No disinfectants
- No application of oils or lotions
- No nail polish remover
- No wet grinding!

Points to note after treatment

- Seal the affixed B/S brace with our B/S glue.
- The sealant is resistant to acetone or nail polish remover.
- Final treatments can only be carried out after sealing, e.g., applying lotion, oil, nail polish, tamponade or wound care.

Further important notes

- The batch number of the glue can be found on the cap and is the expiration date. **Example: batch number 1025, expiration date October 2025.**
- Do not wipe the brush of the glue on the neck of the bottle as this will result in the cap getting glued. This also makes sure you have the perfect amount of glue.
- The glue and the activator must not come into contact with each other to avoid negative consequences:
 - › The shelf life of the glue will be shorter
 - › The glue will get stringy
 - › The bottle content will dry up
- Freeze the glue! The shelf life of the glue is extended by the time it is kept in the freezer. **CAUTION:** The glue must be applied at room temperature, approx. 21 to 24° C. Please also ensure that the B/S adhesive is applied at an air humidity of at least 40%, as it no longer polymerizes below 40%.
- Test the activator: Run the pen over a sheet of paper. If it releases liquid, it is still working. **Important:** Always close the activator pen tightly!
- If the tensile force of the brace is too strong after affixing it, just carefully grind both ends of the brace thinner. If necessary, grind the entire brace thinner.

Training video

If you prefer to watch a video, you are welcome to visit our youtube channel:

www.youtube.com/@BSSpange/videos



How do I scan a QR code?

Open the camera or a QR code scanner app on the smartphone. Hold the camera over the QR code. A link will appear, please tap it and now you can select the appropriate training video.



1. Preparing the nail

For the firm and secure fixing of the B/S brace, it is necessary to grind and roughen the nail plate. Thick nails should be filed down to make them thinner.

Use a diamond or corundum grinder to roughen the surface of the nails.



2. Determining the brace size

The B/S brace must be applied in such a way that it spans the visible nail plate, but it should not come through the ends or be too short, if possible. The rough side of the B/S brace is the side which should be measured and affixed to the nail.

Placing


The best result with the B/S brace is achieved when the B/S brace is affixed 1 to 2 mm behind the area which is causing pain.



3. Degreasing

After determining the size of the B/S brace, **degrease the entire nail surface with one drop of B/S cleaner on a swab.**

Wipe the nail dry with the dry part of the swab.

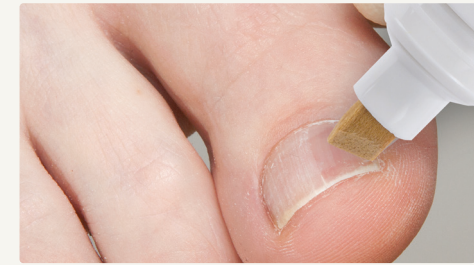

4. Applying B/S activator

Apply activator **only to the side edges of the nail** with the B/S activator pen.

The activator is necessary for the secure and immediate bonding of the adhesive.



5. Preparing the B/S brace

Using a diamond or corundum grinder, grind the ends (approx. 2 to 4 mm) of the B/S brace thinner on the rough side. Through this, the ends can be pressed to the nail more easily and do not stick out.



6. Applying the B/S glue

Hold the chosen B/S brace as shown and, taking into consideration the size of the brace, apply a small drop of B/S glue to the centre of the brace.

Hint:

Our instruction is based on the standard size 18. The larger the B/S brace, the more glue is needed.

Possible alternative: Apply the glue directly to the nail and the activator to the brace.

You can find more about the alternative methode on our website www.bs-spange.de



6.1. Fixing the centre part correctly

Place the B/S brace 1 to 2 mm behind the area which is causing pain and gently press it onto the centre of the nail **for approx. 10 seconds** without slipping.

Our tip:

Keep calm! If possible, avoid slipping with the tweezers during application and do not work without supporting yourself. Support yourself on the patient's foot.



6.2. Fixing the right and left end

Now press the B/S brace onto the right lateral nail plate with the non-slip and blunt part of the applicator - for approx. 10 seconds.

Caution:

Make sure that the ends lie flat against the nail and do not stick out

Repeat this procedure on the left side.



7. Individually adjusting the tensile forces

If your patient feels an **unpleasant pulling or tension pain** after affixing the B/S brace, the tension is too strong.

In this case, you have to gently grind the ends of the affixed B/S brace thinner or, if this does not prove sufficient, the entire B/S brace.

Your advantage:

Corrective adjustments to the B/S brace can be made on the nail!

The patient feels the correction immediately.

Brace ends sticking out can be ground off if less than 1 mm.



8. Cleaning the nail with B/S cleaner

Please remove the activator before sealing the nail. Put a **drop of B/S cleaner** onto a swab and wipe over the nail plate. This removes the activator previously applied.

Caution:

The activator must not come into contact with the glue, as otherwise the glue will dry out.



9. Final sealing

Seal the nail with the B/S brace already affixed with the B/S adhesive to make it particularly secure and water-resistant.

This sealant is resistant to acetone or nail polish remover!

Brush the B/S glue like nail polish over the entire nail surface.



10. Further treatment

Once the sealant is dry (after approx. 5 minutes), standard nail care can be continued:

- Placing cotton into the nail fold
- Applying ointment dressings
- Applying nail polish
- Footbaths
- Application of lotion and oils



11. Removing the B/S brace

It is best to use straight blade nail pliers to remove the B/S brace. Place the pliers at one edge ,between the brace and the nail’.

This way you can easily burst the sealing coat and carefully remove the B/S brace from the nail (see pictures above). The glue, which acts as a sealant, can now be milled off.

Optimisation of the corrective treatment

The best correction of an ingrown nail is achieved when a new B/S brace is applied **every 4 weeks**.


It can also stay affixed to the nail for a longer time - until it comes off by itself. However, it has **little to no lifting force after 4 weeks**.

If necessary and if there is enough space between the nail root and the first B/S brace, a second B/S brace can be affixed to the space after one month or later.

Caution:



Pay attention to the lifting forces when using two B/S braces.

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

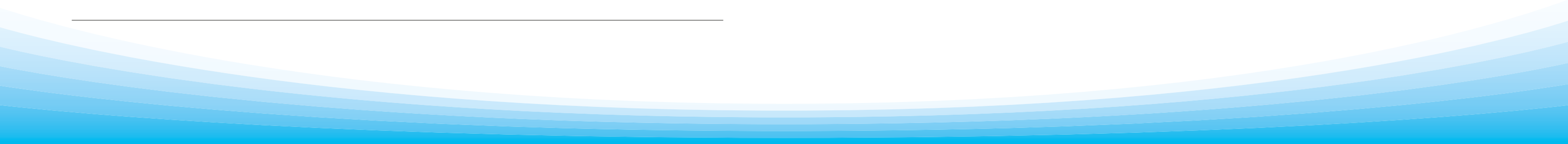

Absolutely original

for over 35 years.

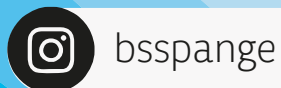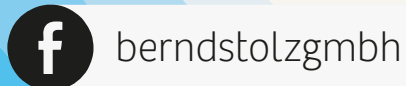

Bernd Stolz GmbH

Fuggerstraße 36 a
D-92224 Amberg

+49 96 21 / 22 1 88
info@bs-spange.de

www.bs-spange.de
www.shop.bs-spange.de