



Clinical Case Cards:
From Protocol to Results

- Proven by Practice -





Alma **Harmony**

The background of the image is a dark, high-contrast photograph of a mechanical device, likely a scale or a precision instrument. The device has various components, including a central vertical column, horizontal arms, and a base. The lighting is dramatic, highlighting the metallic surfaces and creating deep shadows. On the right side, there is a small, partially visible label that reads "Alma o".



Skin Rejuvenation

ClearLift Pro
Iris VL/PL
NIR

Before



After 3 Tx



Treatment Protocol

ClearLift Pro 5x5: -2, -1 and 0 | 1800-2400mJ, Total: 5000J

Iris VL/PL: 16-20J | 12-15msec | 2 passes

NIR: 17-20W | 8 kJ

Sessions intervals: 4 weeks

Courtesy of Dr. Pablo Naranjo, Spain

Alma
Harmony

Alma o



Skin Rejuvenation

ClearSkin Pro
ClearLift Pro
Iris Dye VL/SR
SupErb

Before



1 month post 2 Tx



Treatment Protocol

ClearSkin Pro 7x7: 900-1500 mJ/pulse | 3Hz In-motion | total energy 1500J/face

ClearLift PRO 5x5: -2, -1, 0 | 1800mJ/pulse | 5kJ per face

Iris Dye VL/SR: In-motion | 3.4J/cm² | 3.6kJ total energy

SupErb 7x7: 1500-2500mJ / pulse | 1 pass

Courtesy of Dr. Eyal Kramer, Israel

Alma
Harmony

Alma o




Skin Rejuvenation & Skin Tags

ClearLift Pro
Iris VL/PL
NIR
SupErb

Before



After 3 Tx

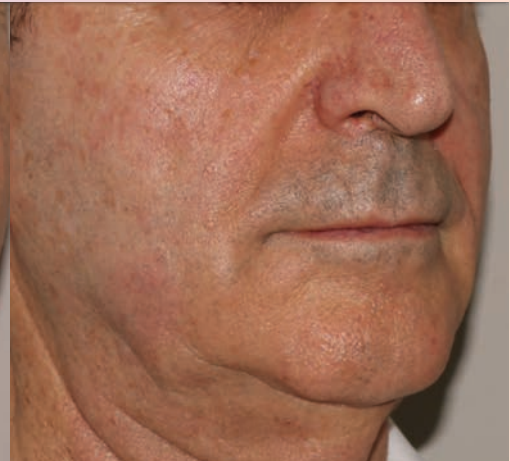

Treatment Protocol

ClearLift Pro 5x5: -2, -1 and 0 | 1800-2400mJ | Total: 5000J


Iris VL/PL: 16-20J | 12-15 ms | 2 passes

NIR: 17-20 W | 8 kJ

SupErb 7x7: 2500-3500 mJ/pulse | PW: 1,5 ms

Sessions intervals: 4 weeks

Courtesy of Dr. Pablo Naranjo, Spain

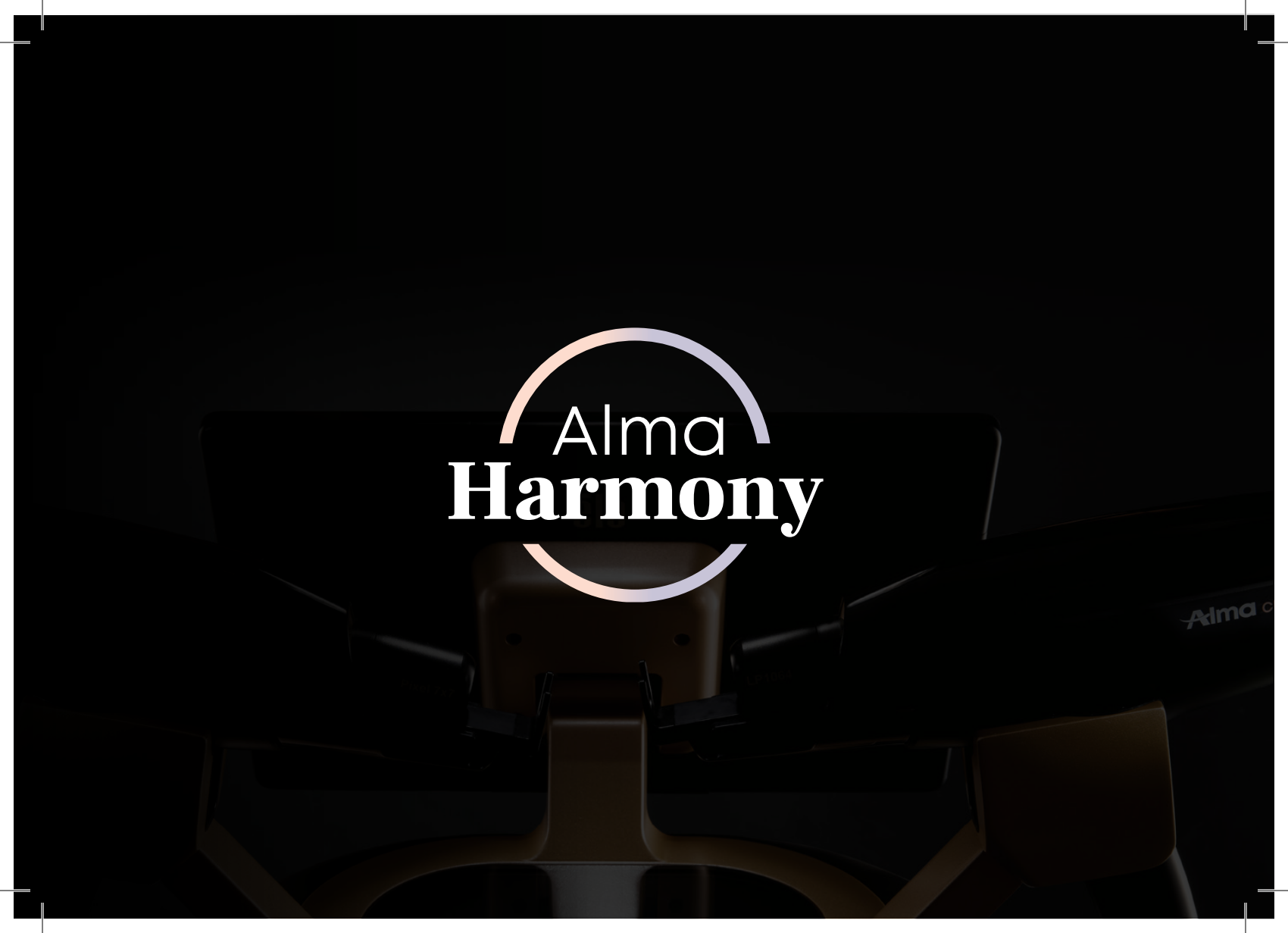The background is a dark, moody photograph of a violin, focusing on the body and f-hole. A semi-transparent logo is centered over the violin. The logo consists of a circular arc with a color gradient from light orange to light purple. Inside the arc, the word "Alma" is written in a white, clean sans-serif font, and below it, the word "Harmony" is written in a larger, white, elegant serif font.

Alma Harmony

Alma o



Texture & Tone

ClearSkin Pro
ClearLift Pro
Iris Dye VL/SR

Before

1 month post 2 Tx



Treatment Protocol

ClearSkin Pro 7x7: 900-1500 mJ/pulse | 3Hz In-motion | total energy 1500J/face

ClearLift PRO 5x5: -2, -1, 0 | 1800mJ/pulse | 5kJ per face

Iris Dye VL/SR: In-motion | 3J/cm² | 3.6kJ total energy

Courtesy of Dr Eyal Kramer and In house clinic



Alma **Harmony**

Alma o



Tone, Texture & Laxity

ClearSkin Pro
Iris Dye-VL
SupErb

Before

1 month post 1 Tx




Treatment Protocol

ClearSkin Pro 7x7: 900-1500 mJ/pulse | 3Hz In-motion | total energy 1500J/face

Iris Dye-VL: In-motion 3J/cm2 | 3.6kJ total energy | Stationary 8-10J/pulse

SupErb 7x7: 1500 - 2500mj/pulse | 1 pass

Courtesy of Dr. Eyal Kramer and In house clinic

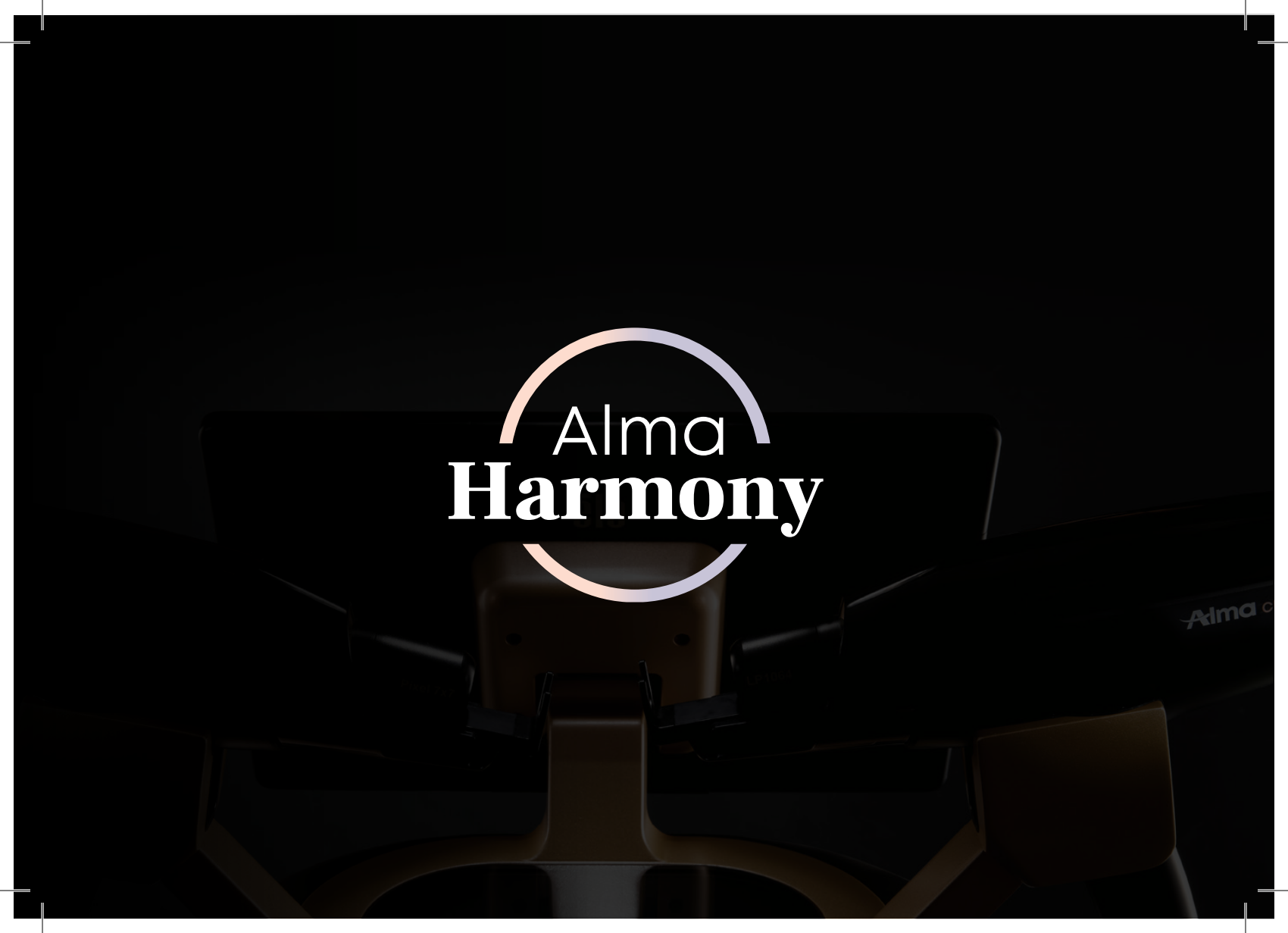

Alma **Harmony**

Alma o



Lip Rejuvenation

ClearLift Pro

Before

Post 6 Tx




Treatment Protocol

ClearLift Pro 5x5:

800-1200 mJ/pulse | 500J (-2), 500J (-1), 500J (0) | 4 Hz

Courtesy of Dr. Hen Ifrach and In house clinic

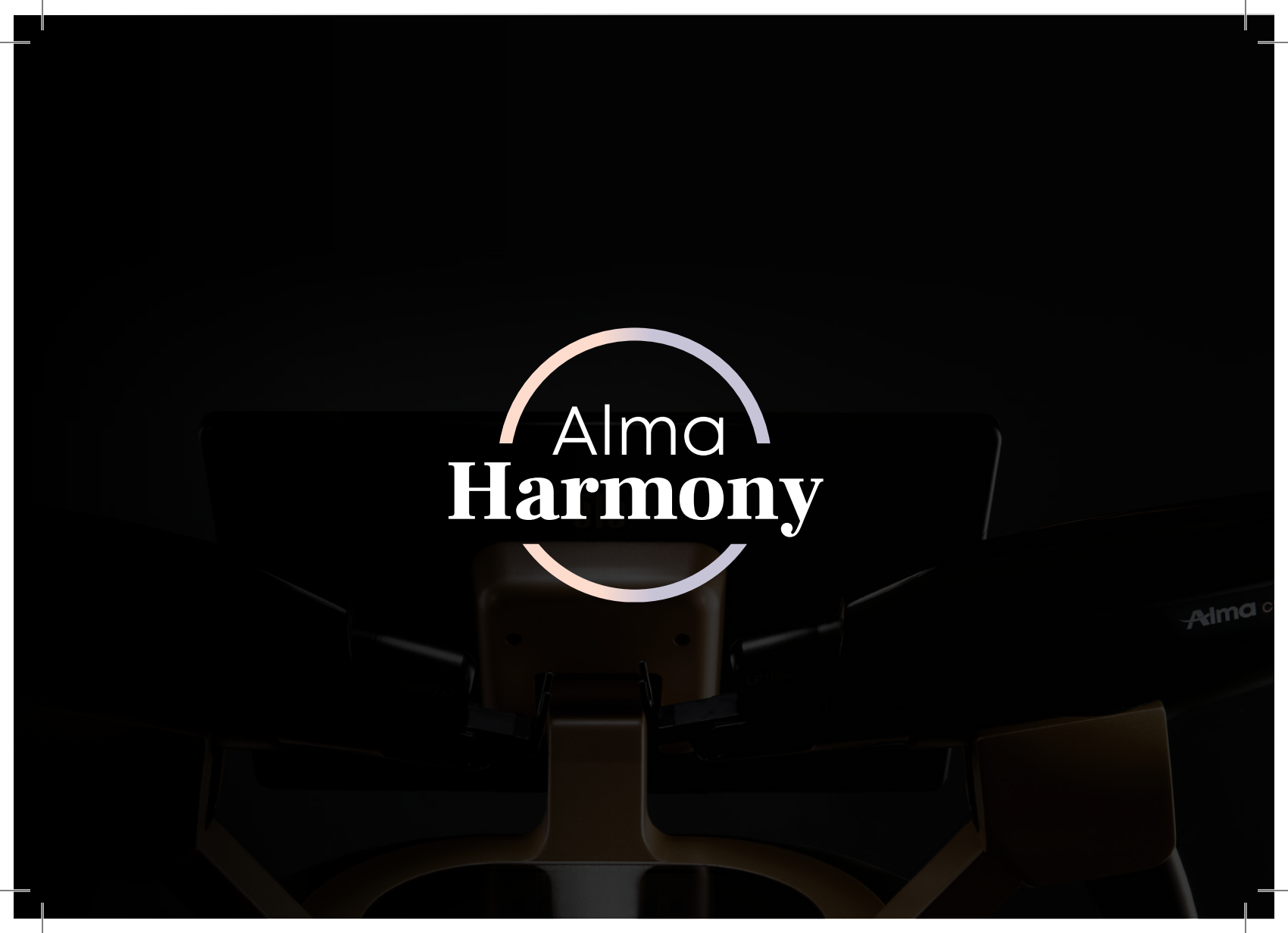

Alma **Harmony**

Alma o



Prejuvenation

ClearSkin Pro
ClearLift Pro


Before

After 6 Tx



Before

After 6 Tx

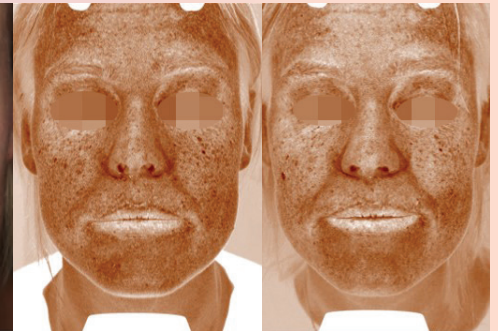

Treatment Protocol

ClearSkin Pro 7x7: 1500 -1800 mJ/pulse | 2 passes

ClearLift Pro 5x5: -2, -1 and 0 and +1 | 5 KJ each side of face

Clearskin Pro: 4mm for active spots

Courtesy of Dr. Munir Somji, UK



Alma **Harmony**

The background of the image is a dark, high-contrast photograph of a mechanical device, likely a scale or a precision instrument. The device has a complex structure with various components, including a central vertical column and several horizontal arms. The lighting is dramatic, highlighting the metallic surfaces and creating deep shadows. On the right side, there is a small, partially visible label that reads "Alma o".



Prejuvenation

ClearSkin Pro
ClearLift Pro
Iris Dye-VL

Before



Post 3 Tx



Treatment protocol

ClearLift Pro 5x5: -2, -1 and 0 | 1200-1800mJ | 2Hz | Total: 5000J

ClearSkin Pro 7x7: 1000-1500J | 3Hz | In-motion | Total: 1500J

Iris Dye-VL: Stationary mode | 8-10J | 12ms | 2 passes

Sessions intervals: 4 weeks

Courtesy of Dr. Pablo Naranjo, Spain

Alma
Harmony

Alma o



Scar

ClearSkin Pro
ClearLift Pro
SupErb

Before

Post 3 Tx



Treatment Protocol

ClearLift Pro 5x5: -2, -1 and 0 | 1200-1500mJ | Total: 200J

ClearSkin Pro 7x7: 600-800J | 2Hz

SupErb 7x7: 2000-2500 mJ/pulse | PW: 1.5 ms

Sessions intervals: 4 weeks

Courtesy of Dr. Pablo Naranjo, Spain

Alma
Harmony

Alma o



Pigmentation

Iris Dye-VL

Before


After 2 Tx



Treatment Protocol

Iris Dye-VL: In-motion and stationary | 12msec | 10 fluence

Courtesy of Dr. Eyal Kramer and In house clinic

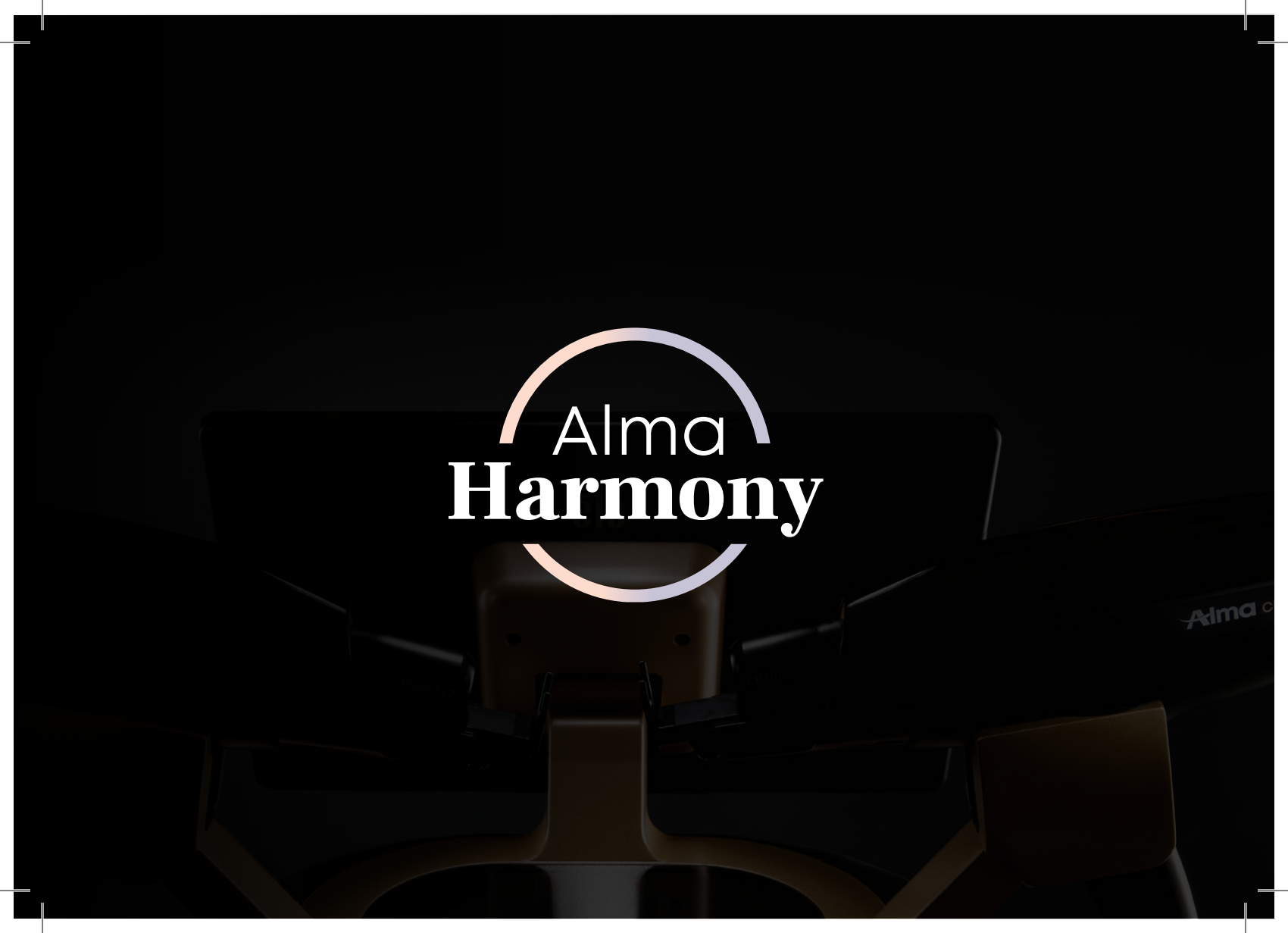

Alma **Harmony**

Alma o



Pigmentation

ClearSkin Pro
Iris Dye-VL
SupErb

Before



After 3 Tx



Treatment Protocol


ClearSkin Pro 7x7: Stationary | 2000-2500J/cm² | 3Hz | Total: 1500J

Iris Dye-VL: Stationary mode | 8-10J | PW: 15 ms | 2 passes

SupErb Pixel 7x7: 2000-2500 mJ/pulse | PW: 1 ms

Sessions intervals: 4 weeks

Courtesy of Dr. Pablo Naranjo, Spain

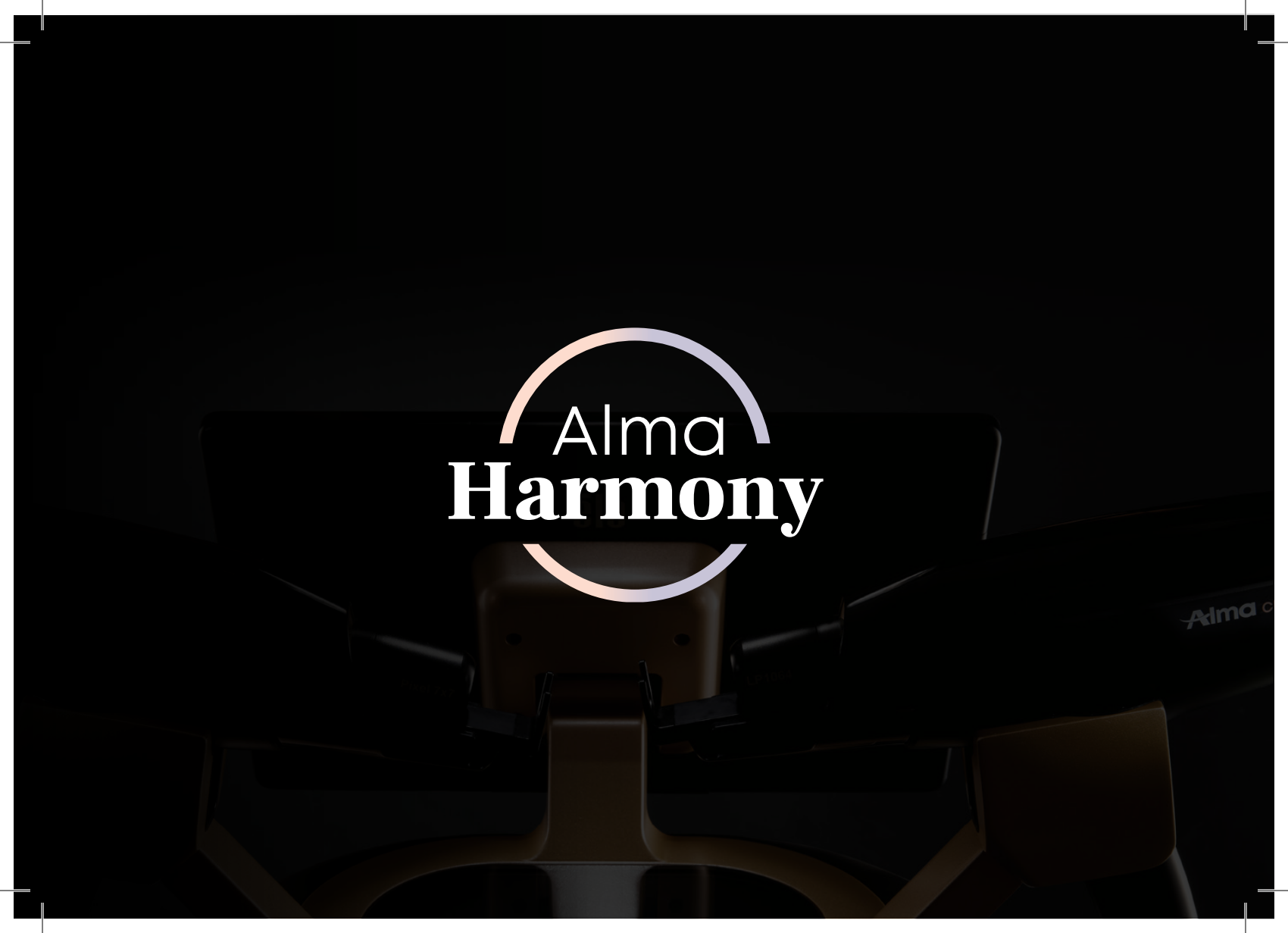

Alma **Harmony**

Alma o



Pigmentations

ClearLift Pro 532nm

Before



After 3 Tx



Treatment Protocol

ClearLift Pro 532nm: 4mm Spot | 400-800mJ/pulse | 1 pass

Courtesy of Dr. Eyal Kramer and In house Clinic

Alma
Harmony

Alma o



Pigmentation

ClearSkin Pro
Iris Dye-VL
SupErb

Before



After 3 Tx



Treatment Protocol

ClearSkin Pro 7x7: Stationary | 1200-1500J/cm² | 2Hz | Total: 800J

Iris Dye-VL: Stationary | 8-10J | PW: 15 ms | 2 passes

SupErb 7x7: 1500-2000 mJ/pulse | PW: 0.5 ms

Sessions intervals: 4 weeks

Courtesy of Dr. Pablo Naranjo, Spain

Alma
Harmony

Alma o



Pigmentation and Rejuvenation

ClearLift Pro
ClearSkin Pro
Iris Dye-VL
SupErb

Before

After 3 Tx



Treatment Protocol


Tx 1 - ClearLift Pro 532nm: 4mm Spot, 400-800mJ/pulse

Tx2 and Tx3 ClearSkin Pro 7x7: 900-1500 mJ/pulse | 3Hz In-motion |
total energy 1500J / face

Iris Dye-VL: In-motion 3J/cm² | 3.6kJ total energy | Stationary: 8-10J/pulse

SupErb 7x7: 1500 - 2500mJ/pulse | 1 pass

Courtesy of Dr. Eyal Kramer and In house clinic

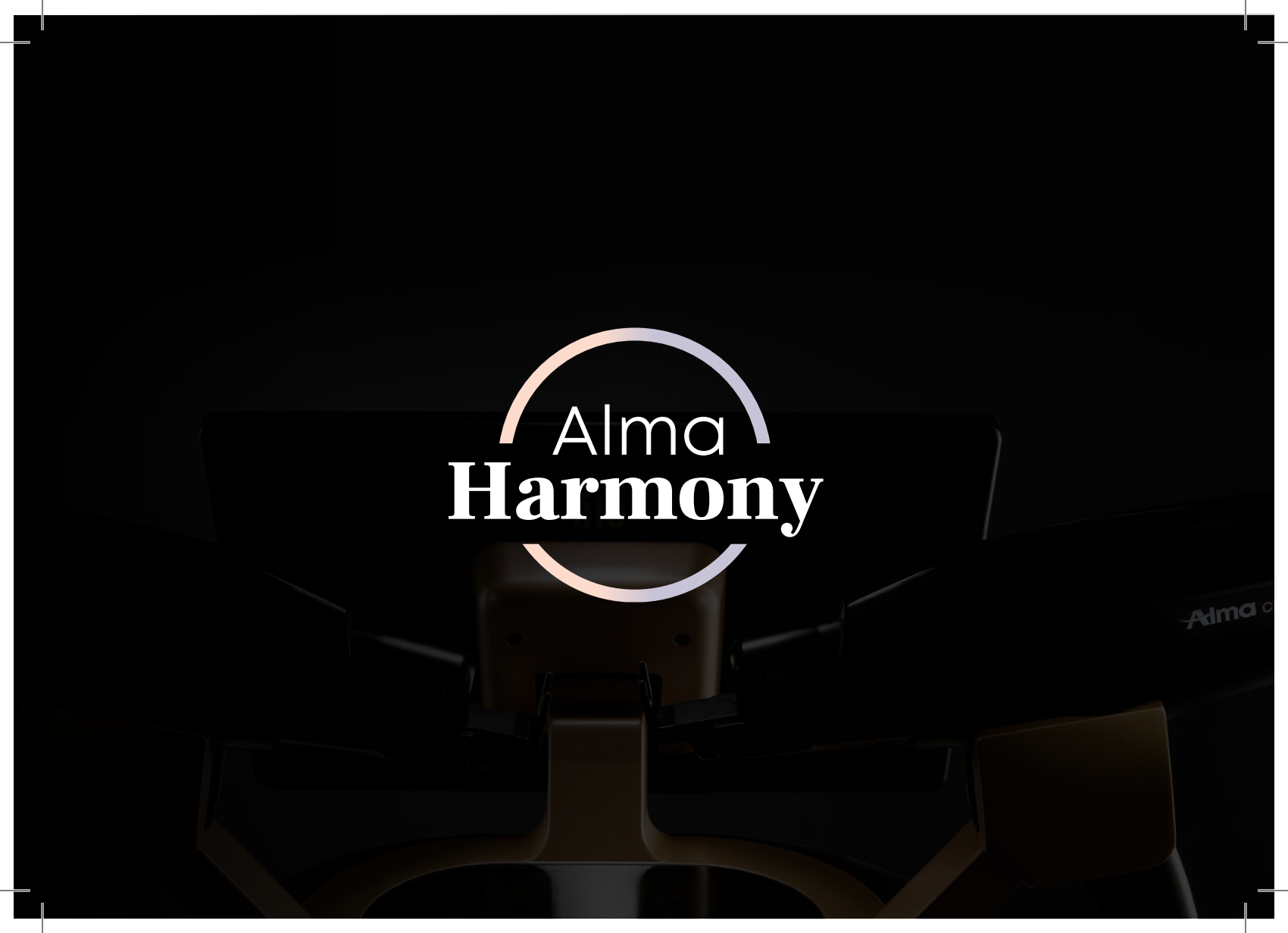

Alma **Harmony**

Alma o



Pigmentation and Rejuvenation

ClearLift Pro
ClearSkin Pro
Iris Dye-VL

Before

1 month post 3 Tx



Before

1 month post 3 Tx



Treatment Protocol

Tx 1 - ClearLift Pro 532nm: 5 mm spot | 800 mJ | 1 pass


ClearSkin Pro 7x7: 3 Hz | 1500 mJ, 2-3 passes, 1285 pulses

Tx 2 - ClearLift Pro 532nm: 8 mm spot | 800 mJ | 4 Hz | 1 pass

ClearSkin Pro 7x7: 3 Hz | 1500 mJ | 31 mJ/Px | 2-3 passes

Tx 3 - Iris Dye-VL: In-motion 3 Hz | 4 J/cm² | 30 sec | 510 pulses | 5.24kj followed by **Iris Dye-VL:** 9 J/cm², 12 ms | 66 pulses | 0.675kj and **ClearSkin Pro 7x7:** (eyes and mouth) | 2 Hz | 2000 mJ

Courtesy of Dr. Jeffrey Hsu, MD, FAAD. USA

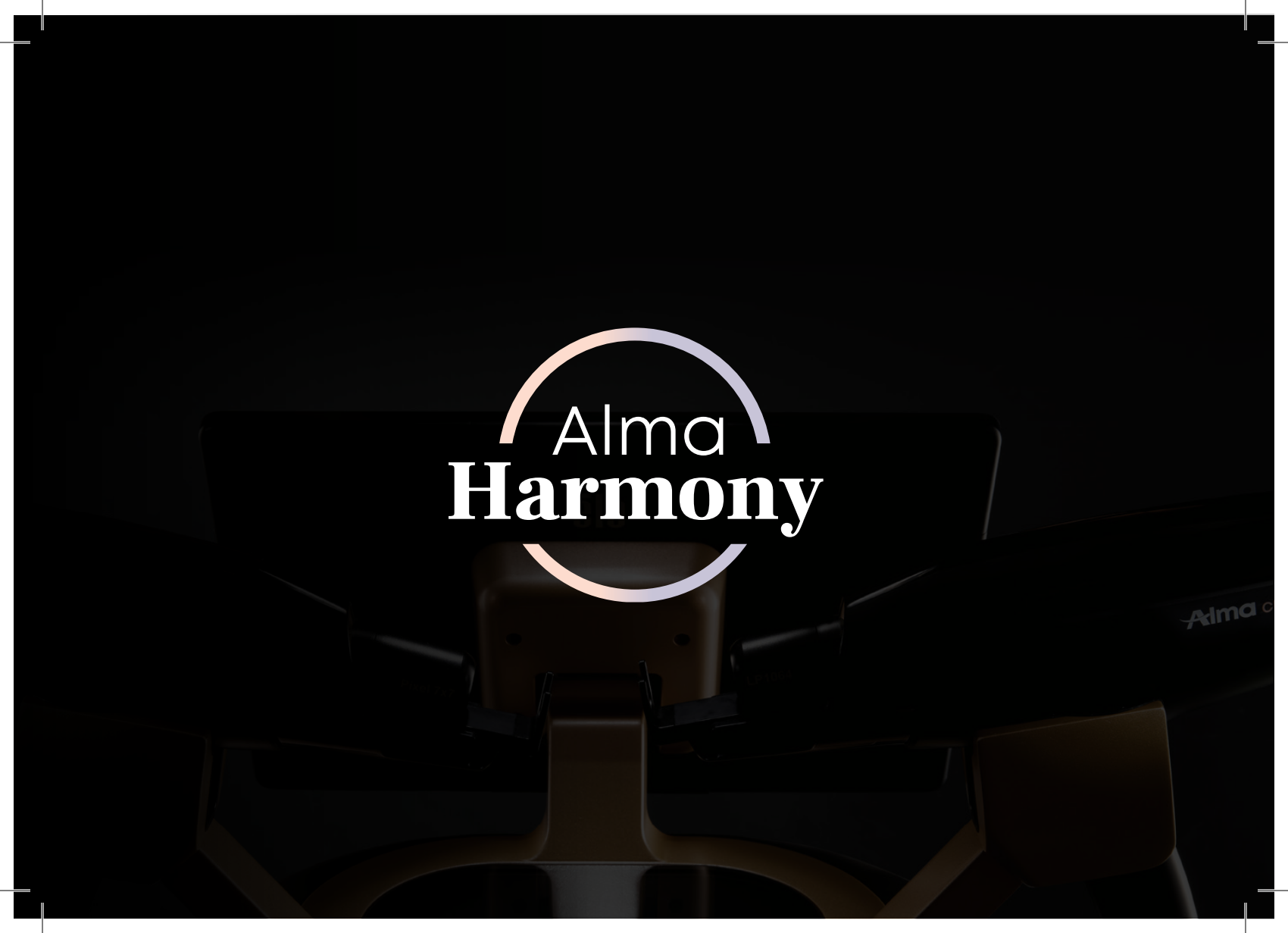

Alma **Harmony**

Alma o



Pigmentation and Rejuvenation

ClearLift Pro
ClearSkin Pro
Iris Dye-VL


Before

1 month post 3 Tx



Before

1 month post 3 Tx

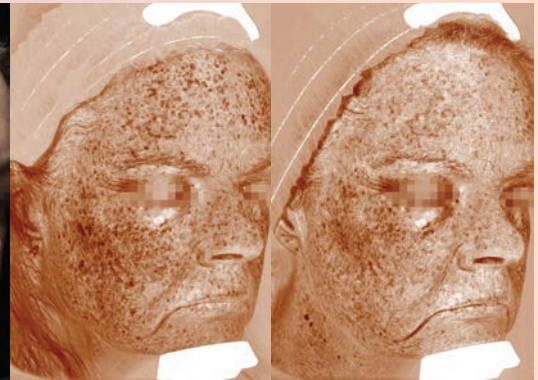

Treatment Protocol

Tx 1 - ClearLift Pro 532nm: 5 mm spot | 800 mJ | 1 pass

ClearSkin Pro 7x7: 3 Hz | 1500 mJ | 2-3 passes | 1285 pulses

Tx 2 - ClearLift Pro 532nm: 8 mm spot, 800 mJ, 4 Hz, 1 pass

ClearSkin Pro 7x7: 3 Hz | 1500 mJ | 31 mJ/Px | 2-3 passes

Tx 3 - Iris Dye-VL: In-motion 3 Hz | 4 J/cm² | 30 sec | 510 pulses | 5.24kJ followed by **Iris Dye-VL:** 9 J/cm² | 12 ms | 66 pulses | 0.675kJ and **ClearSkin Pro 7x7:** (eyes and mouth) | 2 Hz | 2000 mJ

Courtesy of Dr. Jeffrey Hsu, MD, FAAD. USA

Alma
Harmony

Alma o



Melasma

ClearLift Pro

Before

After 3 Tx



Treatment Protocol

ClearLift Pro 5x5: -2, -1, 0 +1 | 1200mJ/pulse | 6-8kJ acc.

Courtesy of Dr. Munir Somji, UK

Alma
Harmony

Alma o



Melasma

ClearLift Pro
ClearSkin Pro
Iris Dye-SR

Before



After 4 Tx



Treatment Protocol

ClearLift Pro 5x5: -2, -1 and 0 | 1500-1800mJ/pulse | 4Hz | Total: 5 KJ

ClearSkin Pro 7x7: In-motion | 1000-1500J | 3Hz | Total: 1500J

Iris Dye-SR: In-motion | 3 -4J | 2 passes | Total: 4 kJ

Sessions intervals: 4 weeks

Courtesy of Dr. Pablo Naranjo, Spain



Alma **Harmony**

The background of the image is a dark, high-contrast photograph of a mechanical device, likely a scale or a precision instrument. The device features various components such as a central vertical column, horizontal arms, and a base. The lighting is dramatic, highlighting the metallic textures and the complex geometry of the machine. The overall aesthetic is professional and technical.

Alma o



Melasma

ClearLift Pro
ClearSkin Pro
Iris Dye-SR

Before



After 3 Tx



Treatment Protocol


ClearLift Pro 5x5: -2, -1 and 0 | 1500-1800mJ/pulse | 4Hz | Total: 5 KJ

ClearSkin Pro 7x7: In-motion | 1000-1500J | 3Hz | Total: 1500J

Iris Dye-SR: In-motion | 3 -4J | 2 passes | Total: 4 kJ

Sessions intervals: 4 weeks

Courtesy of Dr. Pablo Naranjo, Spain

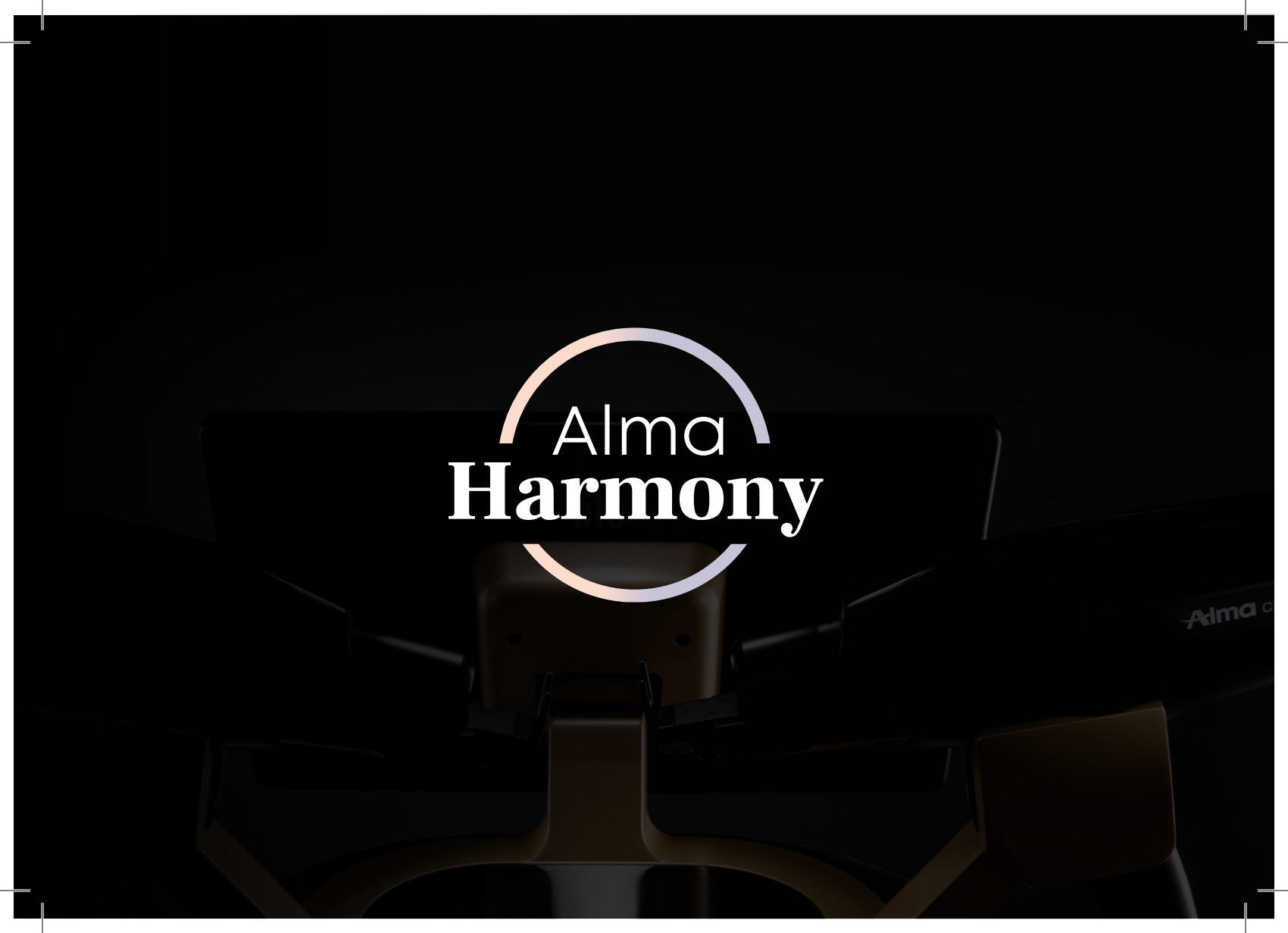

Alma **Harmony**

Alma o



Melasma Skin Type V

ClearLift Pro

Before

1 month post 6 Tx




Treatment Protocol

ClearLift Pro QS 1064nm:

10mm | 1200mJ | 6 Hz, 450-1200 pulses | 4-6 passes

Courtesy of Dr. Jeffrey Hsu, MD, FAAD. USA

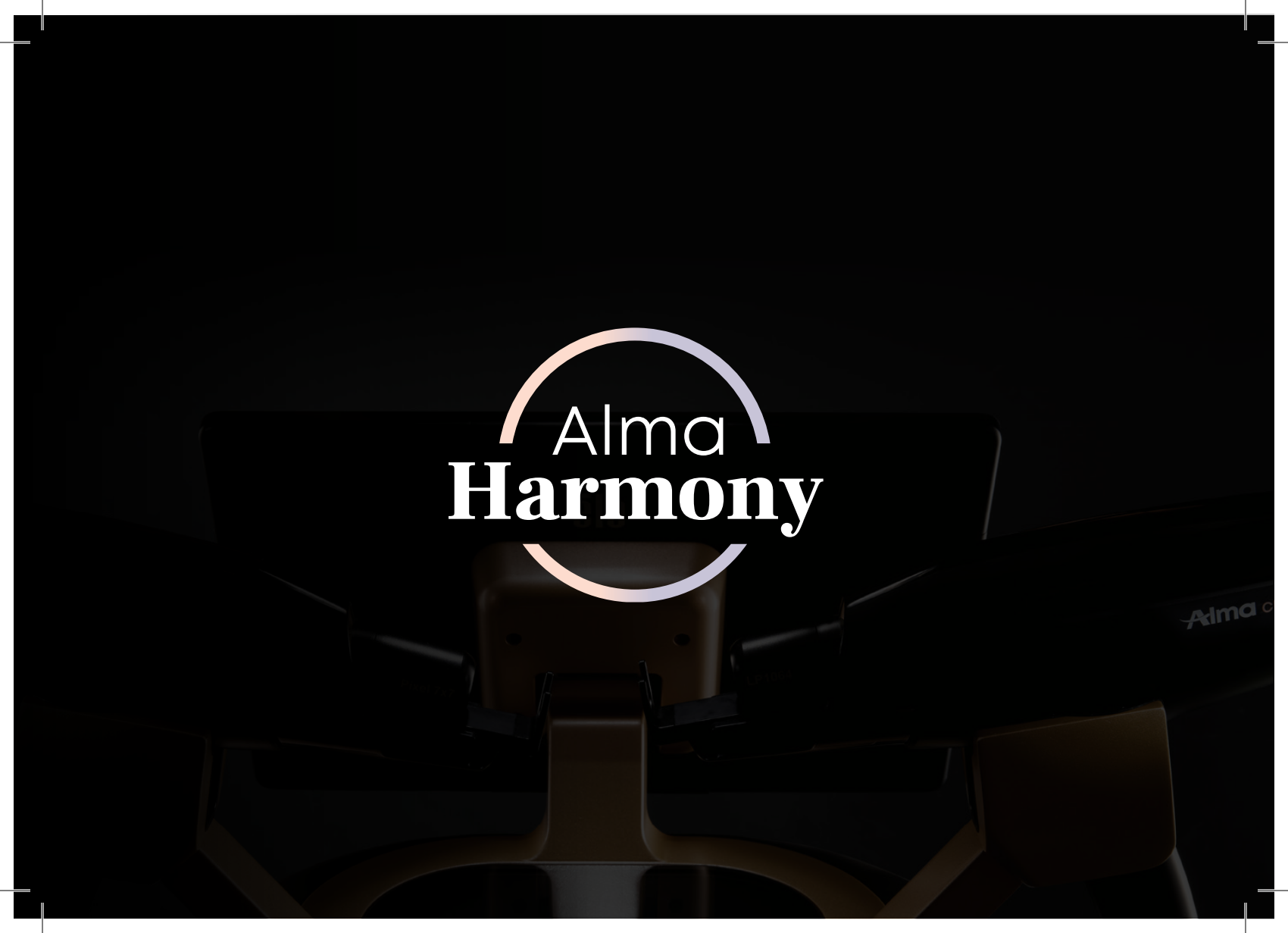

Alma **Harmony**

Alma o



Pigmentation, Texture & Laxity

ClearLift Pro
ClearSkin Pro
Iris Dye VL/SR

Before



1 month post 2 Tx




Treatment protocol

ClearSkin Pro 7x7: 900-1500 mJ/pulse | 3Hz, In-motion, total energy 1500J/face

ClearLift Pro 5x5: -2, -1, 0 | 1800mJ/pulse | 5kJ per face

Iris Dye VL/SR: In-motion | 3J/cm² | 3.6kJ total energy.

Courtesy of Dr. Eyal Kramer and In house clinic

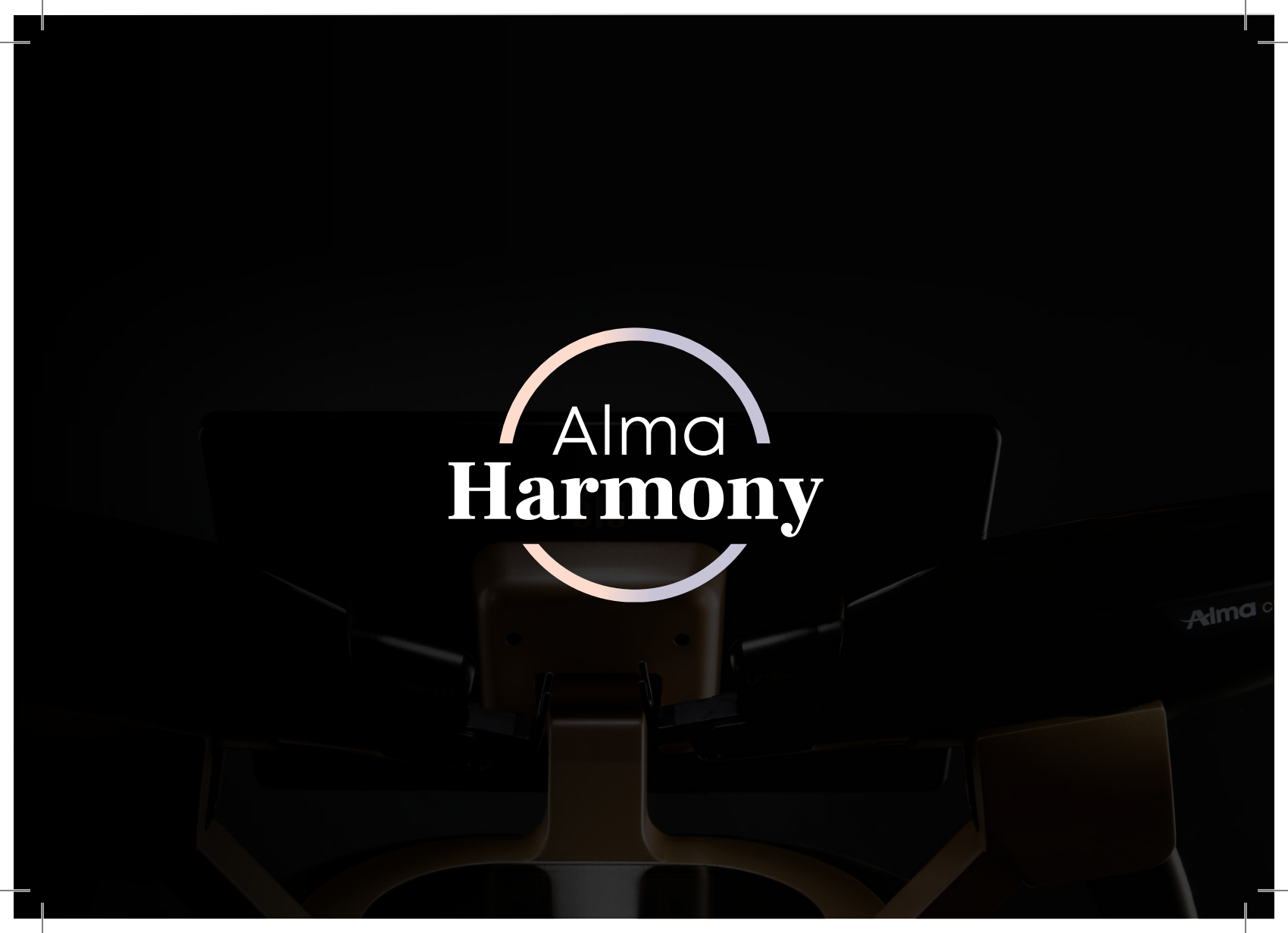

Alma **Harmony**

Alma o



Pigmentation, Texture & Angioma

ClearSkin Pro
ClearLift Pro
SupErb
VascuPen

Before



Post 4 Tx



Treatment protocol

Tx1 - 2 and 4 ClearSkin Pro: 300mJ/pulse | 2 passes.

ClearLift Pro 5x5: 1200mJ/pulse 2 passes -2 | 1200 mJ/pulse 2 passes (-1) | 600 mJ/pulse 1 pass (0)

SupErb 7x7: 1500mJ/pulse 1.5msec | 3 Hz | 1 pass

Tx 3 - VascuPen: 1.0 mm 30 fluence | 21 pw

Courtesy of Dr. Munir Somji, UK

Alma
Harmony

Alma o



Pigmentation & Texture

ClearSkin Pro
ClearLift Pro
Iris Dye VL/SR

Before

1 month post 2 Tx

Before

1 month post 2 Tx



Treatment Protocol

ClearSkin Pro 7x7: 900-1500 mJ/pulse | 3Hz In-motion | total energy 1500J/face

ClearLift Pro QS 7x7: -2, -1, 0 | 1800mJ/pulse | 5kJ per face

Iris Dye VL/SR: In-motion | 3J/cm² | 3.6kJ total energy.

Courtesy of Dr. Eyal Kramer and In house clinic

Alma
Harmony

Alma o



Pigmentation and Texture

ClearSkin Pro
Iris Dye-VL
SupErb

Before



1 month post 2 Tx




Treatment Protocol

ClearSkin Pro 7x7: 900-1500 mJ/pulse | 3Hz, In-motion | total energy 1500J/face

Iris Dye-VL: In-motion | 3J/cm² | 3.6kJ total energy. Stationary: 8-10J/pulse

SupErb 7x7: 1500- 2500mj/pulse | 1 pass

Courtesy of Dr. Eyal Kramer and In house clinic

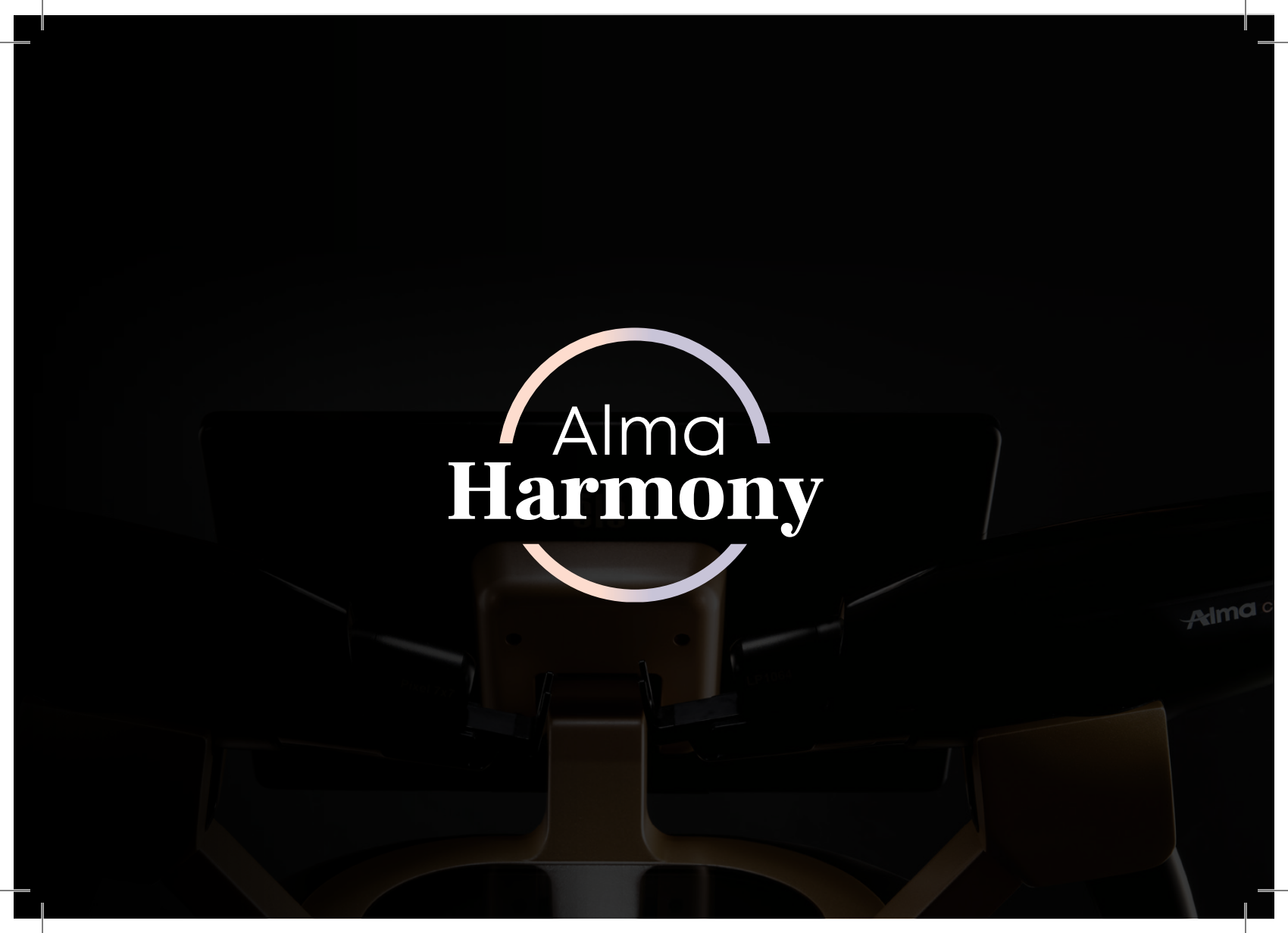

Alma **Harmony**

Alma o



Pigmentation and Texture

ClearSkin Pro
ClearLift Pro
Iris Dye VL/SR

Before

1 month post 2 Tx



Treatment Protocol

ClearSkin Pro 7x7: 900-1500 mJ/pulse | 3Hz, In-motion | total energy 1500J/face

ClearLift Pro 5x5: -2, -1, 0 | 1800mJ/pulse | 5kJ per face

Iris Dye VL/SR: In-motion | 3J/cm² | 3.6kJ total energy

Courtesy of Dr. Eyal Kramer and In house clinic



Alma **Harmony**

Alma o



Pigmented Decollete

ClearLift Pro

Before

Post 2 Tx



Treatment Protocol

ClearLift Pro QS 532nm

10 mm | 800mJ/pulse | 4 Hz | accumulate 400-600 J

Courtesy of Dr. Eyal Kramer and In house clinic



Alma **Harmony**

Alma o



Pigmented Decollete

ClearLift Pro

Before



Post 2 Tx




Treatment Protocol

ClearLift Pro QS 532nm

10mm | 400mJ/pulse | 5 Hz | 4-6 passes | 3313 pulses / 600cm² grid

Courtesy of Dr. Hen Ifrach and In house Clinic

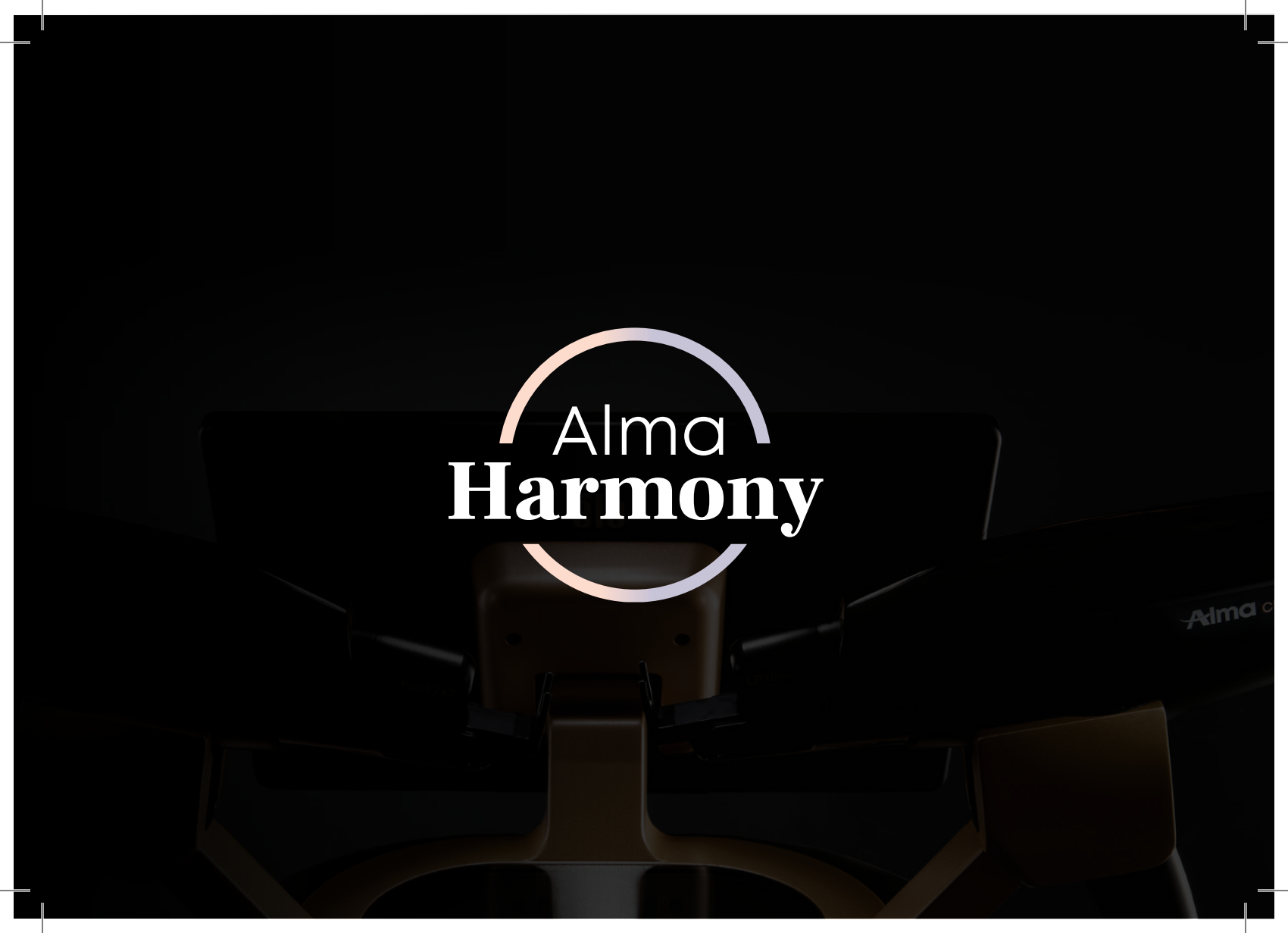

Alma **Harmony**

Alma o



Photodamaged Skin

Iris Dye-VL

Before

Post 2 Tx



Treatment Protocol

Iris Dye-VL Stationary

9J/cm² | 12msec | 1 pass

11J/cm² | 12 msec | 1 pass

Courtesy of Dr. Pablo Naranjo and In house clinic



Alma **Harmony**

Alma o



Solar Lentigo

ClearLift Pro
ClearSkin Pro

Before

1 month after 2Tx



Treatment Protocol

ClearLift Pro QS 532nm: 3 mm | 400 mJ/P | 5.66 J/cm² | 2-4 Hz | 350 Pulses | 1 Pass | Total Energy 0.14 kJ (face)

ClearSkin Pro 7x7: 1,500 mJ/P | 3 Hz | 1,285/1,288 Pulses | 2-3 Passes | Total Energy 1.962/1.971 kJ (face)

Courtesy of Dr. Jeffrey Hsu, USA



Alma **Harmony**

Alma o



Pigmented Leg Scar

SupErb
ClearLift Pro
Iris Dye-VL

Before



After 4 Tx




Treatment Protocol

SupErb 7x7: 0.5msec | 400mJ | one pass

ClearLift Pro: 5mm | 400mJ | one pass

2nd Tx Iris Dye-VL: stationary | 7J/cm² | 15msec

Courtesy of Dr. Eyal Kramer and In house clinic

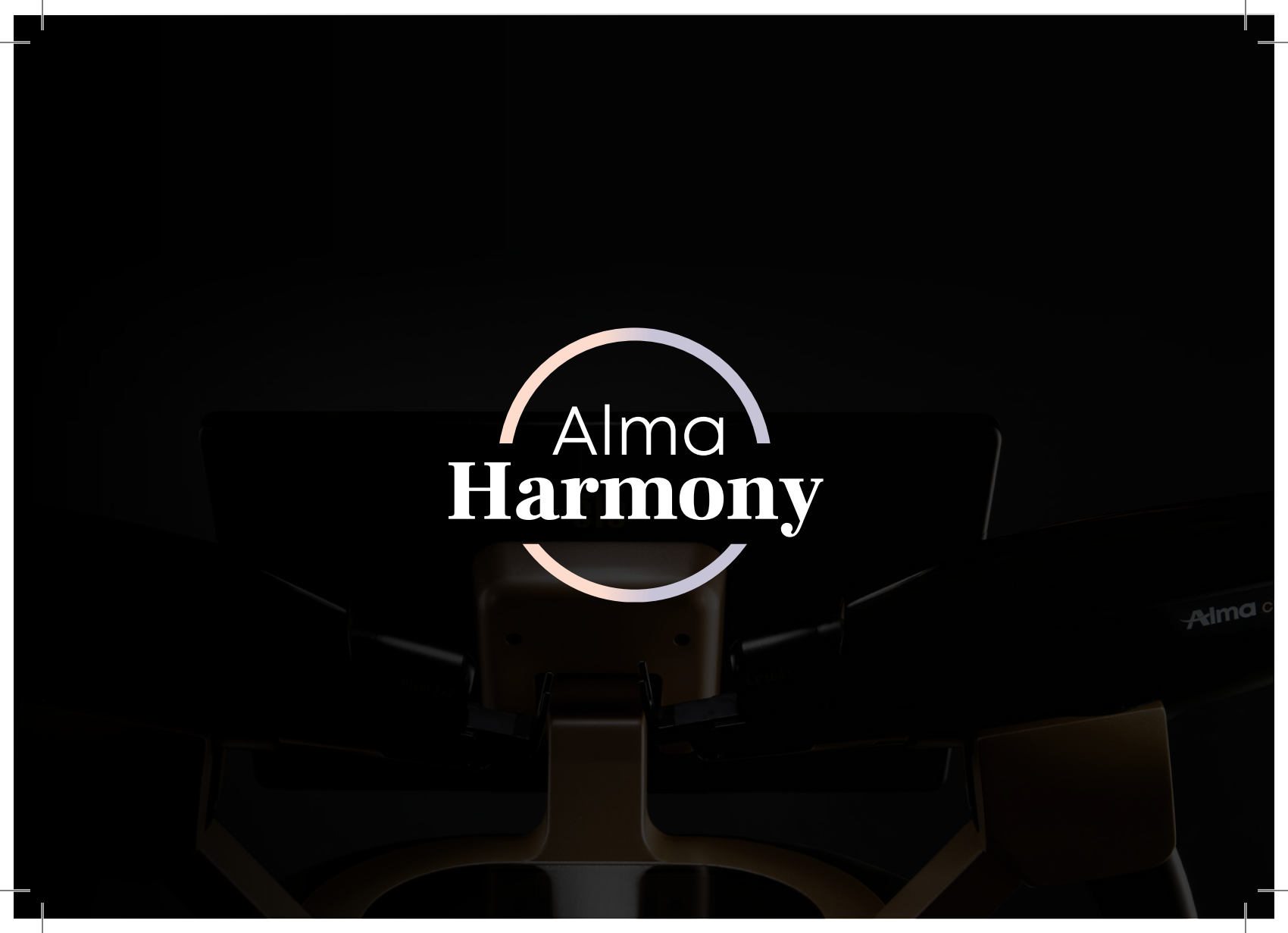

Alma **Harmony**

Alma o



Seborrheic keratosis

SupErb

Before

post 2 Tx



Treatment Protocol

SupErb CVD: 5/6mm | 2000-2500 mJ/pulse | PW 1.5msec

Courtesy of Dr. Eyal Kramer and In house clinic

Alma
Harmony

Alma o



Seborrheic keratosis

SupErb

Before

2 weeks after 1 Tx



Treatment Protocol

SupErb CVD: 4-6mm | 1000 mJ/pulse | PW | 1 msec | 2 Hz

Courtesy of Dr. Eyal Kramer and In house clinic

Alma
Harmony

Alma o



Active Acne

ClearLift Pro
ClearSkin Pro
Iris Acne

Before



After 3 Tx



Treatment Protocol


ClearLift Pro 5x5: -2, -1 and 0 | 1500-2000mJ | Total: 6000J

ClearSkin Pro: 8-10mm with vacuum | 400-600mJ | 2 stacks; 2 passes | Pixel 7x7
In-motion mode | 1500-2000J

Iris Acne: Stationary mode | 8-10J | 40ms

Sessions intervals: 4 weeks

Courtesy of Dr. Pablo Naranjo, Madrid

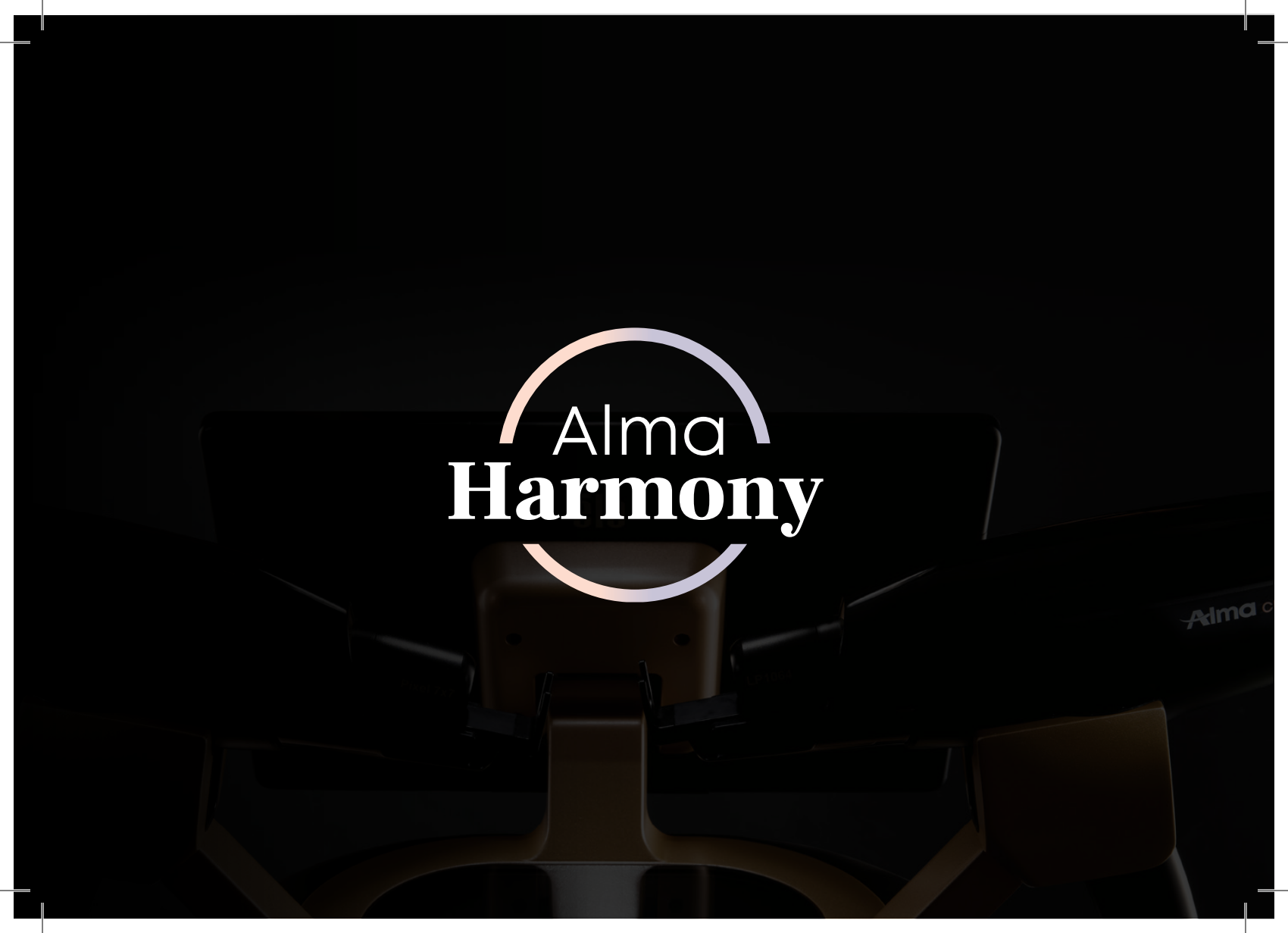

Alma **Harmony**

Alma o



Active Acne

ClearLift Pro
ClearSkin Pro
Iris Acne

Before



1 month after 3 Tx



Treatment Protocol


ClearLift Pro 5x5: -2, -1 and 0 | 1500-2000mJ | Total: 6000J

ClearSkin: 1. CVD tip in papules | 6-10mm with vacuum | 500-700mJ | 2 stacks. 2. Pixel 7x7 | 1500-2000J with 2Hz | 2 passes

Iris Acne: Stationary mode, 10-12J | PW: 40ms

Sessions intervals: 4 weeks

Courtesy of Dr. Pablo Naranjo, Spain

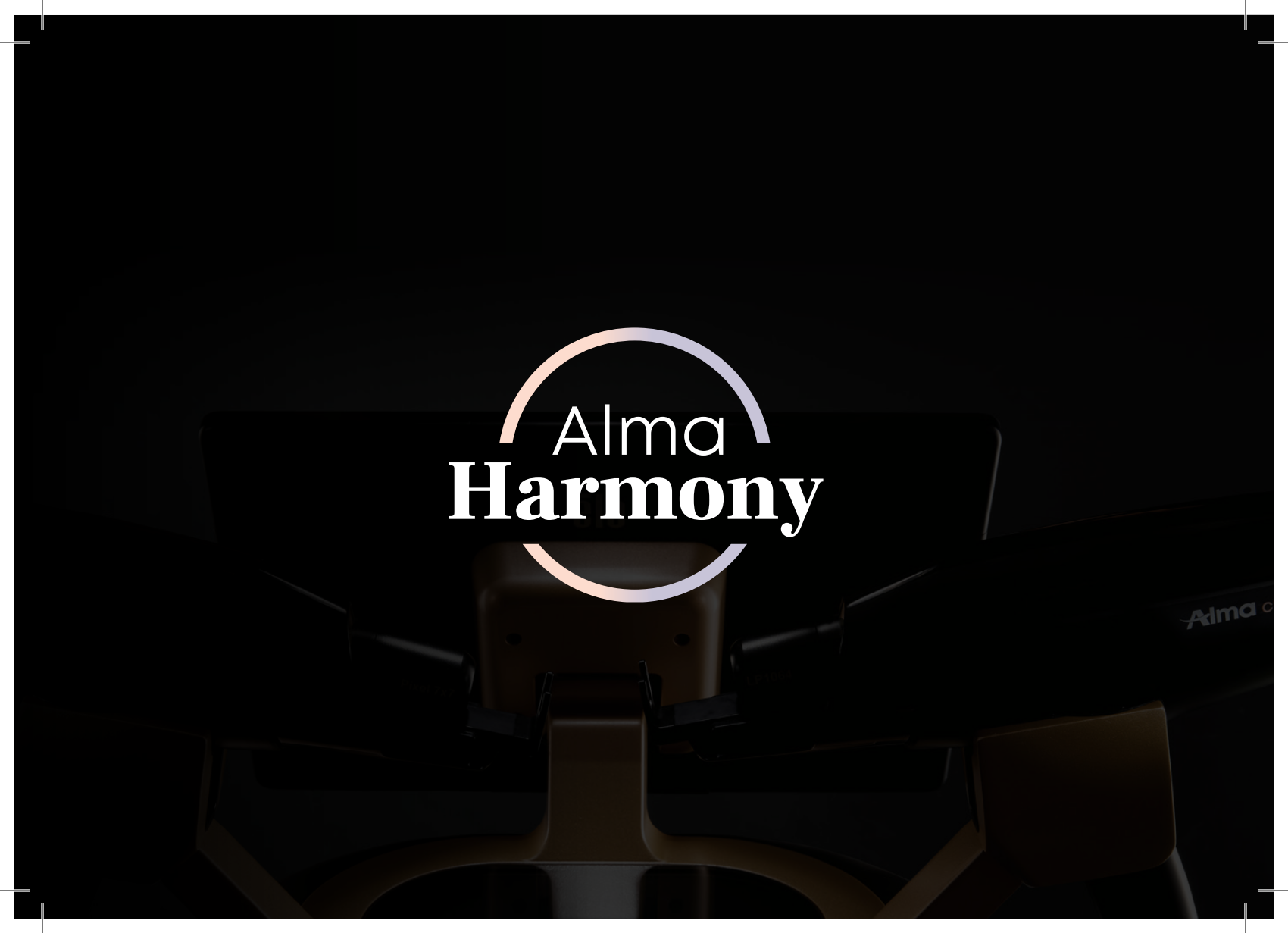

Alma **Harmony**

Alma o



Active Acne

ClearLift Pro
ClearSkin Pro
Iris Acne

Before



After 3 Tx



Treatment Protocol

ClearLift Pro: -2, -1 and 0 | 1500-2000mJ | Total: 6000J

ClearSkin Pro: 1. Full beam in papules | 6-10mm with vacuum

500-700mJ | 2 stacks 2. Pixel 7x7 | 1000-1500J | 2Hz | 2 passes | Total: 1500J (In-motion).

Iris Acne: Stationary mode | 10-12J | PW: 40ms

Sessions intervals: 4 weeks

Courtesy of Dr. Pablo Naranjo, Spain

Alma
Harmony

Alma o



Active Acne

ClearLift Pro
ClearSkin Pro
Iris Acne

Before



After 3 Tx



Treatment Protocol


ClearLift Pro 5x5: -2, -1 and 0 | 1500-2000mJ | Total: 6000J

ClearSkin Pro: 1. CVD Tip in papules | 6-10mm with vacuum, 500-700mJ | 2 stacks. 2. Pixel 7x7 | 1500-2000J with 2Hz | 2 passes

Iris Acne: Stationary mode | 10-12J | PW: 40ms

Sessions intervals: 4 weeks

Courtesy of Dr. Pablo Naranjo, Spain

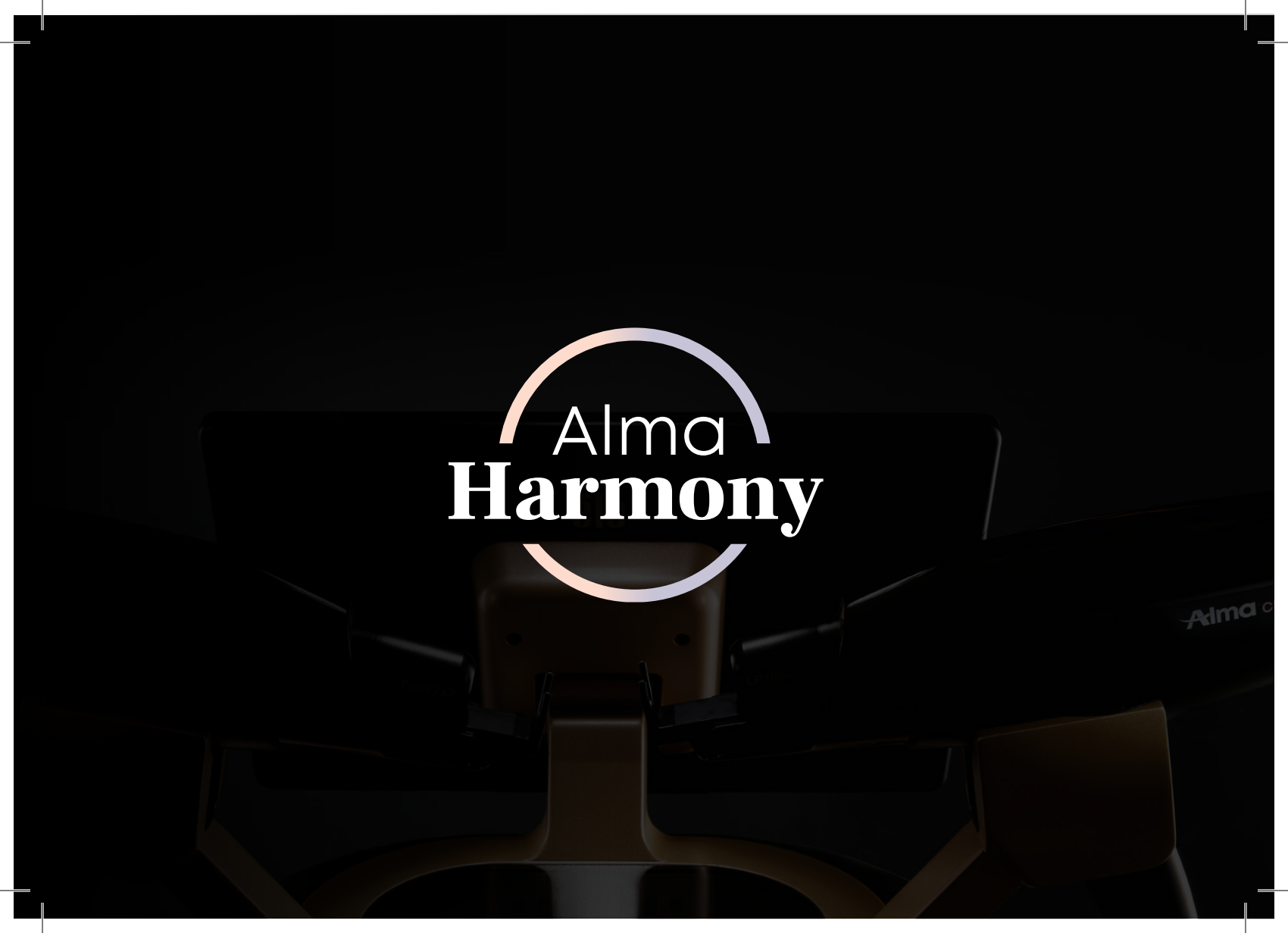

Alma **Harmony**



Acne and Acne Scars

ClearSkin Pro
SupErb

Before

After 4 Tx



Treatment Protocol

ClearSkin Pro 7x7: 1500mJ | 2 Stacks | 2 passes

ClearSkin Pro: 5mm 400mJ | 1 pulse

SupErb 7x7 (Tx#4): 1.5 msec | 1200mJ | one pass

Sessions intervals: 4 weeks

Courtesy of Dr. Chen Ifrach and In House Clinic

Alma
Harmony

Alma o



Acne and Acne Scars

ClearSkin Pro
ClearLift
Iris SHR

Before



After 4 Tx




Treatment Protocol

1st & 2nd Tx - ClearSkin Pro 7x7: 600-800mJ/P | 3Hz | 2 passes. **Iris SHR:** 6J/cm² | 10KJ

3rd Tx ClearSkin Pro 7x7: 800-1000mJ/P | 3Hz | 3 passes. **ClearLift Pro 5x5:** 1200mJ/p | 6Hz | depth 0 | 2 passes **Iris SHR:** 7J/cm² | 10KJ

4th Tx - ClearSkin Pro 7x7: 800-1000mJ/P | 3Hz | 3 passes. **ClearLift Pro 5x5:** 1500mJ/p | 5Hz | depths 0, +1 | 2 passes **Iris SHR:** 7J/cm² | 10KJ

Courtesy of Dr. Hen Ifrach, Italy

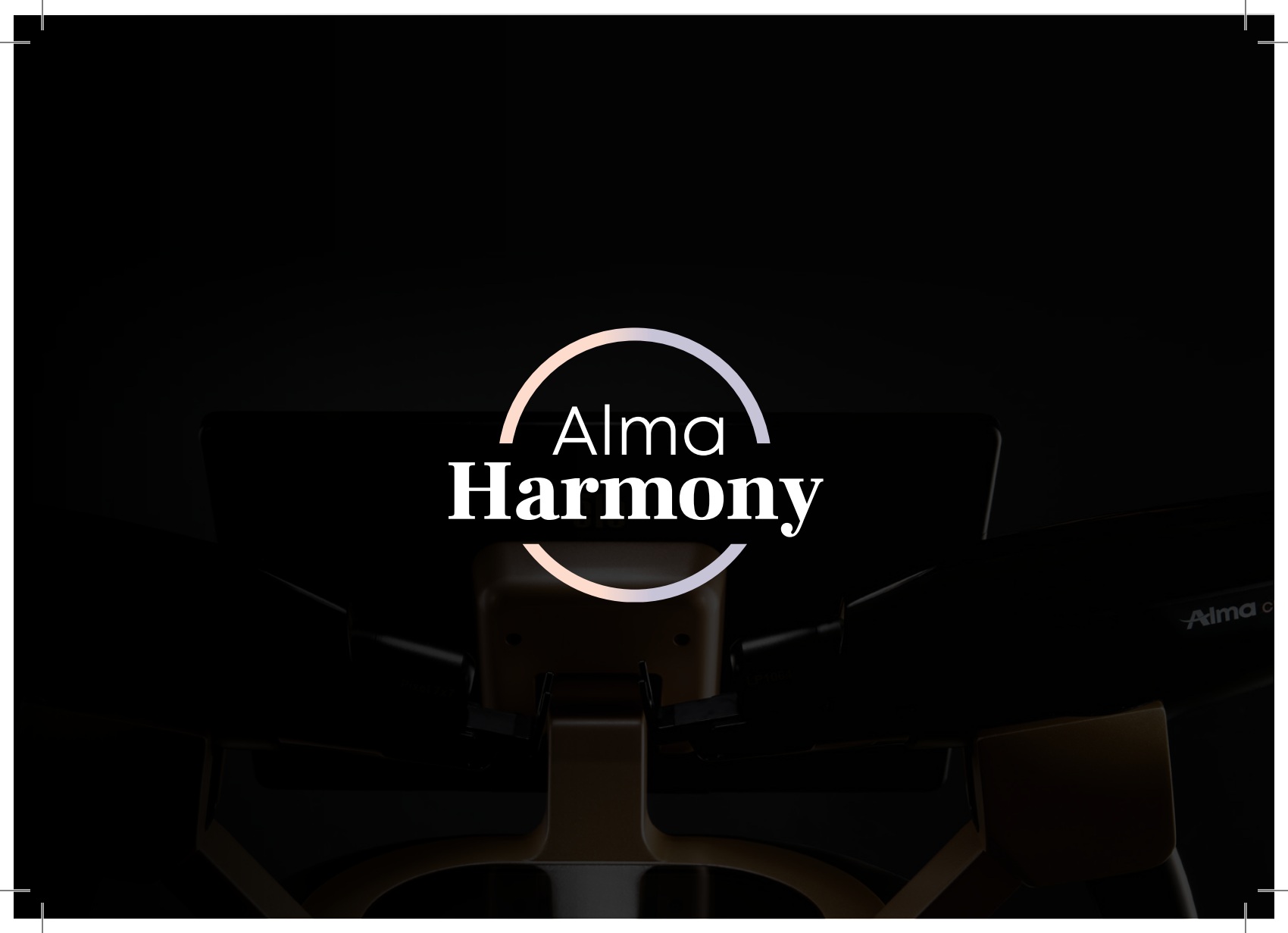

Alma **Harmony**

Alma o



Acne Scars

ClearLift Pro
ClearSkin Pro
SupErb

Before



After 3 Tx



Treatment Protocol

ClearLift Pro 5x5: -2, -1 and 0 | 1500-2000mJ | Total: 6000J

ClearSkin Pro 7x7: 1500-2000J | 2Hz

SupErb: 2500-3500 mJ/pulse | PW: 1.5 ms

Sessions intervals: 4 weeks

Courtesy of Dr. Pablo Naranjo, Spain

Alma
Harmony

Alma o



Acne & Micro Cysts

Iris Acne
ClearSkin Pro

Before



After 4 Tx




Treatment Protocol

Iris Acne: 8J/cm² | 40 msec | 2 passes

ClearSkin Pro 7x7: 800-1100mJ/pulse | 2 passes

Courtesy of Dr. Hen Ifrach, Italy

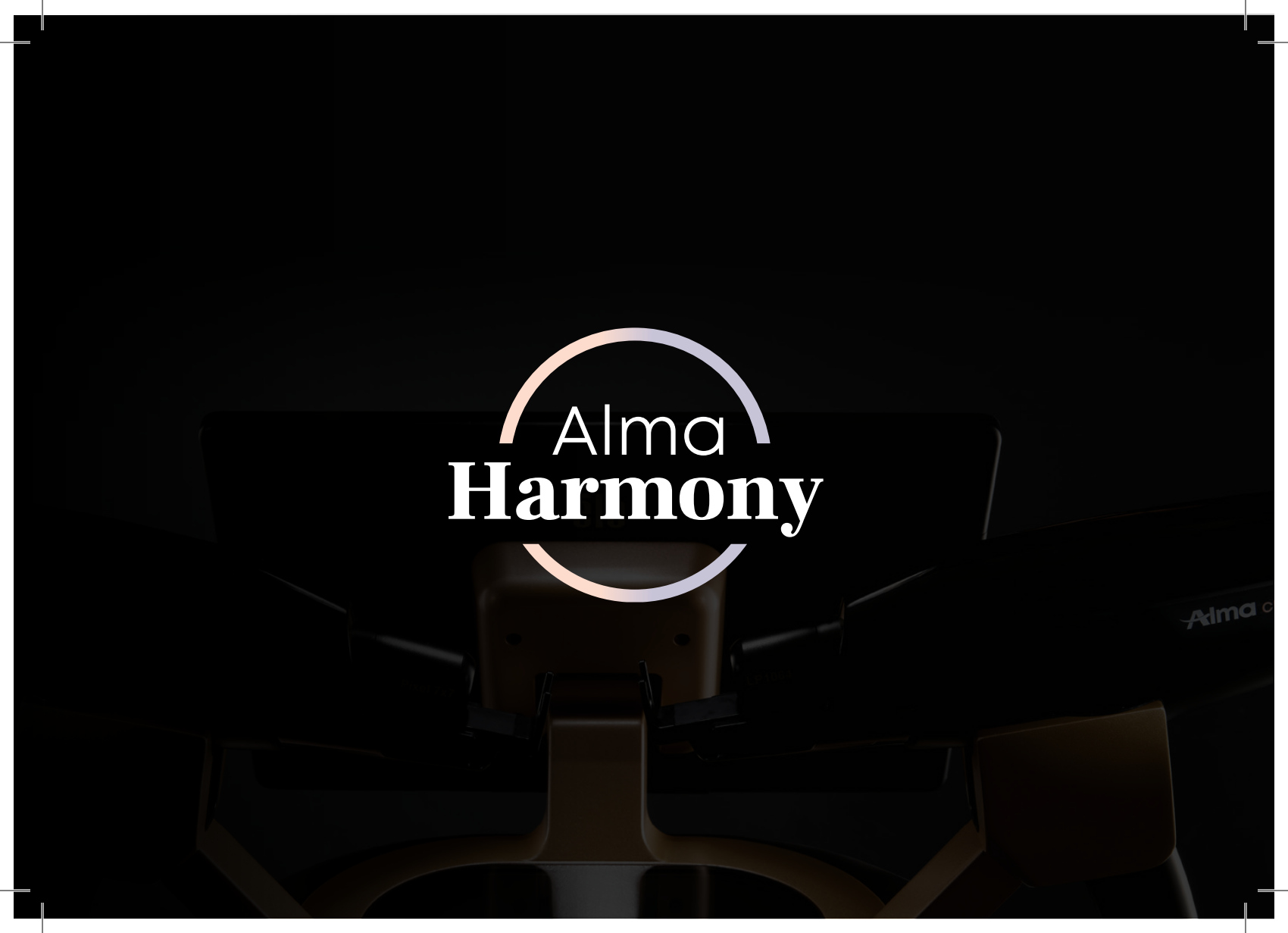

Alma **Harmony**

Alma o



Vascular

ClearVas
Iris Dye-VL
VascuPen

Before



Post 3 Tx



Treatment Protocol


ClearVas: 8mm Spot | 6-8J | 10Hz | 6 passes | Total: 12KJ

Iris Dye-VL: In-motion mode | 2 passes 6-8J

Vascupen: 2 mm Spot (telangiectasia) | 16-20J | PW: 30ms

Sessions intervals: 4 weeks

Courtesy of Dr. Pablo Naranjo, Spain

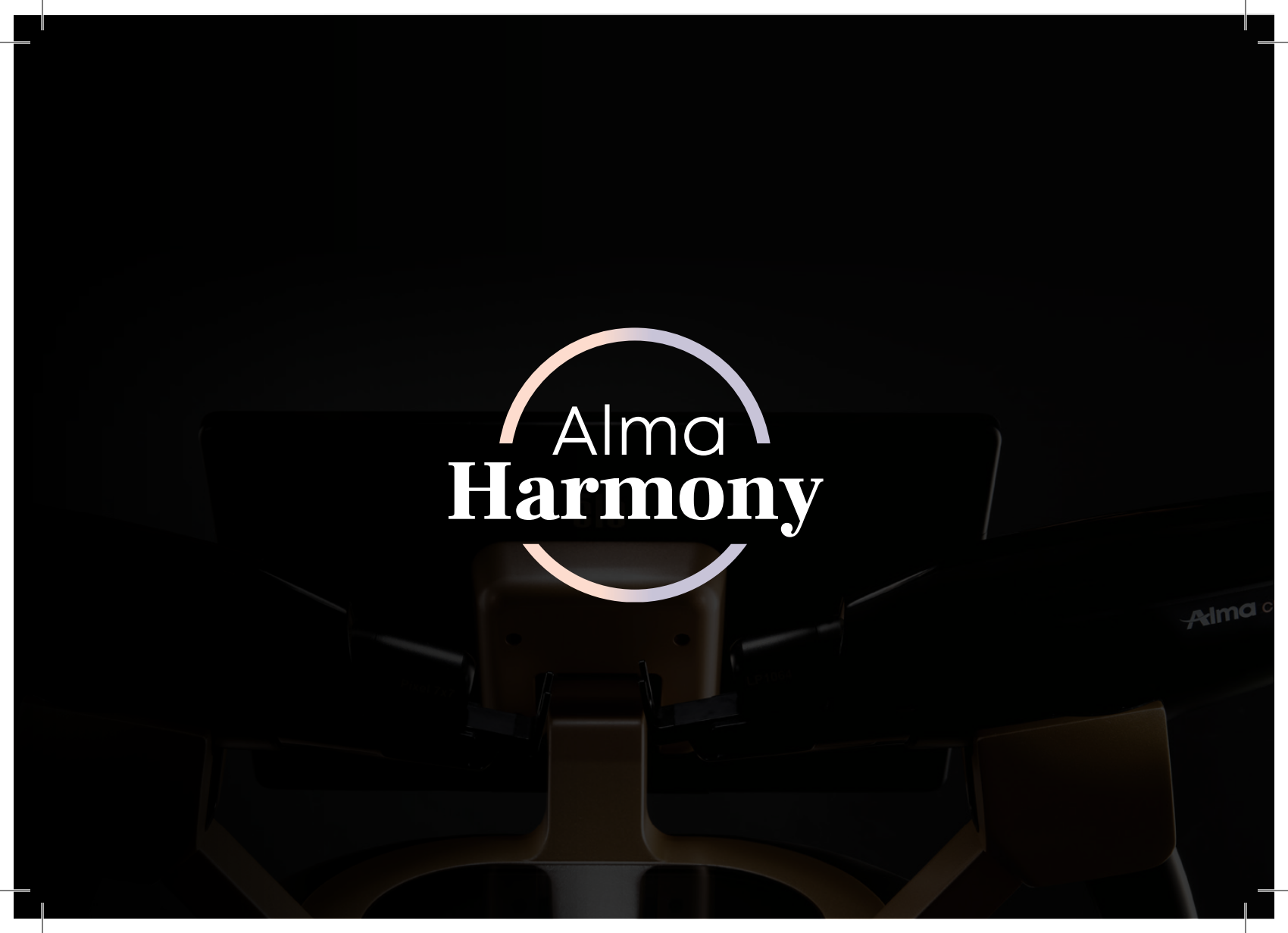A dark, close-up photograph of a mechanical keyboard, likely a Cherry MX series, with a semi-transparent logo overlay. The logo consists of a circular arc with a gradient from orange to purple, framing the text "Alma Harmony". The keyboard's keys and mechanical components are visible in the background, with some keys having a light-colored top. The overall lighting is low, emphasizing the keyboard's texture and the logo's glow.

Alma **Harmony**

Alma o



Vascular

ClearVas
Iris Dye-VL
VascuPen

Before



Post 3 Tx



Treatment Protocol

ClearVas: 8mm Spot | 6-8J | 10Hz | 6 passes | Total: 12KJ

Iris Dye-VL: Stationary mode | 2 passes | 8-10J | PW: 12ms

Vascupen: 2 mm Spot (papule) | 16-20J | PW: 30ms

Sessions intervals: 4 weeks

Courtesy of Dr. Pablo Naranjo, Spain



Alma **Harmony**

The background of the image is a dark, high-contrast photograph of a mechanical device, likely a scale or a precision instrument. The device has a complex, symmetrical structure with various components, including what appears to be a central weighing pan or platform. The lighting is dramatic, highlighting the metallic surfaces and the intricate details of the machinery. The overall aesthetic is professional and technical.

Alma o



Vascular

ClearVas
Iris Dye-VL
VascuPen

Before

Post 3 Tx



Treatment Protocol

ClearVas: 8mm Spot | 6-8J | 10Hz | 6 passes | Total: 12KJ

Iris Dye-VL: In-motion mode (2 passes) | 6-8J

VascuPen: 2 mm Spot (spider angioma) | 16-20J | PW: 30ms

Sessions intervals: 4 weeks

Courtesy of Dr. Pablo Naranjo, Spain

Alma
Harmony

Alma o



Facial Telangiectasia

VascuPen

Before

4 months post 1 Tx



Treatment Protocol

VascuPen 2mm:

16J/cm² | 37msec | 1Hz | Total energy 119mJ

Courtesy of Dr. Jeffrey Hsu, MD, FAAD. USA

Alma
Harmony

Alma o




Leg Veins

ClearVas

Before



Post 8 Tx

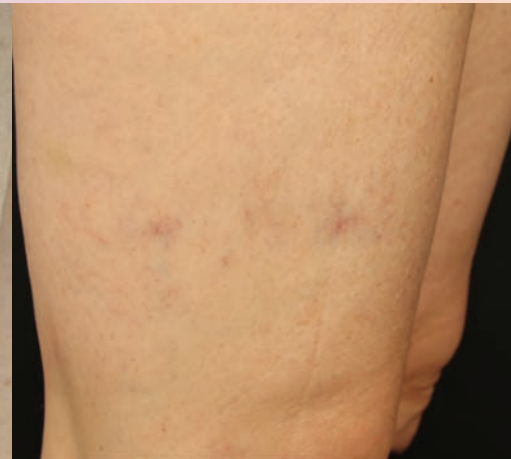

Treatment Protocol

ClearVas: CVD 6mm | 15/45 msec | 110/120 J/cm²

Treatment Intervals: 6 weeks

Courtesy of Dr. Eyal Kramer and In house clinic



Alma **Harmony**

Alma o



Pigmentation and Vascular

ClearLift Pro
ClearSkin Pro
Iris Dye-VL

Before

1 month post 3 Tx

Before

1 month post 3 Tx



Treatment Protocol

Tx 1 - ClearLift Pro 532nm: 5 mm | 800 mJ | 1 pass. **ClearSkin 7x7:** 3 Hz | 1500 mJ | 2-3 passes | 1285 pulses

Tx 2 - ClearLift Pro 532nm: 8 m | 800 mJ | 4 Hz | 1 pass
ClearSkin 7x7: 3 Hz | 1500 mJ | 31 mJ/Px | 2-3 passes.

Tx 3 - Iris Dye SVL: 3 Hz | 4 J/cm² | 30 sec | 510 pulses | 5.24kJ. **Iris Dye-VL:** 9 J/cm² | 12 ms | 66 pulses | 0.675kJ. **ClearSkin Pro 7x7:** (eyes and mouth) 2 Hz | 2000 mJ

Courtesy of Dr. Jeffrey Hsu, MD, FAAD. USA

Alma
Harmony

Alma o




Pigmentation, Rosacea, Vascular lesions

Iris VL/PL
Iris Dye-VL
ClearLift Pro

Before

After 2 Tx

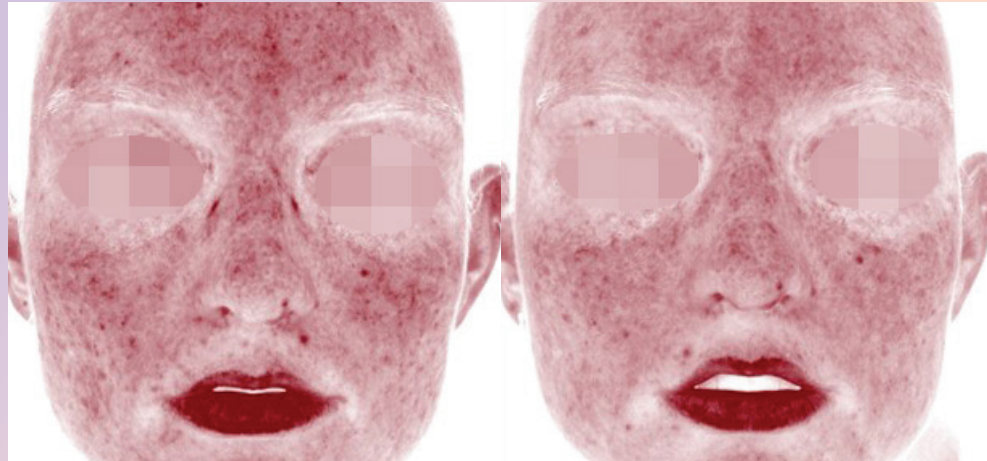

Treatment Protocol

Tx1 - Iris VL/PL: 14J/cm² | 12 msec | 2 passes. **ClearLift Pro 5x5:** 0
2400 mJ/pulse | 3 Hz | 5 kJ each side of face

Tx2 - Iris Dye-VL: 15J/cm² | 12 msec | 2 passes. **ClearLift Pro 5x5:** -1 and 0 |
2100mJ/pulse | 5 kJ total energy each side

Courtesy of Dr. Munir Somji, UK

Alma
Harmony

Alma o



Rejuvenation and Vascular

ClearLift Pro

Before

After 4 Tx



Before

After 4 Tx



Treatment Protocol

ClearLift Pro 5x5:

2100 mJ/pulse | -1 and 0 | 3 Hz | 6 kJ to total energy

Sessions intervals: 2 weeks

Courtesy of Dr. Munir Somji, UK

Alma
Harmony

Alma o



Colored Tattoo Removal

ClearLift Pro

Before



After 8 Tx




Treatment Protocol

Black ink: **ClearLift Pro 1064nm:** 4/5mm | 600 - 1400mJ/pulse

Red Ink: **ClearLift Pro 532nm:** 5mm | 400mJ

Courtesy of Dr. Eyal Kramer and In house clinic

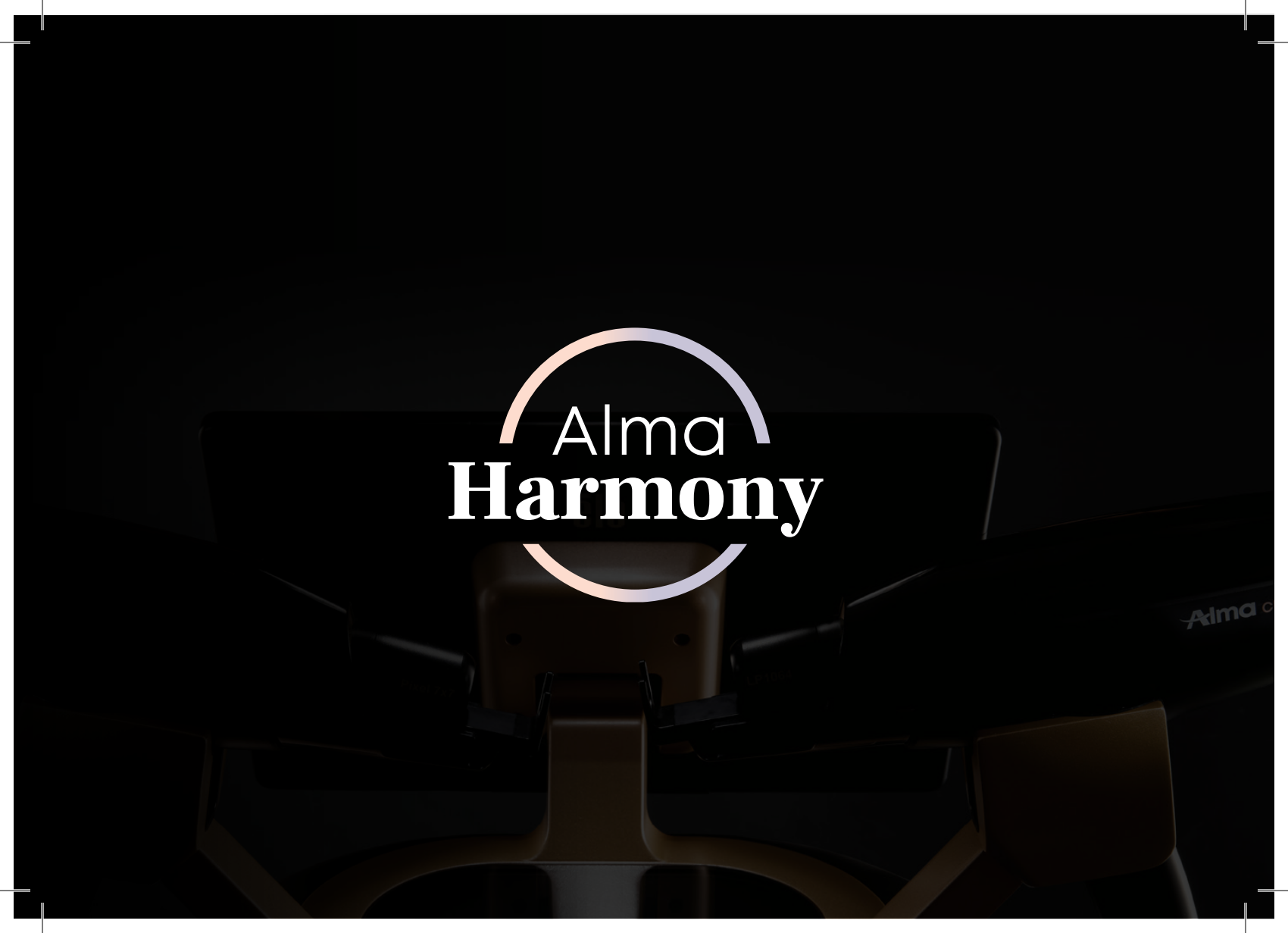

Alma **Harmony**

Alma o



Tattoo Removal

ClearLift Pro

Before

after 3 Tx

after 4 Tx


after 5 Tx



Treatment Protocol

ClearLift Pro 1064nm: 4/5mm | 600 - 2400mJ/pulse

Courtesy of Dr. Eyal Kramer and In house clinic

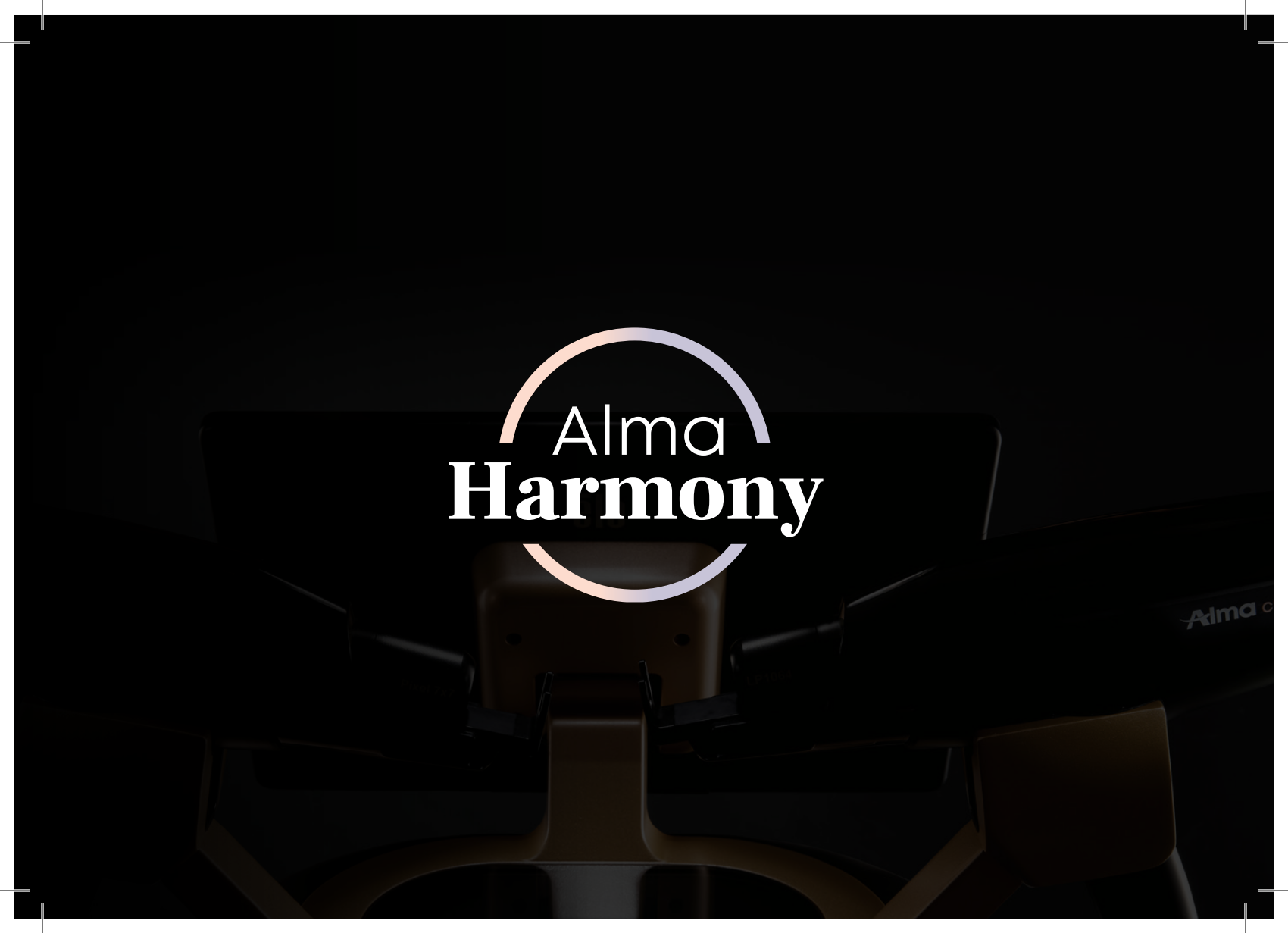

Alma **Harmony**

Alma o



Tattoo Removal Skin Type V

ClearLift Pro

Before



After 6 Tx



Treatment Protocol

ClearLift Pro 1064nm: 4-6mm | 400-800 mJ/pulse | 1-3 Hz |
1.41-6.37 J/cm²

Courtesy of Dr. Eyal Kramer and In house clinic

Alma
Harmony

Alma o



www.almalasers.com