


Acryfold[®] SERIES

Hydrophilic Single Piece Foldable
Intraocular Lens

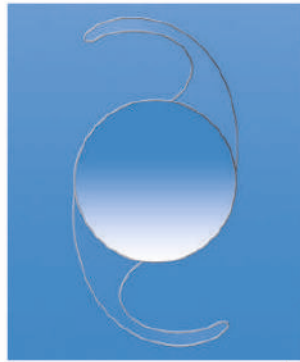

More than One Million Acryfold IOL Implanted Worldwide


Acryfold is manufactured using the latest technology, ensuring good memory and high accuracy. The combination of the best material and choice of design, makes the hydrophilic single piece a complete range of product. The advantage of a foldable lens is clearly that the size of the incision can be made very much smaller.

Model : 601(Clear IOL) & 601Y (Yellow IOL)



Optical Diameter	6.00mm
Overall Length	12.50mm
Positioning Holes	0
Optic Design	Equi-Biconvex
A Constant	118.0
Angulation	0°

Model : 502 (Clear IOL)

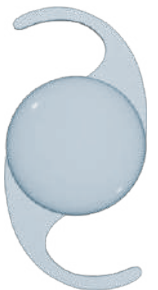

Optical Diameter	6.00mm
Overall Length	12.50mm
Positioning Holes	0
Optic Design	Equi-Biconvex
A Constant	118.0
Angulation	0°

Model : 701 (Clear IOL)



Optical Diameter	6.00mm
Overall Length	11.00mm
Positioning Holes	0
Optic Design	Equi-Biconvex
A Constant	118.0
Angulation	0°

SPECIAL FEATURES & BENEFITS

- Acryfold series of lense are made with high quality bio-compatible UV absorbable hydrophilic material and sterilized with moist heat sterilization.
- The fine polish optic surface with easy un folding of haptics.
- The lens is manufactured by lathe cut and tumble polish with excellent dimensional accuracy.
- All the models are more stable in the capsular bag during capsular contraction.
- Model 701 is more stable and proven for less than 2% PCO in our Multi Center Clinical Studies.
- The flexibility of the material and the design of the lens, facilitates implantation through regular Phaco (2.8mm) incision and above, with proven delivery system.
- Diopter range : -15.0D to 53.0D with 0.5 steps

Manufactured by

 **AI OPTICS PRIVATE LIMITED,**
An ISO 13485:2016 Certified Company

84/1, Aavin Dairy Farm Road,
SIDCO Industrial Estate,
Ambattur, Chennai-600098, INDIA

Marketed by

 **APPASAMY ASSOCIATES**
Empowering Vision*

APPASAMY ASSOICATES (P) LTD.,
20, SBI Officers' Colony, First Street, Arumbakkam, Chennai - 600 106, INDIA
Tel : (91-44) 30101400, 30101401
Email: info@appasamy.com Website: www.appasamy.com