

Maximizing penetration synergy
of Microjetting, Electroporation
and Plasma

1 SYNERGY JETTING | SJ

- Ultrafast Microjetting
- Jet Subcision

The SJ handpiece employs a solenoid mechanism to generate ultra-fast fluid jets. The jet pressure applied to the skin's surface creates micro-openings, enabling drugs to penetrate deeply into the skin.

2 ELECTROPORATION | EP

- Gentle Micro-Channelling
- Enhanced Penetration



The electroporation tip uses electrical microcurrents to temporarily create tiny pores on the skin's surface, significantly enhancing the penetration of active agents.

3 PLASMA | PS

- Modifying Skin Barrier
- Enhanced Skin Permeability



TRIPLE SYNERGY:
BETTER ABSORPTION,
BETTER EFFECT*

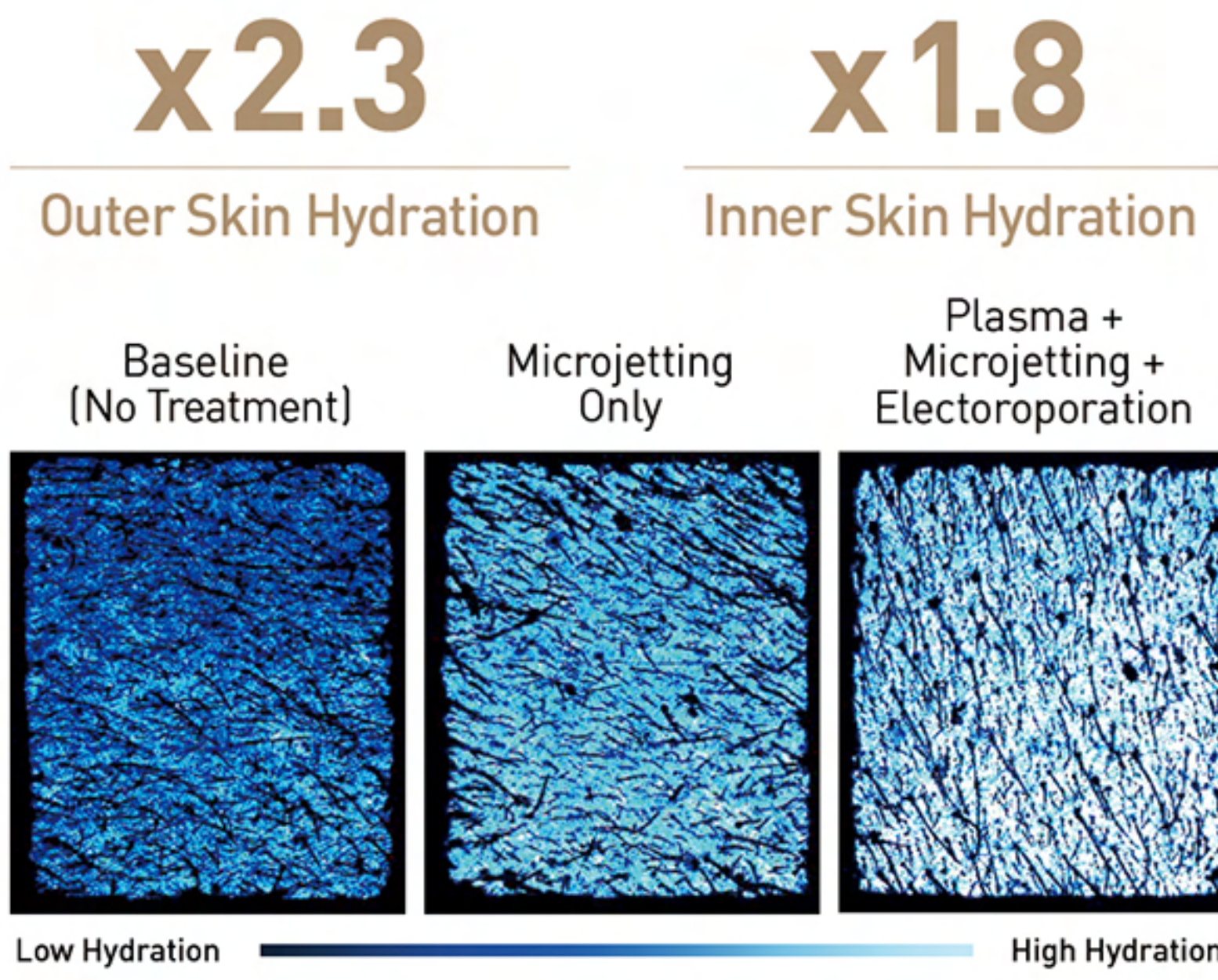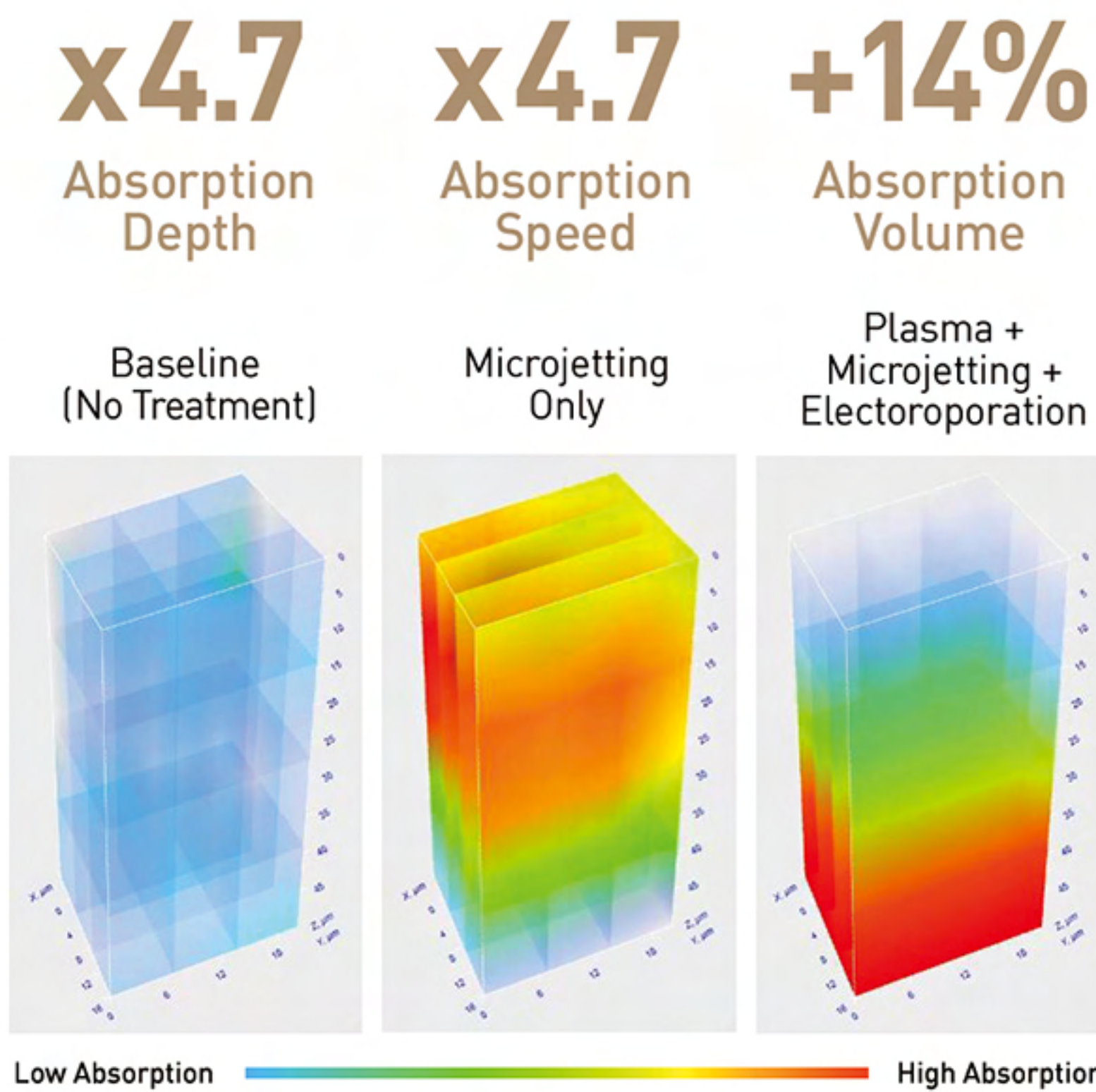

Specifications


Convenience system	RFID, TIS, RMS, Cleaning mode , Jet set-up mode	
Consumables	Nozzle tip, ELPORA tip, PS tip, Brush tip	
Nozzle	Diameter	180μm, 230μm
	Shots	9,000
EL-PORA tip	Distance control	5mm, 10mm, 15mm
Energy generation method	Solenoid	
Energy source	Microjet, Electroporation, Plasma	
LCD screen	12.1"	
Power consumption	700VA	
Dimensions	510(W) x 530(L) x 1,200(H) mm	
Weight	43KG	
EL-PORA Handpiece	Jet speed	150~580m/s
	Injection pressure	100 bar (101.971621kgf/cm ²)
	Injection speed	1~25Hz
EL-PORA Handpiece	EP intensity	1~4Lv
	Power Level (Lv1~5) + Brush tip Level (Lv1~4)	

- This brochure is exclusively crafted for export purposes. For information pertaining to the South Korean market, please consult the localized version of our official website.
- This brochure is intended solely for international application. The device specifications and indications for the usage of the present device may vary based on regional certifications.
- To procure information pertinent to your market, please consult with our regional managers or authorized distributor representatives in your jurisdiction.



HIRONIC Co., Ltd. 19F U-Tower, 767 Sinsu-ro, Suji-gu, Yongin-si, Gyeonggi-do, 16827, South Korea
Sales Support | Tel. +82-31-525-7000 Fax. +82-31-525-7010 E-mail. info@hironic.com
Customer Service | Tel. +82-31-525-7600 E-mail. ocs@hironic.com

SP_2507_V005

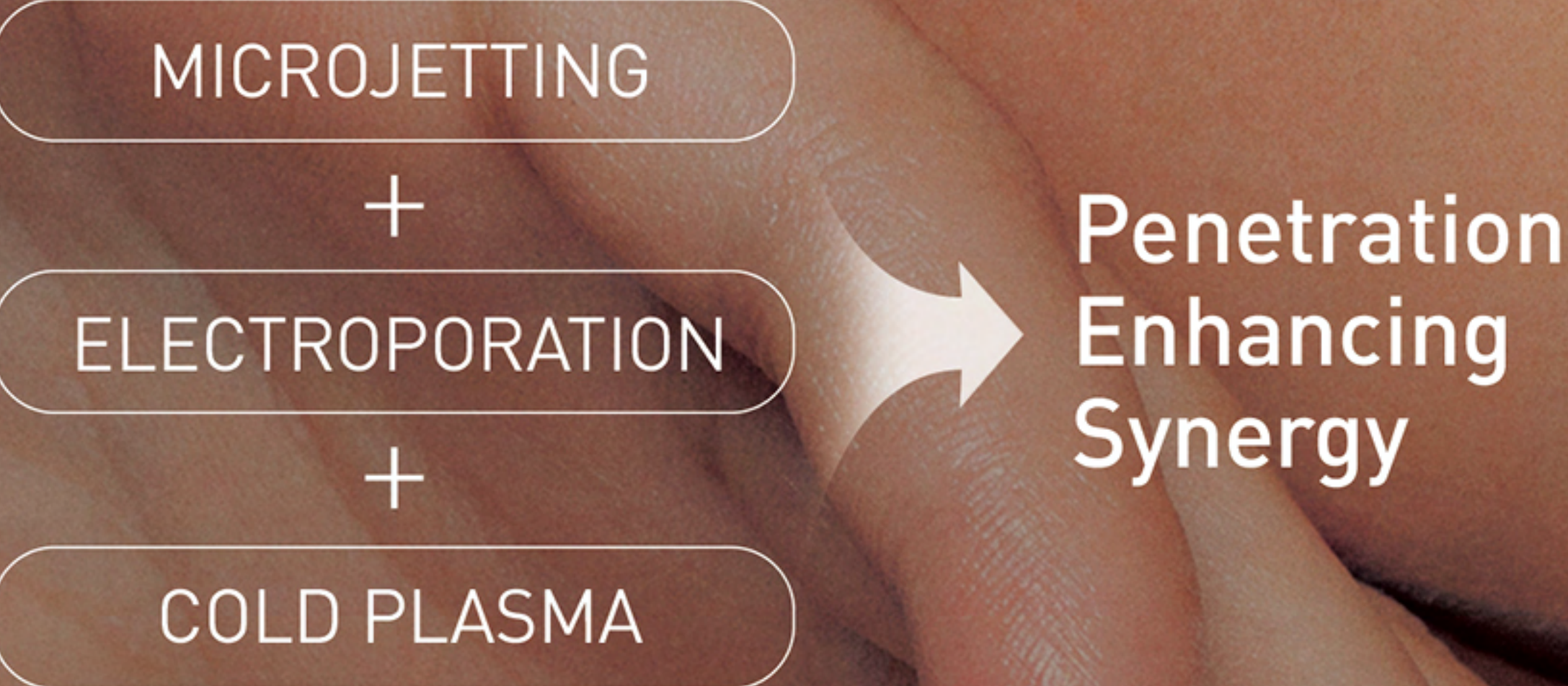

The World's First
Triple Synergy Injector



For Export



3rd Generation Needle-Free Triple Synergy Injector



100% Needle-Free.
Faster. Less Pain.
Less Risks




Uniform Delivery.
Precise Depth.
Volume Control



Needle-Free Subcision.
Scar Treatment.
Enhanced Regeneration

SYNERJET Versatility


Jetting Treatment

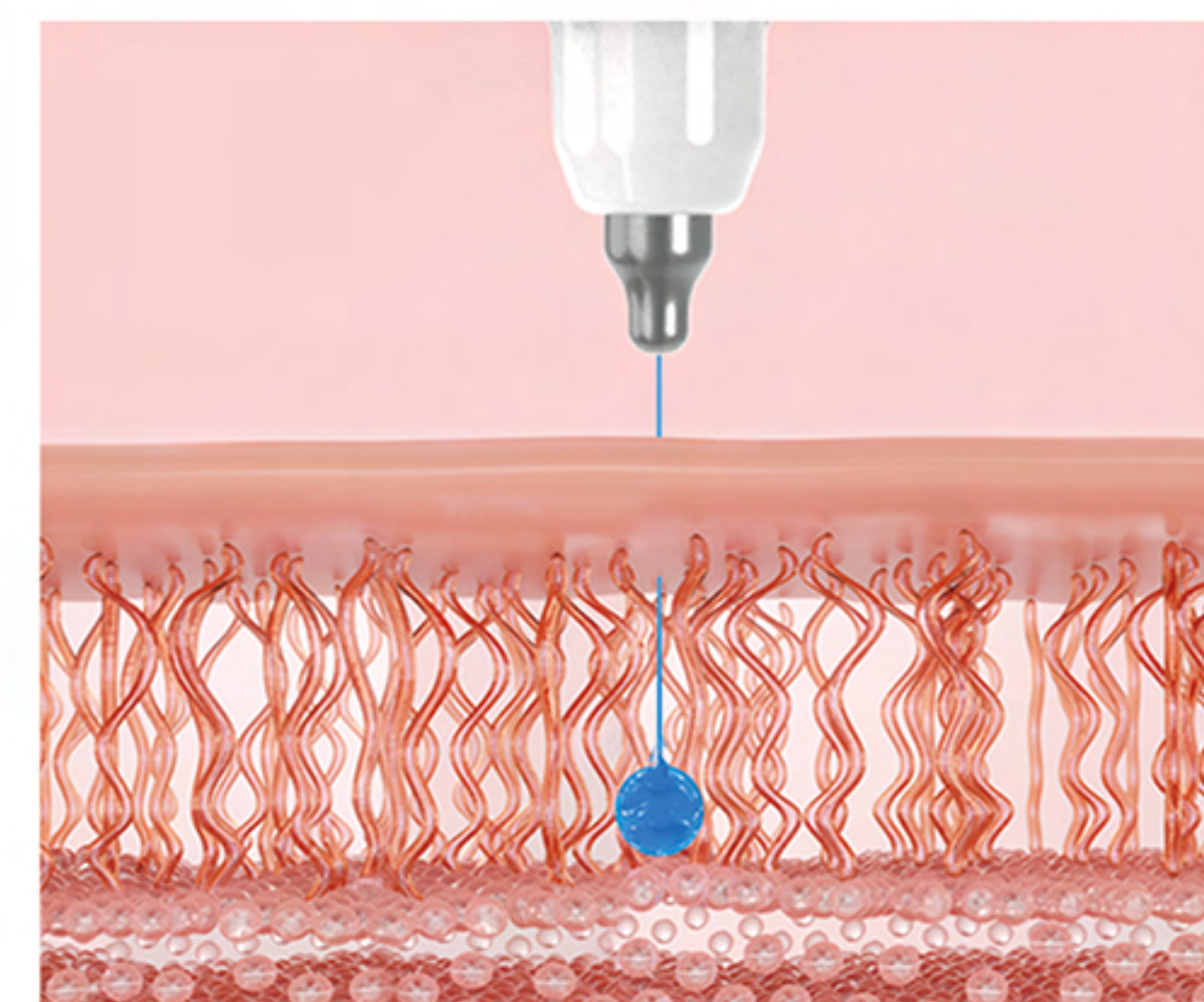

Jetting rapidly and precisely delivers active agents to the upper and mid dermis, with depth controlled by jet power.

- Fast, precise, uniform delivery
- Depth and volume control
- Less pain, trauma, no thermal damage


Jet Subcision Treatment

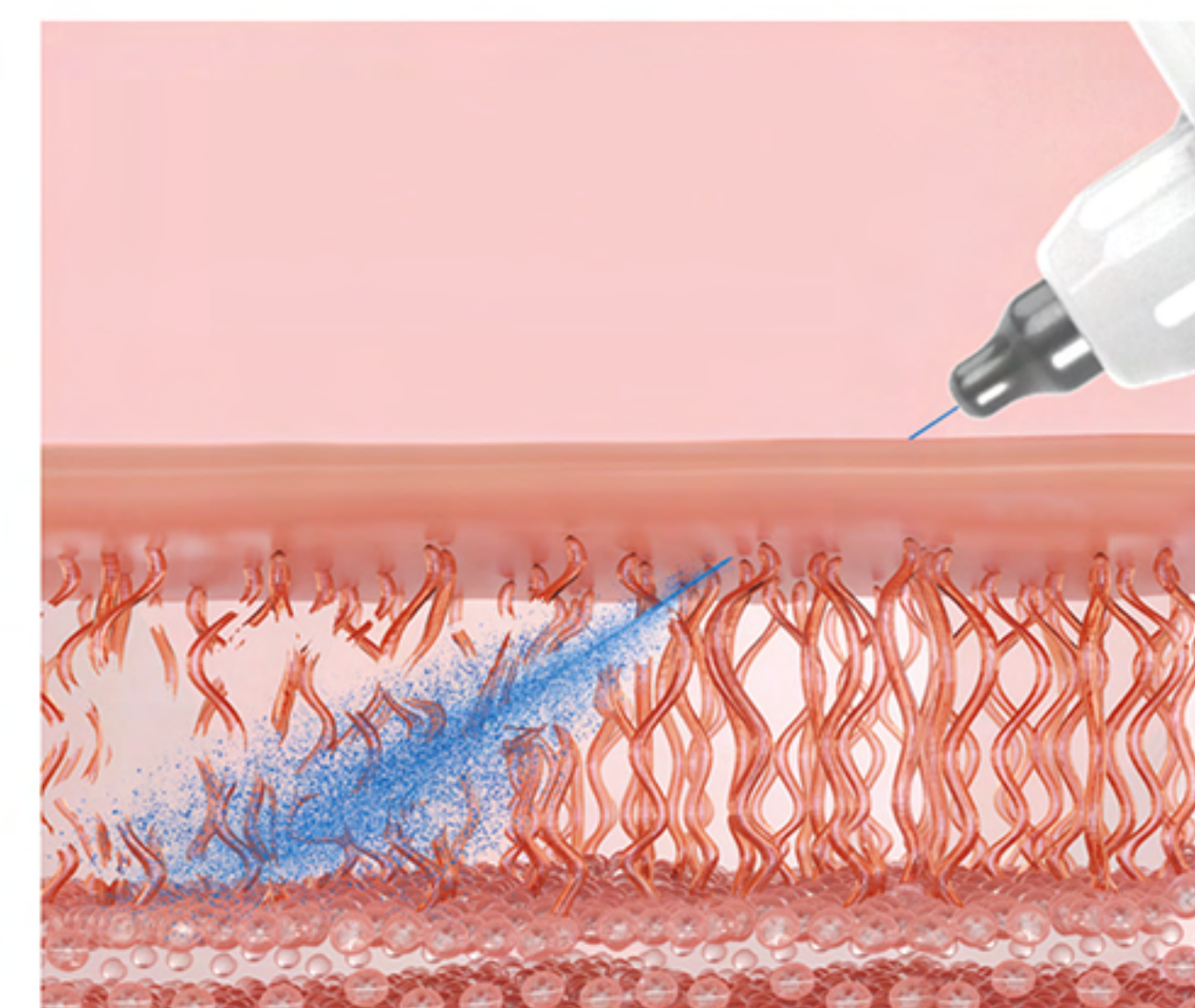

Jet subcision treatment severs tethered fibrous bands, delivering active agents, and promoting volume restoration.

- Gentle and less traumatic micro-subcision
- Simultaneous delivery of therapeutic agent
- Immediate bands release and volume fill-up

Synerjetting™ Treatment

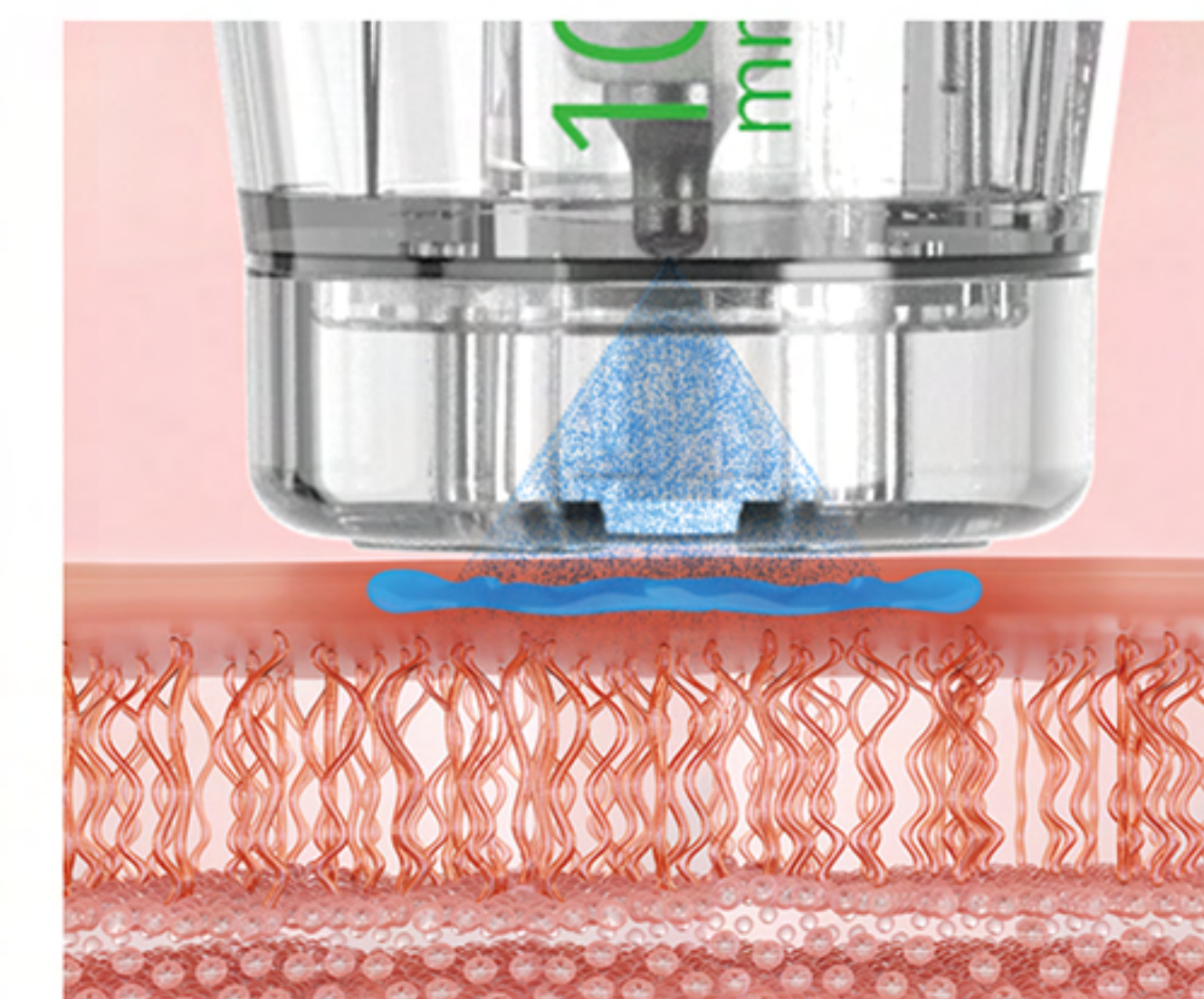

Utilizing the permeation enhancing power of electroporation, Synerjetting™ allows efficient delivery of topical agents.


- More efficient delivery of non-injectables
- Electroporation-enhanced penetration
- Better treatment efficiency

Special EP Tip Design:


- Distance Control
- Prevention of Drug Splattering
- Practitioner and Patient Protection

Optimal Subcision Angle


- 33° angle is advised as optimal for micro-jet subcision
- Special nozzle design features optimal angle guidance for practitioner convenience

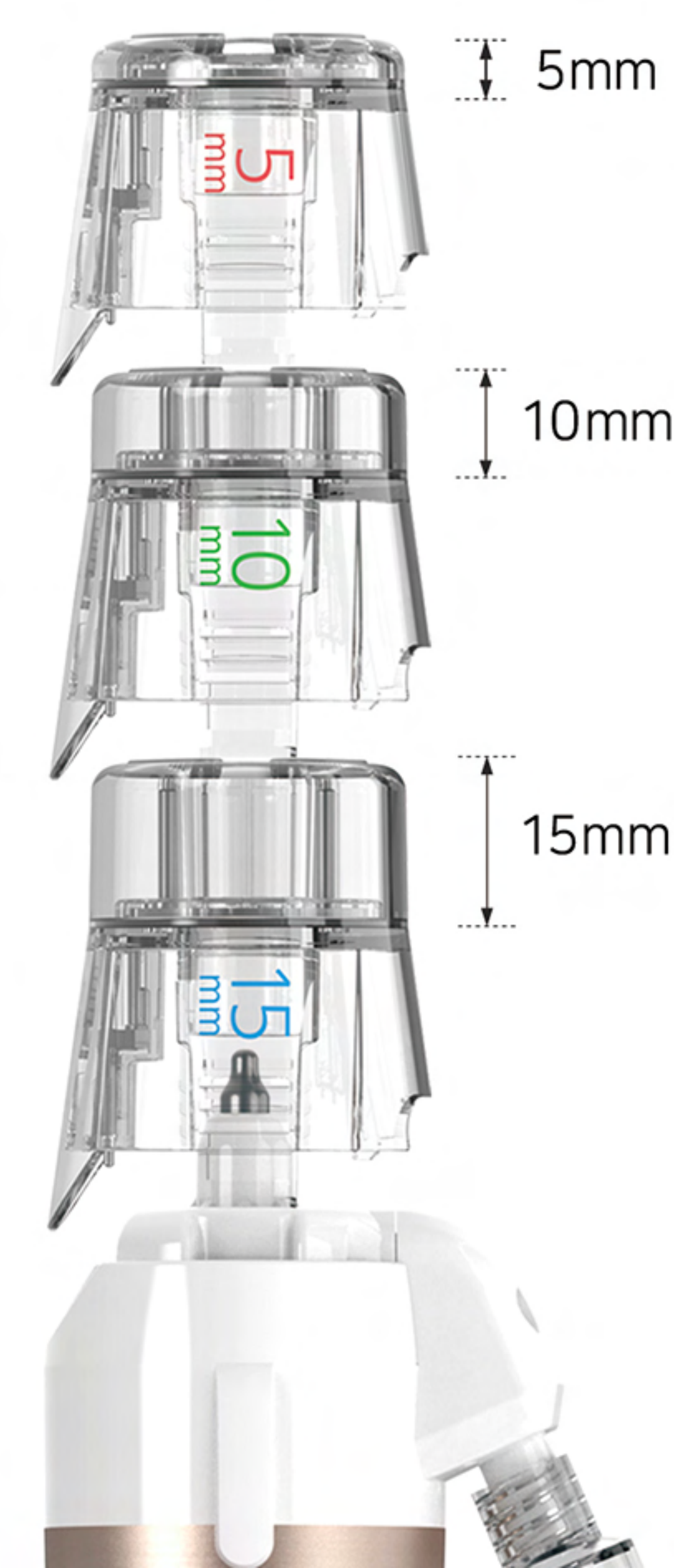

Expanded Treatment Indications

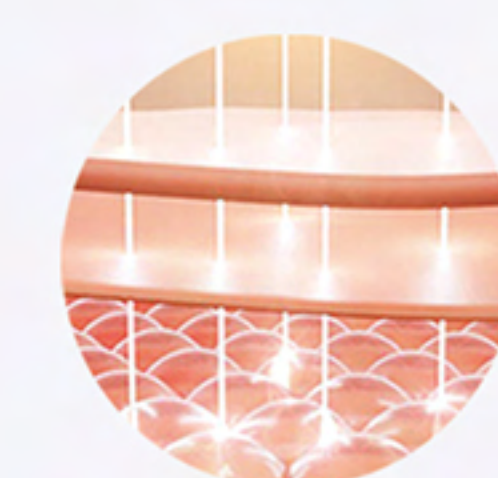

Rejuvenation
Fine lines and wrinkles,
skin texture, hydration,
toning and other

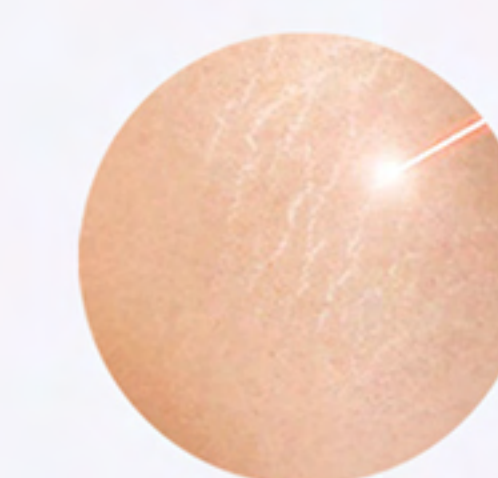

Depressed Skin Conditions
Scars, acne scars,
stretch marks and other



Scalp Treatments
Alopecia

Convenience

1 Treatment Information System (TIS)

Providing the best protocol for each indication and treatment area for the utmost guidance to users. The saved parameters can be managed through the history section for future treatments.

Remote Maintenance System (RMS)

- Reminder to replace the consumable parts
- Storage of protocols
- Device self-diagnosis and inspection
- Malfunction prevention and assistance functions

Big Data System

Accumulating data for enhancing the quality of user experience.




Real-Time
Device Monitoring

Historical
Treatment Data

Pre-A/S

2 Smart GUI

- Jet set up mode
- 9 pre-set modes, 27 slots for saving customized settings
- Nozzle cleaning mode
- Shot count alarm

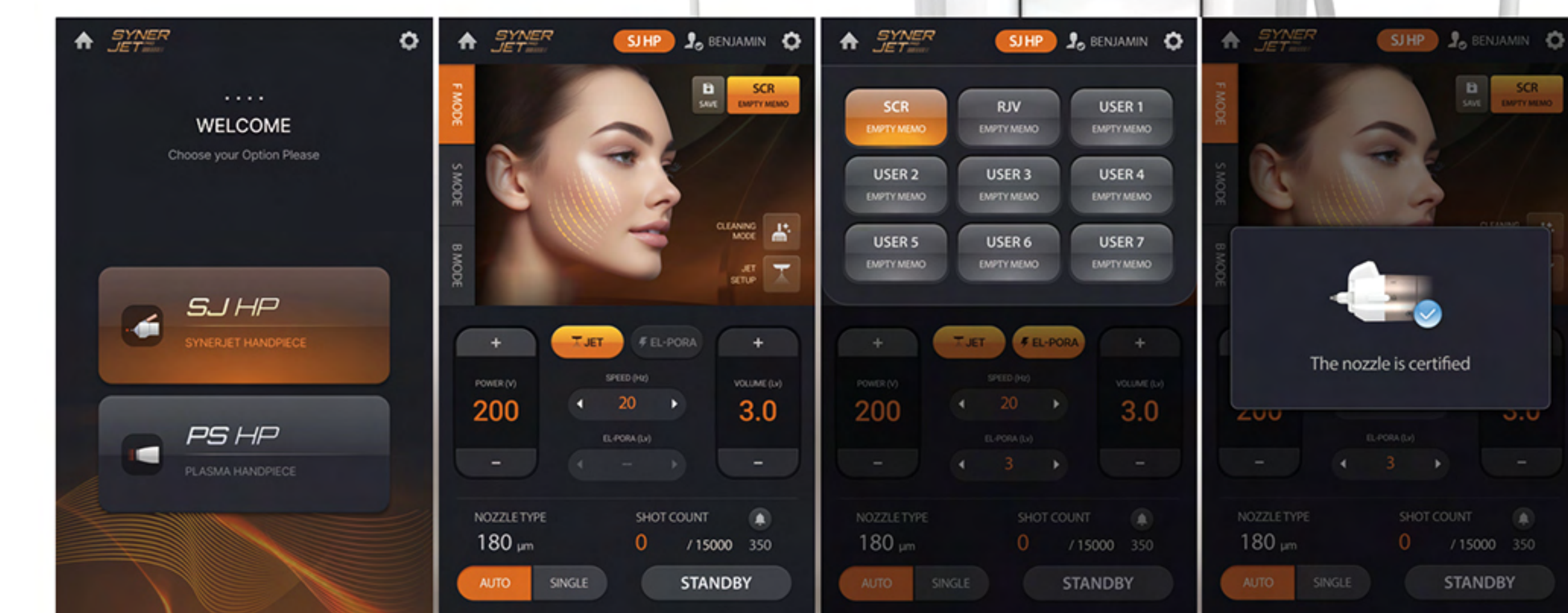

12.1" Wide LCD Display

Max 75° Angle



Two Handpieces

- SJ handpiece
- PS handpiece



RFID Sensor Technology

- Consumable authenticity verification
- Security enhancement
- Contact-based tag activation: automatic consumable recognition



Back Handle

Convenient mobility
with advanced design



Nozzle Management

Clean-up mode for
convenient management



Consumable Nozzles & Tips

9,000 shots with enhanced
nozzle durability

9,000

