



PRECISION-DRIVEN SKIN TIGHTENING

01

2.45 GHz UHF WITH POWER COOLING

This technology efficiently delivers energy to the desired depth, effectively transferring heat to the targeted area while protecting the skin through cooling while efficiently transmitting energy.

02

SUPPORTS SKIN CONTOURING AND TIGHTENING

Targets fat cells and stimulates collagen production, leading to fat reduction and collagen regeneration.

03

SMART HANDPIECES BY PENETRATION DEPTHS

MIGLOU™ offers customized treatment for each patient and enhanced convenience of operation. Choose our various handpieces for different indications for use.

Specifications

Energy	Ultra High Frequency
Frequency	2.45GHz
Power	Max. 200w
Electrical Specifications	100~240Vac / 50~60Hz / 1500VA
Dimensions	700(W) x 435(L) x 1090(H)
Weight	63Kg
Convenience System	Graphical User Interface (GUI)
	12.1" Color Display Touch Screen
Smart Handpiece	- 3mm, 7mm, 12mm - LED indicator - 1.54" LCD screen - Four buttons · Energy Up / Down · Ready / Emission 

WARRANTY

2 YEAR

WARRANTY

All products are covered by a 2-year warranty

Implementation of a Long-Term Warranty Policy

- In accordance with the executed contract
- Consumable items, including cartridges and tips, are covered by a 1-year warranty

Emergency Support

Quick response to urgent repair cases through conference calls when dealers need additional support.

Emergency Conference Call Center [WhatsApp] 010-5013-5343

Emergency Conference Call Schedule [KST] 07:30AM ~ 19:30PM

Copyright(C) HIRONIC. All rights reserved



✓ 10-15 min Procedure



✓ Immediate Return to Daily Life



✓ No Pain, No Anesthesia



HIRONIC Co., Ltd. | [HQ] 19F U-Tower, 767 Sinsu-ro, Suji-gu, Yongin-si, Gyeonggi-do, 16827, South Korea
| 5F Akusan Bldg., 1-5-8 Iidabashi, Chiyoda-ku, Tokyo, 102-0072, Japan

Sales Support | +82-31-525-7400 [FAX] +82-31-525-7010 [E-mail] info@hironic.com

Customer Service | +82-31-525-7600 [E-mail] ocs@hironic.com

MU_2608_V004





Let Me Glow,

MIGLOU™





HIRONIC Ambassador
YANO SHIHO

COMFORTABLE NON-INVASIVE

FAT REDUCTION & CONTOURING & FACE LIFTING


For Internal Use Only

PRECISION-CONTROLLED ENERGY⁺

2.45 GHz UHF^{*} (*UHF : Microwave)


What is Ultra-High Frequency?

- Electromagnetic current is emitted at a frequency of 2.45 GHz, gently vibrating cell tissues within the skin to generate heat
- This energy is precisely delivered to activate fibroblast cells in the dermal layer, without causing irritation to the epidermis
- Consequently, it effectively enhances facial or body contours with minimal discomfort, preserving surrounding tissues from damage

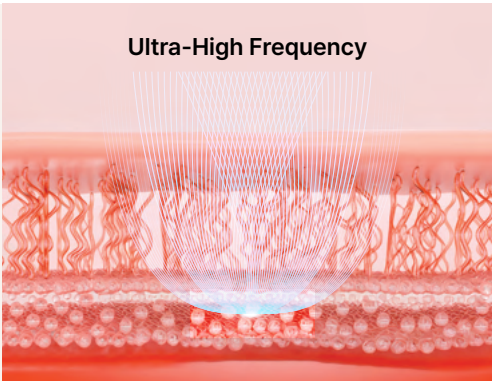

FIM^{*} Smart Matching (*FIM: Focused Intensive Microwave)

- ✓ Focused Intensive Microwave technology enhances **the concentration and efficiency of delivered energy**
- ✓ MIGLOU's proprietary, patented technology precisely **delivers energy to targeted focal zones with high accuracy**

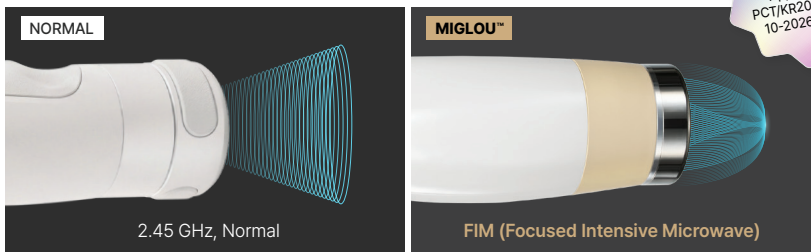

Contact & Moving Smart Sensor



Safer Treatment with Smart Sensor Technology


- Built-in handpiece safety system helps prevent overheating
- Detects both contact and movement for controlled use
- Contact cooling system helps protect the skin and provide comfort treatment

Electromagnetic Shielding Structure

- Reduced airborne energy output for stable operation
- Safe to use with everyday-level electromagnetic exposure

Delimited	MIGLOU™	Daily life estimates
Radiated power (At 200W)	0.10W	At the home Wi-Fi router level
Energy efficient	FIM focusing	Eliminate unnecessary energy dispersion

*IoT Technical Support Center Test result / August 12, 2025 / Antenna characteristics Measurement (Patient absence condition)

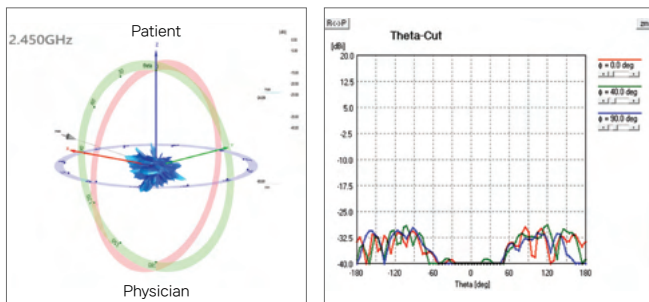

TREATMENT VERSATILITY⁺

Three Handpieces Designed for Different Treatment Areas

Smart Handpieces for Targeted Face to Body Contouring

Customized handpieces selection allows tailored treatments based on individual needs and target areas.

3mm 7mm FACIAL ENHANCEMENT

- Jawline
- Double chin
- Overall facial
- Contouring effects



7mm 12mm BODY SHAPING & FAT REDUCTION

- Thighs
- Abdomen
- Upper arms
- Lovehandles

Facial parts

Body parts

Key Features of the Smart Handpiece

Four Buttons Control

- Energy Up / Down
- Ready / Emission



1.54" LCD Screen

- Standby
- Ready
- Progress mode



LED Indicator

Wavelength control display



Contact Cooling System

Seamlessly built into the handpiece for maximum comfort




DESIGNED FOR CLINICAL EFFICIENCY⁺

User-Friendly System

Intuitive Graphical User Interface (GUI)

The built-in tutorials and treatment protocols simplify the overall workflow, making treatment management more efficient, consistent, and easy to operate.

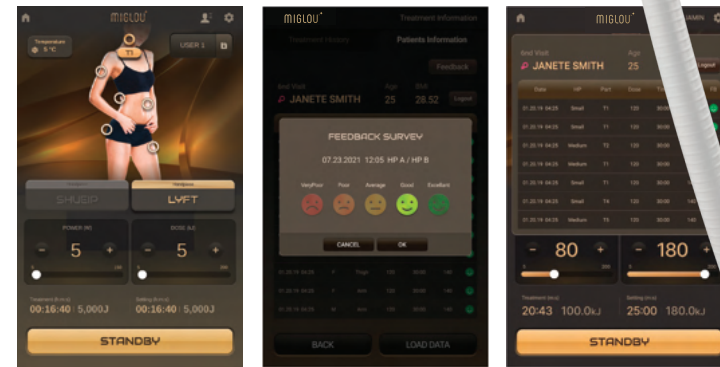

Treatment Information System (TIS)

- Provides optimized treatment protocols for efficient clinical operation
- Enables systematic management of patient information and treatment history


Remote Maintenance System (RMS)

- Supports built-in device diagnostics and performance monitoring
- Helps prevent errors and enables timely maintenance support

Convenience

12.1" Wide LCD Display

Adjust the viewing angle up to 60° vertically and horizontally for optimal visibility during treatment.

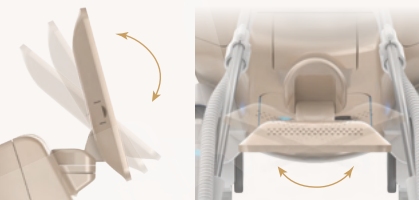

A 2-year warranty

Implementation of a Long-Term Warranty Policy



The Lightest Handpiece

Minimizes wrist strain during procedures

200g

Cable Holder

For more comfortable and lightweight handpiece use

