

# QTekLaser™ 770 nm Fiber Laser

#### **Product Description**

The QTekLaser™ 770 nm Fiber Laser is a high-performance source tailored for atomic physics and quantum research, particularly applications involving potassium transitions, laser cooling, and interferometry. Using advanced Second-Harmonic Generation (SHG), it delivers up to 2 W of narrow-linewidth, polarization-maintaining output at 770 nm. With a compact 3U, 19″



rack-mount chassis, robust beam quality, and an IoT-enabled control system, this laser combines scientific precision with industrial reliability for seamless laboratory integration.

#### **Features**

- High output power (2 W)
- High reliability with all-fiber design
- Narrow linewidth (<20 kHz)
- Excellent power stability (<1%)
- User-friendly interface via IoT technology
- 3U 19" Rack mount
- Certified to IEC 60825-1:2014 safety standards

#### **Applications**

- Quantum computing
- · Quantum sensing
- Atomic interferometry
- Laser cooling and trapping
- Research

# **Second-Harmonic Generation (SHG) Technology**

QTekLaser™ offers laser systems with extended wavelength range by combining nonlinear frequency conversion technologies. Through second-harmonic generation (SHG) we achieve significant laser power at the visible and NIR regime (figure 2; red-shaded cells of table 1). With the development of waveguide technology, nonlinear frequency conversion has become a powerful tool to extend the application scope of fiber lasers.

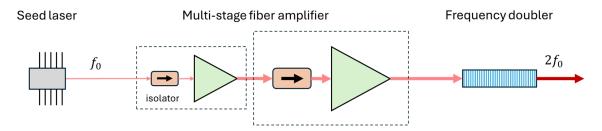


Figure 2. Laser system with integrated frequency doubler.



Table 1. Wavelength extension of QTekLaser™ products using nonlinear frequency conversion.

	SF-AMP-Yb	SF-AMP-Er	SF-AMP-Tm
	(1012 – 1080 nm)	(1530 – 1610 nm)	(1776 – 2040 nm)
SF-AMP-Yb	506 – 540 nm		
(1012 – 1080 nm)			
SF-AMP-Er	609 – 646 nm	765 – 805 nm	
(1530 – 1610 nm)			
SF-AMP-Tm	645 – 706 nm	822 – 900 nm	888 – 1020 nm
(1776 – 2040 nm)			

#### **Specifications**

Parameter	Value	
Wavelength	770 nm	
Max. Output Power	2 W	
Polarization Extinction Ratio (PER)	17 dB	
Power Stability	< 1%	
Laser Linewidth	< 20 kHz	
Operating Temperature	17-25 °C	
Cooling Mode	Air or Water cooling	
Beam Quality	TEM00; $M^2 \sim 1.0$	
Output Mode	PM780-HP with FC/APC connector	
Side-Mode Suppression Ratio (SMSR)	70 dB	
Chassis	3U 19" rack mount	

## Safety & Retro-Reflection Advisory

The 770 nm laser system complies with 21 CFR Subchapter J, Part 1040 (as enforced by FDA) and IEC 60825-1:2014 safety standards.

Users must ensure that significant light is not retroreflected into the system—doing so may harm performance or damage the laser. Use of an external optical isolator is strongly recommended. Damage from retroreflected light is not covered under warranty.

## **Ordering Information**

Part Number: QT-LASR-SHG-770-2-A-2-2-0

Configuration: Seed laser + Er-doped fiber amplifier + SHG



#### **Performance Figures**

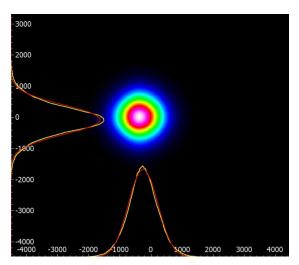


Figure 3. Beam Profile. Collimated with Thorlabs F220APC-780 collimator (not included). Excellent round Gaussian profile with a  $1/e^2$  beam diameter of 1.9 mm.

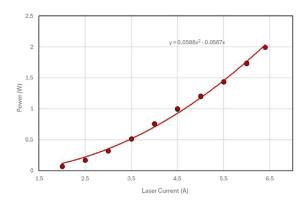


Figure 4. Laser power vs. current.

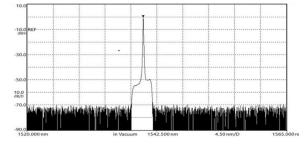


Figure 5. Optical Spectrum: SMSR 70 dB.

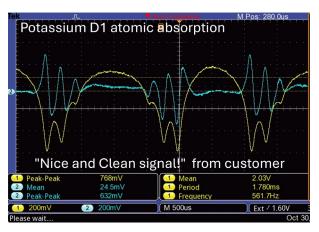


Figure 6. AMO Application: Potassium D1 Atomic Absorption.



#### **Mechanical Details**

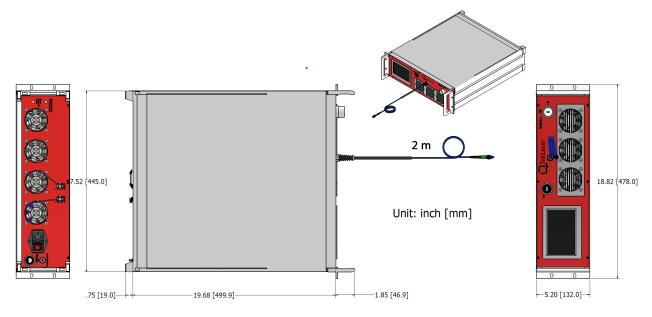


Figure 6. Mechanical dimensions of the 770 nm fiber laser.

# **Product Photographs**



Figure 7. 770 nm fiber laser 3U chassis air cooled.



Figure 8. 770 nm fiber laser 3U chassis back panel air cooled.