

RODIN[®] LUTING CEMENT

RESIN-MODIFIED GLASS IONOMER CEMENT

PROVEN STRENGTH,
EFFORTLESS VERSATILITY



Developed to tackle the demanding challenges of modern dentistry, **Rodin[®] Luting Cement** is the most advanced resin-modified glass ionomer cement.

RODIN[®]

LUTING CEMENT

RESIN-MODIFIED GLASS IONOMER CEMENT



Veneer, Adhesive Bridge



Crown & Bridge



Post, Bonding for Build-up, Crown



Inlay/Onlay

Product Features



Reinforced Strength & Minimal Shrinkage:

Proprietary chemistry strengthens the cement matrix for durable bonds across printed resins, lithium disilicate, zirconia, alloys, and titanium.



Sustained Fluoride Release:

Continuously releases fluoride to support remineralization, reduce secondary caries risk, and enhance long-term clinical success.



Superior Protection Against Microbial Challenges:

Low shrinkage reduces marginal gaps to prevent bacterial infiltration and maintain long-term restoration integrity.



Enhanced Bond Strength & Unmatched Reliability:

Engineered for strong, lasting adhesion to natural abutments and implants.




Optimized Aesthetics for Seamless Integration:

With natural translucency and high radiopacity, four shades blend seamlessly and mask dark preps or metal abutments.



Comparison of Bond Strength for Different Materials and Cements

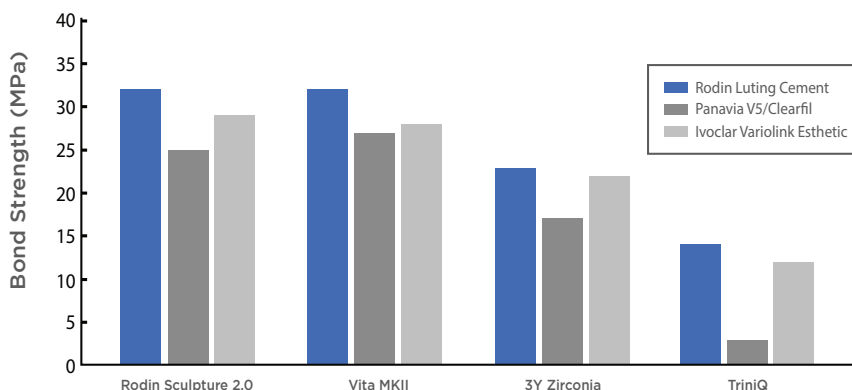

Data provided by Dr. Russell Giordano of Boston University

Best in Class Bond Strength on:

- ✓ 3D Resin based
- ✓ Glass-ceramic
- ✓ Zirconia Restorations

RODIN[®]

LUTING CEMENT

RESIN-MODIFIED GLASS IONOMER CEMENT

WORKFLOW



STEP 1

Etch enamel and dentin with **EtchPro™**



STEP 2

Apply **Rodin® Bond** Adhesive (optional)



STEP 3

Apply **Rodin® Silane Primer™** (Crown and Veneer)



STEP 4

Apply **Rodin® Luting Cement** to Restoration



STEP 5

Seat Restoration and Remove Excess Cement



STEP 6

Light Cure Surface from All Angles.

COMPLETE RESTORATIVE WORKFLOW

Product Name					
	Luting Cement	Veneer Cement	Rodin Bond	Glaze N2-Free	SensiGuard
Product Type	Resin Modified Glass Ionomer Cement	Light-Cured Resin Cement	Universal Adhesive/Bonding Agent	Restorative Glaze	Dentin Desensitizer
Indications	Cementation of crowns, bridges, inlays, onlays, posts, some ceramics	Bonding of porcelain/composite veneers (thin restorations)	Bonding to enamel, dentin, & various restorative surfaces	Glazing & finishing of all surfaces including porcelain, zirconia, & ceramic restorations	Reduction of dentinal hypersensitivity
Setting Mechanism	Dual-cure (light + chemical)	Light-cure only	Light-cure	Light-cure	Self-Drying (No Light Cure)
Available Shades	Multiple Esthetic Shades (Clear, White, White Opaque, Yellow)	Multiple Esthetic Shades (Clear, White, White Opaque, Yellow)	Clear/Universal	Clear/Universal	N/A
Delivery Method	Automix Syringe	Syringe with Tips	Bottle (Brush Application)	Bottle (Brush Application)	Bottle (Brush Application)
Key Features	High bond strength to multiple substrates, easy cleanup, low film thickness	Ultra-thin film thickness, excellent esthetics, color stability	One-step universal application, strong adhesion to tooth structure & restorations	High-gloss finish, N2-free formulation, easy application, enhances esthetics & longevity	Seals dentinal tubules, immediate & long-lasting sensitivity relief





Item#	Product Description
-------	---------------------

451101	Rodin Luting Cement Starter Kit
---------------	--

5x	Rodin Luting Cement Syringe, 5 g Shades: 2 x Clear, 1 x White, 1 x White Opaque, 1 x Yellow	10x	Tapered Mixing Tips
4x	Rodin Try-in Gel Clear, 3 g Shades: 1 x Clear, 1 x White, 1 x White Opaque, 1 x Yellow	10x	Flat Mixing Tips
1x	Rodin Bond Adhesive Bottle, 5 mL	10x	Intra-Oral Tips
1x	Rodin Silane Primer Bottle, 5 mL	10x	Dispensing Tips
2x	Rodin EtchPro Etching Gel Syringe, 1.2 mL	10x	MicroPro Micro Applicators Ultra Fine
		10x	4-Well Mixing Wells



451102	Rodin Luting Cement Refill Clear
---------------	---

2x 5g Clear Syringes

451103	Rodin Luting Cement Refill White
---------------	---

2x 5g White Syringes

451104	Rodin Luting Cement Refill White Opaque
---------------	--

2x 5g White Opaque Syringes

451105	Rodin Luting Cement Refill Yellow
---------------	--

2x 5g Yellow Syringes

All refills include 5x Intraoral Tips, 5x Tapered Mixing Tips, 5x Flat Mixing Tips



451107	Rodin Try-in Gel Refill: 1x 3g Try-in Gel Syringe, Clear
---------------	---

451108	Rodin Try-in Gel Refill: 1x 3g Try-in Gel Syringe, White
---------------	---

451109	Rodin Try-in Gel Refill: 1x 3g Try-in Gel Syringe, White Opaque
---------------	--

451110	Rodin Try-in Gel Refill: 1x 3g Try-in Gel Syringe, Yellow
---------------	--



471102	Rodin Bond Dental Adhesive Refill, 5mL Bottle
---------------	--



471201	Silane Primer Refill, 5ml Bottle
---------------	---



135	Rodin EtchPro 38% Phosphoric Acid Etch
------------	---

4x 1.2 mL Syringes, 8x Dispensing Tips

135E	Econo pack: 36x 1.2 mL Syringes
-------------	--

135C	Clinic pack: 144x 1.2 mL Syringes
-------------	--

231	25 Ga tips, Etch, Blue, 100/pk
------------	---------------------------------------



901	Disposable MicroPro Applicators Standard 2.0mm, 400/pk
------------	---

902	Disposable MicroPro Applicators Fine 1.5mm, 400/pk
------------	---

903	Disposable MicroPro Applicators Ultra fine 1.0mm, 400/pk
------------	---

904	Disposable MicroPro Applicators Tapered, 400/pk
------------	--



906	Mixing Well, 4-Well, 200/pk
------------	------------------------------------