



AlloMend[®] DUO

ACELLULAR DERMAL MATRIX

The Natural Solution for Soft Tissue Reconstruction

AlloMend Duo Acellular Dermal Matrix (ADM) is a human-derived, non-directional allograft that incorporates into surrounding tissue and provides a scaffold for cell repopulation and regeneration of the natural host tissues.¹ This deep reticular dermal allograft is used for soft tissue reconstruction, helping surgeons repair and reconstruct soft tissue defects for long-term structural support.




⊕ THE BENEFITS OF ALLOMEND DUO

Flexible and Pliable	Provides optimal handling characteristics enabling precision placement
Packaged Pre-Hydrated	Packaged in sterile water, immediately ready to use with room temperature storage
High Suture Retention Strength	Average suture retention strength of >39N/mm ² for assurance during application ²
Terminally Sterilized	Irradiated using e-beam technology to minimize the risk of infection
Precisely Processed	Consistent thickness and size through proprietary splitting and cutting technologies

⊕ Dual-sided Reticular Dermal Matrix for Non-Directional Implantation

Nearly equivalent tissue structures on both sides of the allograft for ease of use and increased confidence.

ALLOMEND DUO SUPERFICIAL DERMAL SURFACE

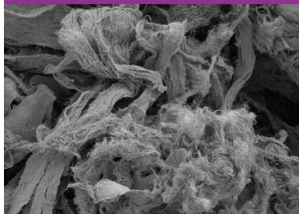


200X Magnification



1000X Magnification

ALLOMEND DUO DEEP DERMAL SURFACE

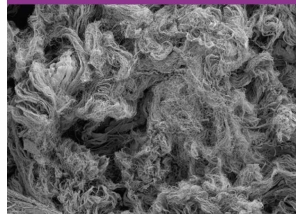

200X Magnification



1000X Magnification

AlloMend® Duo



AlloMend® Duo ACELLULAR DERMAL MATRIX


Deep Reticular Dermal Allograft

MESH	THICKNESS	WIDTH	LENGTH	AREA	REF/PRODUCT #
Non-Meshed	1.0-2.0 mm	2 cm	4 cm	8 cm ²	DLIN3966
Non-Meshed	1.0-2.0 mm	4 cm	8 cm	32 cm ²	DLIN3967
Non-Meshed	1.0-2.0 mm	6 cm	12 cm	72 cm ²	DLIN3968
Non-Meshed	1.0-2.0 mm	8 cm	10 cm	80 cm ²	DLIN3969
Non-Meshed	1.0-2.0 mm	8 cm	12 cm	96 cm ²	DLIN3970
1:1	1.0-2.0 mm	4 cm	16 cm	64 cm ²	DLIN3971
1:1	1.0-2.0 mm	6 cm	12 cm	72 cm ²	DLIN3972
1:1	1.0-2.0 mm	8 cm	10 cm	80 cm ²	DLIN3973
1:1	1.0-2.0 mm	8 cm	12 cm	96 cm ²	DLIN3974
1:1	1.0-2.0 mm	8 cm	16 cm	128 cm ²	DLIN3975
1:1	1.0-2.0 mm	8 cm	20 cm	160 cm ²	DLIN3976
1:1	1.0-2.0 mm	16 cm	20 cm	320 cm ²	DLIN3977

⊕ DermaTrue™ DECELLULARIZATION PROCESS


AlloMend Duo ADM is created using AlloSource's proprietary DermaTrue Decellularization Process to remove cellular debris (including DNA, RNA, proteins and antigens), without the use of harsh detergents or enzymes which can leave residuals in the tissue. The dermal tissue is rendered acellular, contributing to a low immunologic response, while retaining growth factors and maintaining the morphological collagen structure.³

H&E (hematoxylin and eosin) stain review of "before and after" decellularization process.

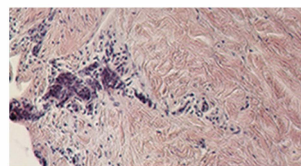

BEFORE

Noticeable large number of well-defined cell nuclei (purple)

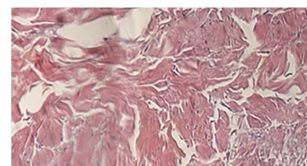

AFTER

Absence of identifiable defined nuclei; no viable cells present

1. Stillwell R., Delaney R. The biomechanics of acellular dermal matrix biocompatibility study. *AlloSource White Paper*. 2016; 00088-LIT(001).
2. Data on File.
3. Delaney R., Stilwell R. The biologic properties of allomend acellular dermal matrix: growth factor study. *AlloSource White Paper*. 2016; 00104-LIT [001].

00177-LIT(001) - 02/2023

LEARN MORE:

Contact Spartan Medical Inc.

(888)240-8091

CS@SPARTANMEDICAL.COM



SPARTAN MEDICAL INC.

1445 Research Blvd, Suite 520,
Rockville, MD 20850

spartanmedical.com

UEI: LCJBD88KQ8G3

DUNS: 807438341 | CAGE:65VT5

