

Leaseir SRF

Cellulite and remodeling
integral platform



Leaseir

www.leaseir.com

Technology



SUPERFICIAL RADIO FREQUENCY
4 modes to treat at different depths



INTELLIGENT PARAMETERS
6 power levels to minimize power losses




COMFORTABLE AND PAINLESS
Suction and cutaneous cooling systems for better results and greater comfort



SIMULTANEOUS HEADPIECES
3 simultaneously connected headpieces allow continuous operation (FC, FS and BS)



WITHOUT PROGRAMMED OBSOLESCENCE
The 3 headpieces have a non-limited duration per use or for number of shots

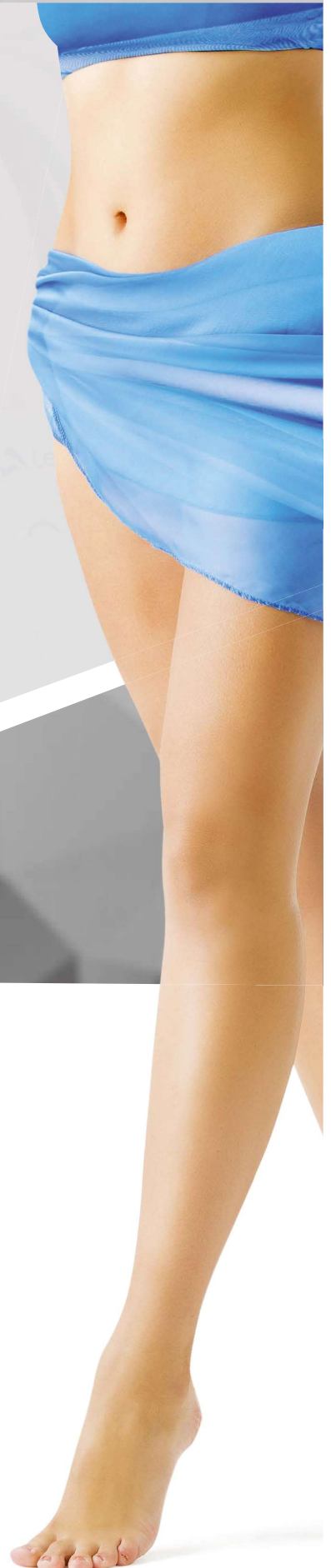

Applications



**PRACTICAL
FUNCTIONAL
NO CONSUMABLES
BEST RESULTS**

**EFFECTIVE TREATMENT
FOR CELLULITE AND
FACIAL AND BODY
REMODELING**

Cellulite

Leaseir SRF is a bipolar radio frequency, with two Facial Suction and Body Suction heads, designed to treat cellulite nuclei, while strengthening existing collagen fibers, giving a more toned and firmer appearance to the skin.

In both headpieces, the radio frequency and the suction favor the natural regeneration processes of the organism.

RF energy is specifically targeted at fatty accumulations, while vacuum manipulation increases blood circulation and lymphatic drainage, increasing treatment efficiency.


Flaccid Skin

Collagen fibers weaken over time due to aging and stress.

Leasier SRF, with its Facial Cold and Facial Suction handpieces, heats the skin between 39 ° and 41 °, stimulating collagen synthesis to strengthen existing fibers and improving laxity, flaccid skin and texture.

The isolation and filling results are immediate, given the contraction and long-term strengthening of fibers in skin tissues, making the skin smoother and firmer.

Treatment

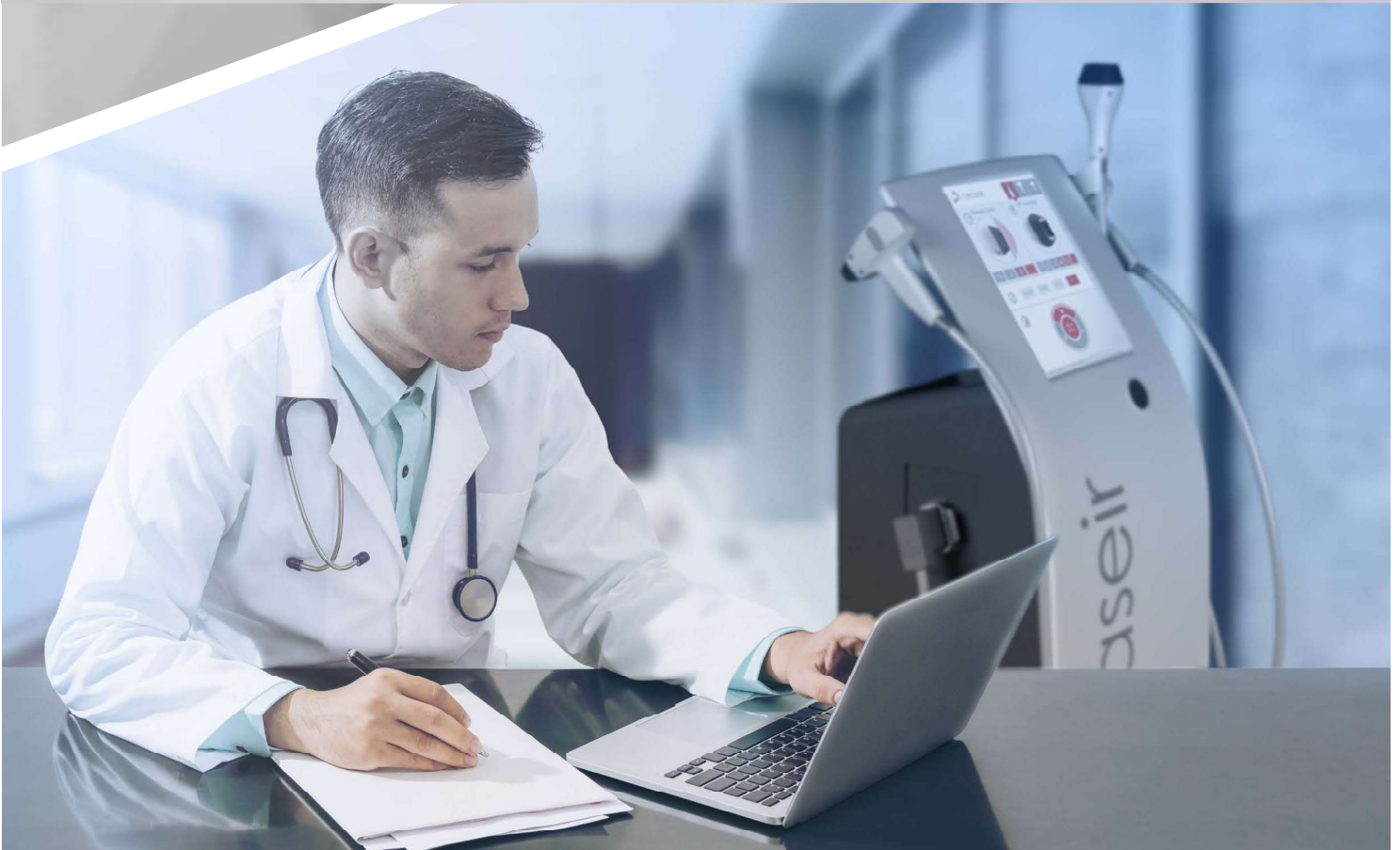

A quality device, easy-to-maintain and work that is controlled by a 15-inch touch interface.

Highly effective treatment with visible results from the first session.

It can be applied to all kinds of people and at any time of the year. Safe, non-invasive treatment without any adverse effects.

CELLULITE

**20 to 30
minute sessions**

**Cycle between 6-8
sessions depending
on the case**

**1 treatment per week
For maintenance
1 every 3-4 months**

FLACCID SKIN

**15 to 20
minute sessions**

**Cycle between 3-5
sessions depending
on the case**

**1 treatment every 3 weeks
For maintenance
1 every 3-5 months**

Features

DEVICE

DIMENSIONS

475mm × 758mm × 1133mm (W × L × H)

WEIGHT

38 Kg

INTERFACE

Touch LCD Display

FEATURES

RF INTENSITY

Up to 50 Watts (BS Applicator)
Up to 20 Watts (FS Applicator)
Up to 130 J/cm³ (FC Applicator)

RF MODE

Frequency Mode I - 0.8 MHz
Frequency Mode II - 1.7 MHz
Frequency Mode III - 2.45 MHz
Frequency Mode IV - 0.8 MHz, 1.7 MHz, 2.45 MHz

SUCTION - COOLING

BS Applicator - Suction when pressing
FS Applicator - Suction when pressing
FC Applicator - Active cooling

TREATMENT AREA SIZES

BS Applicator : 35mm X 30mm
FS Applicator: 15mm X 10mm
FC Applicator : 8mm X 8mm

ELECTRICAL REQUIREMENTS

90-264 VAC; 50/60Hz; Single Phase

VOLTAGE REQUIREMENTS

Maximum 700 Watts

OTHER INFORMATION

Insured with maintenance contract, after the first year

Consumables without programmed obsolescence

Own technical service

Designed and manufactured in Spain

THE MOST ADVANCED
RADIO FREQUENCY
TECHNOLOGY FOR
CELLULITE AND
FACIAL AND BODY
REMODELING



The professionals who integrate Leaseir Medical Light have been working at the forefront of aesthetic equipment technology for more than 18 years.

We are committed to the needs of professionals and are confident of providing the results your customers demand.

Our main objective is quality, which is why our devices have the best features to optimize results and ensure customer satisfaction.



Leaseir

CONTACT DETAILS

OFICINA CENTRAL

Edificio Asturias
Parque Científico y Tecnológico
C/ Jimena Fernandez de la Vega, 140
33203 Gijón (Asturias) - SPAIN
+34 984 390 933
info@leaseir.com