

# SIGNATURE SERIES

## Wood Laboratory Furniture



KEWAUNEE®

# Signature Series

Wood Laboratory Furniture



# Kewaunee Quality Construction

**KEWAUNEE'S SIGNATURE SERIES WOOD FURNITURE** is engineered and fabricated to exacting specifications. From raw material to finished products, our quality assurance program is designed to certify that Kewaunee products exceed established industry standards.

## DOVETAILED DRAWER CONSTRUCTION



**Drawer Sides** - ½" nine-ply birch plywood, dovetailed to drawer sub-front and back.

**Drawer Bottom** - ¼" hardboard set and glued into grooves on all four sides.

**Optional Drawer Bottom:** white melamine on hardboard or plywood.



## PULL STYLE OPTIONS



## DOOR HINGE OPTIONS



Full swing opening allows easier access to the cabinet interior. Hinges are rated for 150 lb. load.

## SIGNATURE SERIES WOOD FURNITURE CABINET FEATURES

- **QUALITY CONSTRUCTION** - Joints are multiple doweled, glued, and screwed. Top rails and solid bottom provide complete case rigidity.
- **RIDGED END PANELS** - ¾" hardwood plywood with front exposed edges faced with solid hardwood.
- **ONE-PIECE CABINET BACK** - ¼" hardboard, dadoed into end panels and securely fastened.
- **FULL DEPTH SHELVES** - Hardwood plywood, with ⅛" hardwood facing, matching cabinet interior. ¾" in Base cabinets - 1" in Wall and Tall cabinets.
- **FULL EXTENSION**, ball bearing, drawer slides on all drawers, pullboards, and pullout trays and pans. Load rated at 100 lb.
- **OPTIONAL ECO-FRIENDLY** - chain of custody FSC® certified plywood, NAUF plywood and particleboards.

## Signature Series Wood Cabinet Styles



### Style-1

#### Conventional Radius Lip


Conventional radius lip doors and drawers. Available in oak, maple, or beech with a wide selection of pull options.



### Style-5

#### Contemporary Full Overlay

Square-edged full overlay doors and drawers with vertical matching grain. Available in oak, maple, or beech with a wide selection of pull options.

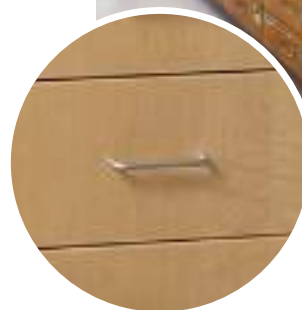



# Educational Assemblies



# Typical Installations





# Typical Installations



# SIGNATURE SERIES

## Wood Laboratory Furniture

Kewaunee Scientific Corporation is dedicated to manufacturing high quality products for the laboratory marketplace. We offer **Total Laboratory Solutions, Innovative Designs, and Technological Expertise** with a worldwide distribution network to provide laboratories that are truly World Class.

### Laboratory Products from Kewaunee®:

- ALPHA® – Adaptable Casework System
- ALPHA Columns – Expandable Workstations
- ALPHA Overhead Carrier Systems
- BasikBench® – Simple 4-Leg Benches or Tables
- Element Bench – Efficient Benches with Plumbed Posts
- Enterprise® – Free-Standing Workstations
- Evolution® – Column Based Workstations
- Kemresin® – Epoxy Resin Counters & Sinks
- Kewaunee Matrix® – Educational Stations
- Research Collection® – Steel Casework
- Signature® Series – Wood Casework
- Supreme Air® – Fume Hoods

For our complete product offering, specifications, catalogs, brochures, case studies, and dealer list, go to [kewaunee.com](http://kewaunee.com), or contact your local sales representative directly.



The mark of responsible forestry



### Kewaunee Scientific Corporation

P.O. Box 1842 • Statesville, NC 28687-1842

704-873-7202 • [kscmarketing@kewaunee.com](mailto:kscmarketing@kewaunee.com) • [kewaunee.com](http://kewaunee.com)

### Kewaunee Labway India Pvt. Ltd.

No. CA-9A, 2nd Floor, Jigani Link Road

Jigani, Bangalore 562 106, India

Phone: 9180-27826725 • Fax: 9180-27826724

Web Site: [www.kewaunee.in](http://www.kewaunee.in)

### Kewaunee Labway Asia Pte. Ltd.

194 Pandan Loop, #06-22

Pantech Business Hub, Singapore 128 383

Phone: 65-6773-0288 • Fax: 65-6773-2322

Kewaunee, Alpha, BasikBench, Enterprise, Evolution, Kemresin, Kewaunee Matrix, Research Collection, Signature, and Supreme Air are registered trademarks of Kewaunee Scientific Corporation

