

Re-Prime Tubing with the twiist™ System

The twiist Automated Insulin Delivery System allows you to re-prime your infusion set tubing with a few taps in your twiist App. You may need to use this function if you experience air bubbles in your tubing.

To start re-priming your infusion set tubing, tap the Pump Status icon on the Home Screen to open the twiist Pump Menu.

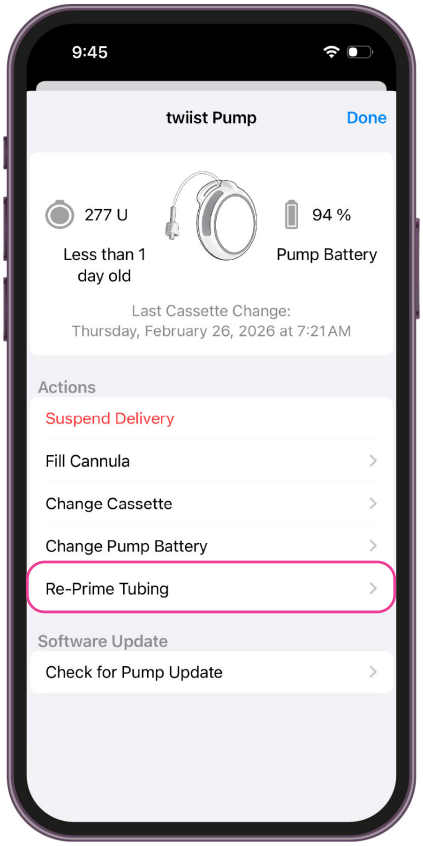

1 Select **Re-Prime Tubing** from the twiist Pump Menu.

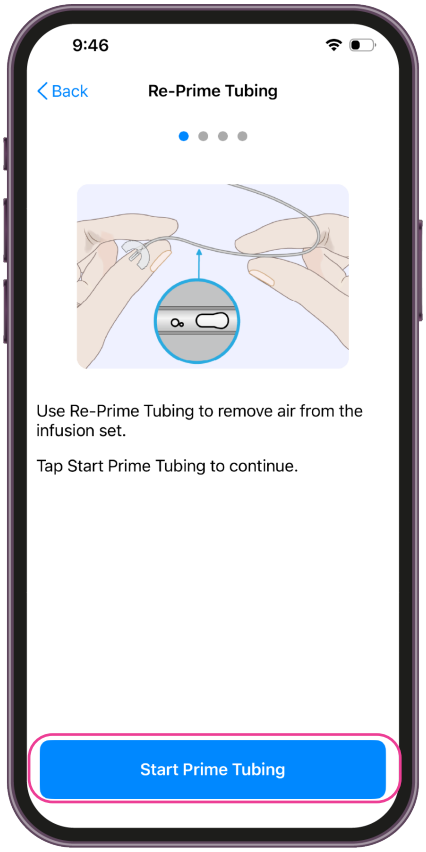

2 Select **Start Prime Tubing**.

Follow the step-by-step instructions in the twiist App to finish re-priming your infusion set tubing.

Do not re-prime infusion set tubing while connected to your infusion site. Connecting the infusion set tubing while re-priming may lead to the unintended delivery of insulin which may result in low blood glucose.

Customer Care: 1 877-4TWIIST (1-877-489-4478)

The twiist™ Automated Insulin Delivery (AID) System is for people with type 1 diabetes, 6 years of age and above, available with a prescription, only in the U.S.

© 2026 Sequel Med Tech, LLC. All rights reserved. All other trademarks are the property of their respective owners. Visit **twiist.com/safety** for important safety information.