twiist™ App Menu Map


The twiist App helps support diabetes management by putting pump controls, therapy settings, and notifications in one place. Use this guide to quickly find options in the **twiist Pump Menu** and **Settings Menu**.

twiist Pump Menu

The **twiist Pump Menu** has options to view information about your twiist System and perform pump-related tasks. It is accessed by tapping the **twiist Pump Status icon** on the twiist App home screen. Use the twiist Pump Menu whenever pump information needs to be reviewed or a pump-related action needs to be completed.

Settings Menu

The **Settings Menu** has options to review and update settings including Therapy Settings, Bolus Settings, Notification Settings, and CGM Settings. Here you can also turn Loop On and Off and use the **Contact Sequel** menu option to access training resources and videos for using the twiist System.

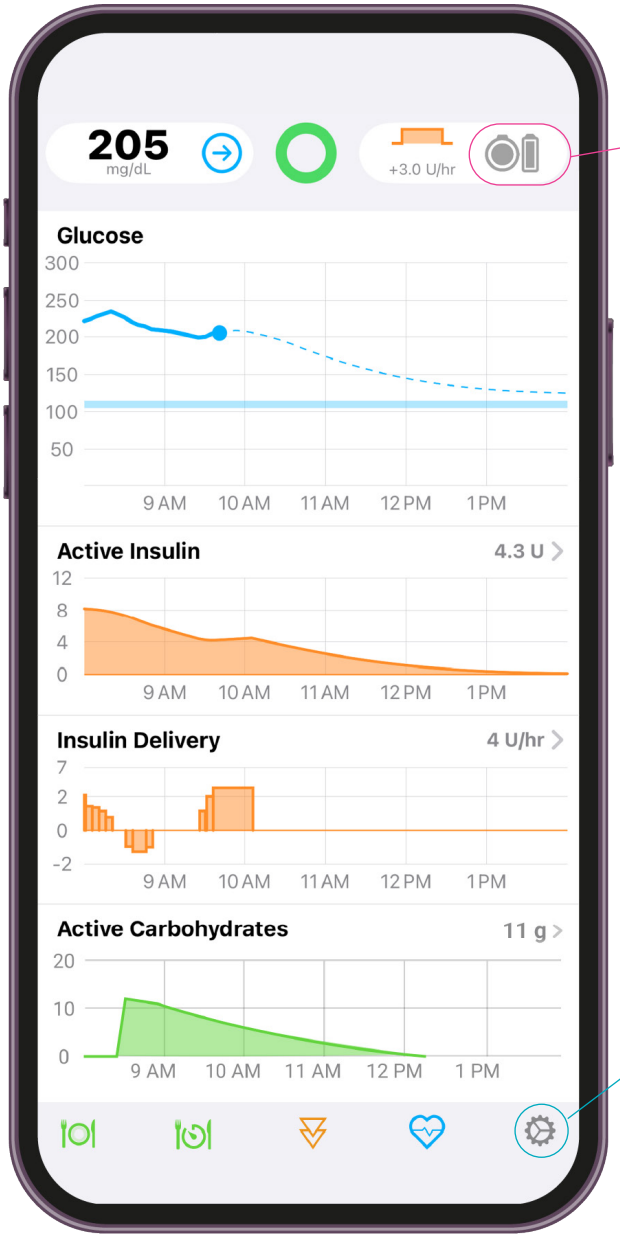

twiist Pump Menu


- Current Cassette Volume & Age
- Pump Battery Percentage
- Suspend Delivery
- Fill Cannula
- Change Cassette
- Change Pump Battery
- Re-Prime Tubing
- Check for Pump Update

Settings Menu

- Loop Toggle (On/Off)
- Therapy Settings 1
- Bolus Settings 2
- Notification Settings 3
- CGM Settings 4
- Help
- Contact Sequel 5
- Configure Devices 6
- About

TRNG-00553 Ver. 1.0

Settings Menu: The Submenus

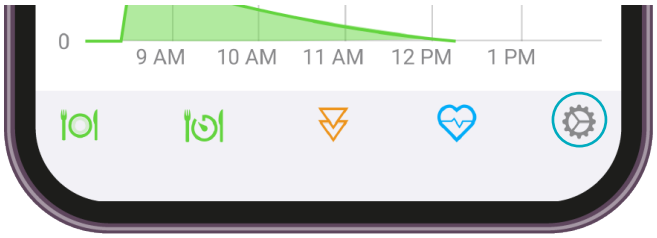

1 Therapy Settings

- Glucose Safety Limit (67-110 mg/dL)
- Android: Glucose Units (mg/dL or mmol/L)
- Correction Ranges (87-180 mg/dL)
- Pre-Meal Range (67-130 mg/dL)
- Workout Range (87-250 mg/dL)
- Carb Ratios (2-150 g/U)
- Basal Rates (0-30 U/hr)
- Delivery Limits (Maximum Basal Rate & Max Bolus)
- Insulin Model (Fast-Acting Adults & Fast-Acting Children)
- Insulin Sensitivities (10-500 mg/dL)
- Temporary Basal (Loop Off*)

2 Bolus Settings

- One-Button Bolus
- Extended Bolus (Loop Off*)
- Dual Bolus (Loop Off*)

3 Notification Settings

- iPhone: Critical Alert Volume
- Android: Do Not Disturb Override
- Pump Sounds and Quiet Mode
- Low Insulin Alert #1
- Low Insulin Alert #2
- Phone Out of Range
- Pump Inactive Alert
- Cassette Change Reminder

4 CGM Settings

- CGM Alerts
- Target Range Setting
- Start New Sensor (FreeStyle Libre 3 Plus)
- Pair Smart Transmitter (Eversense 365)

5 Contact Sequel

- Support and Training Resources & Videos

6 Configure Devices

- Delete Pump
- Delete CGM

*__Loop Off__: The twiist AID System is not automatically adjusting insulin delivery. Your pump will continue delivering your programmed basal insulin.

Customer Care: 1 877-4TWIIST (1-877-489-4478)

The twiist™ Automated Insulin Delivery (AID) System is for people with type 1 diabetes, 6 years of age and above, available with a prescription, only in the U.S.

© 2026 Sequel Med Tech, LLC. All rights reserved. All other trademarks are the property of their respective owners.
Medical disclaimer pending, to be provided by Sequel once approved.

Visit twiist.com/safety for important safety information.

TRNG-00553 Ver. 1.0