
CardioScan

Holter Monitor

What do the lights mean?


When the device is fitted on the patient

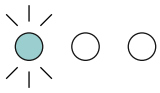


Green

The green light will flash **twice** every 2 seconds to indicate the device is receiving a good signal and the test is underway.

After a few minutes, the green light will flash **once** every 2 seconds when recording starts. It will continue to flash for the test duration.

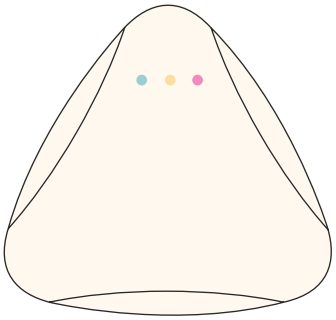

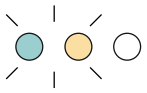


Green and Yellow

The green and yellow lights will flash **once** every 2 seconds to indicate the device battery level is low. The device will continue to record normally in this mode.



Red

When the red light remains permanently on, it indicates a device error. Please contact your local support team.

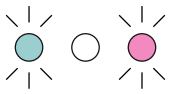

Green and Red

Flashing green and red lights indicate that the heart recording duration has been completed.




Green / Yellow / Red

The green, yellow, and red lights will flash for 2 minutes when the patient double-taps the device to record an event.

All 3 lights will flash for 1 to 2 seconds during a self-test.

When the device is being charged

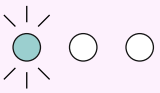

Green

Charging: The green light flashes to indicate that the device is charging.

Fully charged (100%): The green light stops flashing when the device is fully charged.

Need help?

Contact your local support team by visiting **cardioscan.co/contact**