

Certificate of Analysis

BTK (R28H), active

(Recombinant enzyme expressed in Sf21 insect cells)

Item # 14-765, 14-765-K, 14-765M

Parent Lot # 2067978

The data presented in this document apply to the parent lot shown above and to all pack sizes derived from subsequent vialling runs of this parent lot. An alphabetical suffix after the parent lot number is used to denote each vialling run.

Product Description: N-terminal 6His-tagged, recombinant, full-length, human BTK containing the R28H mutation. Expressed by baculovirus in Sf21 insect cells. Purified using Ni²⁺/NTA-agarose.

R28H is one of a number of mutations in BTK shown to cause X-linked agammaglobulinemia (XLA), an immune deficiency disorder characterized by a lack of mature, immunoglobulin-producing, peripheral B cells (Ohta *et al.*, (1994), Proc. Natl. Acad. Sci. USA, **91**, 9062-9066; Hyvonen & Saraste, (1997), EMBO J. **16**, 3396-3404). Purity 57.3% by SDS-PAGE and Coomassie blue staining. MW = 78.4kDa.

Specific Activity (Parent lot# 2067978): 21U/mg, where one unit of BTK (R28H), active activity is defined as 1nmol phosphate incorporated into 300µM LKBtide (LSNLYHQGKFLQTFCGSPLYRRR) per minute at 30°C with a final ATP concentration of 100µM.

Formulation: 0.441mg/ml of enzyme in 50mM Tris/HCl pH7.5, 300mM NaCl, 0.1mM EGTA, 0.03% Brij-35, 270mM sucrose, 0.2mM PMSF, 1mM benzamidine, 0.1% 2-mercaptoethanol. Frozen solution.


Storage and Stability: On receipt of material store at -70°C. Unopened reagent is stable for a minimum of 1 year from date of shipment when stored at recommended storage temperature. Avoid repeat freeze/thaw cycles. For maximum recovery of product, centrifuge original vial prior to removing the cap.

Handling Recommendations: Rapidly thaw the vial under cold water and immediately place on ice. Aliquot unused material into pre-chilled micro-centrifuge tubes and immediately snap-freeze the vials in liquid nitrogen prior to re-storage at -70°C.


**FOR IN VITRO RESEARCH USE ONLY
NOT FOR USE IN HUMANS OR ANIMALS**

Quality Control Testing

Kinase Assay: 25–551ng of this lot of enzyme phosphorylated 300µM LKBtide in the assay described on page two. Assay background was subtracted from the actual counts to yield the results shown below.

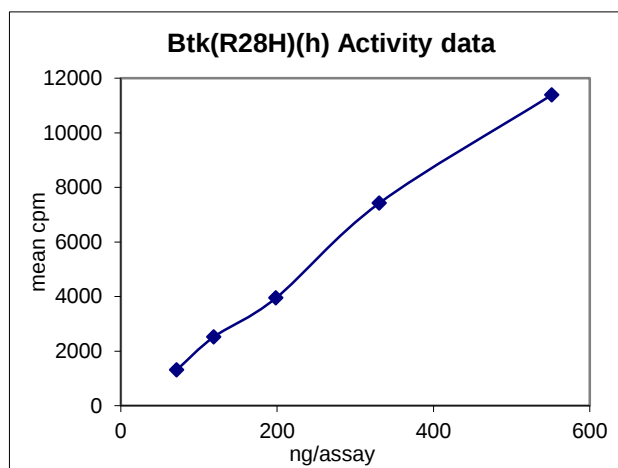

MS Tryptic Fingerprint: Confirmed identity as BTK with the translated sequence listed on page three.

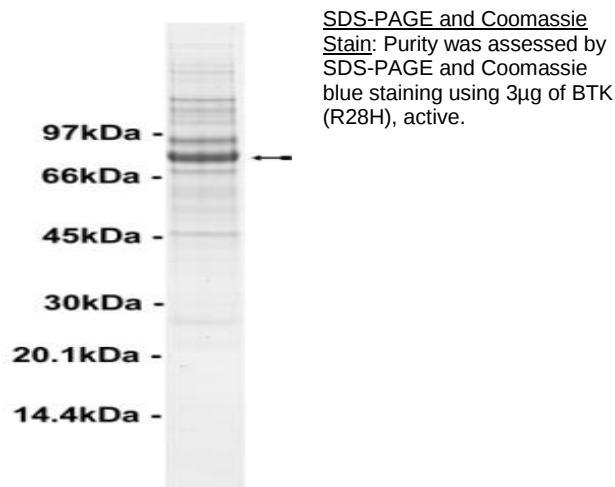

Certificate of Analysis

Kinase Assay Protocol

Stock Solutions:

1. **5 x Reaction Buffer:** 40mM MOPS/NaOH pH7.0, 1mM EDTA.
2. **LKBtide (LSNLYHQGKFLQTFCGSPLYRRR):** Use at a final assay concentration of 300 μ M. Prepare a 3mM stock and add 2.5 μ l of stock per assay point.
3. **NaCl:** Use at a final assay concentration of 50mM. Make a 3M stock. Add 0.42 μ l of stock per assay point.
4. **Triton X-100:** Make up a 10% (w/v) stock. Add 2.5 μ l of stock per assay point.
5. **BTK (R28H), active:** Dilute with 20mM MOPS/NaOH pH7.0, 1mM EDTA, 5% glycerol, 0.01% Brij-35, 0.1% 2-mercaptoethanol, 1mg/ml BSA. Use 25–551ng per assay point.
6. **[γ -³³P]ATP:** 2.5 x magnesium acetate/[γ -³³P]ATP cocktail: 25mM MgAc and 0.25mM ATP to which is added [γ -³³P]ATP (specific activity approximately 500 - 800cpm/pmol as required.)

Assay Procedure (96 well plate format):

1. Add 5 μ l of 5 x reaction buffer per well.
2. Add 2.5 μ l of **LKBtide**.
3. Add **2.5 μ l (25–551ng) BTK (R28H), active**.
4. Add 0.42 μ l of 3M NaCl
5. Add 2.5 μ l of 10% (w/v) Triton X-100.
6. Make up to 15 μ l with dH₂O.
7. Add 10 μ l of diluted [γ -³³P]ATP mixture.
8. Incubate for 10 minutes at 30°C.
9. Stop the reaction by adding 5 μ l of 3% phosphoric acid.
10. Transfer a 10 μ l aliquot onto the appropriate area of a **P30 Filtermat**.
11. Wash the filtermat three times for 5 minutes with 50mM phosphoric acid.
12. Wash the filtermat once for 2 minutes with methanol.
13. Transfer the filtermat to a sealable plastic bag and add 4ml of scintillation cocktail.
14. Read in a scintillation counter. Compare cpm of enzyme samples with cpm of control samples that contain all assay components plus 1 μ l of 30% phosphoric acid.

Certificate of Analysis

BTK (R28H) Sequence Information

<u>Protein</u>	human BTK (R28H)
<u>Tags</u>	N-terminal 6His
<u>Native sequence</u>	M16 of the recombinant protein is equivalent to M1 of human BTK (R28H)
<u>Accession number</u>	GenBank NM_000061

Recombinant BTK amino acid sequence:

```

1  MHHHHHHEFK GLRRRMAAVI LESIFLKRSQ QKKKTSPLNF KKHLLFTLVH KLSYYEYDFE
61  RGRRGSKKGS IDVEKITCVE TVVPEKNPPP ERQIPRRGEE SSEMEQISII ERFYPYFQVV
121 YDEGPLYVFS PTEELRKRWI HQLKNVIRYN SDLVQKYHPC FWIDGQYLCC SQTAKNAMGC
181 QILENRNGSL KPGSSHRKTK KPLPPTPEED QILKKPLPPE PAAAPVSTSE LKKVVALYDY
241 MPMNANDLQL RKGDEYFILE ESNLPWWRAR DKNQGEGYIP SNYVTEAEDS IEMYEWYSKH
301 MTRSQAELL QKEGKEGGFI VRDSSKAGKY TVSVFAKSTG DPQGVIRHYV VCSTPQSQYY
361 LAEKHLFSTI PELINYQHQN SAGLISRLKY PVSQQKNAP STAGLGYSW EIDPKDLTFL
421 KELGTGQFGV VKYGKWRGQY DVAIKMIKEG SMSEDEFIEE AKVMMNLSHE KLVQLYGVCT
481 KQRPIFIITE YMANGCLLNY LREMRHRFQT QQLLEMCKDV CEAMEYLESK QFLHRDLAAR
541 NCLVNDQGVV KVSDFGLSRY VLDDEYTSSV GSKFPVRWSP PEVLMYSKFS SKSDIWAQGV
601 LMWEIYSLGK MPYERFTNSE TAEHIAQGLR LYRPHLASEK VYTIMYSCWH EKADERPTFK
661 ILLSNILDVM DEES

```

Recombinant BTK nucleotide sequence:

```

1  atgcatcatc accatcacca tgaattcaaa ggcctacgtc gacgaatggc cgcagtgatt
61  ctggagagca tctttctgaa gcgatcccaa cagaaaaaga aaacatcacc tctaaacttc
121 aagaagcacc tgtttctctt gaccgtgcac aaactctcct actatgagta tgactttgaa
181 cgtgggagaa gaggcagtaa gaagggttca atagatgttg agaagatcac ttgtgttgaa
241 acagtggttc ctgaaaaaaa tcctcctcca gaaagacaga ttccgagaag aggtgaagag
301 tccagtgaat tggagcaaat ttcaatcatt gaaagggtcc cttatccctt ccaggttgta
361 tatgatgaag ggcctctcta cgtcttctcc ccaactgaag aactaaggaa gcggtggatt
421 caccagctca aaaacgtaat ccggtacaac agtgatctgg ttcagaaata tcacccttgc
481 ttctggatcg atgggcagta tctctgctgc tctcagacag ccaaaaatgc tatgggctgc
541 caaattttgg agaacaggaa tgggaagcta aaacctggga gttctcaccg gaagacaaaa
601 aagcctcttc cccaacgcc tgaggaggac cagatcttga aaaagccact accgcctgag
661 ccagcagcag caccagtctc cacaagttag ctgaaaaagg ttgtggccct ttatgattac
721 atgccaatga atgcaaatga tctacagctg cggaagggtg atgaatattt tatcttggag
781 gaaagcaact taccatggtg gagagcacga gataaaaatg ggcaggaagg ctacattcct
841 agtaactatg tcaatgaagc agaagactcc atagaaatgt atgagtggta ttccaaacac
901 atgactcgga gtcaggctga gcaactgcta aagcaagagg ggaaagaagg aggtttcatt
961 gtcagagact ccagcaaaagc tggcaaatat acagtgtctg tgtttgctaa atccacaggg
1021 gaccctcaag ggggtgatac tcattatgtt gtgtgttcca cacctcagag ccagtattac
1081 ctggctgaga agcacctttt cagcaccatc cctgagctca ttaactacca tcagcacaac
1141 tctgcaggac tcatatccag gctcaaatat ccagtgtctc aacaaaacaa gaatgcacct
1201 tccactgcag gcctgggata cggatcatgg gaaattgatc caaaggacct gaccttcttg
1261 aaggagctgg ggactggaca atttggggta gtgaagtatg ggaaatggag aggccagtac
1321 gacgtggcca tcaagatgat caaagaaggc tccatgtctg aagatgaatt cattgaagaa
1381 gccaaagtca tgatgaatct ttcccatgag aagctggtgc agttgtatgg cgtctgcacc
1441 aagcagcgcc ccatcttcat catcactgag tacatggcca atggctgcct cctgaactac
1501 ctgagggaga tgcgccaccg cttccagact cagcagctgc tagagatgtg caaggatgtc
1561 tgtgaagcca tggaatacct ggagtcaaag cagttccttc accgagacct ggcagctcga
1621 aactgtttgg taaacgatca aggagttgtt aaagtatctg atttcggcct gtccagggtat

```

Certificate of Analysis

```
1681 gtcctggatg atgaatacac aagctcagta ggctccaaat ttccagtccg gtggtcccca
1741 ccggaagtcc tgatgtatag caagttcagc agcaaactctg acatttgggc ttttgggggtt
1801 ttgatgtggg aaatttactc cctggggaag atgccatatg agagatttac taacagtgag
1861 actgctgaac acattgcca aggcctacgt ctctacaggc ctcatctggc ttcagagaag
1921 gtatatacca tcatgtacag ttgttggcat gagaaagcag atgagcgtcc cacttttaaa
1981 attcttctga gcaatattct agatgtcatg gatgaagaat cctga
```

Reviewed and approved by site quality representative.

Unless otherwise stated in our catalogue or other company documentation accompanying the product(s), our products are intended for research use only and are not to be used for any other purpose, which includes but is not limited to, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses or any type of consumption or application to humans or animals.

© 2014 **Eurofins Pharma Discovery Services UK Limited** is an independent member of Eurofins Discovery Services