Voyager III



The Beauty of Dentistry

FREE 5 Year Warranty

for Voyager Dental Chairs and Dental Units, and Dental Operating Lights. Visit belmontdental.co.uk or refer to the back cover for further details.

Voyager III

Key improvements

Improved vacuum air-flow
Chair base and delivery arm re-designed
Chair ergonomics and comfort improved
700 Series light (720) with touchless sensor available as standard

Key features

Ambidextrous
Ideal console arrangement for Assistant's 'Preparation zone'
Below-the-Patient swing-arm delivery system
Mobile delivery systems options
Seamless upholstery
Standard or Ultrasoft upholstery colour options

Air operating system

Upgradeable Operator's Console functionality

Unit-mounted Operating Light options

720SLR or 300 LED

Upgrading and accessorising

Belmont's range of factory built-in options can be used to increase the functionality of a treatment centre and the services that your Practice offers to the Patients. Visit belmontdental.co.uk or see our price list for further details.

Specification

135kg max. Patient weight | 480mm min. seat height | 780mm max. seat height




Very versatile and compact ambidextrous treatment centres, designed and built for the most demanding working environments.



Voyager III ambidextrous treatment centres

Understated versatility & value

Voyager III's versatile Below-the-Patient delivery system and ultra-slim backrest design eliminate physical stresses associated with some treatment centres.

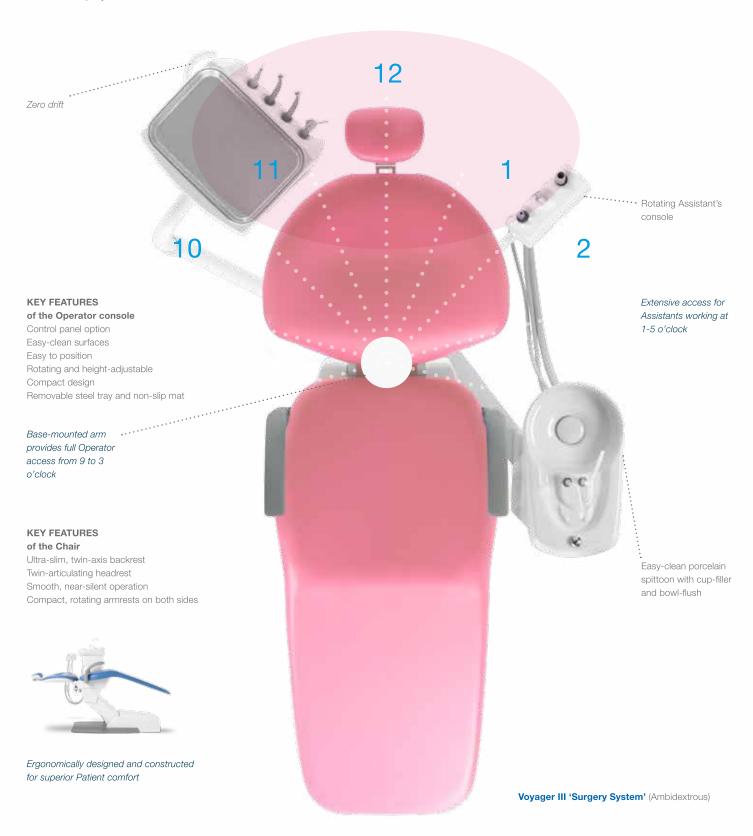



Left- or right- anytime

Switchover in no time at all

The central pivot design enables the very easy changeover procedure and provides 100% equality for left- and right-handed Operators.



Right-handed set up.
To begin the switchover, raise chair height.



Incline chair.

1-2. Swing table round to left side.



3-4. Swing cuspidor round to right side.



5. Secure cuspidor in position.Secure cuspidor locking button.2. Rotate console.Left-handed set up.




Below-the-Patient advantages

Ready to welcome the next Patient, instruments are out of sight and preparation has begun.



Be more efficient

Intelligently designed for a comfortable and compact, easy-clean working environment.

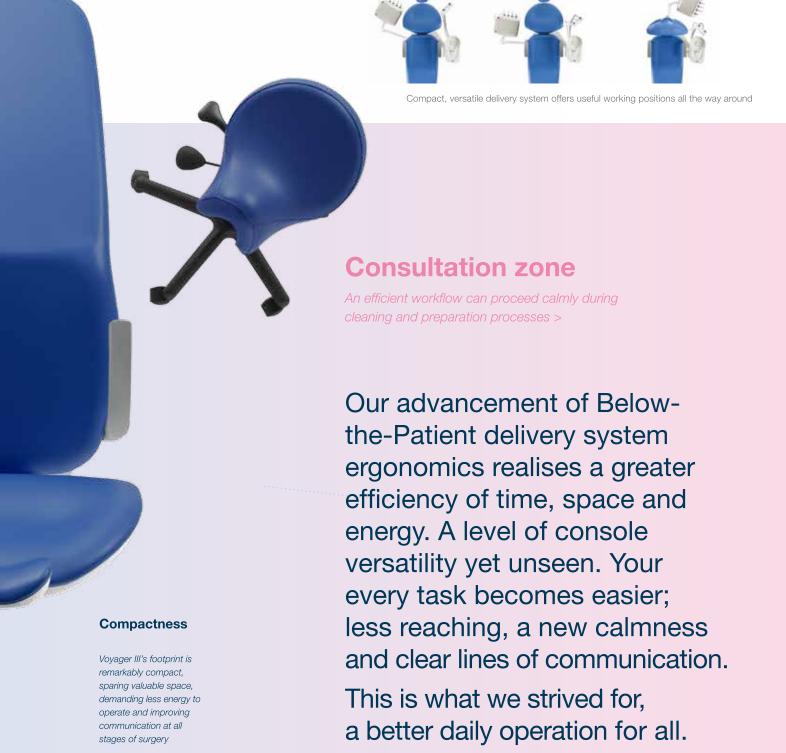

Control is simple

320 LED Operating Light Control chair movements and instruments at the foot-controls. Voyager III is an uncluttered, simple and reliable working partner which allows the Operator full control with the option of mounting a control panel on the Operator's console. Foot control for instruments and chair

Flexible options to suit

Voyager III 'SP' ambidextrous chair and cuspidor



Operator's control panel offers two presets, Last-Position Memory, Auto-Return and manual control for chair movements (optional)





Voyager 'MC' Mobile Cart

Voyager 'M' Module



720SLR for Voyager III

KEY FEATURES

Unit-mounted for Voyager III White natural beam (4,200 Kelvin) 3 light settings: Low (18,000 Lux), High (28,000 Lux) and Composite Curing Safe Mode (at 8,000 Lux) 650mm focal distance Touchless sensor switch Patient mirror Auto ON/OFF with preset chair positions Halogen bulb Twin-axis head movement





650mm focal distance and the very compact light head allow extensive working space for the dental team

Light intensity button and indicator: soft-push-up button is easily and hygienically activated with the back of a hand



320 LED for Voyager III

KEY FEATURES

Unit-mounted for Voyager III
White natural beam (5000 Kelvin)
Stepless intensity control:
3,100 to 28,000 Lux
650mm focal distance
Composite Curing Safe Mode (at 4,300 Lux)
Touchless sensor switch
Auto ON/OFF with preset chair positions*
5 LEDs
Triple-axis head movement







Unique stepless control allows adjustment to any light intensity between 3,100 and 28,000 Lux



5 Our warranties are extended and FREE

On condition that the purchaser returns the Extended Warranty form to Takara Belmont (UK) Ltd, the Standard Warranty will be extended to ~ **5 Years** for Voyager Dental Chairs and Dental Units, and Dental Operating Lights. **2 Years** for Intra-Oral X-Ray equipment.

Standard Warranty Policy – All Belmont products are guaranteed against mechanical, electrical and workmanship faults, under normal working conditions, for a period of one year from date of delivery. Paintwork, plating, natural wear/tear of all tubings are not warranted unless they are proved to be defective at time of delivery. OEM factory built-in options (i.e. scalers, motors and CVS kits) carry a warranty as provided by the individual manufacturer.

(]Belmont

belmontdental.co.uk
@BelmontDentalUK
London showroom 020 7515 0333
Manchester showroom 0161 743 9992
All showroom visits by appointment only

The manufacturers reserve the right to alter the specifications without prior notice.





