



**Remote Dispense Option Shown**

D

C

B

A

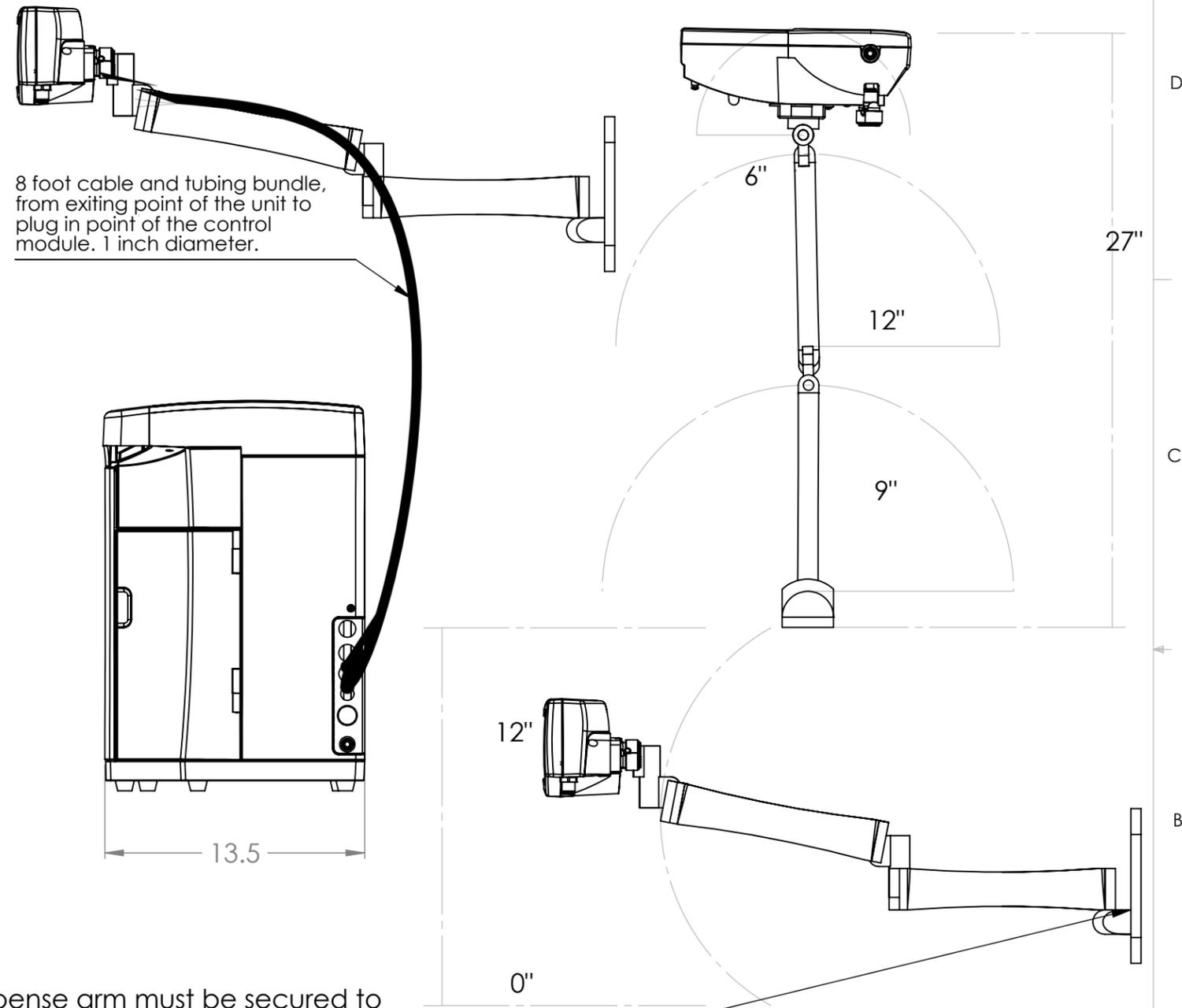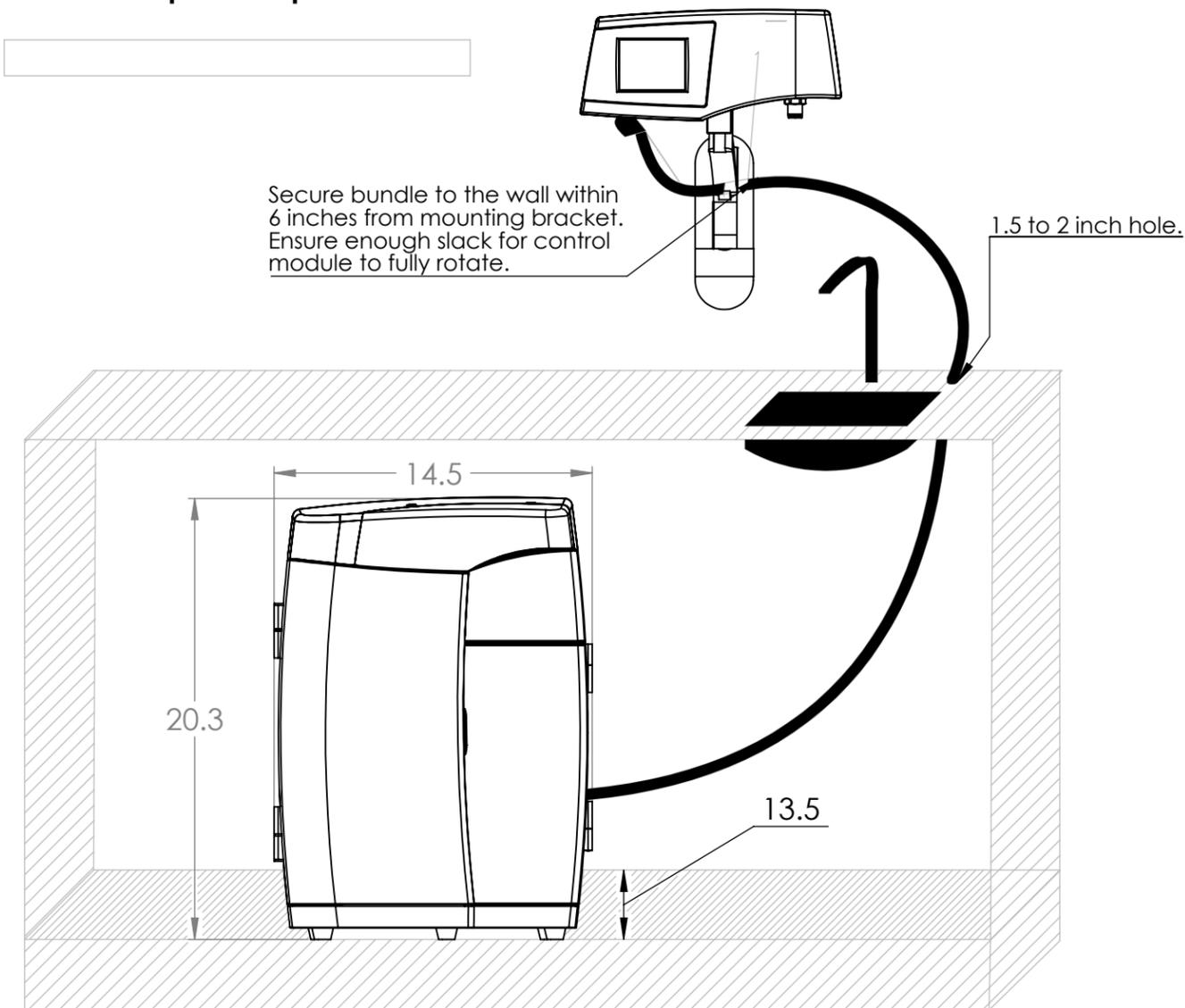

Remote dispense arm must be secured to a structural stud. If no stud available, a proper backing should be mounted prior to, such as a plywood sheet.

**Requirements:**

- GFI within 6 feet
- Feed water within 6 feet
- Drain within 6 feet

**NOTES:**

- **WEIGHT** - 42 LBS. DRY / 47 LBS. OPERATING
- **INLET CONNECTION** - 3/8" O.D. TUBING
- **POWER REQUIREMENTS:**
  - "-1" - 120 VAC / 60 HZ @ 1.0 AMP
  - "-2" - 220VAC / 50 HZ @ 0.5 AMP
- **PRESSURE** - 90 PSIG MAXIMUM / 30 PSIG MINIMUM
- **TEMPERATURE (MAX)** - 100 DEGREES F.
- **FLOW RATE** - 2.5 LPM (< 2.0 LPM WITH ULTRAFILTER AND/OR ENDOTOXIN FILTER INSTALLED)

	NAME	DATE	RESINTECH INC CAMDEN, NJ (856) 768-9600	
DRAWN	R DURBANO	12/06/16	TITLE: CLiR - 5000 SERIES ULTRA HIGH-PURITY WATER SYSTEM DIMENSIONS	
DIMENSIONS ARE IN INCHES UNLESS OTHERWISE SPECIFIED			DWG. NO. REMOTE CONFIGURATION	REV 2
SCALE: 1:7			SHEET 2 OF 2	

