




MEDICAL FACE MASK

TIES • HIGH FILTRATION • GREY

ASTM LEVEL 1

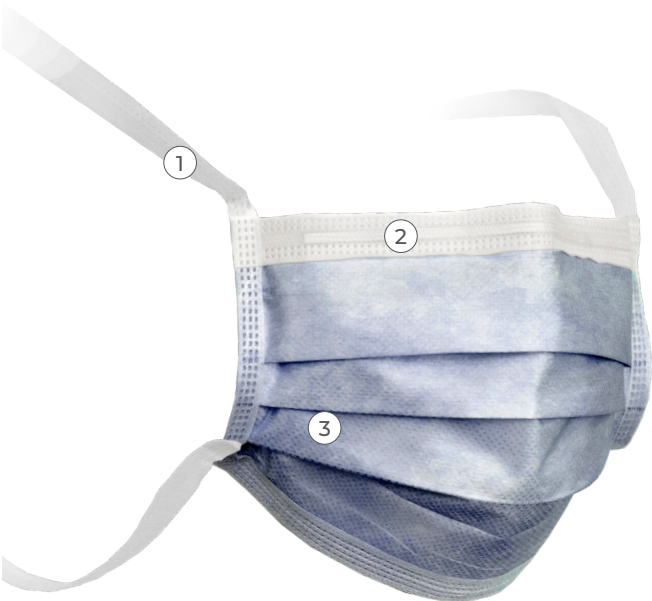

10317TM

For elevated protection during surgical cases with exposure to laser plume, this ASTM Level 1 high filtration mask offers fluid resistance with enhanced filtration efficiency. Crafted with premium filter material and a breathable inner layer, clinicians can feel both protected and comfortable through even the most rigorous procedures.

BFE	PFE	DELTA P $\Delta P(\text{mmH}_2\text{O}/\text{cm}^2)$	FLUID RESISTANCE	FLAME SPREAD
> 95%	> 80%	< 5.0	80 mm Hg	Class 1
MASKS/BOX		BOXES/CASE		MASKS/CASE
50		6		300

① Premium Tie Material and Design

- Tie material that is comfortable on the skin throughout extended wear.
- Long ties create a flexible fit while tying over top of a bouffant or surgical cap.

② High Quality Nosepiece with Ultra Soft Foam

- Aluminum nosepiece forms smoothly across the nose and cheekbones.
- Remains securely in place and reduces fogging of eyewear.

③ High-Performing Filter Material

- Greater filtration efficiency to reduce the risk of inhaling potentially hazardous laser plume during procedures that use surgical lasers and electrosurgical units.

