
4 Return

Place the tube(s) into the biohazard bag and seal.

There is an absorbent blotter inside of the biohazard bag in case of leakage—please do not remove it.

Insert the completed sample information card into the specimen collection kit along with the sealed biohazard bag containing the tube. Close the box.

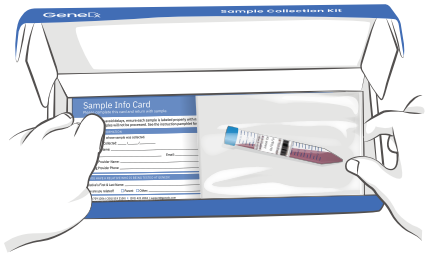

Place the Skin Punch Biopsy Kit inside the FedEx Pak, remove the adhesive strip, and seal the FedEx Pak. Attach the FedEx shipping label to the outside of the Pak.

Drop off the FedEx Pak at a local FedEx location. You can find the nearest FedEx location to you by calling FedEx or visiting FedEx.com/locate.



If you are not able to access a FedEx drop-off location, please contact **GeneDx** for help scheduling a pickup.

Important

Before sealing and shipping the completed kit, make sure to include:

- ☐ Sample labeled with name, date of birth, and date of collection
  Unlabeled samples will not be processed.
- ☐ Completed sample info card

GeneDx

205 & 207 Perry Parkway, Gaithersburg, MD 20877
T 1 (888) 729–1206 or 1 (301) 519–2100
F 1 (201) 421–2100
support@genedx.com • www.genedx.com

©2026 GeneDx, LLC. All Rights Reserved 18988

GeneDx.com | MKT-00301-v1i Jan 2026

GeneDx

Skin Punch Biopsy Kit

Before you begin, please read through all instructions.

Important

The media is light sensitive; please leave all contents inside the box and minimize light exposure to the media until use (e.g. prolonged exposure to fluorescent light or direct sun light). The media does not need to be refrigerated.

Once a tissue specimen is collected, if it is not being shipped immediately, place it in the refrigerator until shipment.

Your kit includes the following:

- FedEx Clinical Pak with UN3373 label
- FedEx return shipping label (addressed to 7331 Calhoun Place, Rockville, MD)
- Instruction brochure
- Biohazard bag with absorbent sheet
- One 4”x4” sheet of Parafilm
- GeneDx specimen label sheet with 5 labels
- GeneDx sample information card
- One 15 mL conical tube containing 10 mL of BenchStable RPMI media


*Test requisition form not included

1 Collection

One or two 4 mm skin punch biopsies should be collected. The tissue should be placed into tube media and submerged.


Place the lid(s) back on the tube(s) and close tightly. Please ensure that each cap is screwed on evenly and completely to avoid any leaking. Wrap each tube cap tightly with parafilm to prevent leaking during transit.

Samples in the media should be sent at room temperature the same day as the collection.

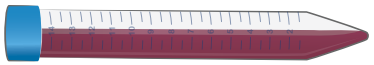

2 Label

Fill out the remaining provided labels with the full name and date of birth of the person whose sample is contained in each tube, as well as the date each sample was collected. Place the labels onto each of the tubes.
Reminder: All tubes must be labeled. Unlabeled samples will not be processed.

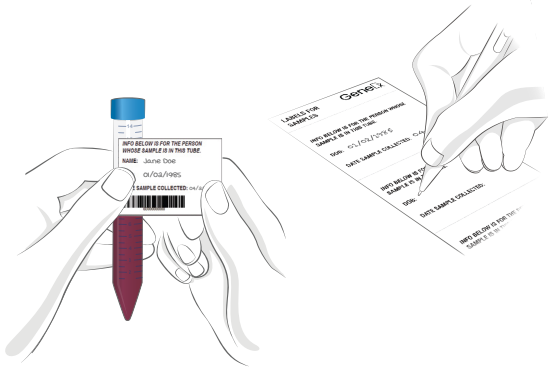

3 Sample Information

Reminder: Every sample MUST be accompanied by a completed sample info card. Samples submitted without this card will not be processed. Please indicate the tissue type on the sample information card.

