

Ophthalmic Res , DOI: 10.1159/000534098

Received: February 17, 2023

Accepted: September 8, 2023

Published online: September 27, 2023

Diagnostic accuracy of Automated Diabetic Retinopathy Image Assessment Softwares: IDx-DR and MediosAI

Grzybowski A, Rao DP, Brona P, Negiloni K, Krzywicki T, Savoy FM

ISSN: 0030-3747 (Print), eISSN: 1423-0259 (Online)

<https://www.karger.com/ORE>

Ophthalmic Research

Disclaimer:

Accepted, unedited article not yet assigned to an issue. The statements, opinions and data contained in this publication are solely those of the individual authors and contributors and not of the publisher and the editor(s). The publisher and the editor(s) disclaim responsibility for any injury to persons or property resulting from any ideas, methods, instructions or products referred to the content.

Copyright:

This article is licensed under the Creative Commons Attribution-NonCommercial 4.0 International License (CC BY-NC) (<http://www.karger.com/Services/OpenAccessLicense>). Usage and distribution for commercial purposes requires written permission.

© 2023 The Author(s). Published by S. Karger AG, Basel

Research Article

Diagnostic accuracy of Automated Diabetic Retinopathy Image Assessment Softwares: IDx-DR and MediosAI
Andrzej Grzybowski¹, Divya Parthasarathy Rao², Piotr Brona³, Kalpa Negiloni⁴, Tomasz Krzywicki⁵, Florian M Savoy⁶

¹Institute for Research in Ophthalmology, Foundation for Ophthalmology Development, Poznan, Poland

²Department of AI R & D, Remidio Inc, Glen Allen, Virginia, USA

³Department of Ophthalmology, Poznan City Hospital, Szwajcarska 3, Poznan 60-285, Poland

⁴Department of Clinical Research, Remidio Innovative Solutions Pvt Ltd, Bangalore, India

⁵Department of Mathematical Methods of Informatics, University of Warmia and Mazury, Olsztyn, Poland

⁶Department of AI R & D, Medios Technologies, Remidio Innovative Solutions, Singapore

Short Title: Performance of two DR screening AI softwares

Corresponding author:

Divya Parthasarathy Rao

Department of AI R & D, Remidio Innovative Solutions Inc.,

11357 Nuckols Rd, #102, Glen Allen, Virginia, 23059, USA

Tel +1 (855) 513-3335

Email drdivya@remidio.com

Number of Tables: 2

Number of Figures: 1

Word count: 3018

Keywords: Diabetic Retinopathy, AI, Medios, IDx-DR

Abstract

Introduction: Numerous studies have demonstrated the use of Artificial Intelligence for early detection of referable diabetic retinopathy (RDR). A direct comparison of these multiple Automated DR Image Assessment Softwares (ARIA) is however challenging. We retrospectively compared the performance of two modern ARIAs, IDx-DR and Medios AI.

Methods: In this retrospective-comparative study, retinal images with sufficient image quality were run on both ARIAs. They were captured in 811 consecutive patients with Diabetes visiting diabetic clinics in Poland. For each patient, four non-mydratic images, 45-degree field of view i.e two sets of one optic disc and one macula-centered image using Topcon NW400 were captured. Images were manually graded for severity of DR as no DR, any DR (mild NPDR or more severe disease), RDR (moderate NPDR or more severe disease and/or clinically significant diabetic macular edema (CSDME)) or sight-threatening DR (severe NPDR or more severe disease and/or CSDME) by certified graders. The ARIAs output was compared to manual consensus image grading (reference standard).

Results: On 807 patients, based on consensus grading, there was no evidence of DR in 543 patients (67). Any DR was seen in 264 (33%) patients, of which 174 (22%) were referable DR, and 41 (5%) sight-threatening DR. The sensitivity

of detecting RDR against reference standard grading was 95% (95%CI 91, 98%) and the specificity was 80% (95%CI 77, 83%) for Medios AI. They were 99% (95%CI 96, 100%) and 68% (95%CI 64, 72%) for IDx-DR respectively.

Conclusion: Both the ARIAs achieved satisfactory accuracy, with few false negatives. Although false-positive results generate additional costs and workload, missed cases raise the most concern whenever automated screening is debated.

Introduction

Diabetes is a global epidemic and one of the world's fastest-growing diseases. The number of patients with diabetic retinopathy (DR) and sight-threatening DR is also expected to rise. There are only a few established nationwide DR screening programmes and overall DR screening services remain inadequate in most of the developing world and even some developed countries [1]. This is further compounded by the increasing resources needed for the implementation and maintenance of comprehensive DR screening programs [2].

One of the proposed solutions to this global issue is the use of automated diabetic retinopathy image assessment software (ARIA) to grade fundus images instead or alongside human graders. There are multiple ARIAs currently available with many more being developed worldwide [1]. Although there is an abundance of studies looking into the performance of a single ARIA, studies comparing multiple ARIAs are currently rare, as direct comparison is often difficult [3]. Based on previous studies it is clear that the performance of even state-of-the-art algorithms may vary considerably [3,4]. We set out to analyze the performance of two modern ARIAs, IDx-DR and MediosAI.

Materials and Methods

Study Design: In this retrospective comparative study, the performance of two different ARIAs in screening for DR were compared to human graders (reference standard). The screening for DR was conducted and retinal images were obtained from diabetic clinics in Poznan, Poland between March 2020 and April 2021. The Institutional Review Board (Ophthalmology 21, Foundation for the Advancement of Ophthalmology) waived the need for IRB approval and written informed consent from participants for this retrospective study. The study was in adherence to the tenets of the Declaration of Helsinki. All the extracted images were anonymized, and no change in the clinical pathway was anticipated.

The primary outcome of the study was to assess the sensitivity and specificity of ARIAS in detecting referable Diabetic Retinopathy (RDR). The secondary outcomes were to assess the positive & negative predictive values of ARIAS to detect RDR and to assess the sensitivity of ARIAS in detecting sight-threatening diabetic retinopathy.

Sample size: Using an alpha error of 0.05, a precision rate of 10% (two sided), an estimated sensitivity of 85%, and an estimated incidence of RDR (International Clinical Diabetic Retinopathy (ICDR) -Moderate Non-proliferative DR (NPDR) and/or presence of clinically significant diabetic macular edema (CSDME)) to be 7%, the sample size calculated was 700 participants. Given these assumptions and expecting that 10% of subjects may be qualified as insufficient quality, a sample size of 800 subjects was chosen.

Inclusion & Exclusion criteria: The retinal images of subjects with established diabetes mellitus that were captured at the time of DR screening were included. Those that did not have at least one disc and one macula-centred image of sufficient quality were excluded from the study. Additionally, subjects who received treatment for DR (lasers or intraocular injections) were excluded.

Retinal Image Acquisition: The screening process involved undilated fundus images captured using a Topcon camera Nw-400 by trained operators who followed a specific imaging protocol. For each patient, a total of four images (45 degrees field of view each) were captured. They included one image centred on the optic disc and one centred on the macula for each eye. Additional images were taken to ensure sufficient quality. Retinal images were obtained from 811 consecutive patients with established diabetes mellitus who underwent screening for DR. Images deemed of sufficient quality graded by the IDx-DR AI software were selected. A total of 3200 sufficient quality images from 811 patients were used for the study.

Reference standard grading: The patients with images of sufficient quality were split into two sub-datasets. 362 patients were graded by three Polish retina specialists and 491 by three certified graders in India. All the graders are masked to the output of the AI and to each other's grading. Images were graded for severity of diabetic retinopathy based on International Clinical Diabetic Retinopathy (ICDR) severity classification as no DR, mild NPDR, moderate NPDR, severe NPDR and Proliferative Diabetic retinopathy (PDR). Macular edema was determined by the presence of surrogate markers like hard exudates. If hard exudates were found within 1 DD of the fovea, macular edema was determined as significant and labelled as clinically significant diabetic macular edema (CSDME) present. Image grading was done on a per eye basis. The final diagnosis for each patient was determined by the stage of DR of the more affected eye. Consensus image grading was regarded as the final reference standard based on Polish and Indian graders for the comparison of both AI systems. All the analysis was performed at the patient level.

Definitions: Referable diabetic retinopathy was defined as moderate NPDR and more severe disease (moderate NPDR, severe NPDR, PDR) and/or the presence of CSDME. Sight-threatening diabetic retinopathy (STDR) was defined as severe NPDR and more severe disease (severe NPDR, PDR), and/or the presence of CSDME

Artificial Intelligence (AI) Analysis using automated grading systems: We used two different Automated Diabetic Retinopathy Image Assessment Softwares (ARIAs) i.e., Medios AI for DR (Medios Technologies, Remidio Innovative Solutions, Singapore) and IDx-DR (Digital Diagnostics, Iowa, USA). The retinal images were run on both the ARIAs to screen for DR. Both the systems processed images deemed as sufficient quality by the IDx-DR system. IDx-DR results were recorded during live screening and all images captured for the patient were analyzed on a per patient basis. Two images per eye that passed the AI quality check were submitted to the AI for DR analysis. For Medios AI analysis, anonymized images for each patient were securely transferred to a cloud platform and the images were analyzed on an automated script version of the AI on a server instead of a manual analysis through the standard iPhone app deployment.

Both the AI systems are based on Convolutional Neural Networks (CNN), with the Medios system being based on the Inception-V3 architecture. Detailed description of the model is provided in the literature [5]. In brief, the Medios AI algorithm evaluated two possible outputs: “no signs of DR detected” (non-referable DR), “signs of DR detected” (referable DR). Report was generated on a per patient basis. The IDx-DR system also has an image quality and a diagnostic algorithm. The IDx-DR system outputs the stage of DR and generates a per patient report.

Statistical analysis:

All data was stored in Microsoft Excel sheets and Apache Parquet files and was analyzed using R and Python programming languages along with Numpy, Pandas, Scikit learn and Scipy libraries. The diagnosis of the AI using Medios and IDx-DR AI systems were tabulated against the consensus image diagnosis (reference standard) by constructing 2×2 tables. The sensitivity, specificity, positive predictive value (PPV) and negative predictive value (NPV) with 95% CIs were calculated. Inter-rater agreement for Polish and Indian graders was measured by calculating the kappa statistic.

Results

The study included the images of 811 patients. Image quality analysis was evaluated as part of the clinical workflow using the IDx-DR AI system. Four patients deemed ungradable by the graders were excluded. In total, 807 patients were included for further analysis. An additional two patients were removed from STDR analysis as they did not have a consensus for a STDR diagnosis despite being labelled as RDR by consensus. Figure 1 presents the STARD diagram of retinal image selection in the study.

Grades from the Polish and Indian graders were converted to No DR, any DR, RDR and STDR before computing consensus. Based on consensus grading, there was no evidence of DR in 543 patients (67%). Any DR was seen in 264 (33%), of which 174 (22%) were referable DR, and 41 (5%) sight-threatening DR. The inter-rater agreement (Cohen's kappa) for Poland graders was 0.679 (Ophthalmologist 1), 0.904 (Ophthalmologist 2) and 0.848 (Ophthalmologist 3). For the Indian graders, kappa was 0.632 (Ophthalmologist 1), 0.916 (Ophthalmologist 2) and 0.87 (Ophthalmologist 3).

IDx-DR AI system gives an output at a stage level. 567 patients (567/807; 70.3%) were flagged positive for DR (mild, moderate or threaten). 374 of them were also categorized as RDR (46.3%, moderate or threaten) and 189 as STDR (23.4%, threaten). Out of 174 patients with ground truth labeled as RDR, the IDx-DR system detected 172. This translates to a sensitivity of 99% (95% CI 96%, 100%). The Medios AI gives a binary output. It detected the presence of Referable Diabetic Retinopathy in 291 patients (36.1%). It correctly identified 166 of the 174 patients with a ground truth diagnosis of RDR. This translates to a sensitivity of 95% (95% CI 91%, 98%). The diagnostic ability of both the AI systems including sensitivity, specificity, positive and negative predictive values are tabulated in Table 1 and 2.

Discussion

With the rising burden of DR, the importance of early detection and screening cannot be overstated. To address this glaring need, advanced technologies like AI software have emerged as promising tools for DR screening. These AI systems are meticulously developed and optimized using diverse datasets. Before implementing such AI software in real-world scenarios, it is crucial to conduct comparisons among different solutions available. In our evaluation, we examined the performance of two automated DR Image Assessment Softwares: Medios AI and IDx DR. The results exhibited comparable performance in terms of sensitivity, with Medios AI achieving 95% and IDx-DR achieving 99% in identifying referable DR respectively. These findings underscore the potential of these software solutions in facilitating early detection and screening of DR.

Overall, the prevalence of DR in the sample analyzed was 33% for any DR and 22% for RDR, significantly higher than commonly reported in other studies. Scottish National Diabetic Retinopathy Screening Programme reported rates of RDR between 4.3% and 7%, large primary-care-based screening in California reported RDR rate of 8.2% and a hospital-based study in Ethiopia found any DR rate of 18.9% [6,7,8]. It is also much higher than previous estimates for DR prevalence in Poland [9]. This is likely a side-effect of the original screening set-up. The screening is based around diabetic clinics and diabetes medical centres, therefore selecting for a higher-risk population with other diabetic complications or difficult to control disease. A similarly high prevalence of DR was found in a study of 297 patients attending a tertiary center for diabetes care in India with DR prevalence of 40.8% [5].

Only patients who initially had images of sufficient quality for IDx-DR during the initial screening were included in this study. The IDx-DR image quality assistant process was used as part of the original screening program and is not evaluated herein. Out of 811 patients deemed gradable by IDx-DR only 4 (less than 0.5%) were excluded by the manual graders indicating that overall, the images selected for this study have good image quality.

The accuracy measures for Medios AI are in line with previously published studies. Natarajan et al reported accuracy of the Medios AI offline, smartphone-based version, with sensitivity and specificity pairs of 100% and 88.4% for RDR and 85.2% and 92.0% for any DR [10]. In the aforementioned study based in a tertiary diabetes center, Medios AI achieved 98.8% and 86.7% sensitivity and specificity for any DR [5]. In another India based study of 900 prospectively included patients Medios AI achieved 83.3%, 95.5% sensitivity and specificity for any DR and 93% and 92.5% respectively for RDR [11].

Crucially, all of the above-mentioned studies were done using images gathered with the Remidio FOP mobile smartphone-based camera in contrast to using a stationary, full-size automatic fundus camera for this study. Previous studies describing Medios AI were smartphone-based, with the algorithm app being run on a smartphone, also used to take the fundus pictures. This is the first study outside of India to investigate using Medios AI with images from a stationary fundus camera in a real-world screening scenario. Images captured with different cameras may differ in resolution, level of detail, contrast, noise and other parameters that may influence the accuracy of an algorithm. It is unclear whether the software or human graders may benefit from higher resolution images and provide a more robust golden-standard, and if so, how significant the difference is. For this study the previous smartphone-based results obtained by Medios AI seem to translate into comparable accuracy when using dedicated stationary fundus camera. These results are similar to another study where Medios AI was evaluated on Topcon images in an Indian population. This demonstrates generalizability of the model performance on a desktop system. This device agnostic

approach is particularly useful in screening programs that have already invested in camera systems and would want to move towards an AI-based approach without having to replace expensive cameras. IDx-DR exceeded the sensitivity measures of Medios AI at the cost of lower specificity. We have previously reported sensitivity and specificity of IDx-DR of 94% and 95% when compared to a single reader [3,12]. In this study IDx-DR retained excellent sensitivity at 99%, with a significantly lower specificity. This was more pronounced for any DR, with IDx-DR specificity of only 44%. This may be explained in part by the fact that IDx-DR has been specifically marketed for the detection of more than mild DR, and the specificity for detection of RDR is much higher at 68% with a 99% sensitivity. For comparison, in the pivotal trial that led to IDx-DR receiving FDA approval, where IDx-DR was compared against a diagnosis based on a 7-field ETDRS study with stereoscopic images and OCT, it achieved 87% sensitivity and 90% specificity for more than mild DR [13]. Both systems over-referred mild cases (False positives included 55% milds, 45% no DR by Medios AI; 38% mild, 62% no DR cases by IDx-DR). Another possibility of lower specificity could be referral of patients with similar lesions and concurrent pathologies that were not evaluated as part of this grading. Overall, both systems achieved satisfactory accuracy, particularly when patient safety is concerned with excellent negative predictive values, meaning very few patients receiving a false-negative result. Although false-positive results generate additional costs and workload due to the increase in referrals, it is the patients with missed disease that raise the most concern whenever automated screening is debated.

Many of the studies regarding automated analysis of DR from fundus images are sponsored or even performed directly by the respective software's owner company, which raises questions regarding bias. As previously mentioned Medios AI does not offer a dedicated desktop application at this point, therefore it was necessary to submit the images to Remidio, owner of the Medios AI algorithm, for a remote analysis on their system. As the authors of this study collaborated remotely, we could not directly oversee or verify the Medios AI output on site. Upon reviewing of the study methodology, we considered this to be a source of potential bias and asked Remidio for a way to independently verify some of the software's results. We submitted the subset of images assessed by Polish graders, through a dedicated API (application programming interface) provided by Remidio with live results. Images were anonymized without changing the image content. The results were in line with those previously submitted by Remidio for all but 3 patients.

Out of those three patients, for whom the initial MediosAI output differed from the verification, all three decisions changed from no RDR to RDR. For two of those patients the new MediosAI result now matched the grader decision of RDR, for the remaining patient the new MediosAI now disagreed with the grader consensus. All three of those patients had very subtle retinal signs. This is a study looking into MediosAI outside of the smartphone application which involves a custom-made deployment for the study. The discrepancies are likely due to challenges surrounding the implementation of the algorithm on different hardware or inconsistencies in image compression parameters between the version of the images submitted to Remidio for the first analysis and the version of the images sent for verification through the API.

Inter-grader and intra-grader variability is known amongst DR graders. Previous studies have shown that inter-grader kappa scores typically range from 0.40 to 0.65 in DR grading [13-19]. The kappa values for both Polish (0.68-0.90) and Indian graders (0.63-0.91) in the study showed similar variability but were well within the limits showing overall good agreement. This ensured reliability while having the data split and graded by both groups separately. The possible reasons for variability amongst graders could be identification and differentiation of subtle DR features (retinal hemorrhages, microaneurysms, hard exudates, new vessels, intraretinal microvascular abnormalities, neovascularization, and surrogate markers of macular edema), variation in image quality due to artifacts, brightness or contrast of images. This has been found in other studies and in other fields of medical imaging as well [13-15]. Gold standard grading in the current study was done based on a majority decision by the graders. As the individual grades were converted to a binary decision for each grader before computing consensus grading there was a majority decision for any DR and RDR for each of the patients. The reliability of the human grading could be improved with an adjudication process for patients without a full consensus [14].

This study included only patients with good quality, non-mydriatic images, which may not be representative of the whole screening cohort. Using a non-mydriatic protocol may underrepresent patients with smaller pupils or media opacities, particularly the elderly. In a previous study about DR grader reliability, based on the same screening programme in Poland from which images for this study were taken, out of 495 patients only 335 were deemed as sufficient quality by IDx-DR and all three human graders [15]. How many of those low-quality screening encounters

could we image and diagnose after mydriasis remains to be seen, as is the comparative performance of both systems in those patients.

In conclusion, our study compared the performance of two AI screening software, Medios AI and IDx DR, in detecting RDR. Both software systems demonstrated robust performance, with high accuracy and sensitivity, highlighting their potential as reliable tools for screening DR in real-world settings. Continued research and validation in larger and diverse patient populations will be essential to strengthen the evidence base and ensure the widespread adoption of these AI screening tools. Our study underscores the promise of these AI systems for DR screening, facilitating early detection and timely intervention for improved patient outcomes.

Statement of Ethics: The Institutional Review Board (Ophthalmology 21, Foundation for the Advancement of Ophthalmology, Application No. 2/2022) waived the need for IRB approval and written informed consent from participants for this retrospective study. The study was in adherence to the tenets of the Declaration of Helsinki.

Conflict of Interest: AG has Grants/Contracts from Alcon, Bausch & Lomb, Zeiss, Hoya, Thea, Viatri, Teleon, J&J, Cooper Vision, Essilor and Polpharma. AG has consulting fees / honoraria from Thea, Polpharma, Viatri and stock with GoCheckKids. DRP, FMS and KN are employees of Remidio Innovative Solutions. Remidio Innovative Solutions, Inc, USA and Medios Technologies are wholly owned subsidiary of Remidio Innovative Solutions Pvt Ltd, India. FMS has patents (mentioned in ICMJE) and stock (ESOP and stock, Remidio Innovative Solutions Pvt Ltd). Other authors declare no financial disclosures.

Funding/Support: There was no funding for this study

Author contribution

AG – Conceptualization; Data curation; Formal analysis; Investigation; Methodology; Project administration; Resources; Software; Supervision; Validation; Visualization; Roles/Writing - original draft; Writing - review & editing
DRP – Conceptualization; Formal analysis; Investigation; Methodology; Project administration; Resources;

Supervision; Visualization; Validation; Writing - review & editing

PB – Conceptualization; Data curation; Formal analysis; Investigation; Methodology; Roles/Writing - original draft; Writing - review & editing

KN – Formal analysis; Writing - original draft; Writing - review & editing

TK - Data curation; Formal analysis; Investigation; Methodology; Software; Validation; Visualization; Roles/Writing - original draft; Writing - review & editing

FMS - Conceptualization; Data curation; Formal analysis; Investigation; Methodology; Project administration; Resources; Software; Supervision; Validation; Visualization; Writing - review & editing

Data Availability Statement

The data that support the findings of this study is not publicly available due to ethical reasons and are available on request from the Dr Andrzej Grzybowski (email: ae.grzybowski@gmail.com).

References

1. Grzybowski A, Brona P, Lim G, Ruamviboonsuk P, Tan GSW, Abramoff M, Ting DSW. Artificial intelligence for diabetic retinopathy screening: a review. *Eye* 2020;34(3):451-460.
2. Ting DS, Cheung GC, Wong TY. Diabetic retinopathy: global prevalence, major risk factors, screening practices and public health challenges: a review. *Clin Exp Ophthalmol*. 2016;44(4):260-77.
3. Grzybowski A, Brona P. Analysis and Comparison of Two Artificial Intelligence Diabetic Retinopathy Screening Algorithms in a Pilot Study: IDx-DR and Retinalyze. *J Clin Med*. 2021;10(11):2352.
4. Lee AY, Yanagihara RT, Lee CS, Blazes M, Jung HC, Chee YE, et al. Multicenter, Head-to-Head, Real-World Validation Study of Seven Automated Artificial Intelligence Diabetic Retinopathy Screening Systems. *Diabetes Care*. 2021;44(5):1168-1175.
5. Sosale B, Sosale AR, Murthy H, Sengupta S, Naveenam M. Medios- An offline, smartphone-based artificial intelligence algorithm for the diagnosis of diabetic retinopathy. *Indian J Ophthalmol*. 2020;68(2):391-395.
6. Looker HC, Nyangoma SO, Cromie DT, Olson JA, Leese GP, Black MW, et al. Scottish Diabetes Research Network Epidemiology Group; Scottish Diabetic Retinopathy Collaborative. Rates of referable eye disease in the Scottish National Diabetic Retinopathy Screening Programme. *Br J Ophthalmol*. 2014;98(6):790-5.

7. Cuadros J, Bresnick G. EyePACS: an adaptable telemedicine system for diabetic retinopathy screening. *J Diabetes Sci Technol*. 2009;3(3):509-16.
8. Tilahun M, Gobena T, Dereje D, Welde M, Yideg G. Prevalence of Diabetic Retinopathy and Its Associated Factors among Diabetic Patients at Debre Markos Referral Hospital, Northwest Ethiopia, 2019: Hospital-Based Cross-Sectional Study. *Diabetes Metab Syndr Obes*. 2020;13:2179-2187.
9. Kozioł M, Nowak MS, Udziela M, Piątkiewicz P, Grabska-Liberek I, Szaflik JP. First nation-wide study of diabetic retinopathy in Poland in the years 2013-2017. *Acta Diabetol*. 2020;57(10):1255-1264.
10. Natarajan S, Jain A, Krishnan R, Rogye A, Sivaprasad S. Diagnostic Accuracy of Community-Based Diabetic Retinopathy Screening With an Offline Artificial Intelligence System on a Smartphone. *JAMA Ophthalmol*. 2019;137(10):1182-1188.
11. Sosale B, Aravind SR, Murthy H, Narayana S, Sharma U, Gowda SGV, Naveenam M. Simple, Mobile-based Artificial Intelligence Algorithm in the detection of Diabetic Retinopathy (SMART) study. *BMJ Open Diabetes Res Care*. 2020;8(1):e000892.
12. Grzybowski A, Brona P. A pilot study of autonomous artificial intelligence-based diabetic retinopathy screening in Poland. *Acta Ophthalmol*. 2019;97(8):e1149-e1150.
13. Abramoff MD, Lavin PT, Birch M, Shah N, Folk JC. Pivotal trial of an autonomous AI-based diagnostic system for detection of diabetic retinopathy in primary care offices. *NPJ Digit Med*. 2018;1:39.
14. Krause J, Gulshan V, Rahimy E, Karth P, Widner K, Corrado GS, et al. Grader Variability and the Importance of Reference Standards for Evaluating Machine Learning Models for Diabetic Retinopathy. *Ophthalmology*. 2018;125(8):1264-1272.
15. Grzybowski A, Brona P, Krzywicki T, Gaca-Wysocka M, Berlińska A, Świąch A. Variability of Grading DR Screening Images among Non-Trained Retina Specialists. *J Clin Med*. 2022;11(11):3125.
16. Guan MY, Gulshan V, Dai AM, Hinton GE. Who said what: Modeling individual labelers improves classification. In: *Thirty-Second AAAI Conference on Artificial Intelligence*. 2018. arXiv:1703.08774v2.
17. Sedova A, Hajdu D, Datlinger F, Steiner I, Neschi M, Aschauer J, Gerendas BS, Schmidt-Erfurth U, Pollreisz A. Comparison of early diabetic retinopathy staging in asymptomatic patients between autonomous AI-based screening and human-graded ultra-widefield colour fundus images. *Eye (Lond)*. 2022 Mar;36(3):510-516.
18. Gangaputra S, Lovato JF, Hubbard L, Davis MD, Esser BA, Ambrosius WT, Chew EY, Greven C, Perdue LH, Wong WT, Condren A, Wilkinson CP, Agrón E, Adler S, Danis RP; ACCORD Eye Research Group. Comparison of standardized clinical classification with fundus photograph grading for the assessment of diabetic retinopathy and diabetic macular edema severity. *Retina*. 2013 Jul-Aug;33(7):1393-9.
19. Raumviboonsuk P, Krause J, Chotcomwongse P, Sayres R, Raman R, Widner K, Campana BJ, Phene S, Hemarat K, Tadarati M, Silpa-Acha S. Deep Learning vs. Human Graders for Classifying Severity Levels of Diabetic Retinopathy in a Real-World Nationwide Screening Program. *arXiv preprint arXiv:1810.08290*. 2018 Oct 18.

Figure Legends:

Fig. 1. STARD flow diagram showing the patient breakdown for the Referable Diabetic Retinopathy (RDR) analysis

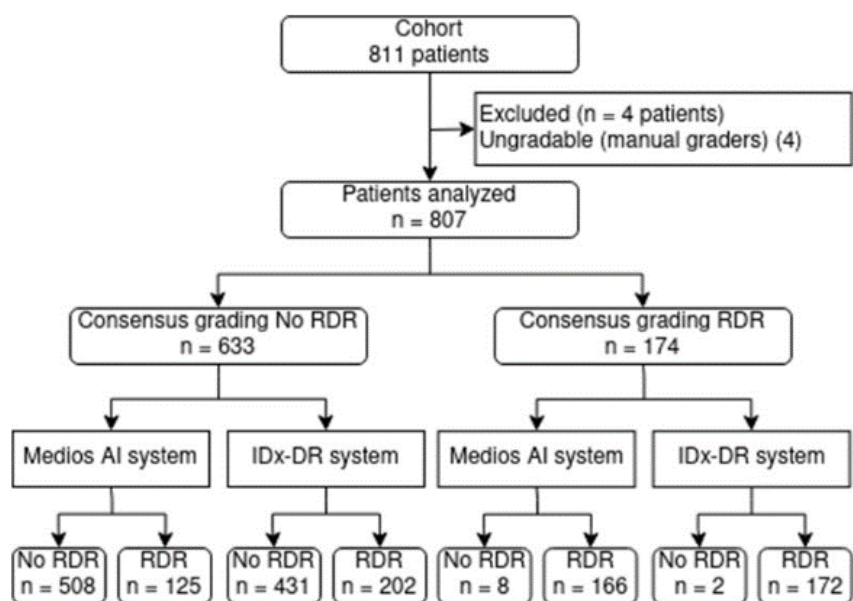


Table 1a) Performance of Medios and IDx DR with ground truth cutoff at any DR (n=264)

		Automated DR Image Assessment Softwares				
		Medios AI		IDx DR		
		Positive	Negative	Positive	Negative	TOTAL
Reference Standard (Consensus grading)	Positive (Mild NPDR and above and/or CSDME)	235 (29%)	29 (4%)	262 (32%)	2 (0.2%)	264
	Negative (No DR)	56 (7%)	487 (60%)	305 (38%)	238 (29%)	543
	TOTAL	291	516	567	240	807

Table 1b) Performance of Medios and IDx DR with ground truth cutoff at RDR (n=174)

		Automated DR Image Assessment Softwares				
		Medios AI		IDx DR		
		Positive	Negative	Positive	Negative	TOTAL
Reference Standard (Consensus grading)	Positive (Moderate NPDR and above and/or CSDME)	166 (21%)	8 (1%)	172 (21%)	2 (0.2%)	174
	Negative (No DR and mild NPDR)	125 (15%)	508 (63%)	202 (25%)	431 (53%)	633
	TOTAL	291	516	374	433	807

Table 1c) Performance of Medios and IDx DR with ground truth cutoff at STDR (n=41)

		Automated DR Image Assessment Softwares				
		Medios AI		IDx DR		
		Positive	Negative	Positive	Negative	TOTAL
Reference Standard (Consensus grading)	Positive (Severe, PDR or CSDME)	40 (5%)	1 (0.1%)	39 (5%)	2 (0.2%)	41

*

DR – Diabetic Retinopathy, NPDR- Non-Proliferative Diabetic Retinopathy, PDR – Proliferative Diabetic Retinopathy, CSDME – Clinically significant Diabetic Macular Edema, RDR - Referable Diabetic Retinopathy, STDR – Sight-threatening Diabetic Retinopathy

Table 2: Performance Analysis of Artificial Intelligence System Compared to Reference Standard

Values in % (95% CI)	For any DR		For referable DR		For sight-threatening DR	
	Medios AI	IDx-DR	Medios AI	IDx-DR AI	Medios AI	IDx-DR AI
Sensitivity	89 (85, 93)	99 (97, 100)	95 (91, 98)	99 (96, 100)	98 (87, 100)	95 (83, 99)
Specificity	90 (87,92)	44 (40, 48)	80 (77, 83)	68 (64, 72)	NA	80 (77, 83)
PPV	81 (76, 85)	46 (42, 50)	57 (51, 63)	46 (41, 51)	NA	21 (15, 27)
NPV	94 (92, 96)	99 (97, 100)	98 (97, 99)	100 (98, 100)	NA	100 (99, 100)