
How to use SensaCath® Flex Compact

1



Wash your hands with soap and water then dry them.

2

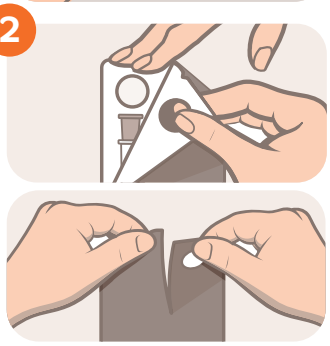

Using the pull tabs, or tearing down the middle, open the pouch until the catheter funnel is visible.

3




You can hang the pouch using the adhesive spot or rest it on a clean surface.

4




If your healthcare professional has advised you, try to urinate before catheterising. If you have a foreskin, gently pull it back and clean the tip of the penis.

5

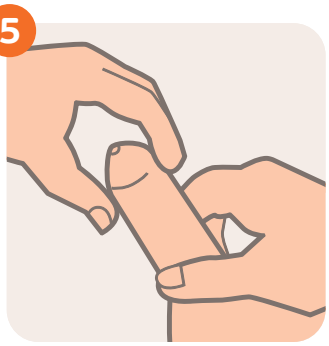

To ensure the urethra is straight, hold your penis pointing upwards close to your stomach. Be careful not to squeeze it.

6

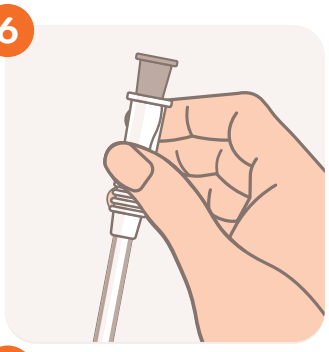


Remove the catheter packaging. There is a white non-touch insertion slider included so it is simple for you to hold and insert the catheter.

7




Position the non-touch insertion slider approximately 10 to 15 centimetres from the catheter tip so that you feel you can control the tip when inserting the catheter.

8

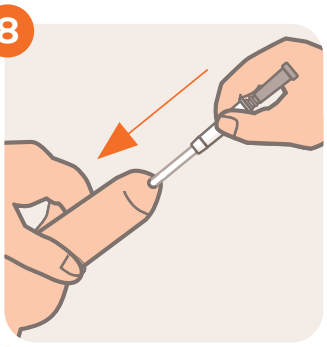

Gently insert the catheter into your urethra until it reaches your bladder (you may feel a little resistance but that's normal). Continue to insert until the urine starts to flow.

9

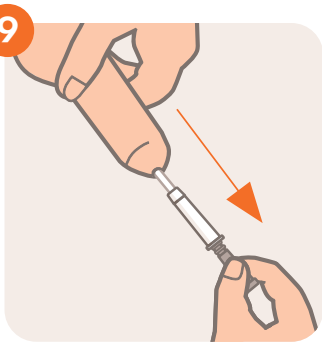

When the flow stops, slowly remove the catheter. If urine flows again, pause until the urine stops flowing.

10



Dispose of **SensaCath® Flex Compact** by placing it back into the pouch before disposing of it in a waste bin. Do NOT flush it down the toilet. Wash your hands.

Clinisupplies

Clinisupplies Ltd, 1 Blackmoor Lane, Croxley Park, Watford, Hertfordshire, WD18 8GA
Telephone: +44 (0) 20 8863 4168 Email: info@clinisupplies.co.uk
www.clinisupplies.co.uk

CSF01469_02/26_v2

A quick-look guide

SensaCath® Flex Compact Intermittent Catheter

Clinisupplies

Nightingale
by Clinisupplies

Ordering is simple and convenient.

Nightingale Home Delivery Service.
For all your bladder and bowel products.

Download the app today!

Download on the App Store

GET IT ON Google Play

Try it now!

Clinisupplies | Nightingale

Nightingale Home Delivery Service is part of the **Clinisupplies Group** of companies.
To find out more, visit: **www.clinisupplies.com**

You can also order online, or by phone:
www.nightingaledelivery.co.uk
0800 304 7434