Clinisupplies

All-Silicone Foley Catheter

User guide



Designed for you.

The **Prosys**® range of continence products have been designed to be comfortable to use and to provide you with confidence in managing your bladder health.

Our comprehensive range of **Prosys**® **All-Silicone Foley Catheters** are available in multiple sizes and two different tips and will sit comfortably in situ for up to 90 days.

Prosys® All Silicone Foley Catheter


COMFORT

- Made with 100% silicone
- Flexible tubing
- Compatible with all urine drainage bags

CONFIDENCE

- Large rounded eyelets for optimal drainage
- Open ended tip for reduced risk of blockages
- Nelaton tip to facilitate catheter insertion

Any questions? Call 0800 304 7434 and we will be happy to help.

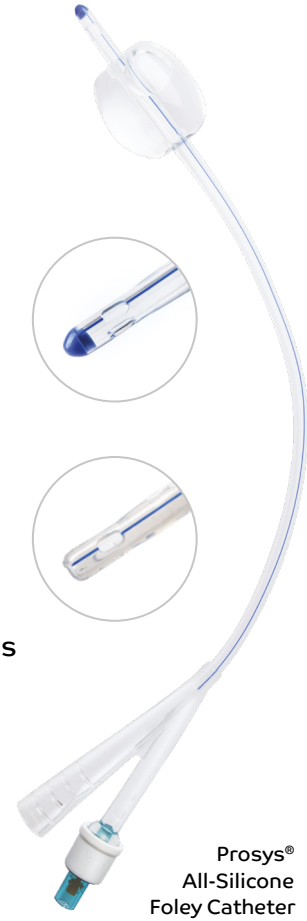

Prosys®
All-Silicone
Foley Catheter

LATEX

Connecting to **Prosys®** Leg Bag

1. To reduce the risk of infection, it is important you wash your hands prior to connecting the **Prosys®** Leg Bag to the **Prosys®** All-Silicone Foley Catheter.



2. Ensure the lever on the outlet tap is in the closed position.
3. Remove the grey protective cap from the universal connector.
4. Insert the connector into the end of the indwelling catheter using a twist motion to ensure a secure connection.

Make sure you do not touch the catheter or leg bag connectors as that will increase chances of cross-contamination and infection.

How to secure

1. The **Prosys®** Leg Bag should be positioned on the thigh or calf, depending on the tube length, and secured appropriately with the leg bag straps provided.
2. Feed the **Prosys®** Leg Bag Straps through the eyelets on the bag and ensure the silicone beaded lines are positioned against the leg for optimal comfort and security.



Caring for an indwelling catheter

1. Ensure that you wash your hands before and after handling your catheter, leg/night bag or catheter valve.
2. Wash the area where your catheter enters the body with mild unperfumed soap and warm water at least once a day.
3. Do not use oil-based creams around the catheter.
4. Avoid unnecessary disconnection of the leg bag and do not touch the connector when replacing. This will help to prevent infection.
5. Ensure the leg/night bag is kept below the level of your bladder to promote urine drainage.
6. To maintain urine flow, ensure your catheter tubing is not kinked.

The National Institute for Health and Care Excellence recommend you secure your catheter with a catheter strap to prevent the risk of urethral trauma.



Ordering information

Prosys® All-Silicone Foley Catheter

Long-term, 2-way with 10ml balloon and pre-filled syringe (sterile water) Box of 1

Reference	Description	Size	Length	PIP code
PCF12M10	Standard length	12Ch	40cm	386-0616
PCF14M10	Standard length	14Ch	40cm	386-0582
PCF16M10	Standard length	16Ch	40cm	386-0574
PCF18M10	Standard length	18Ch	40cm	386-0541
PCF12F10	Female length	12Ch	26cm	386-0624
PCF14F10	Female length	14Ch	26cm	386-0558
PCF16F10	Female length	16Ch	26cm	386-0608
PCF18F10	Female length	18Ch	26cm	386-0590

Prosys® Open Ended Foley Catheter

Long-term, 2-way with 10ml balloon and pre-filled syringe (sterile water) Box of 1

Reference	Description	Size	Length	PIP code
PCF12M10E	Standard length	12Ch	40cm	419-0674
PCF14M10E	Standard length	14Ch	40cm	419-0682
PCF16M10E	Standard length	16Ch	40cm	419-0690
PCF18M10E	Standard length	18Ch	40cm	419-1086

Liberty® Hydro-Caine® Lubricating Gel – 2.00% Lidocaine Hydrochloride, 0.05% Chlorhexidine Gluconate Solution

Reference	Description	Size	PIP code
HYD006	Small	6ml	366-5528
HYD011	Large	11ml	366-5510

Prosys® Catheter Strap

Reference	Description	Length	Quantity	PIP code
PCS40	Short	40cm	5	408-5536
PCS50	Adult	50cm	5	408-5544
PCS85	Abdominal	85cm	5	408-5551

Nightingale by Clinisupplies

NHS
Providing NHS services

Welcome to Nightingale Home Delivery Service.

For all your bladder and bowel products.

Fast, accurate ordering with **discreet** home delivery.



People love us

Trustpilot

Excellent*

*July 2025

Clinisupplies | Nightingale

Nightingale Home Delivery Service is part of **Clinisupplies**, helping you and your healthcare professional to find and order the bladder and bowel products you may need.

To order, register or login online or via the app:
www.nightingaledelivery.co.uk

Download on the App Store

GET IT ON Google Play

or phone: **0800 304 7434**

CS23360_v7