Clinisupplies

Catheter Strap

User guide



Designed for you.

The **Prosys®** range of continence products have been designed to be comfortable to use and to provide you with confidence in managing your bladder health.

Our simple **Prosys®** Catheter Strap provides all the support you need to maintain an active lifestyle and the freedom to move around safely.

Prosys® Catheter Strap

COMFORT

- Elasticated design allows for free movement
- Soft, breathable fabric maximises comfort
- Cut-to-size to remove any excess material

CONFIDENCE

- Clasp ensures the catheter stays fixed in position
- Silicone beading gently secures the strap without irritation
- Compatible with all indwelling catheters and leg bags



Any questions? Call 0800 304 7434 and we will be happy to help.


How to use

The **Prosys® Catheter Strap** is most commonly worn on the leg in the thigh position.

For those with a suprapubic catheter, the strap can be positioned on the lower abdomen.

Securing the strap to your body

Place the strap around either the leg or lower abdomen. Fasten the strap in the desired position on your leg or abdomen. Ensure the strap is snug and the clasp is easily accessible.

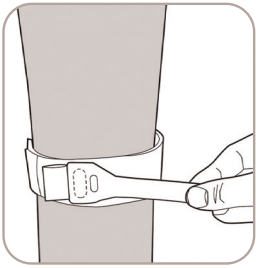


Securing the catheter tube

Position the clasp over the tubing. Wrap the clasp around and under the tubing. Push the tapered end of the clasp through the eyelet provided.




Securing the clasp

Gently pull the clasp to tighten around the tube. Pull the clasp up and over, pressing the tapered end into place and securing it on the main strap.

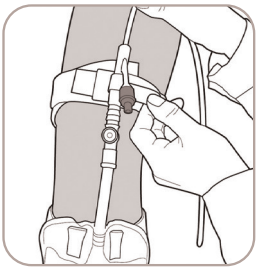



Use the catheter strap to secure either the catheter or the leg bag

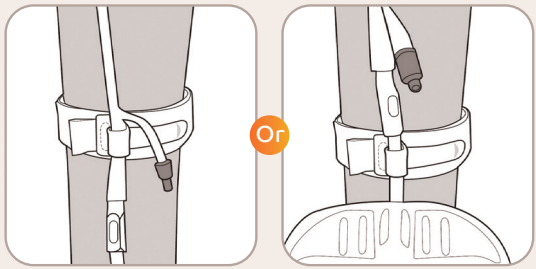


Selecting the appropriate strap length

Use a measuring tape to measure the circumference of the patient's thigh at the widest point, or around the abdomen at the area the strap is to be worn. Choose a **Prosys® Catheter Strap** that is longer than the measured circumference to ensure proper fit and secure placement.

Caring for your strap

To support performance and good hygiene, we recommend you wash the strap when needed. To do so, secure the clasp in place onto the strap and wash according to the instructions on the label.

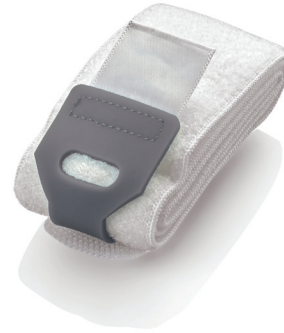

Ordering information

Prosys® Catheter Strap

Reference	Size	Length	Quantity	PIP code
PCS40	Short	40cm	5	4085536
PCS50	Adult	50cm	5	4085544
PCS85	Abdominal	85cm	5	4085551



Nightingale
by Clinisupplies

NHS
Providing NHS services

Welcome to Nightingale Home Delivery Service.

For all your bladder and bowel products.

Fast, accurate ordering with **discreet** home delivery.



People love us
Trustpilot
Excellent*

*July 2025

Clinisupplies | Nightingale

Nightingale Home Delivery Service is part of **Clinisupplies**, helping you and your healthcare professional to find and order the bladder and bowel products you may need.

To order, register or login online or via the app:
www.nightingaledelivery.co.uk

Download on the App Store
GET IT ON Google Play

or phone: **0800 304 7434**

CS2118_0818_v7