Clinisupplies

Flofit® Sheath

User guide



Designed for you.

The **Prosys**® range of continence products have been designed to be comfortable to use and to provide confidence in managing bladder health.

Our comprehensive range of **Prosys**® **Flofit Sheaths** are available in three different styles and five sizes, to provide choice.

Prosys® Flofit Sheath

COMFORT

- Made with 100% silicone
- Includes hair protectors
- Compatible with our drainage bags

CONFIDENCE

- Skin friendly adhesives providing security
- Easy to apply and remove
- Kink resistant spout

Any questions?
Call 0800 304 7434



LATEX

Suitable for 24hrs

Selecting the correct sheath style.

STANDARD



A standard length sheath with centrally positioned adhesive is suitable for men who need a sheath to cover the entire length of the penis.

POP-ON



A pop on sheath with maximum adhesive area is suitable for men who do not need the sheath to cover the entire penis, and is useful for men who experience retraction.

WIDE BAND



A wide band sheath with greater adhesive area is suitable for men who are active and need to cover the entire length of the penis.

Selecting **Prosys**® Flofit Sheath size

- **Prosys**® **Flofit** is available in 5 sizes; 24, 28, 31, 35, 40mm
- For fitting, measure the penis at the widest point by carefully placing the cut-out of the measuring guide around the mid-shaft of the penis
- If you find the penis is in between two sizes, select the larger one



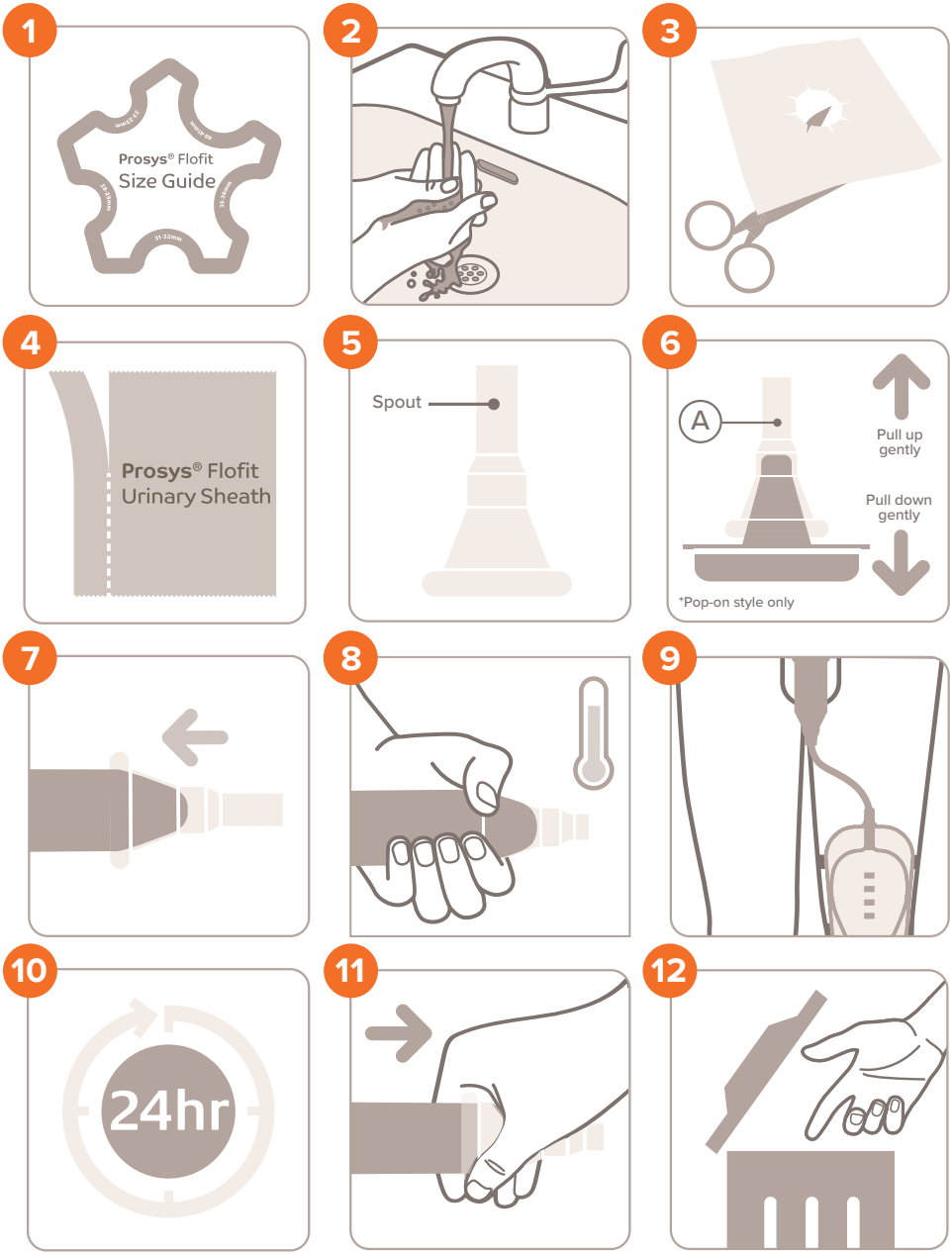
Prosys® **Flofit** Sheath

How to fit

1. Please ensure you are wearing the correct sheath size by utilising the measuring guide.
2. Please ensure you have washed your hands before using the **Prosys® Flofit** Sheath. Gloves should be worn if the sheath is being fitted by a healthcare professional.
3. If necessary, trim any pubic hair or use the hair protectors supplied. Do not shave as this can cause skin irritation. Wash the penis and groin area using a mild unperfumed soap and completely dry before proceeding.
4. Open the sheath pouch – tear along the perforation to access the product.
5. Remove the sheath from the pack by the spout.
6. If using the Pop-On style sheath, take it off the plastic cone by holding the sheath just below the spout (point A on Diagram 6) and pull down gently on the plastic cone.
7. Position the sheath over the penis. Unroll the sheath along the shaft of the penis, keeping it even and straight until fully unrolled.
8. Place your hand around the shaft of the penis and the sheath. Hold the adhesive area for 30 seconds – this will allow the adhesive to warm up and will ensure a secure fit.
9. Fit the leg bag to the leg before connecting it to the sheath to prevent any twisting.

How to remove

1. The sheath can be worn for up to a 24hr period.
2. Removal of the sheath may be easier in the bath or shower. Use warm, non-perfumed, soapy water and a wash-cloth before gently rolling off the sheath.
3. Dispose of your sheath in the household bin and wash your hands.



PLEASE NOTE: if your skin becomes irritated whilst using this product, do not use until you have spoken to your healthcare professional. Always check your skin between applications to ensure there is no irritation. Ensure you always dry your skin properly before applying a new sheath.

Connect to a leg bag.

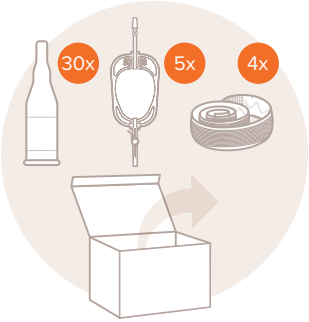
1. It is important to wash your hands prior to connecting the leg bag.
2. Open the leg bag pouch. Inside you will find a pair of gloves. Please put these on before removing the leg bag from its packaging.
3. Ensure the lever on the outlet tap is in closed position.
4. Remove the blue protective cap from the inlet connector.
5. Insert the inlet connector into the end of the sheath.

Reconnect to the bag when the sheath is changed every day, unless using a single use night bag. The reusable leg bag will need to be changed every 5-7 days.

How to secure

1. The leg bag should be positioned should be positioned on the thigh or calf, depending on the tube length, and secured appropriately with the leg bag straps provided.
2. Feed the leg bag straps through the eyelets on the bag and ensure the silicone beaded lines are positioned against the leg for optimal comfort and security.
3. Position the leg bag below the bladder for continuous urine flow and storage.

Prosys® FloPack



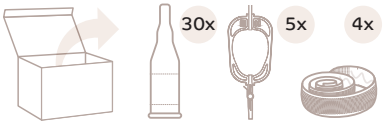
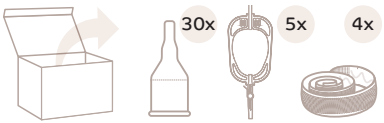
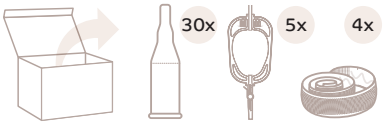
One box,
One month,
One prescription.

Prosys® FloPack is an all-in-one box containing **Prosys® Flofit** Sheaths, **Prosys®** Leg Bags and **Prosys®** Leg Bag Straps. **Prosys® FloPack** contains one months' worth of supply for patients with urinary incontinence.

- Includes sheaths and leg bags
- One regular monthly order
- Less packaging waste
- 100% latex free
- Includes hair protectors

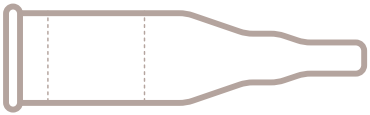
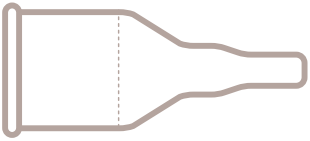
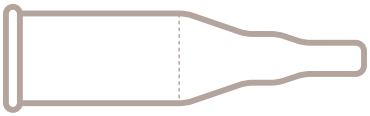
Prosys® FloPack is available in three different diameters 28mm, 31mm and 35mm and three different styles, **Standard**, **Pop-on** and **Wide band**.

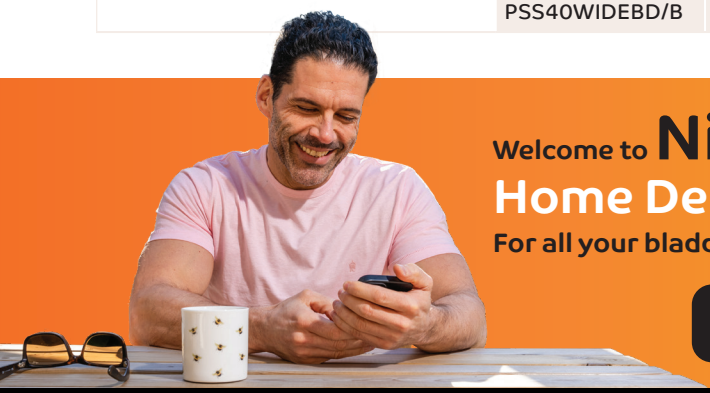
Prosys® FloPack 30 Day Box

Standard style	Reference	Size	PIP code
	PFPS28	28mm	422-3772
	PFPS31	31mm	422-3780
	PFPS35	35mm	422-3798
Pop-On style	**Reference**	**Size**	**PIP code**
	PFPP28	28mm	422-3806
	PFPP31	31mm	422-3814
	PFPP35	35mm	422-3822
Wide Band style	**Reference**	**Size**	**PIP code**
	PFPW28	28mm	422-3749
	PFPW31	31mm	422-3756
	PFPW35	35mm	422-3764

Ordering information

Prosys® Flofit Urinary Sheath — Box of 30

Standard style	Reference	Size	PIP code
	PSS24STD/B	24mm	406-0059
	PSS28STD/B	28mm	405-9762
	PSS31STD/B	31mm	405-9770
	PSS35STD/B	35mm	405-9788
	PSS40STD/B	40mm	405-9796
Pop-On style	**Reference**	**Size**	**PIP code**
	PSS24POPON/B	24mm	405-9804
	PSS28POPON/B	28mm	405-9812
	PSS31POPON/B	31mm	405-9754
	PSS35POPON/B	35mm	405-9820
	PSS40POPON/B	40mm	405-9846
Wide Band style	**Reference**	**Size**	**PIP code**
	PSS24WIDEBD/B	24mm	405-9911
	PSS28WIDEBD/B	28mm	405-9903
	PSS31WIDEBD/B	31mm	405-9887
	PSS35WIDEBD/B	35mm	405-9879
	PSS40WIDEBD/B	40mm	405-9853



Welcome to **Nightingale Home Delivery Service**

For all your bladder and bowel products.

Fast, accurate ordering with discreet home delivery.

Clinisupplies | Nightingale

Nightingale Home Delivery Service is part of **Clinisupplies**, helping you and your healthcare professional to find and order the bladder and bowel products you may need.

To order, register or login online or via the app:
www.nightingaledelivery.co.uk

Download on the App Store | GET IT ON Google Play

or phone: 0800 304 7434

CS23362_1020 v6