Clinisupplies

Catheter Valve

User guide



Designed for you.

The **Clinisupplies** range of continence products have been designed to be both comfortable to use and to provide you with confidence in managing your bladder health.

Our **Prosys®** and **Libra® Catheter Valves** are manufactured with a simple open and closed lever tap function, allowing you to take care of yourself with ease and independence.

Catheter Valve

COMFORT

- Compact size
- Soft outlet tube
- Low profile lever

CONFIDENCE

- Coloured tap easy to identify
- Compatible with all drainage bags
- Ridged connector with open and closed indication on tap

Any questions? Call 0800 304 7434 and we will be happy to help.



How to use

With the catheter valve in place, you can maintain bladder function and capacity by allowing the filling and emptying of the bladder, mimicking normal function.


It is important you remember how to use the valve and must empty your bladder at regular intervals. Failure to do so can result in urinary infection.

Frequency of emptying bladder is recommended by your healthcare professional.


How does the Catheter Valve work?

1. Wash your hands thoroughly before opening the peel pouch, and then remove the valve from its packaging.


2. Attach the ridged connector of **Catheter Valve** to the indwelling catheter.

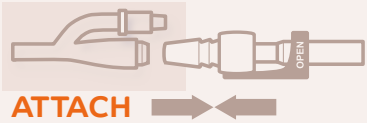

3. To open the **Catheter Valve** and allow urine to drain, push the lever fully down until the OPEN symbol is facing upwards.

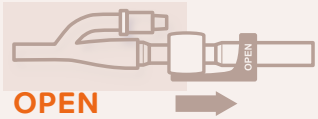

4. To close the **Catheter Valve**, pull the lever up so it sits with the catheter and the CLOSED symbol is facing upwards.

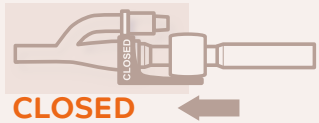

Caring for my **Prosys®** and **Libra®** Catheter Valve


1. Do not touch the open end of the catheter or valve when connecting the two devices, or emptying the bladder, as this can lead to infection.

2. Make sure you have washed your hands thoroughly with soap and warm water. If your carer is changing the catheter valve, gloves must be worn prior to replacing the valve.

3. When wearing a catheter, it is important you have secured your catheter with a securement device to prevent risk of urethral trauma when mobilising. See **Clinispplies range of Catheter Straps** for more information.

4. You can have a bath or shower with the catheter and valve in situ.

5. Cleanse your catheter and valve daily with warm soapy water.

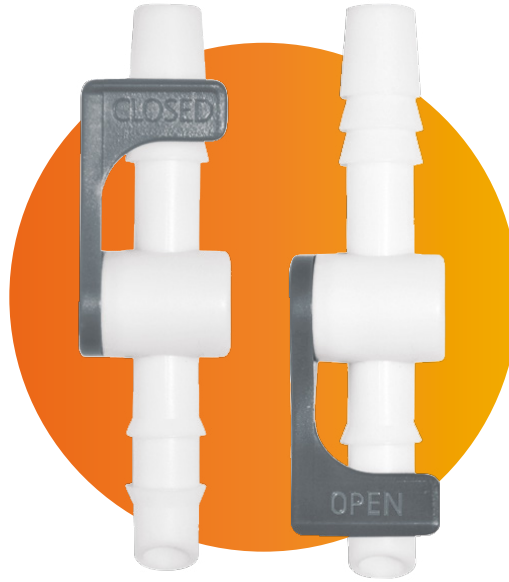

Ordering information

Clinisupplies Catheter Valves are available on prescription

Reference	Description	Size	PIP code
PCV3942	**Prosys**® Catheter Valve	5 in a pack	386-3800
10540A	**Libra**® Catheter Valve	5 in a pack	342-9487



Nightingale
by Clinisupplies

NHS
Providing NHS services

Welcome to Nightingale Home Delivery Service.

For all your bladder and bowel products.



People love us

Trustpilot

Excellent*

Fast, accurate ordering with discreet home delivery.

*Feb 2026

Clinisupplies | Nightingale

Nightingale Home Delivery Service is part of **Clinisupplies**, helping you and your healthcare professional to find and order the bladder and bowel products you may need.

To order, register or login online or via the app:
www.nightingaledelivery.co.uk

Download on the App Store

GET IT ON Google Play

or phone: **0800 304 7434**

CS23361_1020_v7