Clinisupplies

Aquaflush®

Compact and Compact+ Mini and Midi Systems

EN: **Transanal Irrigation Systems:** Instructions for Use
NO: **Transanal Irrigasjonssystem:** Brukerveiledning
DA: **Transanalt irrigationssystem:** Brugsanvisning
SV: **Transanal Irrigationssystem:** Bruksanvisning
FI: **Suolihuuhtelujärjestelmä:** Käyttöohje
NL: **Transanaal irrigatiesysteem:** Gebruiksaanwijzing
ES: **Sistema de irrigación transanal:** Instrucciones de uso
DE: **Transanale Irrigationssysteme:** Gebrauchsanweisung
FR: **Systèmes d'irrigation transanale:** Mode d'emploi
PT: **Sistemas de Irrigação Transanal:** Instruções de utilização
IT: **Sistemi per irrigazione transanale:** Istruzioni per l'uso

EN: Thank you for choosing Aquaflush. Aquaflush Transanal Irrigation (TAI) systems are designed to be a SIMPLE and EFFECTIVE way for you to control your bowel management routine.

Please read this document carefully, including the warnings and precautions. If you have any questions, please contact us at Clinisupplies or discuss any concerns with your Healthcare Professional who should also provide you with some training before commencing use.

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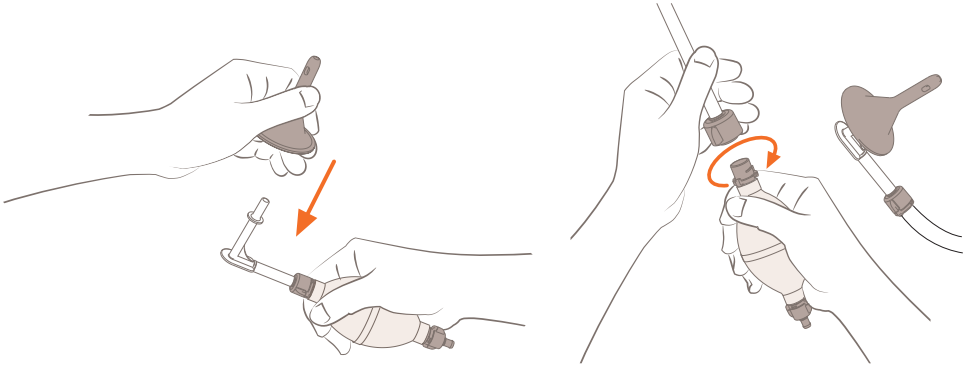
with optional extension tube

2

EN: Preparing the system for use: Fill water between 35°C to 37°C (95°F to 98°F), which is close to body temperature in a clean bowl or wash basin (Do not use cold or hot water). Squeeze the pump to let the air out as it is important you do not pump air into the bowel. Squeeze the pump and then submerge the bottom of the pump into the water and relax your grip to draw the water into the pump and tubing. Repeat this process until you get a steady stream of water from the cone tip. Your Aquaflush system is now ready for use.

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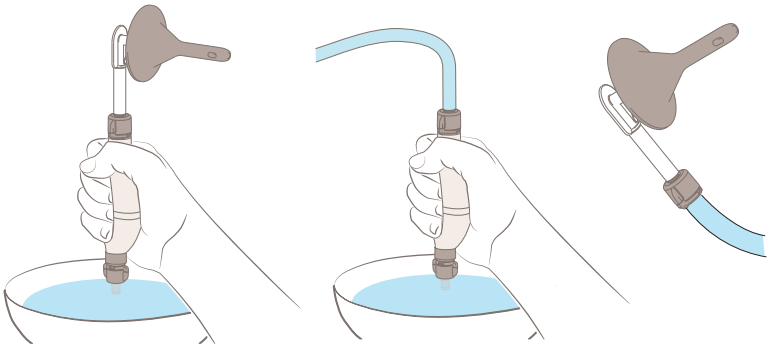
3

EN: Apply the lubricant provided to the tip and shaft of the cone and to the anus if required. While hovering over the toilet or seated, gently push the cone into the rectum at the angle taught by your Healthcare Professional. Tilt the pump upward slightly to ensure air is not pumped into the bowel. Gently squeeze the pump to allow the water to flow into your rectum as directed by your Healthcare Professional. (If any discomfort is felt stop pumping immediately and seek advice from your Healthcare Professional). Once the process is complete, remove the cone and remain seated on the toilet. Your bowel may take up to 30 minutes to empty, and it is recommended to wait an additional 5 minutes after everything appears to have been expelled.

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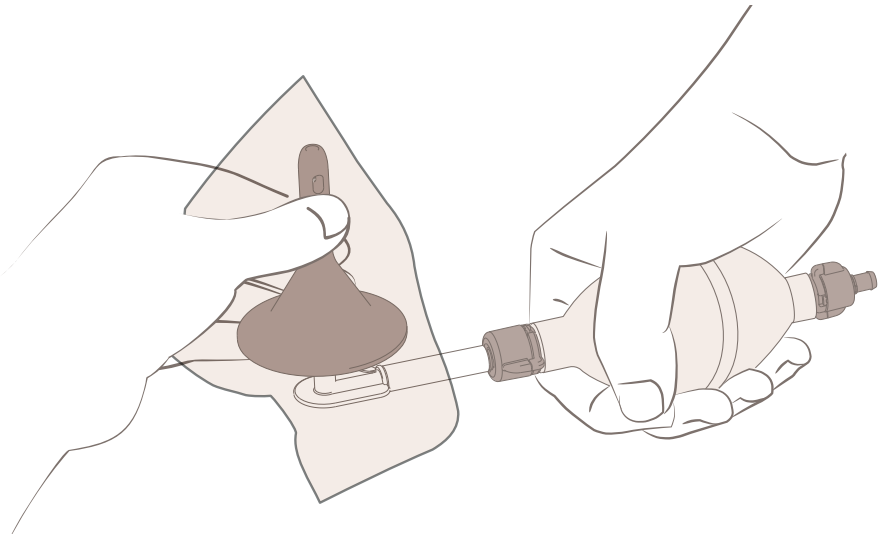
EN: Aquaflush® Transanal Irrigation systems must only be used by the person that it is prescribed for following a thorough assessment by a Healthcare Professional.

Caution: The following list does not rule out transanal irrigation as a bowel management programme, but you should have consulted a Healthcare Professional to ensure that the system is right for you. Children under 5 must be reviewed by a Paediatric Specialist before starting transanal irrigation.

- Inflammatory bowel disease (e.g. Crohn's disease or Ulcerative Colitis and Ischaemic Colitis)
- Rectal pain or bleeding e.g. anal fissure, severe haemorrhoids
- Diverticular Disease
- Autonomic Dysreflexia
- Recent spinal cord shock following spinal trauma
- Radiation therapy in the abdominal or pelvic region
- Abdominal or anal surgery
- Recent colonic biopsy or endoscopic polypectomy
- Cancer in the abdominal or pelvic region
- Long term steroid therapy, Anticoagulant Therapy
- Coronary Disease
- Pregnancy
- Nursing mothers

Intended use: The intended use of the product is for transanal irrigation of the bowel in adults and children 3 years and above with a bowel dysfunction, who have been assessed and trained by an appropriate healthcare professional.

WARNINGS: You must consult a Healthcare Professional if any of the following occur:

- Blood in the faeces, weight loss or abdominal pain
- Significant changes in the frequency, colour and consistency of the stool

Rare occurrence: Bowel perforation is extremely rare but is serious and life threatening if it occurs and will require urgent attention. If you experience any of the following symptoms during or after transanal irrigation, contact your Healthcare Professional immediately:

- Severe and sustained abdominal pain or back pain, especially if combined with fever
- Sustained bleeding

SAFETY: Ensure you store your Aquaflush® systems away from the reach of small children to avoid risk of strangulation from tubing, suffocation from plastic bags or choking in the event of small components coming loose from the product.

Do not use after the date of expiry, which can be found on the carton label.

Any serious incident that has occurred with this device should be reported to the manufacturer or competent authority.

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Clinisupplies Ltd, 1 Blackmoor Lane, Croxley Park, Watford, Hertfordshire WD18 8GA, UK
Clinisupplies (Nanjing) Co., Ltd.,
EC REP **Advena Ltd.** Tower Business Centre, 2nd Flr., Tower Street, Swatar, BKR 4013 Malta
VitalCare Limited, 6-9 Trinity Street, Dublin 2, Ireland