
Clinisupplies


Aquaflush®

Aquaflush® Maxi Systems


EN: **Transanal Irrigation Systems:** Instructions for Use
NO: **Transanal Irrigasjonssystem:** Brukerveiledning
DA: **Transanalt irrigationssystem:** Brugsanvisning
SV: **Transanal Irrigationssystem:** Bruksanvisning
FI: **Suolihuuhtelujärjestelmä:** Käyttöohje
NL: **Transanaal irrigatiesysteem:** Gebruiksaanwijzing
ES: **Sistema de irrigación transanal:** Instrucciones de uso
DE: **Transanale Irrigationssysteme:** Gebrauchsanweisung
FR: **Systèmes d'irrigation transanale:** Mode d'emploi
PT: **Sistemas de Irrigação Transanal:** Instruções de utilização
IT: **Sistemi per irrigazione transanale:** Istruzioni per l'uso

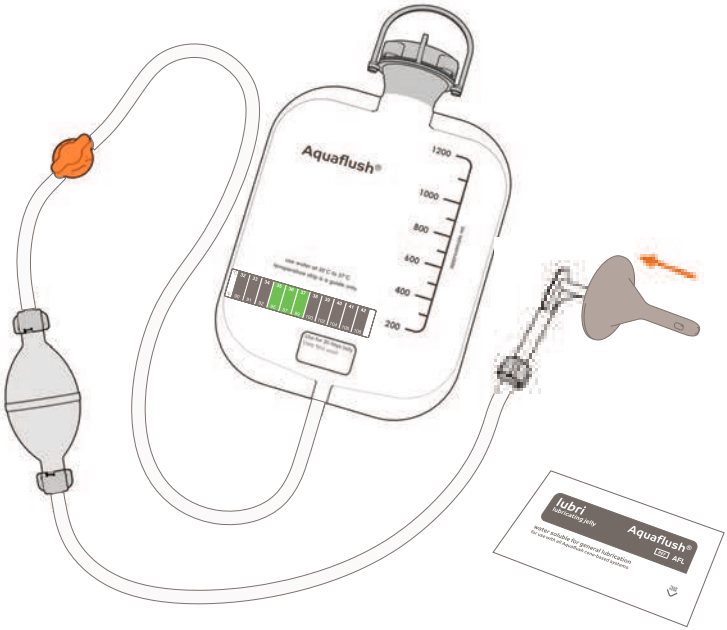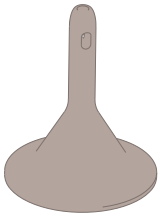
Lite


REF AFLS, AFLM

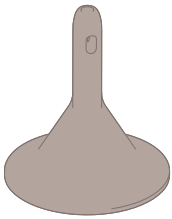Actif

REF AFAS, AFAM

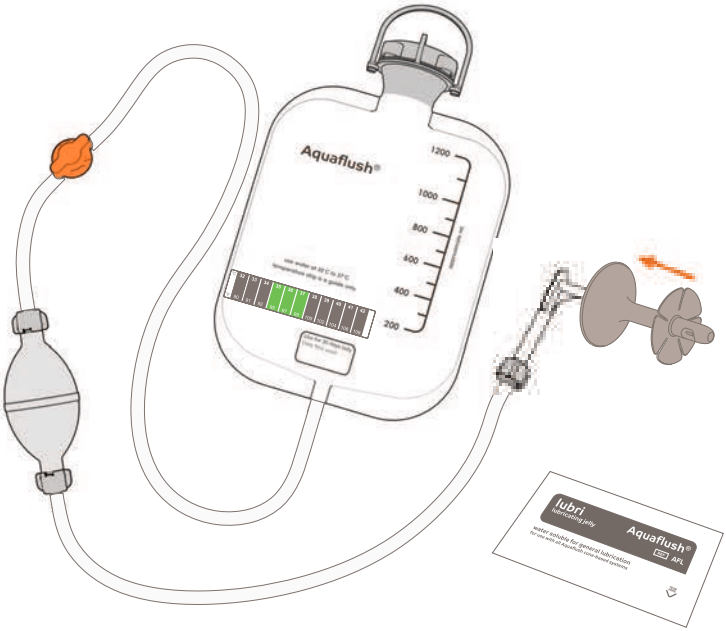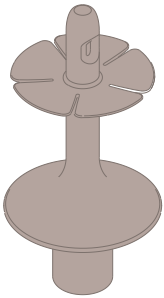Self-Retaining Cone (SRC)

REF AFSRCS, AFSRCM

EN: Thank you for choosing Aquaflush. Aquaflush® Transanal Irrigation (TAI) systems are designed to be a SIMPLE and EFFECTIVE way for you to control your bowel management routine.

Please read this document carefully, including the warnings and precautions. If you have any questions, please contact us at Clinisupplies or discuss any concerns with your Healthcare Professional who should also provide you with some training before commencing use.

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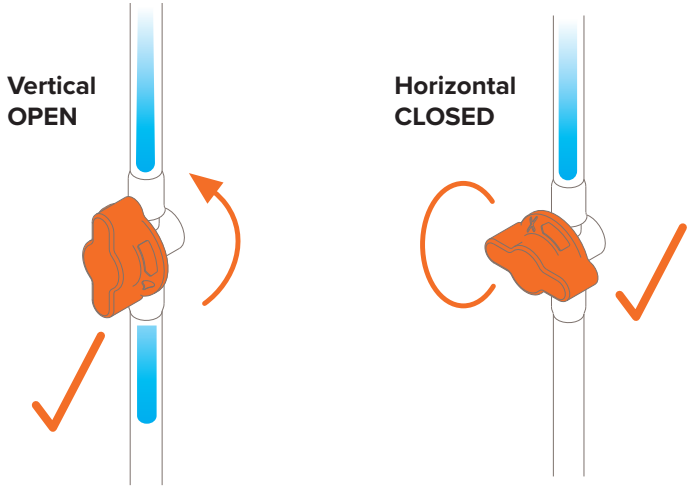
2

EN: Preparing your system for use: Unscrew the waterbag cap and fill the bag with water at a temperature between 35°C to 37°C (95°F to 98°F), which is close to body temperature. (Do not use cold or hot water). Fold the bag at the water line to remove any trapped air before securely closing the cap. This process ensures that no air is trapped in the bag, preventing the risk of pumping air into the system. It is important you do not pump air into the bowel.

Helpful Hint: The temperature gauge on the bag will indicate when the water in the correct range.

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3

EN: Place the waterbag on a flat surface or on the floor using the Aquaflush® Bag Stand provided or hang it above shoulder height using the built-in hanger with an 'S' hook*. The hanging option uses gravity, allowing the system to free flow without the need for a pump and should be hung above the height of the toilet to avoid getting air into the system. You can only free flow if the water bag is above the height of the person doing the irrigation.

*If you need a bag stand or an 'S' hook, contact your delivery company or supplier.

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

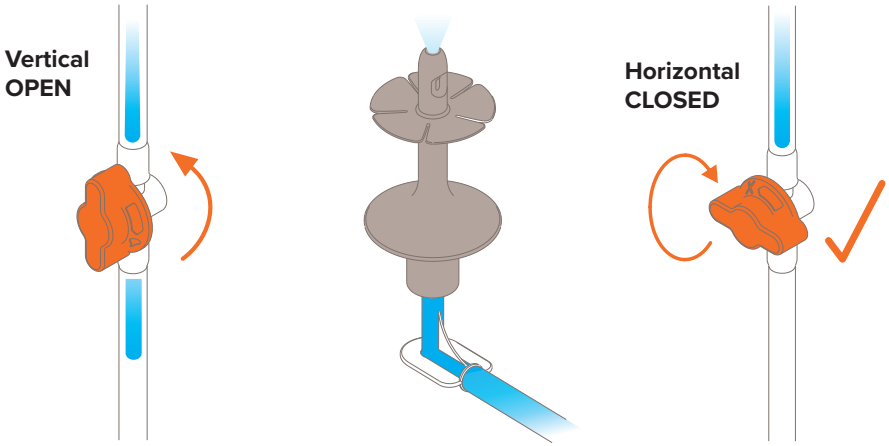
5

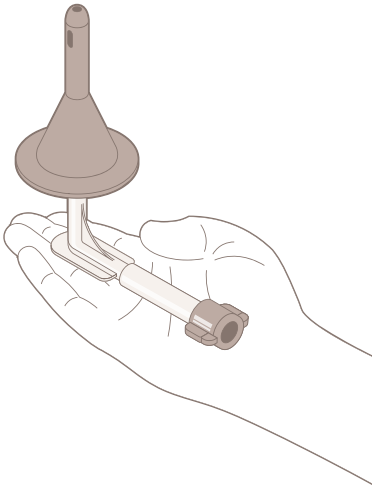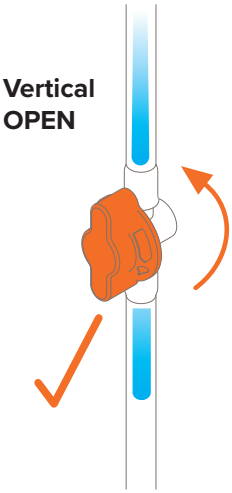EN: Apply the lubricant provided to the tip and shaft of the cone and the anus if required. Sit on the toilet and gently push the cone into the rectum at the correct angle of insertion advised by your Healthcare Professional. [If you are using a Self Retaining Cone (SRC) then pull back the cone slightly so that the petals sit firmly in the rectum]. Turn the orange tap to the vertical open position fully.

Helpful Hint: If you have difficultly correctly inserting the cone, consult your Healthcare Professional.

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
6

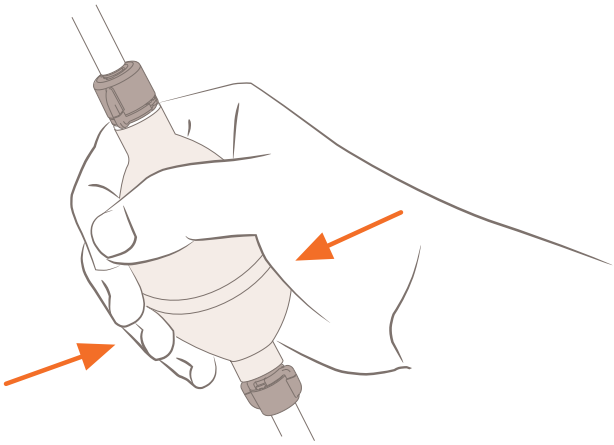EN: Gently pump the volume of water advised by your Healthcare Professional into the rectum, monitoring the volume scale on the water bag (usually between 300ml to 500ml). (If any discomfort is felt stop pump immediately, close the orange tap to the horizontal position and seek advice from your Healthcare Professional). When you have finished pumping, turn the orange tap to the horizontal closed position and remove the cone from the rectum. Remain seated on the toilet. Your bowel may take up to 30 minutes to empty, and it is recommended to wait an additional 5 minutes after everything appears to have been expelled.

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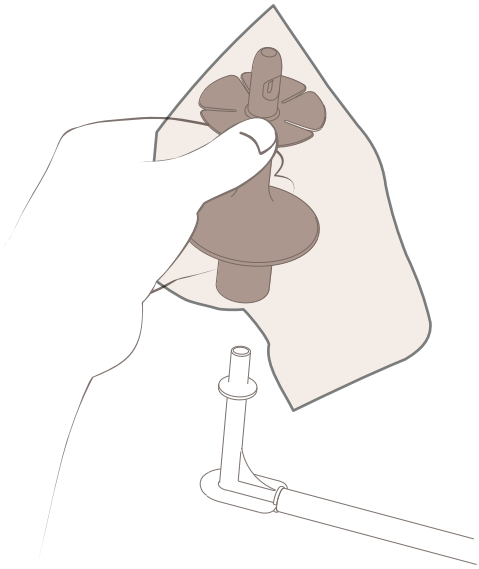
EN: Aquaflush® Transanal Irrigation systems must only be used by the person that it is prescribed for following a thorough assessment by a Healthcare Professional.

Caution: The following list does not rule out transanal irrigation as a bowel management programme, but you should have consulted a Healthcare Professional to ensure that the system is right for you. Children under 5 must be reviewed by a Paediatric Specialist before starting transanal irrigation.

- Inflammatory bowel disease (e.g. Crohn's disease or Ulcerative Colitis and Ischaemic Colitis)
- Rectal pain or bleeding e.g. anal fissure, severe haemorrhoids
- Diverticular Disease
- Autonomic Dysreflexia
- Recent spinal cord shock following spinal trauma
- Radiation therapy in the abdominal or pelvic region
- Abdominal or anal surgery
- Recent colonic biopsy or endoscopic polypectomy
- Cancer in the abdominal or pelvic region
- Long term steroid therapy, Anticoagulant Therapy
- Coronary Disease
- Pregnancy
- Nursing mothers

Intended use: The intended use of the product is for transanal irrigation of the bowel in adults and children 3 years and above with a bowel dysfunction, who have been assessed and trained by an appropriate healthcare professional.

WARNINGS: You must consult a Healthcare Professional if any of the following occur:

- Blood in the faeces, weight loss or abdominal pain
- Significant changes in the frequency, colour and consistency of the stool

Rare occurrence: Bowel perforation is extremely rare but is serious and life threatening if it occurs and will require urgent attention. If you experience any of the following symptoms during or after transanal irrigation, contact your Healthcare Professional immediately:

- Severe and sustained abdominal pain or back pain, especially if combined with fever
- Sustained bleeding

SAFETY: Ensure you store your Aquaflush® systems away from the reach of small children to avoid risk of strangulation from tubing, suffocation from plastic bags or choking in the event of small components coming loose from the product.

Do not use after the date of expiry, which can be found on the carton label.

Any serious incident that has occurred with this device should be reported to the manufacturer or competent authority.

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Clinisupplies Ltd, 1 Blackmoor Lane, Croxley Park, Watford, Hertfordshire WD18 8GA, UK
Clinisupplies (Nanjing) Co., Ltd.,
EC REP **Advena Ltd.** Tower Business Centre, 2nd Flr., Tower Street, Swatar, BKR 4013 Malta
VitalCare Limited, 6-9 Trinity Street, Dublin 2, Ireland