EN: Importer
NO: Importør
DA: Importør
SV: Importör
FI: Maahantuoja
NL: Importeur
ES: Importador
DE: Importeur
FR: Importateur
PT: Importador
IT: Importatore

EN: Distributor
NO: Distributør
DA: Distributør
SV: Distributör
FI: Jakelija
NL: Distributeur
ES: Distribuidor
DE: Lieferant
FR: Distributeur
PT: Distribuidor
IT: Distributore

EN: Manufacturer
NO: Produsent
DA: Producent
SV: Tillverkare
FI: Valmistaja
NL: Fabrikant
ES: Fabricante
PT: Fabricante
FR: Fabricant
DE: Hersteller
IT: Produttore

NON STERILE

EN: Non-sterile
NO: Usteril
DA: Ikke-steril
SV: Icke-steril
FI: Ei-steriili
NL: Niet steriel
ES: No estéril
DE: Unsteril
FR: Non-stérilisé
PT: Não estéril
IT: Non sterile

EN: Keep dry
NO: Må holdes tørt
DA: Holdes tør
SV: Förvaras torrt
FI: Pidä kuivana
NL: Droog bewaren
ES: Mantener seco
DE: Trocken aufbewahren
FR: À conserver au sec
PT: Manter seco
IT: Tenere asciutto

EN: Use-by-date
NO: Holdbarhetsdato
DA: Sidste anvendelsesdato
SV: Utgångsdatum
FI: Viimeinen käyttöpäivä
NL: Te gebruiken tot
ES: Fecha de caducidad
PT: Data de validade
FR: Date limite d'utilisation
DE: Verfallsdatum
IT: Da utilizzarsi entro

MD

EN: Medical device
NO: Medisinsk enhet
DA: Medicinsk enhed
SV: Medicinteknisk produkt
FI: Lääketieteellinen laite
NL: Medisch hulpmiddel
ES: Dispositivo médico
PT: Dispositivo médico
FR: Dispositif médical
DE: Medizinprodukt
IT: Dispositivo medico

LOT

EN: Batch code
NO: Batch-kode
DA: Batchkode
SV: Satskod
FI: Eräkoodi
NL: Batchcode
ES: Código de lote
PT: Código do lote
FR: Numéro de lot
DE: Chargencode
IT: Lotto di produzione

EN: Caution
NO: Forsiktig
DA: Forsigtig
SV: Försiktighet
FI: Varoitus
NL: Let op
ES: Advertencia
DE: Achtung
FR: Attention
PT: Cuidado
IT: Attenzione

EU REP

EN: Authorized representative in the European community/European union
NO: Autorisert representant i EU
DA: Bemyndiget repræsentant i Det Europæiske Fællesskab/Den Europæiske Union
SV: Auktoriserad representant i Europeiska gemenskapen/Europeiska unionen
FI: Valtuutettu edustaja Euroopan yhteisössä/Euroopan unionissa
NL: Bevoegd vertegenwoordiger in de Europese Gemeenschap/Europese Unie
ES: Representante autorizado en la Comunidad Europea/Unión Europea
DE: Bevollmächtigter in der Europäischen Gemeinschaft/Europäischen Union
FR: Mandataire agréé dans la communauté européenne/l'union européenne
PT: Representante autorizado na comunidade Europeia/União Europeia
IT: Rappresentante autorizzato nella comunità Europea/Unione Europea

EN: Do not re-use
NL: Niet hergebruiken
FR: Ne pas réutiliser
DE: Nicht wiederverwenden
DA: Må ikke genbruges
FI: Älä käytä uudelleen
IT: Non riutilizzare
NO: Må ikke brukes om igjen
PT: Não reutilizar
ES: No reutilizar
SV: Får inte återanvändas

EN: Keep away from sunlight
NO: Må ikke utsettes for sollys
DA: Holdes væk fra sollys
SV: Förvara skyddad från direkt solljus
FI: Pidä poissa suorasta auringonvalosta
NL: Beschermen tegen zonlicht
ES: Mantener fuera de la luz del sol
PT: Manter afastado da luz solar
FR: À conserver à l'abri de la lumière directe du soleil
DE: Von Sonnenlicht fernhalten
IT: Tenere lontano dalla luce del sole

Clinisupplies

Clinisupplies Ltd, 1 Blackmoor Lane, Croxley Park, Watford, Hertfordshire, WD18 8GA, UK
Tel: 020 8863 4168 Email: info@clinisupplies.co.uk **www.clinisupplies.co.uk**


Clinisupplies (Nanjing) Co., Ltd.
No. 86-8 Shuanggao Rd, Gaochun District, Nanjing, 211316 Jiangsu, P.R. CHINA

EU REP **Advena Ltd.** Tower Business Centre, 2nd Flr., Tower Street, Swatar, BKR 4013 Malta

VitalCare Limited, 6-9 Trinity Street, Dublin 2, Ireland

CE


CSF00916_V3EU 11/25

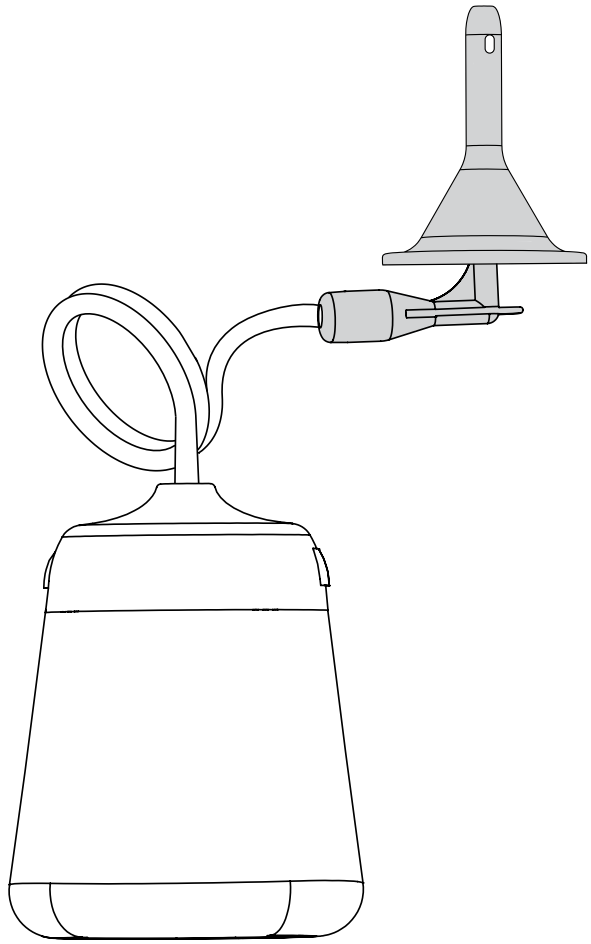AQUAFLUSH® 200X

TRANSANAL IRRIGATION SYSTEM

INSTRUCTIONS FOR USE Aquaflush®

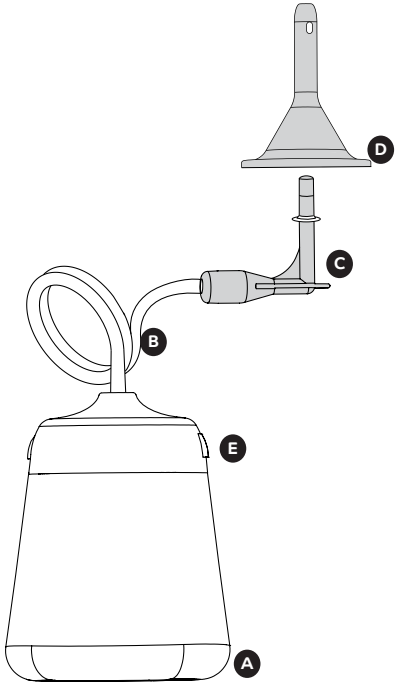Aquaflush®

KIT CONTENTS

A: WATER BOTTLE

B: EXTENSION TUBE

C: ELBOW CONNECTOR

D: STANDARD CONE

E: BOTTLE CAP

DE: PACKUNGSINHALT **NL:** NHOUD VAN HET PAKKET **FR:** CONTENU DU KIT **DA:** INDHOLD AF SÆTTET **FI:** SARJAN OSAT
IT: CONTENUTO DEL KIT **NO:** INNHOLD I SETTET **PT:** CONTEÚDO DO KIT **ES:** CONTENIDO DEL KIT **DE:** KITTETS INNEHÅLL


A: **DE:** WASSERFLASCHE **NL:** WATERFLES **FR:** BOUTEILLE D'EAU **DA:** VANDFLASKE **FI:** VESIPULLO **IT:** CONTENITORE DELL'ACQUA
NO: VANNFLASKE **PT:** GARRAFA DE ÁGUA **ES:** BOTELLA DE AGUA **SV:** VATTENFLASKA

B: **DE:** VERLÄNGERUNGSSCHLAUCH **NL:** SLANG **FR:** TUBE D'EXTENSION **DA:** FORLÆNGERSLANGE **FI:** JATKOLETKU
IT: TUBO DI PROLUNGA **NO:** FORLENGELSESSLANGE **PT:** TUBO DE EXTENSÃO **ES:** TUBO DE EXTENSIÓN **SV:** FÖRLÄNGNINGSSLANG

C: **DE:** WINKELANSCHLUSS **NL:** ELLEBOOGKOPPELING **FR:** CONNECTEUR COUDÉ **DA:** ALBUEFORBINDELSE **FI:** KULMALIITIN
IT: CONNETTORE A GOMITO **NO:** ALBUEKOBLING **PT:** CONECTOR EM COTOVELO **ES:** CONECTOR ACODADO **SV:** KRÖKT ANSLUTNING

D: **DE:** STANDARD-KEGEL **NL:** STANDAARD KEGEL **FR:** CÔNE STANDARD **DA:** STANDARD KAPPE **FI:** VAKIOKOKOINEN KARTIO
IT: CONO STANDARD **NO:** STANDARDKONUS **PT:** CONE PADRÃO **ES:** CONO ESTÁNDAR **SV:** STANDARDKON

E: **DE:** FLASCHENVERSCHLUSS **NL:** DEKSEL VOOR FLES **FR:** BOUCHON DE BOUTEILLE **DA:** FLASKELÅG **FI:** PULLON MUOVIKORKKI
IT: TAPPO DEL CONTENITORE **NO:** FLASKEKORK **PT:** TAMPA DA GARRAFA **ES:** TAPA DE LA BOTELLA **SV:** FLASKLOCK

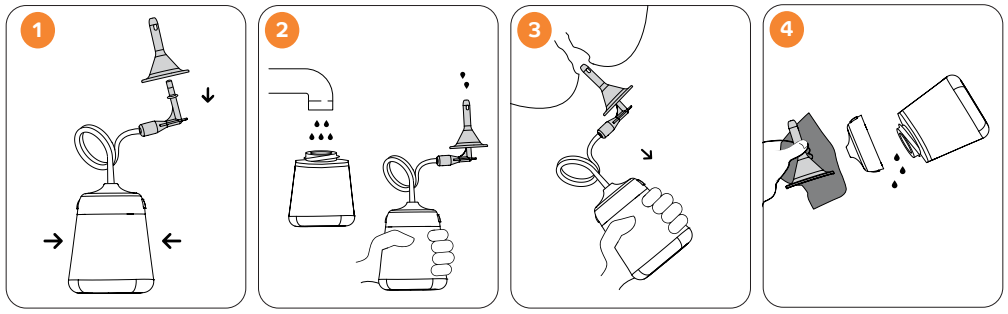EN

DESCRIPTION:

The **Aquaflush®** 200X Transanal Irrigation System (TAI), featuring an integrated flexible extension tube and a convenient elbow connector, is designed as a simple and effective bowel management solution. It's ideal for individuals who require a low water volume (up to 200ml) to seamlessly fit into their daily routine.

INTENDED USE:

The intended use of the product is for transanal irrigation of the bowel in adults and children 3 years and above with a bowel dysfunction, who have been assessed and trained by an appropriate healthcare professional.

INDICATIONS FOR USE:

Transanal irrigation (TAI) is intended for patients with neurological bowel dysfunction (NBD), gastrointestinal disorders of bowel innervation (DGBI), and post-surgical bowel disorders.

TAI helps manage faecal incontinence and/or mild constipation in conditions such as: post-defaecation seepage or soiling, passive faecal incontinence, incomplete evacuation, rectal washout.

The **Aquaflush®** 200X system can be used daily or as advised by your healthcare professional following appropriate assessment and training.

Some individuals may require assistance from a carer to perform the procedure safely and effectively.

USAGE GUIDELINES:

The water bottle, extension tube and elbow connector can be used 30 times (1 month based on daily usage) or up to 60 days, whichever is sooner. Cones and lubricant are single use only.

1. BEFORE YOU START (ASSEMBLY, PREPARATION):

Wash your hands and ensure you have everything you need close by. **Aquaflush®** 200X comes preassembled with the extension tube and elbow connector on the lid. Attach the cone to the elbow connector.

2. PREPARING FOR USE (FILLING & PRIMING):

- Twist to remove the lid from the water bottle. Fill the bottle with water between 35°C to 37°C (95°F to 98°F), which is close to body temperature (Do not use cold or hot water). Replace the lid, ensuring the lid is securely fixed.
- While holding over the sink, gently squeeze the bottle until water comes out of the cone. Your **Aquaflush®** system is now ready for use.
- A one-way valve in the elbow connector avoids water being drawn back in during use.

3. USING AQUAFLUSH®200X:

- Apply the water-based lubricant provided to the tip of the cone and to the anus if required.
- Whilst sitting on the toilet with your feet flat on the floor insert the cone into the anus. Position yourself towards the back of the toilet seat to insert the cone from in between your legs.
- Gently squeeze the bottle until it has been compressed fully to instil the water.
- When you have instilled the water, remove the cone, and remain seated on the toilet.
- You may experience a quick rush of water immediately which is normal. It may take up to 10 minutes to fully evacuate the water after using the **Aquaflush®** 200X system.

4. DISPOSAL OF CONE, CLEANING, DRYING:

Cones are single use and should only be used for one irrigation. Using the disposal bag provided, remove the cone from the elbow connector and dispose according to your local authority guidelines. Open the lid and shake out the excess water and leave to air dry. If you wish to wash the system, use a very mild soap and rinse thoroughly.

CONTRAINDICATIONS:

Do not use if you have certain conditions that may make TAI unsafe, including:

- Known anal stenosis
- Colorectal cancer
- Active inflammatory bowel disease (IBD)
- Ischemic colitis
- Acute diverticulitis and diverticular abscess
- Patients within 3 months of anastomotic colerectal surgery

Always consult your healthcare professional (HCP) before starting TAI.

PRECAUTIONS:

Always consult your HCP before starting TAI.

Regular monitoring during treatment is recommended to make sure the therapy remains safe and effective.

- Painful anorectal conditions—any condition which may cause pain and bleeding, e.g., anal fissure, anal fistula, grade three or grade four haemorrhoids
- Known colorectal stenosis
- Severe diverticulosis
- Patients at risk of uncontrolled autonomic dysreflexia
- Within 3 months of irradiation therapy in the abdominal or pelvic region and cessation of proctitis
- Cognitive impairment
- Change of stool pattern (e.g., sudden diarrhoea)
- Within 4 weeks of previous endoscopic polypectomy
- Previous colorectal surgery
- Patients within 3 months of anal surgery
- Patients within 3 months of radical prostatectomy
- Following pregnancy in an existing user

WARNINGS: Aquaflush® Transanal Irrigation Systems must only be used by the person that it is prescribed for following a thorough assessment by a Healthcare Professional.

You must consult a Healthcare Professional if any of the following occur:

- Blood in the faeces, weight loss or abdominal pain
- Significant changes in the frequency, colour and consistency of the stool

CAUTION: The following list does not rule out transanal irrigation as a bowel management programme, but you should have consulted a Healthcare Professional to ensure that the system is right for you. Children under 5 must be reviewed by a Paediatric Specialist before starting transanal irrigation.

- Inflammatory bowel disease (e.g. Crohn's disease or Ulcerative Colitis and Ischaemic Colitis)
- Rectal pain or bleeding e.g. anal fissure, severe haemorrhoids
- Recent spinal cord shock following spinal trauma

- Long term steroid therapy, anticoagulant therapy
- Coronary disease
- Nursing Mothers

RARE OCCURRENCE: Bowel perforation is extremely rare but is serious and life threatening if it occurs and will require urgent attention. If you experience any of the following symptoms during or after transanal irrigation, contact your Healthcare Professional immediately:

- Severe and sustained abdominal pain or back pain, especially if combined with fever
- Sustained bleeding

SAFETY: Ensure you store your **Aquaflush®** systems away from the reach of small children to avoid risk of strangulation from tubing or choking in the event of small components coming loose from the product.

Do not use after the date of expiry, which can be found on the carton label.
For water-based lubricant composition please check lubricant packaging.

Any serious incident that has occurred with this device should be reported to the manufacturer or competent authority. Details of your local Competent Authority can be found on your local government website.

STORAGE:

- Store the product in a cool, dry environment.
- Avoid direct sunlight and exposure to sources of heat
- Keep the product in a wash bag or its original packaging
- Keep tubing loosely coiled in its protective packaging — do not bend or fold sharply to prevent kinking
- Always inspect tubing for signs of flattening or kinking before use

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