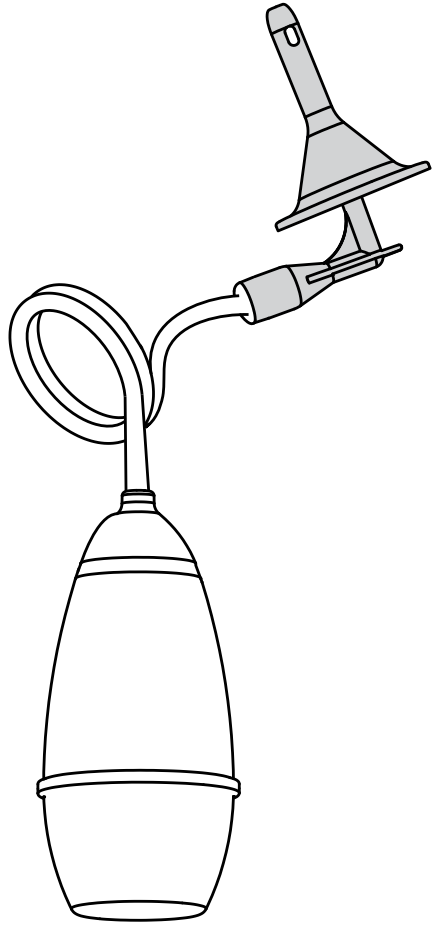
AQUAFLUSH® 500X

(STANDARD CONE SET)

TRANSANAL IRRIGATION SYSTEM

INSTRUCTIONS FOR USE Aquaflush®

Clinisupplies

Clinisupplies Ltd, 1 Blackmoor Lane, Croxley Park, Watford, Hertfordshire, WD18 8GA, UK
Tel: 020 8863 4168 Email: info@clinisupplies.co.uk **www.clinisupplies.co.uk**


Clinisupplies (Nanjing) Co., Ltd.
No. 86-8 Shuanggao Rd, Gaochun District, Nanjing, 211316 Jiangsu, P.R. CHINA

EU REP **Advena Ltd.** Tower Business Centre, 2nd Flr., Tower Street, Swatar, BKR 4013 Malta

VitalCare Limited, 6-9 Trinity Street, Dublin 2, Ireland

CE

CSF00918_V3EU 11/25

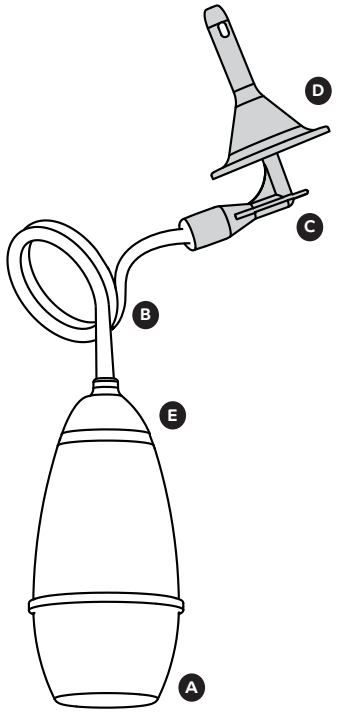Aquaflush®

KIT CONTENTS

A: WATER BOTTLE
B: EXTENSION TUBE
C: ELBOW CONNECTOR
D: STANDARD CONE
E: BOTTLE CAP

DE: PACKUNGSINHALT **NL:** NHOUD VAN HET PAKKET **FR:** CONTENU DU KIT **DA:** INDHOLD AF SÆTTET **FI:** SARJAN OSAT
IT: CONTENUTO DEL KIT **NO:** INNHOLD I SETTET **PT:** CONTEÚDO DO KIT **ES:** CONTENIDO DEL KIT **DE:** KITTETS INNEHÅLL


A: **DE:** WASSERFLASCHE **NL:** WATERFLES **FR:** BOUTEILLE D'EAU **DA:** VANDFLASKE **FI:** VESIPULLO **IT:** CONTENITORE DELL'ACQUA
NO: VANNFLASKE **PT:** GARRAFA DE ÁGUA **ES:** BOTELLA DE AGUA **SV:** VATTENFLASKA

B: **DE:** VERLÄNGERUNGSSCHLAUCH **NL:** SLANG **FR:** TUBE D'EXTENSION **DA:** FORLÆNGERSLANGE **FI:** JATKOLETKU
IT: TUBO DI PROLUNGA **NO:** FORLENGELSESSLANGE **PT:** TUBO DE EXTENSÃO **ES:** TUBO DE EXTENSIÓN **SV:** FÖRLÄNGNINGSSLANG

C: **DE:** WINKELANSCHLUSS **NL:** ELLEBOOGKOPPELING **FR:** CONNECTEUR COUDÉ **DA:** ALBUEFORBINDELSE **FI:** KULMALIITIN
IT: CONNETTORE A GOMITO **NO:** ALBUEKOBLING **PT:** CONECTOR EM COTOVELO **ES:** CONECTOR ACODADO **SV:** KRÖKT ANSLUTNING

D: **DE:** STANDARD-KEGEL **NL:** STANDAARD KEGEL **FR:** CÔNE STANDARD **DA:** STANDARD KAPPE **FI:** VAKIOKOKOINEN KARTIO
IT: CONO STANDARD **NO:** STANDARDKONUS **PT:** CONE PADRÃO **ES:** CONO ESTÁNDAR **SV:** STANDARDKON

E: **DE:** FLASCHENVERSCHLUSS **NL:** DEKSEL VOOR FLES **FR:** BOUCHON DE BOUTEILLE **DA:** FLASKELÅG **FI:** PULLON MUOVIKORKKI
IT: TAPPO DEL CONTENITORE **NO:** FLASKEKORK **PT:** TAMPA DA GARRAFA **ES:** TAPA DE LA BOTELLA **SV:** FLASKLOCK

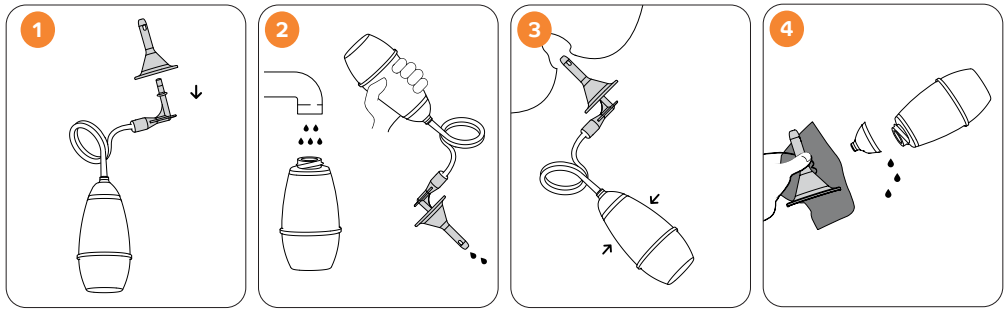DESCRIPTION:

The Aquaflush® 500X Transanal Irrigation System (TAI), featuring an integrated flexible extension tube and a convenient elbow connector, is designed as a simple and effective bowel management solution. It's ideal for individuals who require water volume of up to 500ml to achieve a regular bowel routine.

INTENDED USE:

The intended use of the product is for transanal irrigation of the bowel in adults and children 3 years and above with a bowel dysfunction, who have been assessed and trained by an appropriate healthcare professional.

INDICATIONS FOR USE:

Transanal irrigation (TAI) is intended for patients with neurological bowel dysfunction (NBD), gastrointestinal disorders of bowel innervation (DGBI), and post-surgical bowel disorders.

TAI helps manage faecal incontinence and/or mild constipation in conditions such as: post-defaecation seepage or soiling, passive faecal incontinence, incomplete evacuation, rectal washout.

The **Aquaflush®** 500X system can be used daily or as advised by your healthcare professional following appropriate assessment and training.

Some individuals may require assistance from a carer to perform the procedure safely and effectively.

USAGE GUIDELINES:

The water bottle, extension tube and elbow connector can be used 30 times (1 month based on daily usage) or up to 60 days, whichever is sooner. Cones and lubricant are single use only.

1. BEFORE YOU START (ASSEMBLY, PREPARATION):

Wash your hands and ensure you have everything you need close by. **Aquaflush®** 500X comes preassembled with the extension tube and elbow connector on the lid. Attach the cone to the elbow connector.

2. PREPARING FOR USE (FILLING & PRIMING):

- Twist to remove the lid from the water bottle. Fill the bottle with water between 35°C to 37°C (95°F to 98°F), which is close to body temperature (Do not use cold or hot water). Replace the lid, ensuring the lid is securely fixed.
- While holding over the sink, gently squeeze the bottle until water comes out of the cone. Your **Aquaflush®** system is now ready for use.
- A one-way valve in the elbow connector avoids water being drawn back in during use.

3. USING AQUAFLUSH®500X:

- Apply the water-based lubricant provided to the tip of the cone and to the anus if required.
- Whilst sitting on the toilet with your feet flat on the floor insert the cone into the anus. Position yourself towards the back of the toilet seat to insert the cone from in between your legs.
- With the bottle in the required position, gently squeeze and release the bottle to pump the water into the bowel. Keep pumping until the required amount of water has been dispensed into the bowel. This can be monitored using the graduations found on the side of the bottle.
- If you are using Aquaflush® 500X as a gravity fed system, keep the water bottle upright until you have inserted the cone and are ready to irrigate. Turn the water bottle upside down and raise it higher than the cone (above your waist) to allow the water to flow into your rectum. To stop the water flow, lower the water bottle and turn it upright.
- When you have instilled the water, remove the cone and remain seated on the toilet.
- You may experience a quick rush of water immediately which is normal. Time may vary to fully evacuate the water after using the Aquaflush® 500X system.

4. DISPOSAL OF CONE, CLEANING, DRYING:

Cones are single use and should only be used for one irrigation. Using the disposal bag provided, remove the cone from the elbow connector and dispose according to your local authority guidelines. Open the lid and shake out the excess water and leave to air dry. If you wish to wash the system, use a very mild soap and rinse thoroughly.

CONTRAINDICATIONS:

Do not use if you have certain conditions that may make TAI unsafe, including

- Known anal stenosis
- Colorectal cancer
- Active inflammatory bowel disease (IBD)
- Ischemic colitis
- Acute diverticulitis and diverticular abscess
- Known colorectal stenosis
- Chronic renal disease
- Congestive cardiac failure
- Patients within 3 months of anastomotic colorectal surgery
- Within 4 weeks of previous endoscopic polypectomy
- Following pregnancy in an existing user

Always consult your healthcare professional (HCP) before starting TAI.

PRECAUTIONS:

Always consult your HCP before starting TAI.

Regular monitoring during treatment is recommended to make sure the therapy remains safe and effective.

- Painful anorectal conditions - any condition which may cause pain and bleeding, e.g., anal fissure, anal fistula, grade three or grade four haemorrhoids
- Severe diverticulosis
- Patients at risk of haemorrhage from a systemic condition
- Patients at risk of uncontrolled autonomic dysreflexia
- Within 3 months of irradiation therapy in the abdominal or pelvic region and cessation of proctitis
- Cognitive impairment—appropriate support is required
- Change of stool pattern (e.g., sudden diarrhoea)
- Faecal impaction—need to disimpact prior TAI
- Previous major pelvic surgery
- Previous colorectal surgery
- Patients within 3 months of anal surgery
- Patients within 3 months of radical prostatectomy
- Anticoagulant therapy (not including aspirin or clopidogrel)

WARNINGS: Aquaflush® Transanal Irrigation Systems must only be used by the person that it is prescribed for following a thorough assessment by a Healthcare Professional.

You must consult a Healthcare Professional if any of the following occur:

- Blood in the faeces, weight loss or abdominal pain
- Significant changes in the frequency, colour and consistency of the stool

CAUTION: The following list does not rule out transanal irrigation as a bowel management programme, but you should have consulted a Healthcare Professional to ensure that the system is right for you. Children under 5 must be reviewed by a Paediatric Specialist before starting transanal irrigation.

- Inflammatory bowel disease (e.g. Crohn's disease or Ulcerative Colitis and Ischaemic Colitis)
- Rectal pain or bleeding e.g. anal fissure, severe haemorrhoids
- Recent spinal cord shock following spinal trauma
- Long term steroid therapy, anticoagulant therapy
- Coronary disease
- Nursing Mothers

RARE OCCURRENCE: Bowel perforation is extremely rare but is serious and life threatening if it occurs and will require urgent attention. If you experience any of the following symptoms during or after transanal irrigation, contact your Healthcare Professional immediately:

- Severe and sustained abdominal pain or back pain, especially if combined with fever
- Sustained bleeding

SAFETY: Ensure you store your **Aquaflush®** systems away from the reach of small children to avoid risk of strangulation from tubing or choking in the event of small components coming loose from the product.

Do not use after the date of expiry, which can be found on the carton label.
For water-based lubricant composition please check lubricant packaging.

Any serious incident that has occurred with this device should be reported to the manufacturer or competent authority. Details of your local Competent Authority can be found on your local government website.

STORAGE:

- Store the product in a cool, dry environment.
- Avoid direct sunlight and exposure to sources of heat
- Keep the product in a wash bag or its original packaging
- Keep tubing loosely coiled in its protective packaging — do not bend or fold sharply to prevent kinking
- Always inspect tubing for signs of flattening or kinking before use

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