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SensaCath® Compact
Intermittent Catheter – Female

EN INSTRUCTIONS FOR USE – PLEASE READ CAREFULLY

GENERAL INFORMATION

SensaCath® Compact Intermittent Catheter is a sterile, single use, latex-free catheter made of phthalate-free PVC (Polyvinyl Chloride) medical device.

SensaCath® has two smooth eyelets for efficient drainage. The catheter is coated with a hydrophilic coating to reduce friction during its insertion and removal from the patient's urethra. The **SensaCath®** is stored in a catheter case with sterile water, keeping its surface continuously smooth and fully activated, ready for use at any time. Once the **SensaCath®** is removed from the catheter case, it can be immediately used.

SensaCath® comes with one catheter handle per box. The catheter handle can be used to grip the catheter funnel and act as a handle.

INTENDED USE

SensaCath® are self-lubricating catheters intended to be inserted one-time through the urethra to the urinary bladder, by the patient/user or a healthcare provider, to provide urine drainage, typically for an individual who is physiologically incapable of voiding. It is intended that each user waits until all the urine is drained completely from the bladder. Catheter should be removed once urine has stopped flowing.

INDICATIONS

SensaCath® is used for bladder drainage with individuals who have bladder dysfunction. **SensaCath®** aids in the effective management of urinary retention and incontinence associated with neurological disorders, post-operative care, and other medical conditions.

WARNINGS

Intermittent self-catheterisation is a common procedure, however, you should always seek guidance from your healthcare professional and carefully read the instructions before use. Intermittent self-catheterisation should only be carried out under medical advice e.g. physician/nurse – and in accordance with the instructions given.

Do not use the catheter if the package is damaged, the seal sticker is broken, or it has passed its expiry date.

Reuse, reprocessing, cleaning, disinfection or sterilisation may compromise product characteristics which in turn may create a potential risk to the user.

In case of symptoms such as fever, discomfort/pain in the urinary tract, the desire to frequently empty bladder or having blood in the urine then consult your professional healthcare provider for advice and guidance. In the event of a serious side effect like bleeding from the urethra, urethral injury or urinary tract infection then please contact your doctor and or specialist nurse for treatment advice.

COMPLICATIONS

- Urinary tract infection
- Bleeding linked to urethral trauma or infection
- Epididymitis, urethral stricture and bladder stones – rare

POSSIBLE SIDE EFFECTS

- Urinary tract infection/Sepsis
- Burning and/or pain in urethra
- Urethral bleeding
- Urinary system stone disease

CONTRAINDICATIONS

- Acute urethritis
- Patients with PVC allergy
- Patients have calcareous urolithiasis (bladder stones)


PRECAUTIONS

- If you feel any resistance during insertion, do not force the catheter. Pull the catheter back slightly, try and relax and then continue to introduce the catheter
- If you cannot advance the catheter into the bladder seek medical advice
- If there is no urine flowing out, please check if the catheter is blocked or bent

INSTRUCTIONS FOR USE

SensaCath® is ready to use as soon as you open the packaging, so ensure you are prepared to proceed with catheterisation immediately after opening it.


1. Begin by washing your hands and the area around the urethral orifice. It is important to avoid transferring bacteria onto the catheter from your hands or surrounding areas.
2.  To access the catheter, twist the lid to break the seal between the lid and the catheter case. Hold the product vertically. You can use the lid as a base to stand the catheter case on a surface while you prepare for catheterisation.
3. Once the lid is removed, you can take out the catheter and use it immediately. The catheter can be used with or without the handle; however, using the handle offers increased reach and a more hygienic catheterisation experience. To remove the catheter from its case, gently grip the funnel using either the catheter handle or your fingers. Please take care not to drop or contaminate the catheter.
4. With one hand, part your labia to expose the urethral orifice.
5. Without touching the clear catheter part, use the Catheter Handle or grip the catheter funnel to gently insert the catheter into the urethra approximately 2cm at a time. Just before the catheter reaches the bladder you may feel a slight resistance due to the sphincter.

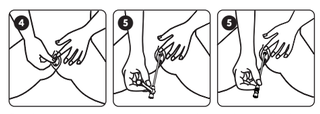6. Continue to insert the catheter into the urethra until urine starts to flow.
7. When the urine stops flowing, ease out the catheter 2 to 3cm. If the urine restarts, wait a few seconds before easing out the catheter another 1 to 2cm. Once the bladder is completely empty, remove the catheter slowly.
8. Place the catheter back in the container. Dispose of it or carry it in your bag until you can dispose of it in accordance with local authority guidelines.

 PLEASE NOTE! The black catheter case and Catheter Handle are recyclable.

SensaCath® can also be used with the re-useable handle included. After use the handle should be washed in mild soapy water and can be re-used up to 30 times.

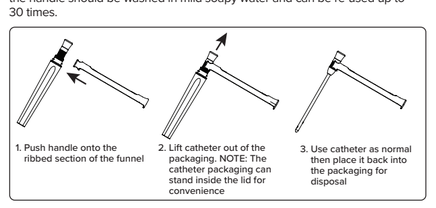Any serious incident that has occurred with this device should be reported to the manufacturer or competent authority. Details of your local Competent Authority can be found on your local government website.

SensaCath® Compact Intermittent Catheter

Female Box of 30

Reference	Description	Size	Length
PSCC10	Female length	10Ch	9cm
PSCC12	Female length	12Ch	9cm
PSCC14	Female length	14Ch	9cm

UK CA 0086

CE 2797

Clinisupplies

Clinisupplies Ltd, 1 Blackmoor Lane, Croxley Park, Watford, Hertfordshire, WD18 8GA, UK
Tel: +44 (0) 20 8863 4168
Email: info@clinisupplies.co.uk
www.clinisupplies.co.uk

EC REP Advena Ltd. Tower Business Centre, 2nd floor, Tower Street, Swatar, BKR 4013. Malta.

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