
SensaCath®

Intermittent Catheter

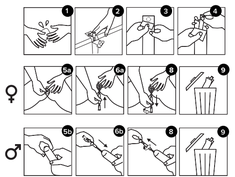The device is intended for people requiring intermittent bladder drainage due to conditions such as:

1. Urinary retention (e.g., bladder outlet obstruction).
2. Neurogenic bladder dysfunction leading to incomplete bladder emptying.
3. Pre/post-operative bladder management.

The catheter may also be indicated for people with anatomical challenges (e.g., urethral strictures, enlarged prostate) where the Tiemann tip/Flex variant provides easier passage.

ESSENTIAL PERFORMANCE

An intermittent catheter is a sterile, flexible tube intended to be inserted via the urethra into the bladder to enable urine drainage. It is used by the patient/user or health care professional, typically for individuals who are physiologically unable to completely empty their bladder.

INDICATIONS FOR USE

Intermittent catheters are used in people with various conditions including but not limited to:

- Bladder outlet obstruction (urinary retention).
- Patients with neurogenic conditions affecting bladder function i.e. Spinal-cord injury, multiple sclerosis, spina bifida, cauda equina syndrome.

WARNINGS

- Intermittent catheterisation should only be performed under medical advice e.g. physician/nurse – and in accordance with the instructions given.
- User should be trained in the proper use of the product.
- In the event of any signs of discomfort, bleeding, or infection, discontinue use and consult a healthcare professional.
- Do not use if the packaging has been deliberately or accidentally unsealed before intended use, product damaged or expired.
- Do not re-sterilise or re-use.
- Keep out of reach of children to prevent choking on small components.

CAUTIONS

- If you feel any resistance during insertion, do not force the catheter. Withdraw it slightly, try to relax and then continue insertion gently.
- If you cannot advance the catheter into the bladder, seek medical advice.
- If no urine flows, please check whether the catheter is blocked or kinked.
- Handle the catheter via non-touch slider; do not touch the catheter directly.
- Consult instructions for use.

PRECAUTIONS

- Catheterise according to the schedule recommended by your Healthcare Professional to maintain bladder health.
- Catheterisation should be performed as soon as possible after the packaging is opened and the catheter should be removed once the bladder is empty.

CONTRAINDICATIONS

- Do not use in patients with known allergy or hypersensitivity to any materials used in the device.
- Active urinary tract infection
- Undiagnosed haematuria
- Suspected urethral trauma
- Bleeding disorders
- Recent radiotherapy

SIDE EFFECTS

- Potential sensitivity to the materials used in the catheter.

COMPLICATIONS

Possible complications associated with intermittent catheterisation include:

- Urinary tract infection which in rare cases may progress to sepsis
- Urethral trauma
- Urethral strictures
- Bladder stones formation associated with long-term use

STORAGE CONDITIONS

- Keep Dry
- Keep away from sunlight

REPORTING OF SERIOUS INCIDENCE

If a device malfunction or failure results in a serious incident, such as significant deterioration in health or death, you should contact Clinisupplies Ltd and your local Competent Authority as soon as possible. Details of your local Competent Authority can be found on your local government website.

Clinisupplies Ltd accepts no liability for any loss or harm that may arise if this product is not used in accordance with the Instructions for Use and the company's recommendations.

DISPOSAL

Dispose of used product in accordance with local policies and guidelines.

Do not flush any product component down the toilet

CLINICAL BENEFITS

The clinical benefits of intermittent catheterisation include:

- Reduced risk of urinary tract infection through regular and effective bladder emptying.
- Improved bladder management and reduced risk of bladder overflow.
- Avoidance of complications associated with long term use of indwelling catheters.

STERILISATION

Sterilised by irradiation

INSTRUCTION FOR USE

1. Wash your hands and the area around the urethral opening.
2. Using the adhesive dot, attach the foil package to an appropriate surface.
3. Peel open the foil packaging from the top - enough to expose the catheter funnel and slider.
4. Take hold of the catheter funnel and slider and remove it from the packaging, taking care not to drop the catheter.
5. **5a. (Female)**: Using one hand part your labia to expose the urethral opening.

 5b. (Male): Straighten your urethra by lifting your penis up. Do not squeeze it, as they may make insertion more difficult.
6. **6a. (Female):** Use the funnel and/or slider to gently insert the catheter into the urethra, 2cm at a time.

 6b. (Male): Use the funnel and/or slider to gently insert the catheter into the urethra, 2cm at a time. As the tip reaches the bladder sphincter, you may feel a slight resistance.
7. Continue to insert the catheter until urine starts to flow.
8. When urine stops flowing, begin to remove the catheter - slowing rotating it as you go. If urine begins to flow, stop and allow it to drain. Repeat until the flow stops completely and the catheter is removed.
9. Dispose of the catheter according to local waste guidelines.

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Clinisupplies Ltd, 1 Blackmoor Lane, Croxley Park, Watford, Hertfordshire WD18 8GA, UK
EU REP **Advena Ltd.** Tower Business Centre, 2nd floor, Tower Street, Swatar, BKR 4013. Malta.