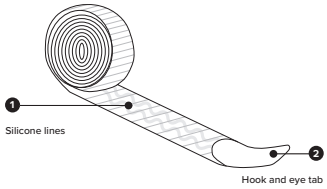
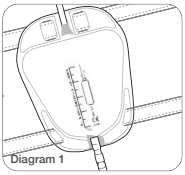
EN Libra® Leg Bag Strap

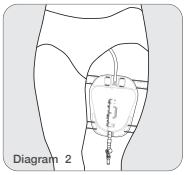
The **Libra®** Leg Bag Strap secures a leg bag to the leg. They feature a hook and eye style fastening for closure.

INSTRUCTIONS FOR USE

PLEASE READ CAREFULLY

Diagram 1

Diagram 2

Important: Do not let your leg bag hang unsupported from your Catheter as this will cause urethral discomfort and trauma.

Before Use

Ensure you are using the correct length strap as indicated in the table overleaf.

1. Feed the longer **Libra® Leg Bag Strap** (60cm or 70cm) through the 4 eyelets at the top of your leg bag **(Diagram 1)**
2. Feed the shorter **Libra® Leg Bag Strap** (40cm or 30cm) through the 2 eyelets at the bottom of your leg bag **(Diagram 1)**
3. Ensure the silicone lines **(Diagram 1)** will be facing your skin
4. Wrap the straps around your thigh or calf in the desired position **(Diagram 2)**
5. Secure using the hook and eye tabs **(Diagram 2)**. Ensure the strap is not overtightened around the leg as this can affect circulation

Cautions

- Use of any medical device carries associated risks. Device failures are rare, however in the event of a device failure of this product, you should contact Clinisupplies GB using the details provided within this leaflet
- If a failure has resulted in a serious incident, such as significant deterioration of health or death, you should contact Clinisupplies GB and your local Competent Authority as soon as possible. Details of your local Competent Authority can be found on your local government website

WASHING INSTRUCTIONS:

1. The **Libra® Leg Bag Strap** can be handwashed. It is not suitable for machine wash.
2. Do not Iron or Tumble dry

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Clinisupplies GB Ltd, 1 Blackmoor Lane, Croxley Park, Watford, Hertfordshire WD18 8GA,
EC REP **MD Reps Ltd**