Nightingale

NHS
Providing NHS services

Welcome to Nightingale Home Delivery Service.

For all your bladder and bowel products.



Clinisupplies | Nightingale

Nightingale Home Delivery Service is part of **Clinisupplies**, helping you and your healthcare professional to find and order the bladder or bowel products you may need.

Place an order online:
www.nightingaledelivery.co.uk


or phone:
0800 304 7434

CS23647_0224_v8

Clinisupplies


Liberty® Fix-it Retaining Strap

A quick look guide

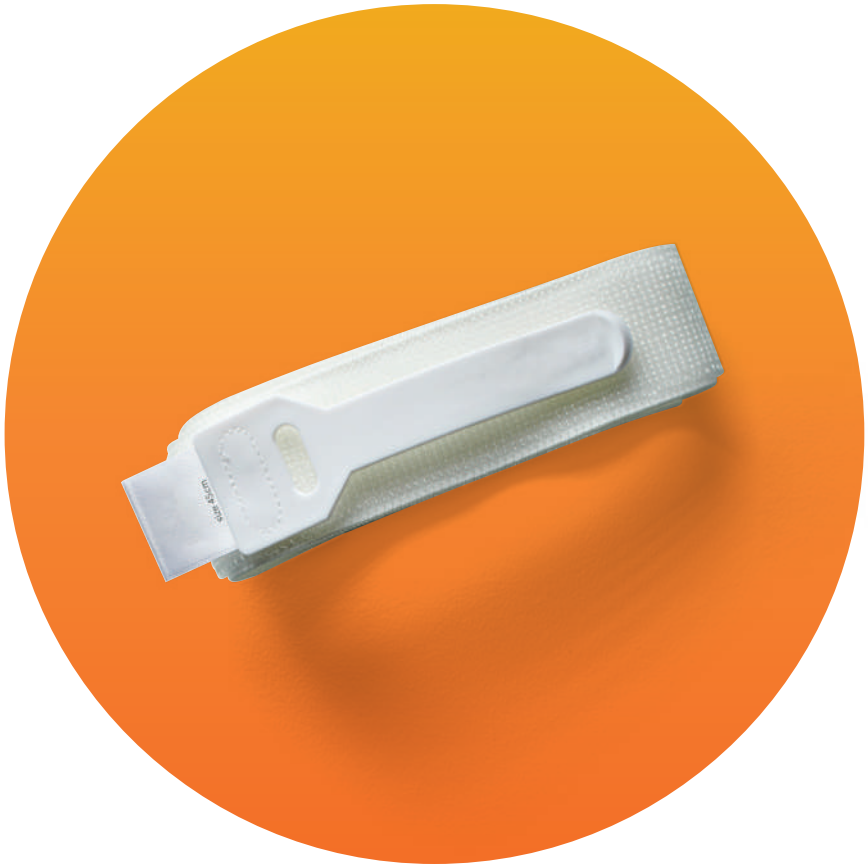

Instructions for use


The Liberty® Fix-it Retaining Strap is designed to secure indwelling and suprapubic catheters, as well as leg bag tubing with a sheath system.

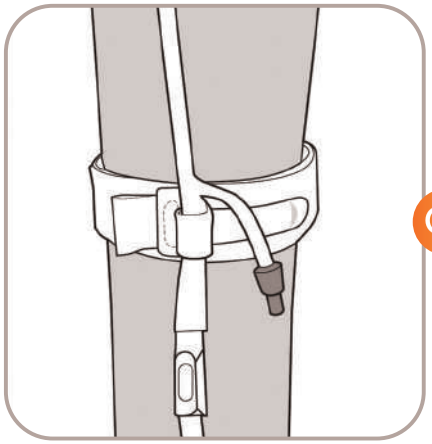


Or




1. Use the Liberty® Fix-it Retaining Strap to secure the catheter at the split point (bifurcation) or around the tubing of the leg bag

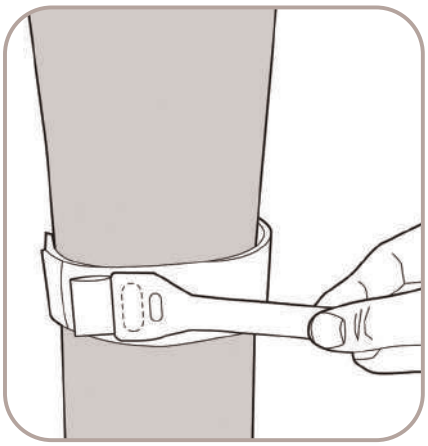

2. Pull Liberty® Fix-it tie away from strap

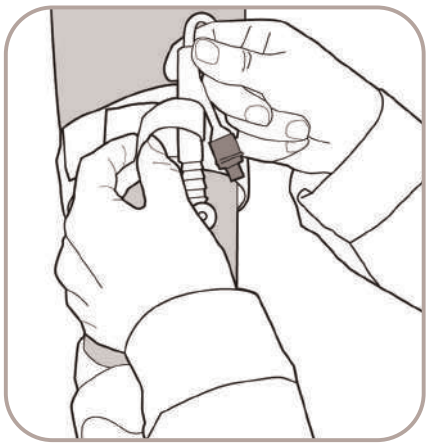

3. Wrap the Liberty® Fix-it tie back on itself around the tubing/ catheter

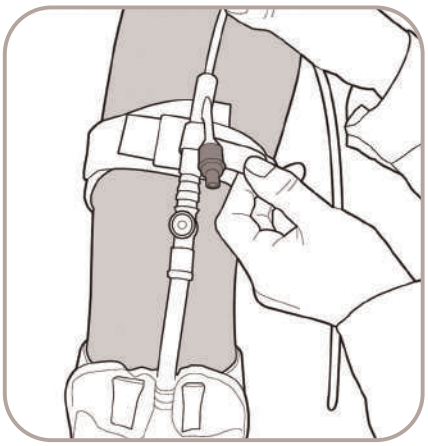

4. Push tab through eyelet and secure in place

Selecting the Appropriate Strap Length

Use a measuring tape to measure the circumference of the patient's thigh at the widest point, or around the abdomen at the area the strap is to be worn. Choose a Liberty® Fix-it strap that is longer than the measured circumference to ensure proper fit and secure placement.

PLEASE NOTE! It is important to secure the catheter after insertion to prevent movement and pulling on urethra[1]

PLEASE NOTE! It is important to secure a catheter to prevent urethral injury, pain, discomfort, and infection[2]

LATEX

Ref. 1. European Association of Urology Nurses (2012) Ref. 2. Fisher, J. (2010) "The importan ce of effective catheter securement" British Journal of Nursing, 19:8