Prosys®
Catheter Retaining Strap


Prosys® Catheter Retaining Strap is designed to keep the catheter or a leg bag tubing in a fixed position and to support the weight of a urinary leg bag.

Prosys® Catheter Retaining Strap can be worn on the lower abdomen or leg, and its position will depend on whether you have a urethral or suprapubic catheter and also on the length of your catheter and/or leg bag tubing.

INSTRUCTIONS FOR USE

Place the strap around either the leg or lower abdomen with the silicone beading facing the skin. Fasten the strap in the desired position on your leg or abdomen. Ensure the strap is snug, but not too tight and the **Prosys®** clasp is positioned where it is easily accessible.

1. Position the **Prosys®** clasp over the tubing.
   Wrap the **Prosys®** clasp around and under the tubing.
2. Push the tapered end of the **Prosys®** clasp through the eyelet provided.
3. Pull the **Prosys®** clasp to gently tighten around the tube.
   Pull the **Prosys®** clasp up and over to the opposite side, then press the tapered end into place, securing it on the main strap.

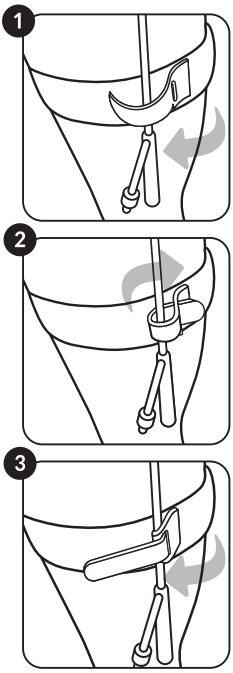

Ensure the **Prosys® Catheter Retaining Strap** fits tightly enough to keep the strap with the tubing secure during movement, but not too tight that it prevents normal blood circulation.

Laundering and care:

Close the catheter retaining strap by securing the **Prosys®** clasp in place and wash according to washing instructions. The strap may shrink up to 5% during washing.

NON STERILE **Non-sterile**	**Keep dry**	LOT **Batch code**	REF **Catalogue number**
Consult instructions for use	**Keep away from sunlight**	**Date of manufacture**	MD **Medical device**
Caution		**Use by date**	UDI **Unique device identifier**
LATEX **Not made with natural rubber latex**	**Do not use if package is damaged**	**Manufacturer**	UK CA **UK conformity mark**
EC REP **Authorized representative in the European Community/European Union**		CE **European conformity mark**	

Clinisupplies

Clinisupplies Ltd, 1 Blackmoor Lane, Croxley Park, Watford, Hertfordshire WD18 8GA, UK
Tel: +44 (0)20 8863 4168 Email: info@clinisupplies.co.uk **www.clinisupplies.co.uk**

EC REP Advena Ltd. Tower Business Centre, 2nd floor, Tower Street, Swatar, BKR 4013. Malta.

CS1090 05/23 v5