CLINIfast®

Tubular Bandage

Elasticated Viscose Stockinette

Clini*fast* for Dressing Retention and Wet and Dry Wrapping

Clini*fast* is a light weight, elasticated tubular bandage suitable for dressing retention and for wet and dry wrapping in the treatment of atopic eczema

Dressing Retention

Instructions for Use

Step 1 Select the correct size and cut to length

Step 2 Stretch Clini*fast* over the affected area

Step 3 Place your hand into the bandage and apply to the limb

Step 4 Clini*fast* may be applied in a double layer for extra security over joints

Wet and Dry Wrapping:

Instructions for Use

Step 1 Select Clini*fast* in the required size to cover the affected area. Measure to ensure full coverage of the area and cut to size

Step 2 Before wet or dry wrapping the person should be bathed, dried thoroughly and have had emollient applied to their skin in a downward motion to prevent the hair follicles becoming blocked

Step 3 Apply Clini*fast* to the affected area (knees, wrists, elbows etc)

Step 4 When wet wrapping soak in warm water, squeeze out excess and apply damp, follow with a dry top layer

Step 5 When dry wrapping apply a single layer only over the emollient which will keep the emollient in place longer

Clini*fast* Sizing Guide

Red Line:	Green Line:	Blue Line:	Yellow Line:	Beige Line:
For small limbs: Unstretched width **3.5cm**, limb circumference **8-15cm**	For small/medium limbs: Unstretched width **5.0cm**, limb circumference **10-25cm**	For large limbs: Unstretched width **7.5cm**, limb circumference **20-45cm**	For children's heads, trunks, extra large limbs and heads: Unstretched width **10.75cm**, limb circumference **35-65cm**	For adult trunks: Unstretched width **17.5cm**, limb circumference **50-120cm**

Washing Instructions



Symptoms may improve by using a non-biological washing powder and avoiding the use of fabric conditioner. Seek advice from your healthcare professional.

Storage

Store in dry conditions out of direct sunlight


Fabric Content

92% Viscose , 5% Elastane , 3% Polyamide This device is free from natural latex

Precautions

- Keep Clini***fast*** bandages away from naked flames
- Only wet or dry wrap if advised to do so and instructed by a healthcare professional
- When dry wrapping, always use an emollient
- Do not use if you think the child's skin is infected or the child is unwell and seek immediate medical advice
- Ensure the correct fit to allow complete freedom of movement

For further information on Clini***fast*** bandages for dressing retention and wet or dry wrapping, contact our Customer Services on **020 8863 4168** or visit our website at **www.clinisupplies.co.uk**.

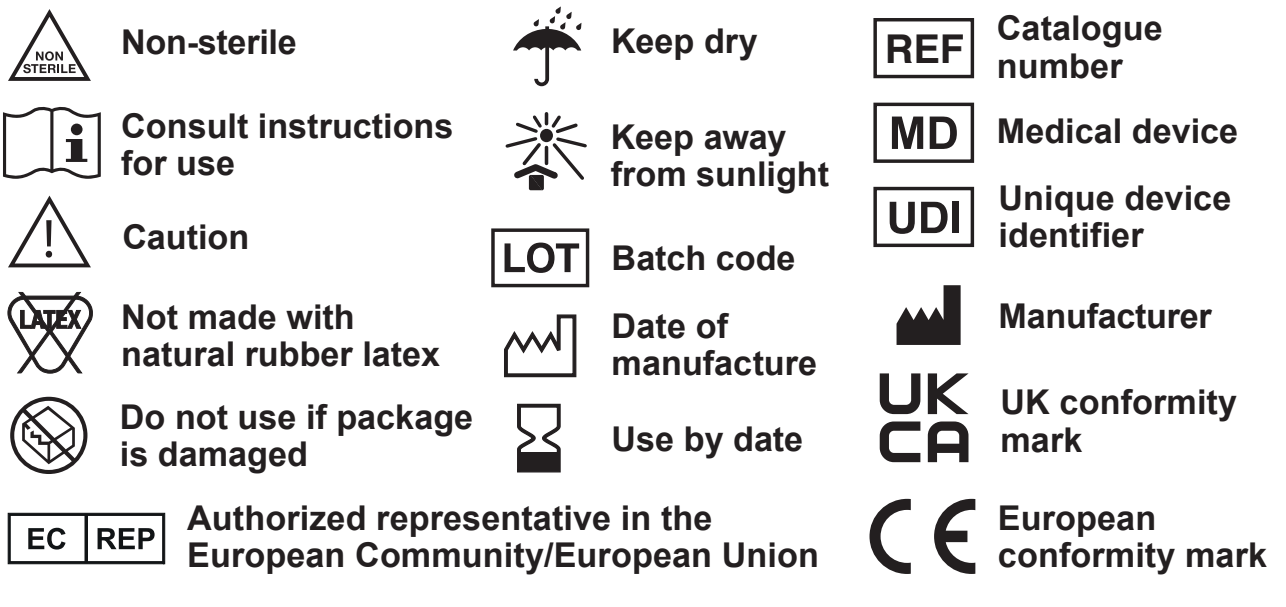

Clinisupplies Ltd,
1 Blackmoor Lane, Croxley Park, Watford, Hertfordshire WD18 8GA, UK
Telephone: 020 8863 4168 Fax: 020 8426 0768
Email: info@clinisupplies.co.uk **www.clinisupplies.co.uk**

EC REP Advena Ltd. Tower Business Centre, 2nd floor, Tower Street, Swatar, BKR 4013. Malta.

Clinisupplies

08/22 v5