SYMBOL EXPLANATION

EN: UK conformity mark
FR: Marque de conformité brittanique
DE: GB-Konformitätszeichen

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Clinisupplies

Clinisupplies GB Ltd, 1 Blackmoor Lane, Croxley Park, Watford, Hertfordshire WD18 8GA, UK
Tel: +44 (0)20 8863 4168 Email: info@clinisupplies.co.uk **www.clinisupplies.co.uk**

EC REP **MD Reps Ltd**
12-16 Bridge Street, Belfast, Northern Ireland, BT1 1LU

CSF00806_1_v3
07/25

UK CA 0120
CE 1639


EN Libra® Lever Catheter Valve


The **Libra®** Lever Catheter Valve is an individually packed, sterile device.


The **Libra®** Lever Catheter Valve is specifically designed to fit all Catheters, and can also be connected to a Night Bag for overnight drainage.

INSTRUCTIONS FOR USE

PLEASE READ CAREFULLY

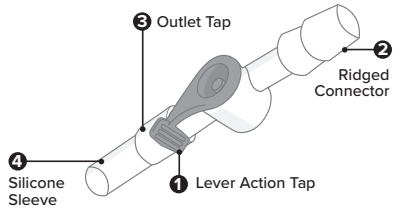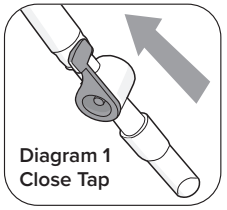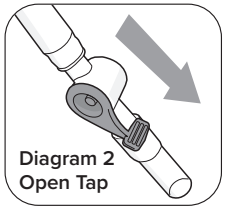Before Use

1. Wash Hands
2. Open peel pouch and remove **Libra® Lever Catheter Valve** from its packaging
3. Ensure the lever on the Lever Action Tap **(1)** is in the closed position **(Diagram 1)**

Direct Connection and Drainage

1. Insert the ridged inlet connector **(2)** into the end of the Catheter. Insert the connector fully to allow a secure connection
2. Drain urine from the bladder at a frequency instructed by your Healthcare Professional
3. To drain the bladder, push the lever action tap **(1)** downwards **(Diagram 2)**
4. After emptying, wrap the outlet tap **(3)** with with a tissue and pat a number of times before gently wiping around the inside of the valve. This will dislodge any residual urine
5. Remember to close the tap after emptying **(Diagram 1)**

Overnight Drainage

1. When you go to bed, do not disconnect your Catheter Valve from your Catheter
2. Insert the connector of your Night Bag to the silicone sleeve **(4)** on your Catheter Valve. Insert the connector fully to allow a secure connection
3. Once connected, turn the Catheter Valve tap to the open position **(Diagram 2)** to allow the overnight drainage into the Night Bag
4. Before disconnecting the Night Bag in the morning, ensure the Catheter Valve tap is in the closed position **(Diagram 1)**
5. Do not reconnect a used Night Bag

Disposal

Dispose with normal household waste – do not flush down the toilet

Cautions

- The **Libra® Lever Catheter Valve** should not be reused after disconnection as this may compromise the integrity of the materials used
- Keep out of the reach of children
- The **Libra® Lever Catheter Valve** should be changed every 5-7 days
- The **Libra® Lever Catheter Valve** can also be used in conjunction with a Sheath, however please seek advice from your Healthcare Professional before use
- Use of any medical device carries associated risks. Device failures are rare, however in the event of a device failure of this product, you should contact Clinisupplies GB using the details provided within this leaflet
- If a failure has resulted in a serious incident, such as significant deterioration of health or death, you should contact Clinisupplies GB and your local Competent Authority as soon as possible. Details of your local Competent Authority can be found on your local government website

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