EN Libra® GB2 Night Bag

Libra® GB2 Night Bag is a 2 litre bag and is for single-use only and cannot be re-used.

INSTRUCTIONS FOR USE

PLEASE READ CAREFULLY

Libra® GB2 is for single-use only and cannot be re-used.

The scale on the bag is an approximation and should not be used as an accurate measuring device.

To reduce the risk of acquiring a catheter associated urinary tract infection do not disconnect your Libra Leg Bag from your Indwelling Catheter. A Libra Leg Bag should only be disconnected when a new Libra Leg Bag is required – normally every 7 days.

A **Libra® GB2 Night Bag** can be connected directly to a sheath.

- Remove grey cap from the **Libra® GB2 Night Bag** connector.
- Insert the ridged connector into the soft silicone tube on the outlet tap of your Libra Leg Bag **(Diagram 1)**. Insert the connector fully to allow a secure connection between the **Libra® GB2 Night Bag** and Libra Leg Bag.
- The **Libra® GB2 Night Bag** should be supported on a night stand using the hanger holes located at the top of the **Libra® GB2 Night Bag**.
- Ensure the **Libra® GB2 Night Bag** is positioned on the night stand at a level below the bladder for optimal drainage.
- Make sure the tubing isn't kinked which can restrict the flow of urine.
- To allow urine drainage from your Libra Leg Bag to the **Libra® GB2 Night Bag** ensure the tap on the Libra Leg Bag is in the open position **(Diagram 2)**.

Emptying the GB2

- Close the Libra Leg Bag tap before disconnecting the **Libra® GB2 Night Bag (Diagram 3)**.
- To empty the **Libra® GB2 Night Bag** push your thumb firmly through the top left hand corner, and continue the tear downwards into the bag **(Diagram 4)**.
- Fully invert the bag, allowing urine to flow into the toilet.


Disposal


- Dispose with normal household waste – do not flush down the toilet.


Caution:


- Use of any medical device carries associated risks. Device failures are rare, however in the event of a device failure, you should contact Clinisupplies GB using the details provided below.
- If a failure has resulted in a serious incident, such as significant deterioration of health or death, you should contact Clinisupplies GB and your local Competent Authority as soon as possible. Details of your local Competent Authority can be found on your local government website.

Warning: Keep out of the reach of children

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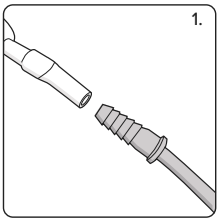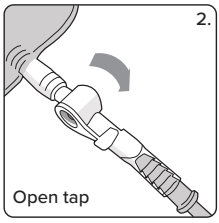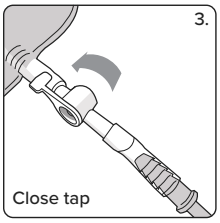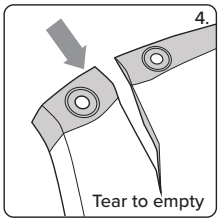A **Libra® GB2 Nachttas** can be connected directly to a sheath.

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Clinisupplies GB Ltd, 1 Blackmoor Lane, Croxley Park, Watford, Hertfordshire WD18 8GA, UK
EC REP **MD Reps Ltd**