



-80°C
-40°C
-20°C
RT


STC
Series


BiO.Li.X!

STC-SERIES - AUTOMATED BIO-LIBRARIES

LiCONiC builds on 25 years of experience in providing climate controlled sample storage automation to offer the STC Bio-Library Storage System. The STC-Series, one of the four Biobanking solutions in the BiOLiX! family of automated Biobanking storage products, offers the broadest range of application and sample capacity requirements. First installed in 2008, the STC is the longest continuously running -80C automated storage system solution on the market. System design features and LiCONiC's continual product development commitment, continue to evolve parallel with changing Biobanking methods and requirements.

Energy Efficiency

The STC-Series delivers an industry leading return on investment by all standards. Through the utilization of a compressible shelf storage system, which is maintained within a chest freezer, our proprietary refrigeration design reduces power consumption to a fraction of the nearest competitor.

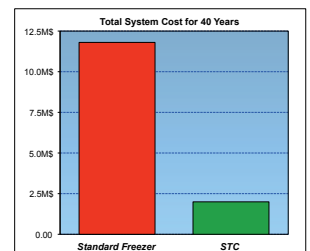

MAIN FEATURES

- Ultra low temperature
- Unique compact shelf design
- Fully automated sample picking & tracking
- Compatible with all SBS rack formatted tubes and vials
- Ultra low energy refrigeration technology
- Multiple levels of sample safety backup
- External transfer station for integration
- Fully featured system control software
- LIMS and robotic workcell integration

STC-SERIES PRODUCT OVERVIEW

The STC-Series storage systems are ideally suited for a wide range of storage applications ranging in size from approximately 100k to millions of samples, and for temperatures ranging from +25°C to -80°C. The STC-ULT series provides all the hallmark LiCONiC advantages such as best in class temperature stability, energy consumption, structural integrity, sample density and simple, effective, software control.

LiCONiC's exclusive compressible, high density shelf system is fully enclosed in an enhanced, insulated chest style compartment. This design exploits the natural property of sinking cold air, providing for an extremely energy efficient enclosure for all temperature ranges.

A multi-function Input/Output Interface with tube picking station is confined to a single system area to prevent moisture ingress and to provide easy service access. A fully redundant refrigeration system is the first of multiple sample safety levels.

Unique frost free chest freezer and compressible shelf designs provide:

- Maximum storage density
- Temperature uniformity and stability
- Extremely energy efficient refrigeration




The STC-Series of automated sample stores address the widest range of sample storage applications at temperatures from +25°C to -80°C.

- Small molecule compound storage
- Diagnostics
- Population based biobanks
- Disease based biobanks
- Research biobanks
- Reagent storage
- Core sample storage facilities


STC STORAGE CONDITIONS

LiCONiC's unique chest freezer design naturally creates numerous benefits. Two of the most important are: Two climate zones, one for active components at -20°C and another for passive elements at -80°C . Proven, stable, horizontal temperature layers within the ULT storage area.

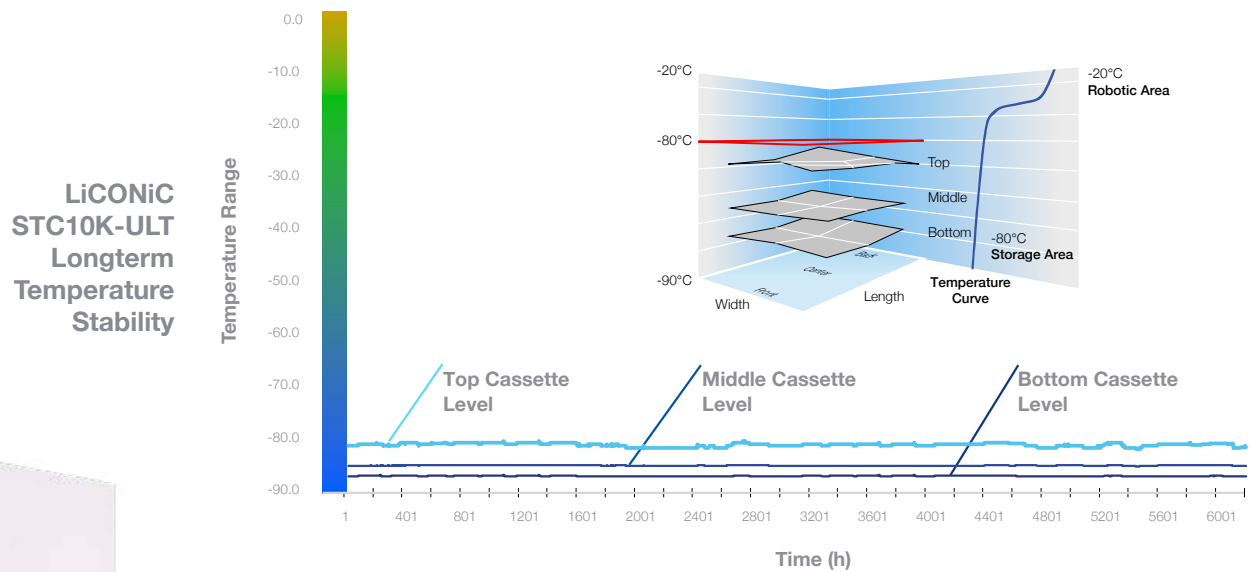

SPACE EFFICIENCY

The STC-Series Automated BioLibraries make best use of given room constraints by its “built to space” modular design. System dimensions and layout are optimised to the individual application’s needs and workflow requirements.

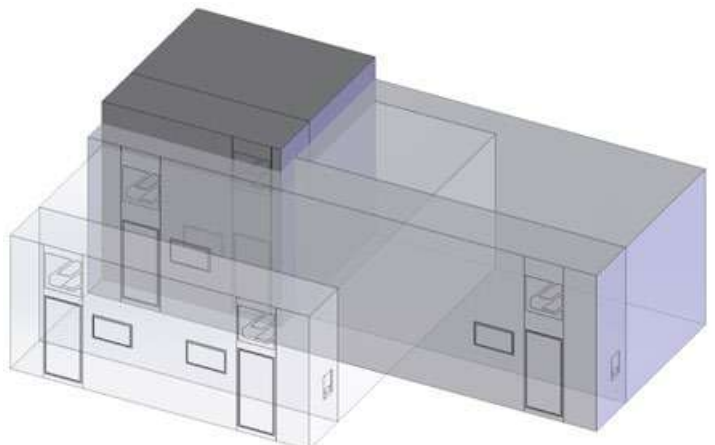

STC CAPACITIES

BiO.Li.X! STC-Series examples			STC6k0-ULT		STC32k5-ULT		STC101k0-ULT	
Labware	Plate height	Pitch	Levels	Samples	Levels	Samples	Levels	Samples
FluidX 0.3 ml Screw Cap	19.1 mm	25 mm	42	524'160	50	2'880'000	78	8'985'600
Micronic 0.5 ml Screw Cap	28.5 mm	35 mm	30	374'400	35	2'016'000	55	6'336'000
LVL LX1000	45.6 mm	52 mm	20	249'600	24	1'382'400	37	4'262'400
Thermo Matrix 1.0ml Screw Cap	54.6 mm	61 mm	17	212'160	20	1'152'000	32	3'686'400
Greiner Bio One 2 ml	48.0 mm	54 mm	19	118'560	23	662'400	36	2'073'600




INTUITIVE SOFTWARE

Graphical Job Creation for Import/Export of Samples
Easy Search of Samples by Barcode, Name, Data etc.
Complete Audit Trail
Automatic Inventory Defragmentation
Automatic Data Backup
Integration to LIMS Database
Remote Support
Mobile Device Control
User Access Security

OPTIONS & ACCESSORIES


Barcode Reader
Flexible Labware Support
Configurable to Room Constraints
Tablet Interface
External Database Integration
Remote Monitoring
Internal Vision System
Multi - Tube / Vial Picker
LN2 - Backup Option
2D - BCR Rack Scanner
Redundant Refrigeration
Dynamically Configurable Storage Cassettes
External Interface Station

STC DIMENSIONS

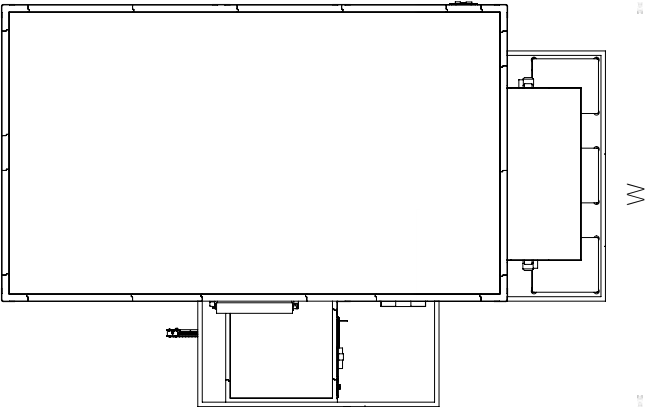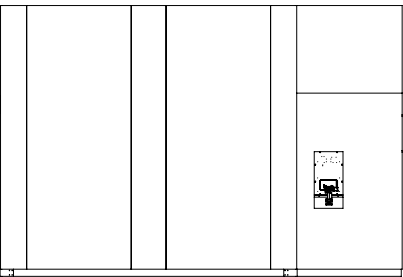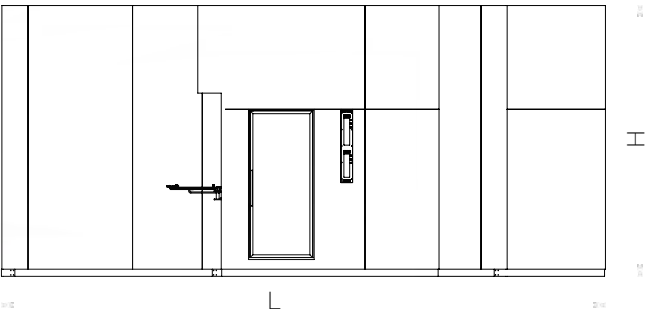

BiO.Li.X! STC Dimensions examples			
Type	L= Length	W= Width	H= Height
STC6k0-ULT	5580 mm	2340 mm	2830 mm
STC32k5-ULT	6755 mm	4495 mm	3030 mm
STC101k0-ULT	11115 mm	4495 mm	4030 mm

WORLDWIDE PRESENCE

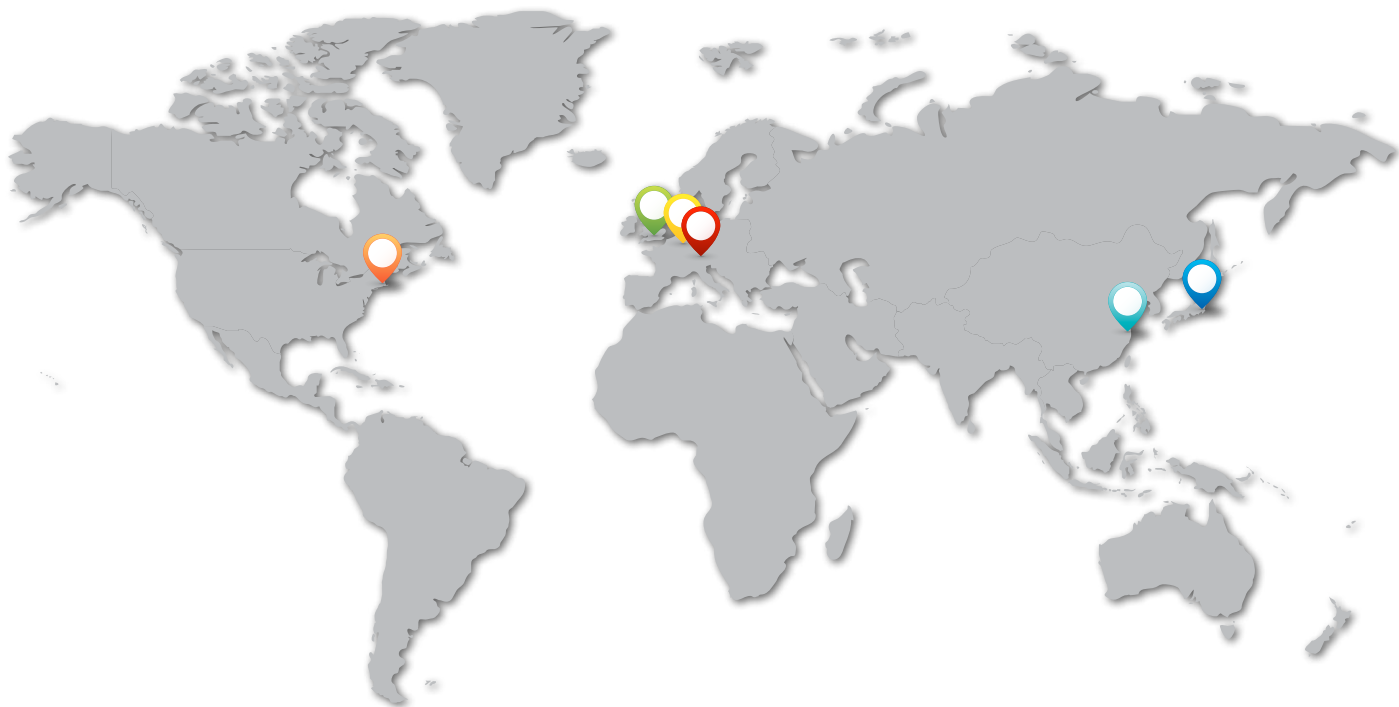

CORPORATE **LiCONiC AG**

LiCONiC AG
Industriestrasse 8-12
9493 Mauren
Principality of Liechtenstein

Phone: +423 340 5000
Fax: +423 37 35 35 9
E-Mail: info@liconic.com
sales@liconic.com
service@liconic.com

LiCONiC Services **Deutschland GmbH**

Bahnhofplatz 3
D-56410 Montabaur
Germany

Phone: +49 2602 997 4557
Fax: +49 2602 997 4559
E-Mail: info.de@liconic.com
sales.de@liconic.com
service.de@liconic.com

LiCONiC UK Ltd

BioHub
Alderley Park
Macclesfield
Cheshire, SK10 4TG
United Kingdom

Phone: +44 (0) 1618 752 542
+44 (0) 1618 752 543
Fax: +44 (0) 1618 752 540
E-Mail: info.uk@liconic.com
sales.uk@liconic.com
service.uk@liconic.com

LiCONiC US, Inc.

21-F Olympia Ave
Woburn, MA 01801
USA

Phone: +1 (781) 309-3340
Fax: +1 (781) 537-6173
E-Mail: info.usa@liconic.com
sales.usa@liconic.com
service.usa@liconic.com

LiCONiC Japan K.K.

2-2-15 Daiya Bldg 2F
Hamamatsu-cho, Minato-ku
Tokyo Japan 105-0013

Phone: +81 3 6841 0707
E-Mail: info.japan@liconic.com
sales.japan@liconic.com
service.japan@liconic.com

LiCONiC China Co. Ltd

1-B102, No.3891, Jinxiu Rd
Shanghai (201204),
R.P.China

Phone: +86 186-2156-0603
E-Mail: info.china@liconic.com
sales.china@liconic.com
service.china@liconic.com

Asia:

Phone: +1 (760) 522-2455
E-Mail: info.asia@liconic.com
sales.asia@liconic.com
service.asia@liconic.com



Information within this document is subject to be changed without prior notice. For details please contact our applications at info@liconic.com.
20190411



www.liconic.com